



CAPSTONE PROJECT REPORT

Report 3 – Software Requirement Specification

CareBridge: A Maternal and Early Childhood
Healthcare Platform

– Hanoi, June 2026 –

Table of Contents

I. Record of Changes	12
II. Software Requirement Specification	14
1. Requirement Overview	14
1.1 Context Diagram	14
1.1.1 Web App Context Diagram	14
1.1.2 Mobile App Context Diagram	17
1.2 Main Workflow	20
1.2.1 User Onboarding, Account & Privacy Setup Workflow	20
1.2.2 Mother, Baby & Health Record Workflow	22
1.2.3 Community Q&A, Expert Answer & Moderation Workflow	24
1.2.4 Expert Verification, Booking & Consultation Workflow	26
1.2.5 Gemini AI Triage & Safe Action Workflow	28
1.2.6 TrackAsia Emergency Map & Nearby Care Workflow	30
1.2.7 Family Sync, Permission & Care Task Workflow	31
1.2.8 Smart Activity Monitoring & Safety Support Workflow	33
1.2.9 Expert App, TrackAsia Nearby Support & Realtime Contact Workflow	36
1.2.10 Pregnancy Exercise & MediaPipe Posture Workflow	39
1.2.11 Smartwatch Data Sync & Health Trend Workflow	41
1.2.12 Partner Clinic & Sponsored Content Governance Workflow	42
1.2.13 Verified Content, Checklist & Gemini RAG Knowledge Workflow	44
1.2.14 Audit, Security Incident & Access Review Workflow	46
2. User Requirements	47
2.1 Actors	47
2.2 Use Cases	49
2.2.1 Diagram(s)	49
2.2.1.1 Use Cases for User	49
2.2.1.2 Use Cases for Mother	50
2.2.1.3 Use Cases for Family Member	51
2.2.1.4 Use Cases for Expert	52
2.2.1.5 Use Cases for Content, Moderation & Operations	53
2.2.1.6 Use Cases for Partner	54
2.2.1.7 Use Cases for Supporting Services	55
2.2.1.8 Use Cases for Smart Activity Monitoring & Safety Support	56
2.2.1.9 Use Cases for Community Q&A Extended Use Cases	57

2.2.1.10 Use Cases for Map, Expert App & Nearby Support	58
2.2.1.11 Use Cases for Pregnancy Exercise & MediaPipe Posture	59
2.2.1.12 Use Cases for Consultation Lifecycle, Dispute & Refund	60
2.2.1.13 Use Cases for Family Sync & Care Task Management	61
2.2.1.14 Use Cases for Verified Content, Checklist & RAG Knowledge	62
2.2.1.15 Use Cases for Vaccination & Growth Tracking	63
2.2.2 Descriptions	63
3. Functional Specifications	73
3.1 Shared	73
3.1.1 Authentication & Profile	73
3.1.1.1 UC-01 Register Account	73
3.1.1.2 UC-02 Verify OTP	74
3.1.1.3 UC-03 Log In	75
3.1.1.4 UC-04 Log Out	76
3.1.1.5 UC-05 Request Password Reset	77
3.1.1.6 UC-06 Reset Password	78
3.1.1.7 UC-07 Change Password	79
3.1.1.8 UC-08 View Account Profile	80
3.1.1.9 UC-09 Update Account Profile	80
3.1.1.10 UC-10 Manage Community Identity	81
3.1.1.11 UC-11 Manage Notification Preferences	82
3.1.1.12 UC-12 View and Mark Notifications	83
3.1.1.13 UC-13 Manage Own Login Sessions	84
3.1.1.14 UC-14 Deactivate or Delete Own Account	85
3.1.2 Consent & Data Permission	86
3.1.2.1 UC-15 Grant Data Permission	86
3.1.2.2 UC-16 Review and Revoke Data Permission	87
3.2 WebApp	88
3.2.1 Account & Security Administration	88
3.2.1.1 UC-17 Administer User Accounts and Role Access	88
3.2.1.2 UC-18 Review Sensitive Access and Security Events	89
3.2.2 Community Moderation & Governance	90
3.2.2.1 UC-56 View Moderation Queue	90
3.2.2.2 UC-57 Moderate Community Content	91
3.2.2.3 UC-58 Resolve Content or Account Report and Apply Enforcement	92
3.2.2.4 UC-59 Manage Community Topics and Visibility	93

3.2.3 Expert Network & Verification	94
3.2.3.1 UC-60 Submit Expert Profile	94
3.2.3.2 UC-61 Update Expert Profile	94
3.2.3.3 UC-62 Submit or Replace Verification Documents	95
3.2.3.4 UC-63 View Verification Status and Renew Submission	96
3.2.3.5 UC-64 Configure Expert Availability and Service Scope	97
3.2.3.6 UC-67 View Expert Question Queue	98
3.2.3.7 UC-68 Post Verified Expert Answer	99
3.2.3.8 UC-69 View Contribution Points and Badges	100
3.2.3.9 UC-70 Review Expert Verification Submission	101
3.2.3.10 UC-71 Restrict, Suspend or Reinstate Expert Trust Status	102
3.2.4 AI Knowledge & Risk Configuration	103
3.2.4.1 UC-75 Manage Approved AI Knowledge Sources	103
3.2.4.2 UC-76 Configure AI Risk and Red-Flag Rules	104
3.2.5 Nearby Expert Support	105
3.2.5.1 UC-82 Manage Expert Nearby Availability and Respond to Nearby Support Request	105
3.2.6 Verified Content Administration	106
3.2.6.1 UC-104 Create Verified Content	106
3.2.6.2 UC-105 Update Verified Content and Sources	107
3.2.6.3 UC-106 Review and Publish Content Version	108
3.2.6.4 UC-107 Unpublish or Archive Content	109
3.2.6.5 UC-108 Manage Content Categories and Stage/Topic Mapping	110
3.3 MobileApp	111
3.3.1 Mother Care Journey	111
3.3.1.1 UC-19 Initialize Mother Care Journey	111
3.3.1.2 UC-20 Update Mother Journey Stage and Dates	112
3.3.1.3 UC-21 View Mother Journey Dashboard	113
3.3.1.4 UC-22 Add Maternal Health Metric	113
3.3.1.5 UC-23 Update or Delete Maternal Health Metric	114
3.3.1.6 UC-24 View Maternal Health Trend	115
3.3.1.7 UC-25 Add Postpartum Recovery Log	116
3.3.1.8 UC-26 Update or Delete Postpartum Recovery Log	117
3.3.1.9 UC-27 Manage Preparation Checklist	118
3.3.1.10 UC-28 Browse Pregnancy Exercise Library	119
3.3.1.11 UC-29 Complete Pre-exercise Safety Check	120

3.3.1.12 UC-30 Conduct Pregnancy Exercise Session with Optional Posture Feedback	121
3.3.1.13 UC-31 View Exercise History and Session Result	122
3.3.2 Baby Care Journey, Growth & Vaccination	123
3.3.2.1 UC-32 Create Baby Profile	123
3.3.2.2 UC-33 Update or Archive Baby Profile	124
3.3.2.3 UC-34 Switch Active Baby Profile	124
3.3.2.4 UC-35 View Baby Care Overview	125
3.3.2.5 UC-36 Add Baby Daily Log	126
3.3.2.6 UC-37 Update or Delete Baby Daily Log	127
3.3.2.7 UC-38 View Baby Log Summary	128
3.3.2.8 UC-39 Record Development Milestone	129
3.3.2.9 UC-40 Update or Delete Development Milestone	130
3.3.2.10 UC-41 Add Growth Measurement	131
3.3.2.11 UC-42 Update or Delete Growth Measurement	132
3.3.2.12 UC-43 View Growth Trend and Measurement History	132
3.3.2.13 UC-44 Manage Vaccination Record	133
3.3.2.14 UC-45 View Vaccination Reference Schedule and Reminder Status	134
3.3.3 Community Q&A	135
3.3.3.1 UC-46 Browse Community Feed and Topics	135
3.3.3.2 UC-47 View Community Question Detail	136
3.3.3.3 UC-48 Create Community Question	137
3.3.3.4 UC-49 Edit or Delete Own Community Post	138
3.3.3.5 UC-50 Post Community Answer	139
3.3.3.6 UC-51 Edit or Delete Own Community Answer	140
3.3.3.7 UC-52 React to Community Answer	141
3.3.3.8 UC-53 Bookmark Community Post	142
3.3.3.9 UC-54 Follow or Unfollow Community Topic	143
3.3.3.10 UC-55 Report Unsafe Community Content or Account	144
3.3.4 Verified Expert Discovery	145
3.3.4.1 UC-65 Browse Verified Expert Directory	145
3.3.4.2 UC-66 View Verified Expert Profile	146
3.3.5 AI Nurse Assistant & Risk Triage	147
3.3.5.1 UC-72 Run AI Symptom Intake	147
3.3.5.2 UC-73 View Risk Triage Result	147
3.3.5.3 UC-74 Open Emergency Support from a Red Risk Result	148
3.3.6 Emergency Map & Nearby Care Support	149

3.3.6.1 UC-77 Open Emergency Map	149
3.3.6.2 UC-78 Find Nearby Care Facilities	150
3.3.6.3 UC-79 View Route, ETA and Quick Call or Navigate	151
3.3.6.4 UC-80 Share Time-Limited Location and Send Family Emergency Alert	152
3.3.6.5 UC-81 Create or Cancel Nearby Support Request	153
3.3.7 Personal Health Records & Source Labeling	154
3.3.7.1 UC-83 Add Personal Health Record and Attachment	154
3.3.7.2 UC-84 Update Health Record Metadata	155
3.3.7.3 UC-85 Archive or Delete User-entered Health Record	156
3.3.7.4 UC-86 View Health Record Timeline and Detail	157
3.3.7.5 UC-87 Generate Health Summary	158
3.3.7.6 UC-88 Share Health Summary or Selected Records Under Consent	159
3.3.8 Reminders, Tasks & Care Plan	160
3.3.8.1 UC-89 Create Appointment Reminder	160
3.3.8.2 UC-90 Create Medicine or Vitamin Reminder	161
3.3.8.3 UC-91 Create Vaccination Reminder	161
3.3.8.4 UC-92 Update, Snooze, Complete, Skip or Delete Care Reminder	162
3.3.8.5 UC-93 View Today Tasks and Care Plan	163
3.3.9 Family Sync & Cooperative Care	164
3.3.9.1 UC-94 Create Care Group	164
3.3.9.2 UC-95 Invite or Revoke Family Member Invitation	165
3.3.9.3 UC-96 Accept or Reject Care Group Invitation	166
3.3.9.4 UC-97 Manage Care Group Membership	167
3.3.9.5 UC-98 Manage Family Permission Scope	168
3.3.9.6 UC-99 Create, Update or Cancel Family Care Task	169
3.3.9.7 UC-100 Update Assigned Task Status	170
3.3.9.8 UC-101 View Shared Care Calendar, Data and Alerts	171
3.3.10 Verified Content & Checklist Hub	172
3.3.10.1 UC-102 Browse Verified Content	172
3.3.10.2 UC-103 View Verified Content Detail	172
3.3.11 Expense & Preparation Planner	173
3.3.11.1 UC-109 Add Expense Entry	173
3.3.11.2 UC-110 Update or Delete Expense Entry	174
3.3.11.3 UC-111 View Expense Summary	175
3.3.12 Connected Device & Health Data Integration	176
3.3.12.1 UC-112 Connect Health Device or Platform	176

3.3.12.2 UC-113 Import or Synchronize Device Observations	177
3.3.12.3 UC-114 View Device Data Trend and Quality	178
3.3.12.4 UC-115 Disconnect Device and Delete Imported Data	179
3.3.13 Smart Activity Monitoring & Safety Support	180
3.3.13.1 UC-116 Manage Emergency Contacts	180
3.3.13.2 UC-117 Configure Smart Activity Monitoring	181
3.3.13.3 UC-118 Enable Smart Activity Monitoring	182
3.3.13.4 UC-119 Disable Smart Activity Monitoring	183
3.3.13.5 UC-120 Handle Suspected Fall or Impact and Safety Check	183
3.3.13.6 UC-121 Review Safety Event History and Report False Positive	184
4. Functional Requirements	185
4.1 System Functional Overview	185
4.1.1 Screens Flow	185
4.1.1.1 Mother Mobile App Screens Flow	187
4.1.1.2 Family Member Mobile App Screens Flow	188
4.1.1.3 Verified Expert App Screens Flow	188
4.1.1.4 Verified Expert Web Portal Screens Flow	189
4.1.1.5 System Admin Web Portal Screens Flow	190
4.1.1.6 Content Admin Web Portal Screens Flow	191
4.1.1.7 Moderator Web Portal Screens Flow	191
4.1.1.8 Partner Representative Web Portal Screens Flow	192
4.1.2 Screen Descriptions	192
4.1.3 Screen Authorization	213
4.1.4 Non-Screen Functions	228
4.1.5 Entity Relationship Diagram	235
4.2 Mobile Screen	310
4.2.1 Common Mobile Account Screens	311
4.2.1.1 Mobile Welcome Screen	311
4.2.1.2 Register Account Screen	312
4.2.1.3 Verify OTP Screen	313
4.2.1.4 Login Screen	314
4.2.1.5 Request Password Reset Screen	315
4.2.1.6 Reset Password Screen	316
4.2.1.7 Account Profile and Settings Screen	317
4.2.1.8 Notification Center Screen	319
4.2.1.9 Login Sessions Screen	320

4.2.1.10 Privacy, Consent and Data Sharing Screen	321
4.2.2 Mother Mobile App	322
4.2.2.1 Mother Journey Setup Screen	322
4.2.2.2 Mother Journey Dashboard Screen	323
4.2.2.3 Maternal Metrics and Trend Screen	324
4.2.2.4 Postpartum Recovery Logs Screen	326
4.2.2.5 Preparation Checklist Screen	327
4.2.2.6 Pregnancy Exercise Library Screen	328
4.2.2.7 Pre-Exercise Safety Check Screen	329
4.2.2.8 Active Exercise Session and Optional Posture Feedback Screen	330
4.2.2.9 Exercise History and Session Result Screen	332
4.2.2.10 Baby Profiles Screen	333
4.2.2.11 Baby Care Overview Screen	334
4.2.2.12 Baby Daily Log Screen	335
4.2.2.13 Baby Log Summary Screen	336
4.2.2.14 Development Milestones Screen	337
4.2.2.15 Growth Measurements and Trend Screen	338
4.2.2.16 Vaccination Records and Reference Schedule Screen	340
4.2.2.17 Community Feed Screen	341
4.2.2.18 Community Question Detail Screen	342
4.2.2.19 Community Post Composer Screen	343
4.2.2.20 Community Answer Composer Screen	345
4.2.2.21 Report Unsafe Content or Account Screen	346
4.2.2.22 Verified Expert Directory Screen	347
4.2.2.23 Verified Expert Profile Screen	348
4.2.2.24 AI Symptom Intake Screen	349
4.2.2.25 Risk Triage Result Screen	350
4.2.2.26 Emergency Map and Nearby Facilities Screen	352
4.2.2.27 Family Emergency Alert and Time-Limited Location Screen	353
4.2.2.28 Nearby Support Request Screen	354
4.2.2.29 Health Record Timeline and Detail Screen	355
4.2.2.30 Add or Edit Health Record and Attachment Screen	357
4.2.2.31 Health Summary and Consent Share Screen	358
4.2.2.32 Today Tasks and Care Plan Screen	359
4.2.2.33 Care Reminder Form Screen	360
4.2.2.34 Care Group Management and Family Permissions Screen	361

4.2.2.35 Verified Content Hub Screen	362
4.2.2.36 Verified Content Detail Screen	364
4.2.2.37 Expense Planner and Summary Screen	365
4.2.2.38 Connected Devices and Synchronisation Screen	366
4.2.2.39 Device Data Trend and Quality Screen	367
4.2.2.40 Emergency Contacts and Smart Activity Monitoring Screen	368
4.2.2.41 Safety Check and Incident Alert Screen	370
4.2.2.42 Safety Event History and False Positive Feedback Screen	371
4.2.3 Family Member Mobile App	372
4.2.3.1 Care Group Invitation Screen	372
4.2.3.2 Family Shared Care Home Screen	373
4.2.4 Verified Expert App	374
4.2.4.1 Expert Profile and Verification Documents Screen	374
4.2.4.2 Expert Verification Status and Renewal Screen	376
4.2.4.3 Expert Availability and Service Scope Screen	377
4.2.4.4 Expert Question Queue Screen	378
4.2.4.5 Expert Answer Composer Screen	379
4.2.4.6 Contribution Points and Badges Screen	380
4.2.4.7 Nearby Availability and Support Requests Screen	381
4.3 Web Screen	383
4.3.1 Common Web Account Screens	383
4.3.1.1 Web Portal Login Screen	383
4.3.1.2 Web Account and Session Settings Screen	384
4.3.2 Verified Expert Web Portal	385
4.3.2.1 Expert Portal Dashboard Screen	385
4.3.2.2 Expert Professional Profile and Documents Screen	386
4.3.2.3 Expert Availability and Nearby Support Screen	388
4.3.2.4 Expert Question Queue and Answer Editor Screen	389
4.3.3 Moderator Web Workspace	390
4.3.3.1 Moderator Dashboard and Moderation Queue Screen	390
4.3.3.2 Moderation Item and Report Resolution Screen	391
4.3.3.3 Community Topics and Visibility Screen	392
4.3.4 Content Admin Web Workspace	393
4.3.4.1 Content Dashboard and Library Screen	393
4.3.4.2 Content Editor and Source Management Screen	395
4.3.4.3 Content Review, Publication and Archive Screen	396

4.3.4.4 Content Categories and Stage Topic Mapping Screen	397
4.3.4.5 Approved AI Knowledge Sources Screen	398
4.3.5 System Admin Web Workspace	399
4.3.5.1 System Administration Dashboard and User Accounts Screen	399
4.3.5.2 Expert Verification Queue and Trust Status Screen	400
4.3.5.3 Sensitive Access and Security Events Screen	402
4.3.5.4 AI Risk and Red-Flag Rules Screen	403
5. Non-Functional Requirements	404
5.1 External Interfaces	404
5.1.1 User Interfaces	404
5.1.2 Software and Service Interfaces	405
5.1.3 Hardware, Device and Sensor Interfaces	405
5.1.4 Communication Interfaces	406
5.2 Quality Attributes	406
5.2.1 Usability	406
5.2.2 Reliability	406
5.2.3 Performance	407
5.2.4 Security and Privacy	408
5.2.5 Safety and Compliance	408
5.2.6 Maintainability and Supportability	409
5.2.7 Scalability and Capacity	409
5.2.8 Compatibility and Portability	409
5.2.9 Data Quality, Integrity and Auditability	410
6. Requirement Appendix	410
6.1 Business Rules	410
6.1.1 Account, Identity, Role, and Access Control	410
6.1.2 Privacy, Consent, Health Data, and	411
6.1.3 Community, Moderation, Content, and Partner Governance	412
6.1.4 Expert Verification, Booking, Consultation, Payment, Pricing, Dispute, Refund, and Realtime Communication	413
6.1.5 AI, Health Guidance, Exercise, Growth, and Reminder Safety	413
6.1.6 Map, Location, Nearby Support, and Connected Devices	414
6.1.7 Smart Activity Monitoring and Safety Support (MF21)	415
6.2 Common Requirements	415
6.3 Application Messages List	417
6.3.1 Common Validation and Access Messages	417

6.3.2 Privacy and Consent Messages	418
6.3.3 Community, Content and Partner Governance Messages	418
6.3.4 Expert, Consultation, Payment and Realtime Messages	419
6.3.5 Community, Content and Partner Governance Messages	419
6.4 Other Requirements	420

I. Record of Changes

Date	A* M, D	In charge	Change Description
5/6/2026	A	PhuongNT	1. Requirement Overview 1.1 Context Diagram 1.1.1 Web App Context Diagram 1.1.2 Mobile App Context Diagram
6/6/2026	A	PhuongNT	1.2 Main Workflow 1.2.1 User Onboarding & Profile Setup Workflow 1.2.2 Mother Care Journey Workflow 1.2.3 Community Q&A and Moderation Workflow 1.2.4 Expert Verification and Consultation Workflow 1.2.5 AI Triage and Safe Action Workflow 1.2.6 Emergency Map and Nearby Care Workflow 1.2.7 Family Sync Workflow
6/6/2026	A	HuyND	4. Functional Requirements 4.1 System Functional Overview 4.1.1 Screens Flow
7/6/2026	A	HuyND	4. Functional Requirements 4.1 System Functional Overview 4.1.1 Screens Flow
8/6/2026	A	HuyND	4. Functional Requirements 4.1 System Functional Overview 4.1.1 Screens Flow
8/6/2026	A	HuyND	4.1.2 Screen Descriptions 4.1.3 Screen Authorization
8/6/2026	M	LamVH	Read and review UCs
9/6/2026	A	PhuongNT	1.2.8 Smart Activity Monitoring & Safety Support Workflow 1.2.9 Expert App & Nearby Support Workflow
9/6/2026	A	ChuongBD	5. Non-Functional Requirements 5.1 External Interfaces 5.1.1 User Interfaces 5.1.2 Software and Service Interfaces 5.1.3 Hardware, Device and Sensor Interfaces 5.1.4 Communication Interfaces
9/6/2026	A	LamVH	Write 6. Requirement Appendix, add 60 BR in 6.1
10/6/2026	A	LamVH	Add BR in 6.1.2; 6.1.3; 6.1.4; 6.1.5; 6.1.6; 6.1.7
10/6/2026	A	LamVH	Docx config footer
10/6/2026	A	ChuongBD	5.2 Quality Attributes 5.2.1 Usability 5.2.2 Reliability 5.2.3 Performance 5.2.4 Security and Privacy 5.2.5 Safety and Compliance 5.2.6 Maintainability and Supportability 5.2.7 Scalability and Capacity 5.2.8 Compatibility and Portability 5.2.9 Data Quality, Integrity and Auditability
10/6/2026	A	LamVH	6.3 Application Messages List

II. Software Requirement Specification

1. Requirement Overview

CareBridge is a community-driven maternal and early childhood health support platform. The system boundary includes Mobile App, Expert Portal, Admin & Moderation Portal, Partner Portal, Backend Services, AI/RAG/Triage Service, Notification Service, Map Service, Health Device Platform and Phone IMU Sensor integration. CareBridge supports mothers, postpartum users, caregivers and families in tracking care journeys, saving personal records, joining a moderated community, reading verified content, receiving reminders, connecting with verified experts and using safety support features.

External entities include human actors such as Guest, User, Care Seeker, Family Member/Caregiver, Verified Expert, Moderator, Content Admin, Admin and Partner Representative; supporting service actors such as Payment Gateway, Notification Service, Map Service, Health Device Platform, AI Service and Phone IMU Sensor; and partner ecosystem entities such as partner clinics and sponsored service providers. Data exchange is controlled by RBAC, consent, scope, expiry, revoke and audit log principles.

1.1 Context Diagram

1.1.1 Web App Context Diagram

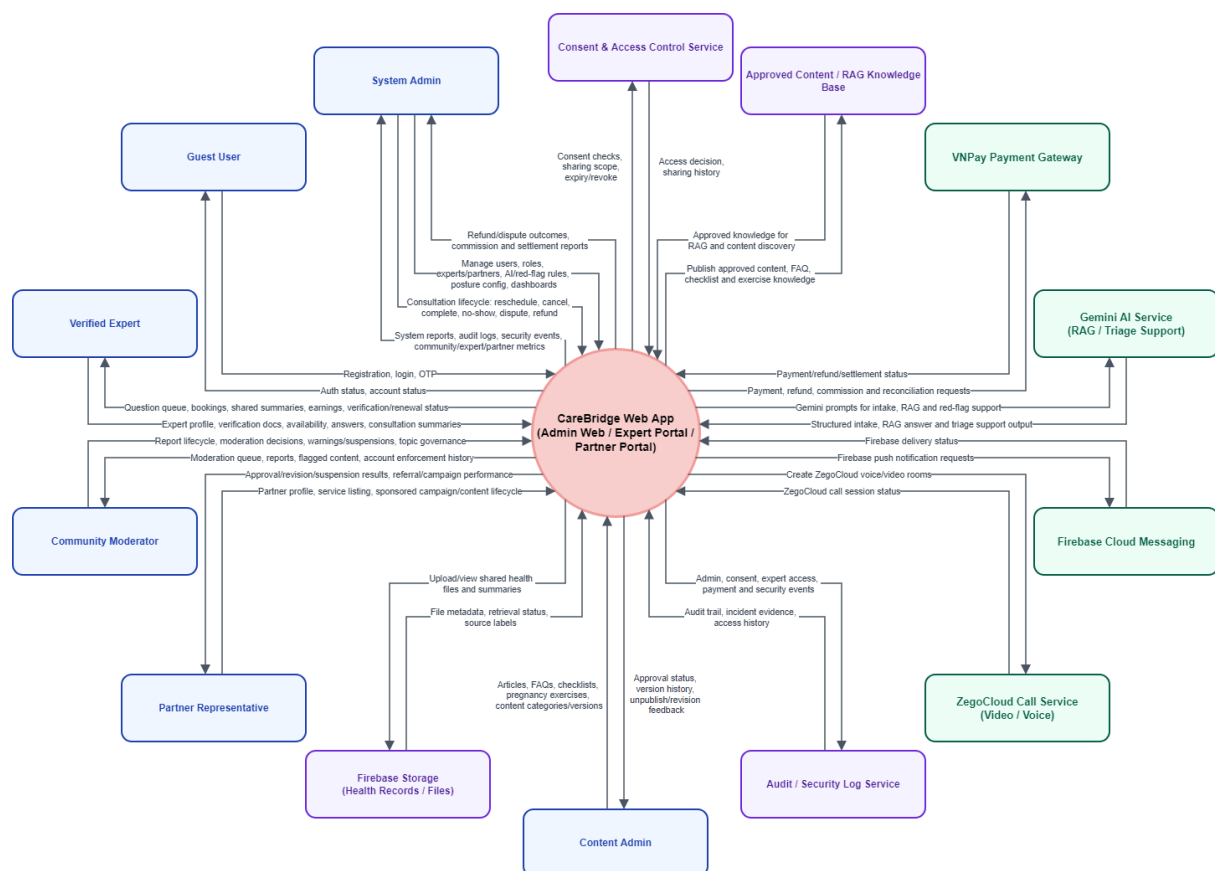


Figure 1: Context Diagram for Web App

Image Detail: [Context Diagram for Web App](#)

Description

External Entity	Input to System Data Flow	Output from System Data Flow	Description
Guest / User	Registration, login, OTP	Auth status, account status	Primary web users for shared account, profile, search, report, and access-controlled functions.
Verified Expert	Expert profile, verification docs, availability, answers, consultation summaries	Question queue, bookings, shared summaries, earnings, verification/renewal status	Uses Expert Portal / Expert App capabilities to provide verified consultation and expert community support.
Community Moderator	Report lifecycle, moderation decisions, warnings/suspensions, topic governance	Moderation queue, reports, flagged content, account enforcement history	Reviews community content and account reports to keep interaction safe and respectful.
Content Admin	Articles, FAQs, checklists, pregnancy exercises, content categories/versions	Approval status, version history, unpublish/revision feedback	Manages verified content and knowledge materials used by users and Gemini RAG support.
System Admin	- Manage users, roles, experts/partners, AI/red-flag rules, posture config, dashboards - Consultation lifecycle: reschedule, cancel, complete, no-show, dispute, refund	- System reports, audit logs, security events, community/expert/partner metrics - Refund/dispute outcomes, commission and settlement reports	Administers platform operation, access, expert governance, pricing configuration, and security review.
Partner Representative	Partner profile, service listing, sponsored campaign/content lifecycle	Approval/revision/suspension results, referral/campaign performance	Represents partner clinics or sponsored services through Partner Portal.
VNPay Payment Gateway	Payment/refund/settlement status	Payment, refund, commission and reconciliation requests	Processes paid consultation transactions and returns payment results to CareBridge.

External Entity	Input to System Data Flow	Output from System Data Flow	Description
Firebase Cloud Messaging	Firebase delivery status	Firebase push notification requests	Delivers push and in-app notifications for consultation, moderation, reminder, and safety events.
Gemini AI Service (RAG / Triage Support)	Structured intake, RAG answer and triage support output	Gemini prompts for intake, RAG and red-flag support	Supports intake, FAQ/RAG response, summarization, and triage assistance without making independent medical decisions.
ZegoCloud Call Service (Video / Voice)	ZegoCloud call session status	Create ZegoCloud voice/video rooms	Provides realtime voice and video call rooms for consultation sessions.
Consent & Access Control Service	Access decision, sharing history	Consent checks, sharing scope, expiry/revoke	Controls whether health records, summaries, location, and shared care data may be accessed.
Firebase Storage (Health Records / Files)	File metadata, retrieval status, source labels	Upload/view shared health files and summaries	Stores health records and supporting files under CareBridge permission and consent rules.
Audit / Security Log Service	Audit trail, incident evidence, access history	Admin, consent, expert access, payment and security events	Preserves traceability for sensitive platform actions and security investigation.
Approved Content / RAG Knowledge Base	Approved knowledge for RAG and content discovery	Publish approved content, FAQ, checklist and exercise knowledge	Holds approved educational and checklist content used for verified content browsing and RAG retrieval.

Table 2: Context Diagram Description For Web App

1.1.2 Mobile App Context Diagram

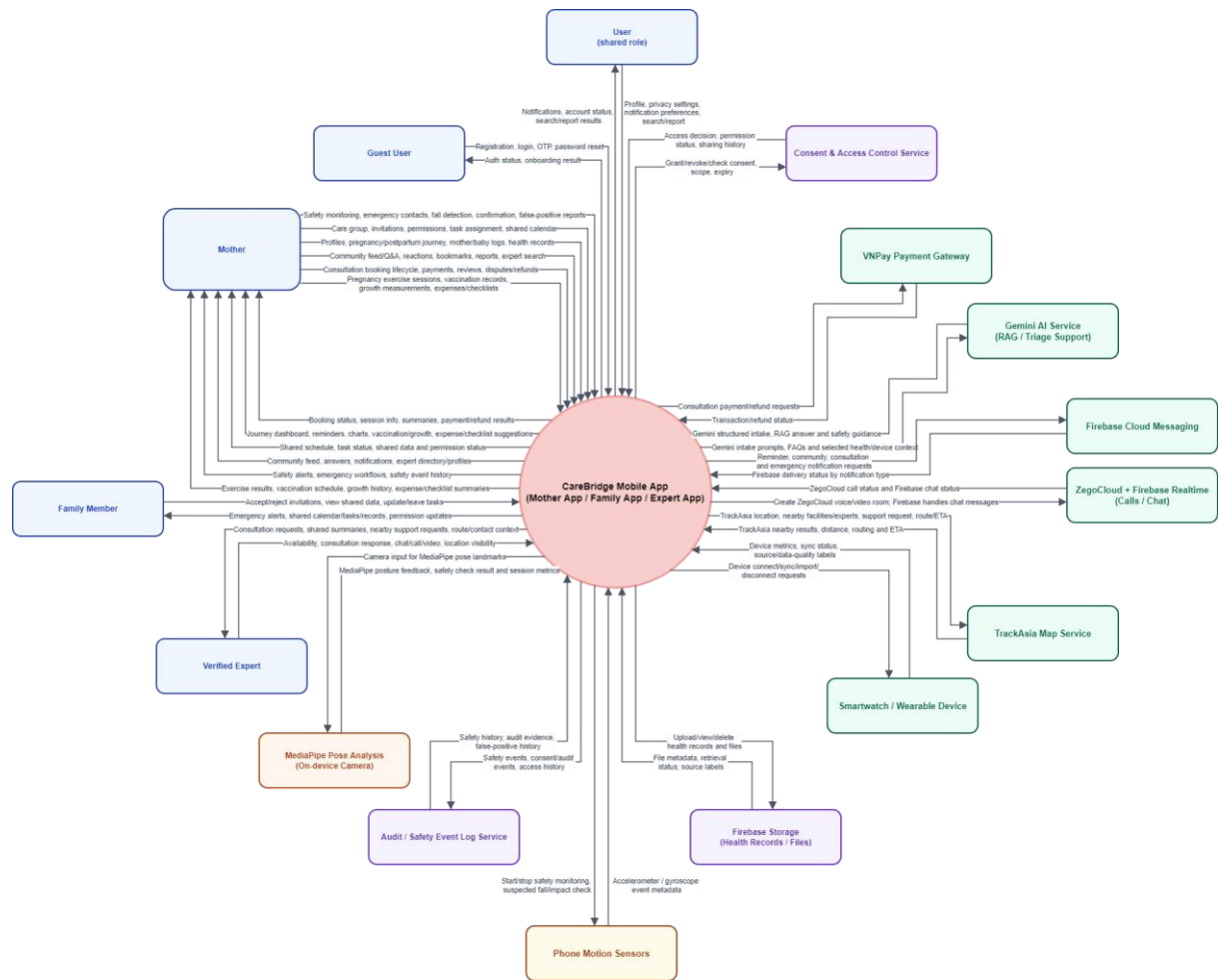


Figure 2: Context Diagram for Mobile App

Image Detail: [Context Diagram for Mobile App](#)

Description

External Entity	Input to System Data Flow	Output from System Data Flow	Description
Mother	- Profiles, pregnancy/postpartum journey, mother/baby logs, health records - Community feed/Q&A, reactions, bookmarks, reports, expert search - Consultation booking lifecycle, payments, reviews, disputes/refunds - Safety monitoring, emergency contacts, fall detection, confirmation, false-positive reports - Care group, invitations, permissions, task assignment, shared calendar - Pregnancy exercise sessions, vaccination records, growth measurements, expenses/checklists	- Journey dashboard, reminders, charts, vaccination/growth, expense/checklist suggestions - Community feed, answers, notifications, expert directory/profiles - Booking status, session info, summaries, payment/refund results - Safety alerts, emergency workflows, safety event history - Shared schedule, task status, shared data and permission status - Exercise results, vaccination schedule, growth history, expense/checklist summaries	Main mobile user managing maternal, baby, health record, consultation, emergency, and safety support features.
Family Member	Accept/reject invitations, view shared data, update/leave tasks	Emergency alerts, shared calendar/tasks/records, permission updates	Supports the mother/care seeker within a consented family care group.
Verified Expert	Availability, consultation response, chat/call/video, location visibility	Consultation requests, shared summaries, nearby support requests, route/contact context	Uses mobile/expert app capabilities for availability, nearby support, and realtime consultation.
Guest / User	- Registration, login, OTP, password reset - Profile, privacy settings, notification preferences, search/report	- Auth status, onboarding result - Notifications, account status, search/report results	Shared mobile account actor for authentication and general user functions.
VNPay Payment Gateway	Transaction/refund status	Consultation payment/refund requests	Processes paid consultation checkout from the mobile booking flow.
TrackAsia Map Service	TrackAsia nearby results, distance, routing and ETA	TrackAsia location, nearby facilities/experts, support request, route/ETA	Supports emergency map, nearby care, and consented expert support location features.

External Entity	Input to System Data Flow	Output from System Data Flow	Description
Smartwatch / Wearable Device	Device metrics, sync status, source/data-quality labels	Device connect/sync/import/disconnect requests	Provides supported wearable health data only after consent and with source labels.
Gemini AI Service (RAG / Triage Support)	Gemini structured intake, RAG answer and safety guidance	Gemini intake prompts, FAQs and selected health/device context	Supports AI nurse assistant and triage flows without replacing professional medical advice.
Firebase Cloud Messaging	Firebase delivery status by notification type	Reminder, community, consultation and emergency notification requests	Delivers mobile notifications and safety alerts.
ZegoCloud + Firebase Realtime (Calls / Chat)	ZegoCloud call status and Firebase chat status	- Create ZegoCloud voice/video room - Firebase handles chat messages	Supports chat through Firebase and voice/video calls through ZegoCloud.

Table 3: Context Diagram Description For Mobile App

1.2 Main Workflow

1.2.1 User Onboarding, Account & Privacy Setup Workflow

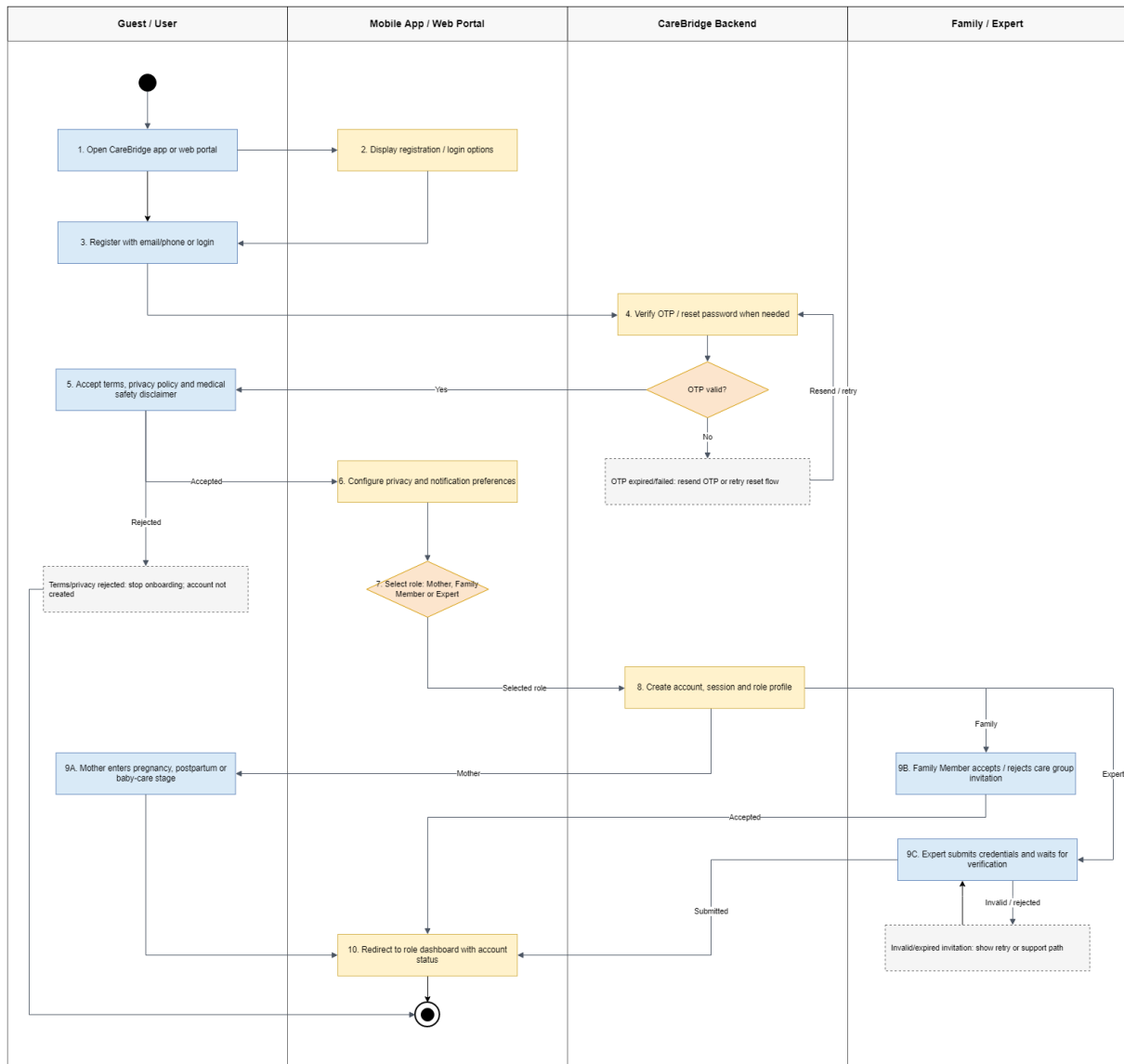


Figure 3: User Onboarding, Account & Privacy Setup Workflow

Image Detail: [User Onboarding, Account & Privacy Setup Workflow](#)

Description

#	Step Name	Detail Description	Role	Note
1	Open CareBridge app or web portal	Guest/User starts onboarding from the mobile app or web portal.	Guest / User	Start node in diagram.
2	Display registration / login options	The system shows sign-up, login and account recovery options.	Mobile App / Web Portal	-
3	Register with email/phone or login	User submits email/phone credentials or existing login data.	Guest / User	-
4	Verify OTP / reset password when needed	Backend validates OTP or supports password reset when required.	CareBridge Backend	Decision follows: OTP valid?
D-1	OTP valid?	If valid, continue to terms/privacy. If failed or expired, user retries OTP or reset flow.	CareBridge Backend	Yes / No; resend or retry path.
5	Accept terms, privacy policy and medical safety disclaimer	User must accept legal, privacy and medical safety terms before account creation.	Guest / User	Rejected path stops onboarding; account not created.
6	Configure privacy and notification preferences	User configures privacy controls and notification preferences.	Mobile App / Web Portal	-
7	Select role: Mother, Family Member or Expert	User selects intended role for profile setup and dashboard routing.	Mobile App / Web Portal	Branches to Mother, Family or Expert.
8	Create account, session and role profile	Backend creates account, session and role-specific profile.	CareBridge Backend	After role selection.
9A	Mother enters pregnancy, postpartum or baby-care stage	Mother completes stage-specific onboarding data.	Guest / User	Mother branch.
9B	Family Member accepts / rejects care group invitation	Family member handles invitation before dashboard access.	Family / Expert	Invalid/expired invitation shows retry or support path.
9C	Expert submits credentials and waits for verification	Expert submits credential documents and remains pending verification.	Family / Expert	Expert branch.
10	Redirect to role dashboard with account status	System opens the correct role dashboard and shows account status.	Mobile App / Web Portal	End of onboarding.

Table 4: User Onboarding, Account & Privacy Setup Workflow Description

1.2.2 Mother, Baby & Health Record Workflow

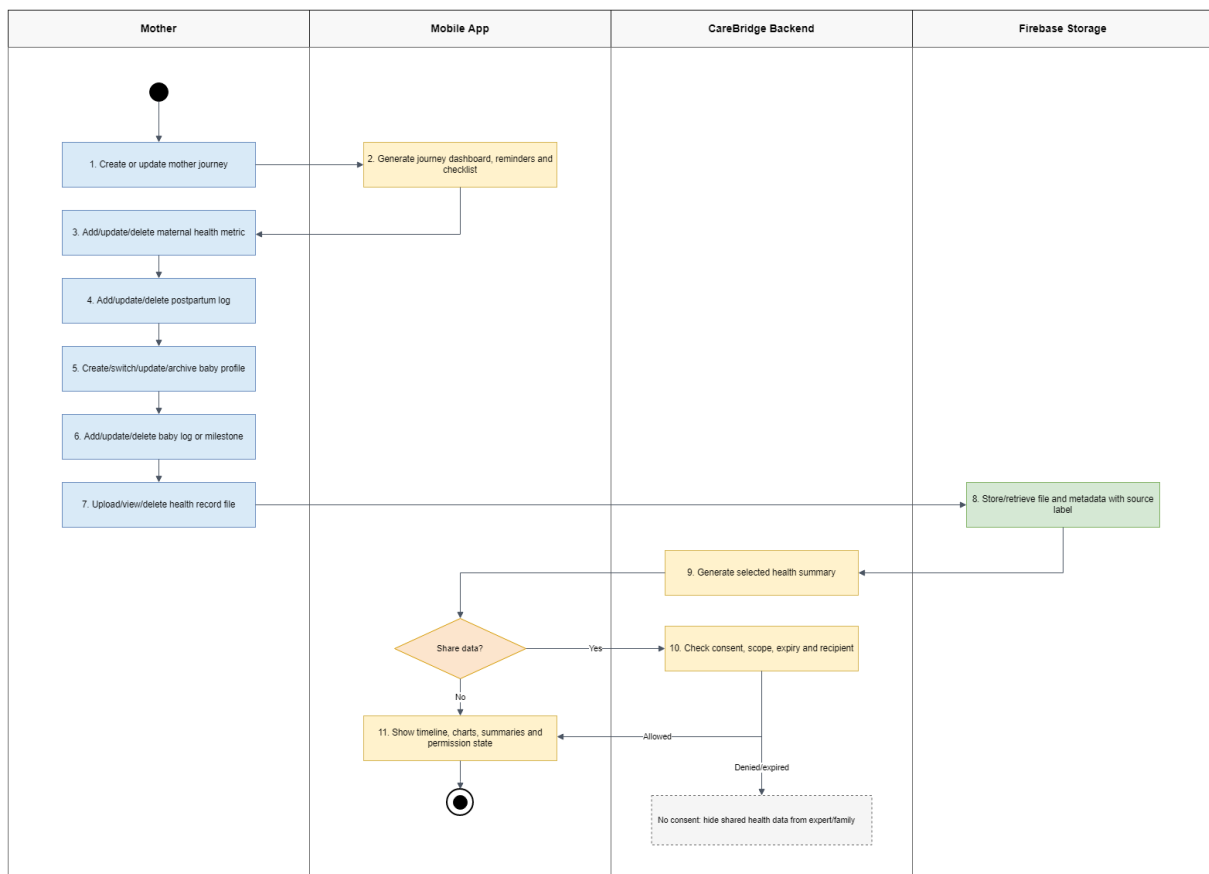


Figure 4: Mother, Baby & Health Record Workflow

Image Detail: [Mother, Baby & Health Record Workflow](#)

Description

#	Step Name	Detail Description	Role	Note
1	Create or update mother journey	Mother creates or updates pregnancy, postpartum or baby-care journey information.	Mother	-
2	Generate journey dashboard, reminders and checklist	Mobile app builds dashboard, reminders and checklist from journey data.	Mobile App	-
3	Add/update/delete maternal health metric	Mother manages maternal health metrics.	Mother	-
4	Add/update/delete postpartum log	Mother manages postpartum logs.	Mother	-
5	Create/switch/update/archive baby profile	Mother manages active baby profile and baby profile lifecycle.	Mother	-
6	Add/update/delete baby log or milestone	Mother records baby daily logs and development milestones.	Mother	-
7	Upload/view/delete health record file	Mother manages health record files from the mobile app.	Mother	-
8	Store/retrieve file and metadata with source label	Firebase Storage stores or retrieves file data and metadata with source labels.	Firebase Storage	-
9	Generate selected health summary	Backend prepares selected health summary for viewing or sharing.	CareBridge Backend	-
D-1	Share data?	If sharing is requested, backend checks consent, scope, expiry and recipient. If not, app shows private timeline and summaries.	Mobile App	Yes / No.
10	Check consent, scope, expiry and recipient	Backend validates whether shared data may be exposed to expert or family.	CareBridge Backend	Allowed / denied or expired.
11	Show timeline, charts, summaries and permission state	Mobile app displays timeline, charts, summaries and current permission state.	Mobile App	-
EX-1	No consent	Shared health data is hidden from expert/family when consent is denied or expired.	CareBridge Backend	Exception path.

Table 5: Mother, Baby & Health Record Workflow Description

1.2.3 Community Q&A, Expert Answer & Moderation Workflow

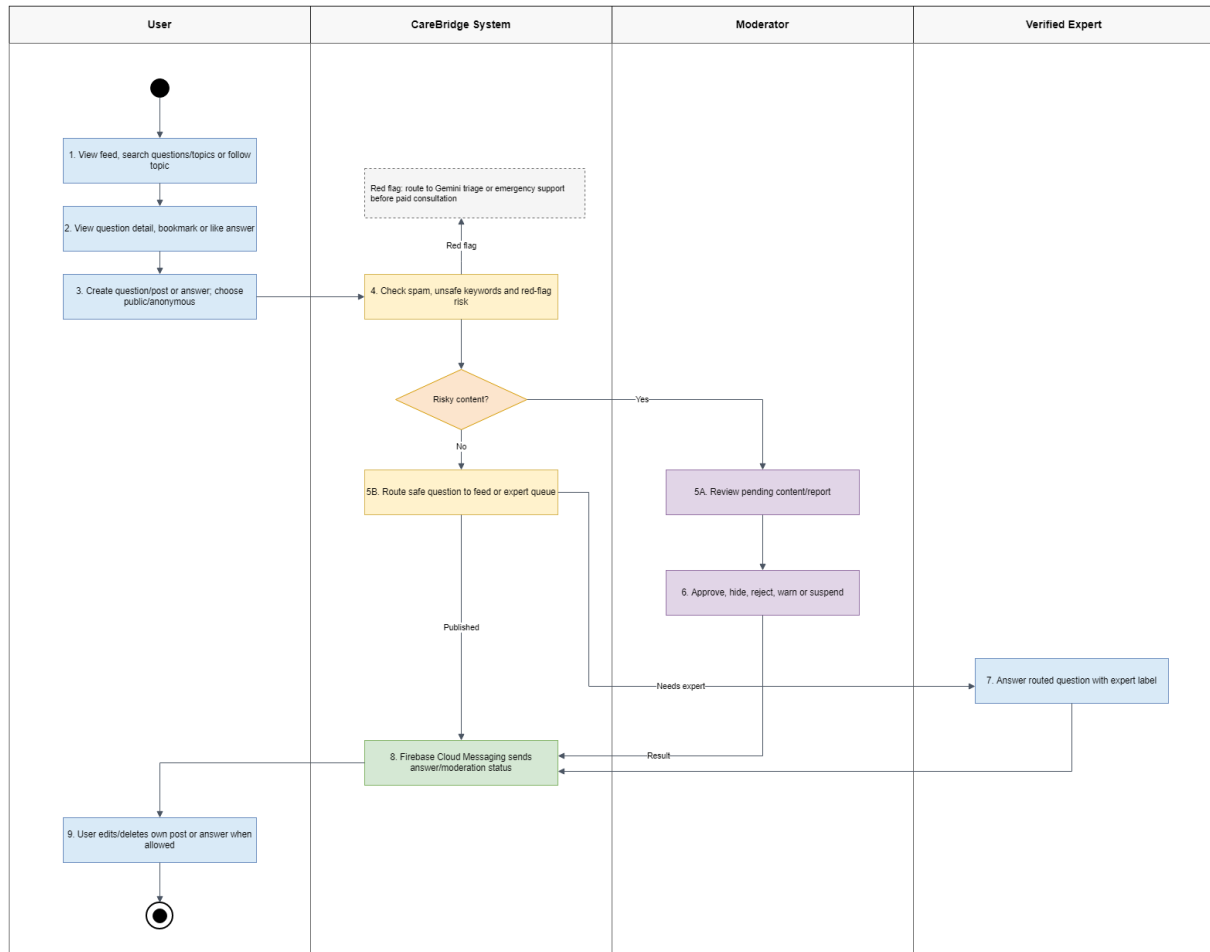


Figure 5: Community Q&A, Expert Answer & Moderation Workflow

Image Detail: [Community Q&A, Expert Answer & Moderation Workflow](#)

Description

#	Step Name	Detail Description	Role	Note
1	View feed, search questions/topics or follow topic	User browses community feed, searches topics/questions or follows a topic.	User	-
2	View question detail, bookmark or like answer	User opens detail and performs allowed engagement actions.	User	-
3	Create question/post or answer; choose public/anonymous	User creates community content and selects visibility identity.	User	-

#	Step Name	Detail Description	Role	Note
4	Check spam, unsafe keywords and red-flag risk	CareBridge System screens the submitted content before publication.	CareBridge System	-
D-1	Risky content?	Risky content goes to moderation; safe content routes to feed or expert queue.	CareBridge System	Yes / No.
5A	Review pending content/report	Moderator reviews reported or pending content.	Moderator	Risky branch.
6	Approve, hide, reject, warn or suspend	Moderator applies content/account enforcement decision.	Moderator	Result returns to notification step.
5B	Route safe question to feed or expert queue	System publishes safe question or routes it for expert answer.	CareBridge System	Needs expert / published.
7	Answer routed question with expert label	Verified Expert answers routed question with expert identity label.	Verified Expert	-
8	Firebase Cloud Messaging sends answer/moderation status	System sends answer, publication or moderation status notifications.	CareBridge System	Uses Firebase Cloud Messaging.
9	User edits/deletes own post or answer when allowed	User updates or deletes own allowed content after publication.	User	-
EX-1	Red flag	Red-flag content is routed to Gemini triage or emergency support before paid consultation.	CareBridge System	Safety branch.

Table 6: Community Q&A, Expert Answer & Moderation Workflow Description

1.2.4 Expert Verification, Booking & Consultation Workflow

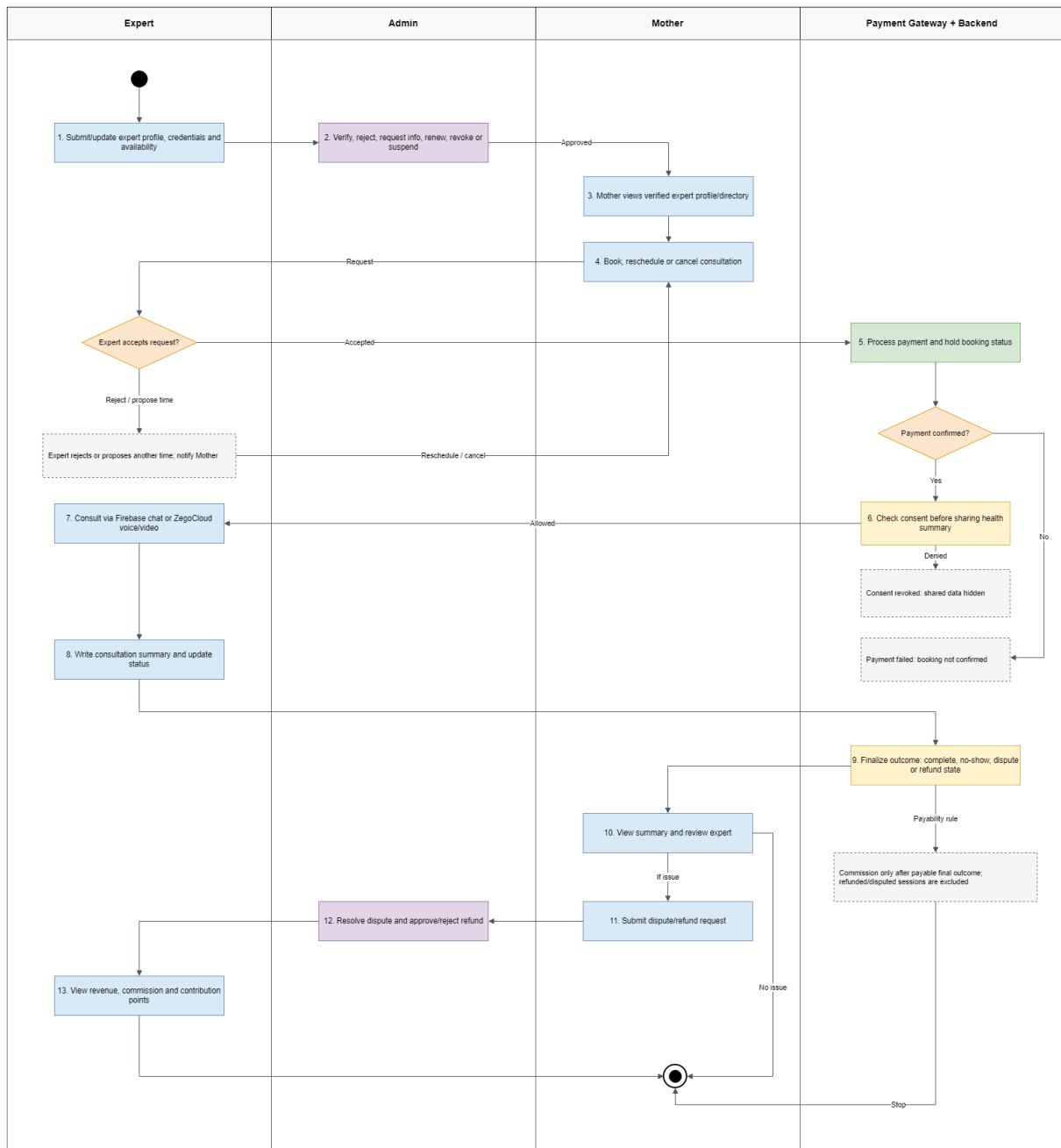


Figure 6: Expert Verification, Booking & Consultation Workflow

Image Detail: [Expert Verification, Booking & Consultation Workflow](#)

Description

#	Step Name	Detail Description	Role	Note
1	Submit/update expert profile, credentials and availability	Expert submits or updates profile, credentials and availability.	Expert	-
2	Verify, reject, request info, renew, revoke or suspend	Admin reviews expert credentials and manages verification lifecycle.	Admin	Approved path exposes expert to directory.
3	Mother views verified expert profile/directory	Mother browses verified expert profiles and directory.	Mother	-
4	Book, reschedule or cancel consultation	Mother requests, reschedules or cancels consultation.	Mother	-
D-1	Expert accepts request?	Experts can accept, reject or propose another time before payment proceeds.	Expert	Reject/propose time notifies Mother.
5	Process payment and hold booking status	Payment gateway and backend process payment while booking remains pending.	Payment Gateway + Backend	-
D-2	Payment confirmed?	Confirmed payment moves to consent check; failed payment stops booking confirmation.	Payment Gateway + Backend	Yes / No.
6	Check consent before sharing health summary	Backend checks permission before exposing selected health summary.	Payment Gateway + Backend	Allowed / denied.
7	Consult via Firebase chat or ZegoCloud voice/video	Expert and Mother conduct consultation through chat, voice or video.	Expert	Only after payment and required consent path.
8	Write consultation summary and update status	Expert writes consultation summary and updates session status.	Expert	-
9	Finalize outcome: complete, no-show, dispute or refund state	Backend/payment layer records final payable or non-payable outcome.	Payment Gateway + Backend	Commission only after payable final outcome.
10	View summary and review expert	Mother views summary and submits expert review.	Mother	-
11	Submit dispute/refund request	Mother submits dispute or refund request when there is an issue.	Mother	If issue.

#	Step Name	Detail Description	Role	Note
12	Resolve dispute and approve/reject refund	Admin resolves dispute and makes refund decision.	Admin	-
13	View revenue, commission and contribution points	Expert views revenue, commission and contribution points after final outcome.	Expert	Refunded/disputed sessions are excluded from commission.
EX-1	Payment failed	Booking is not confirmed when payment fails.	Payment Gateway + Backend	Stop path.
EX-2	Consent revoked	Shared data is hidden when consent is denied or revoked.	Payment Gateway + Backend	Consultation may continue without shared health summary.

Table 7: Expert Verification, Booking & Consultation Workflow

1.2.5 Gemini AI Triage & Safe Action Workflow

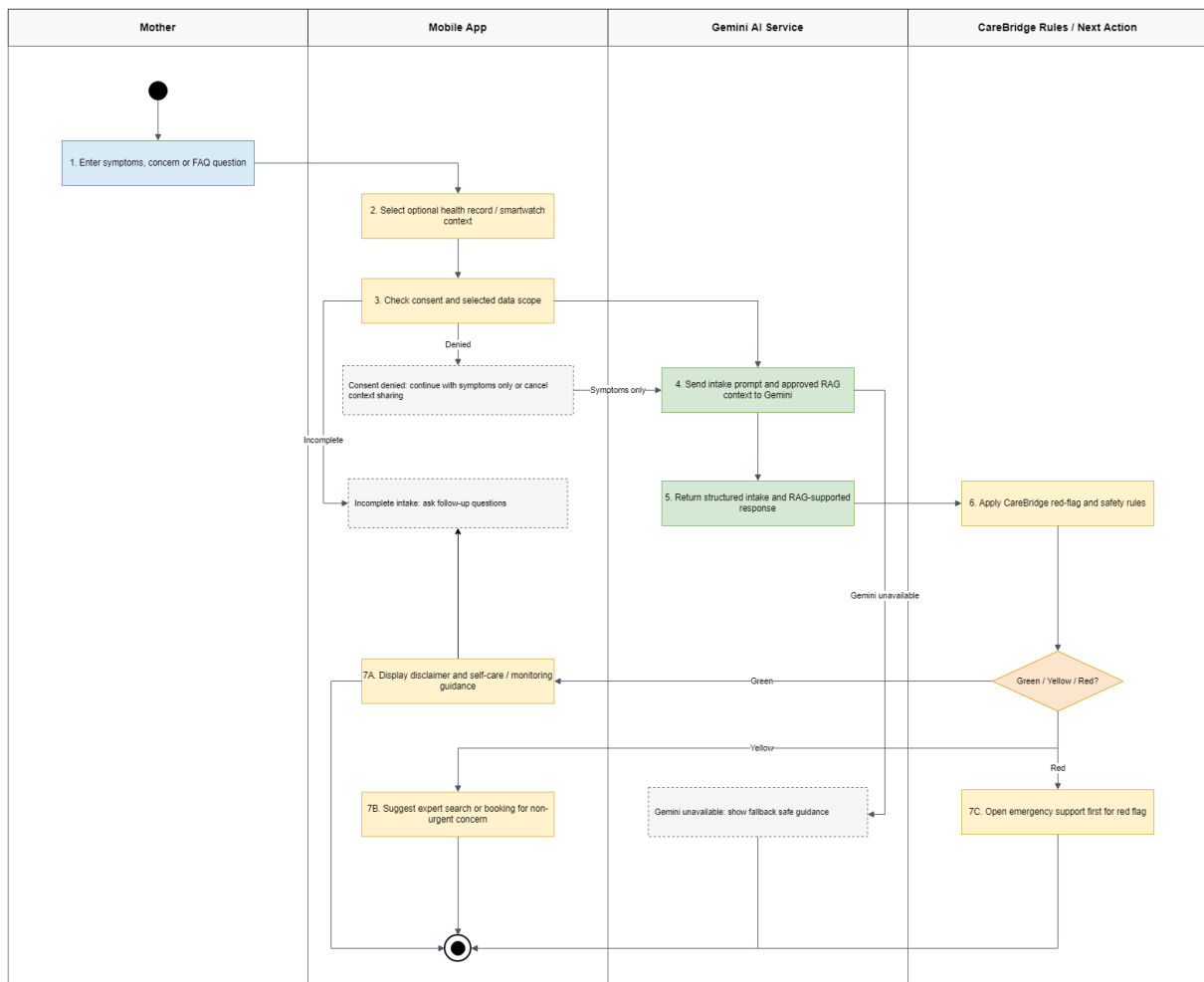


Figure 7: Gemini AI Triage & Safe Action Workflow

Image Detail: [Gemini AI Triage & Safe Action Workflow](#)

Description

#	Step Name	Detail Description	Role	Note
1	Enter symptoms, concern or FAQ question	Mother submits symptoms, concerns or FAQ-style questions.	Mother	-
2	Select optional health record / smartwatch context	Mobile app lets Mother choose optional health or device context.	Mobile App	-
3	Check consent and selected data scope	Mobile app validates consent and selected sharing scope.	Mobile App	Denied path may continue with symptoms only.
4	Send intake prompt and approved RAG context to Gemini	Mobile app sends approved prompt/context to Gemini.	Gemini AI Service	-
5	Return structured intake and RAG-supported response	Gemini returns structured intake and RAG-supported response.	Gemini AI Service	-
6	Apply CareBridge red-flag and safety rules	CareBridge applies safety rules over Gemini output.	CareBridge Rules / Next Action	-
D-1	Green / Yellow / Red?	Safety classification chooses guidance, expert suggestion or emergency-first action.	CareBridge Rules / Next Action	Green / Yellow / Red.
7A	Display disclaimer and self-care / monitoring guidance	App displays safe guidance and disclaimer for green cases.	Mobile App	-
7B	Suggest expert search or booking for non-urgent concern	App suggests expert search or booking for yellow cases.	Mobile App	-
7C	Open emergency support first for red flag	System opens emergency support before other actions for red cases.	CareBridge Rules / Next Action	-
EX-1	Incomplete intake	The app asks follow-up questions before the Gemini response is finalized.	Mobile App	Incomplete branch.
EX-2	Consent denied	Users can continue with symptoms only or cancel context sharing.	Mobile App	Symptoms-only path still reaches Gemini.
EX-3	Gemini unavailable	The app shows fallback safe guidance when Gemini is unavailable.	Gemini AI Service	Fallback branch.

Table 8: Gemini AI Triage & Safe Action Workflow Description

1.2.6 TrackAsia Emergency Map & Nearby Care Workflow

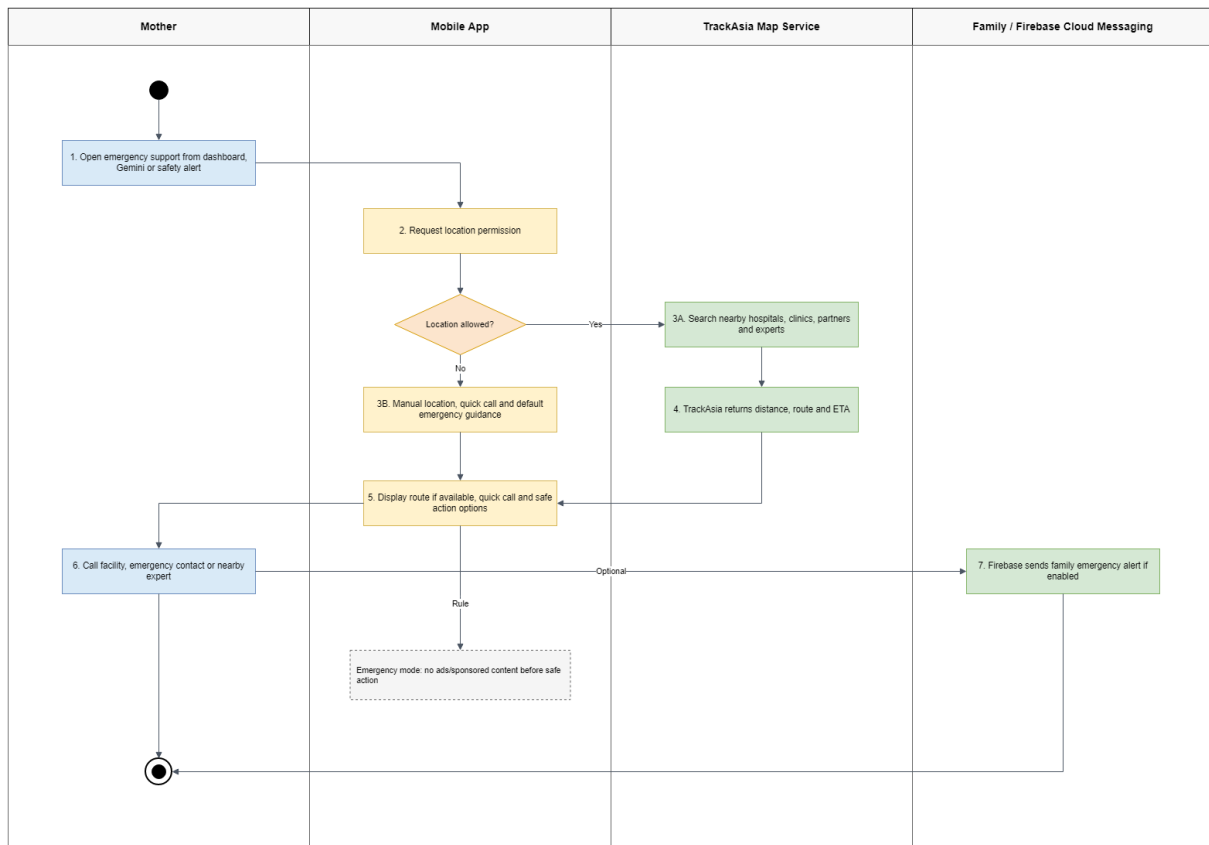


Figure 8: TrackAsia Emergency Map & Nearby Care Workflow

Image Detail: [TrackAsia Emergency Map & Nearby Care Workflow](#)

Description

#	Step Name	Detail Description	Role	Note
1	Open emergency support from dashboard, Gemini or safety alert	Mother opens emergency support from dashboard, AI triage or safety alert.	Mother	-
2	Request location permission	Mobile app requests location permission.	Mobile App	-
D-1	Location allowed?	Allowed location enables nearby search; denied location uses manual/default guidance.	Mobile App	Yes / No.
3A	Search nearby hospitals, clinics, partners and experts	TrackAsia searches nearby care and support options.	TrackAsia Map Service	Location allowed branch.
3B	Manual location, quick call and default emergency guidance	Mobile app provides manual location, quick call and default emergency guidance.	Mobile App	Location denied branch.

#	Step Name	Detail Description	Role	Note
4	TrackAsia returns distance, route and ETA	TrackAsia returns route data, distance and ETA.	TrackAsia Map Service	-
5	Display route if available, quick call and safe action options	Mobile app displays route, quick call and safe action options.	Mobile App	-
6	Call facility, emergency contact or nearby expert	Mother initiates call to facility, emergency contact or nearby expert.	Mother	-
7	Firebase sends family emergency alert if enabled	Firebase Cloud Messaging sends a family emergency alert when enabled.	Family / Firebase Cloud Messaging	Optional.
EX-1	Emergency mode: no ads/sponsored content before safe action	Emergency workflow suppresses ads or sponsored content until safe action is available.	Mobile App	Rule.

Table 9: TrackAsia Emergency Map & Nearby Care Workflow Description

1.2.7 Family Sync, Permission & Care Task Workflow

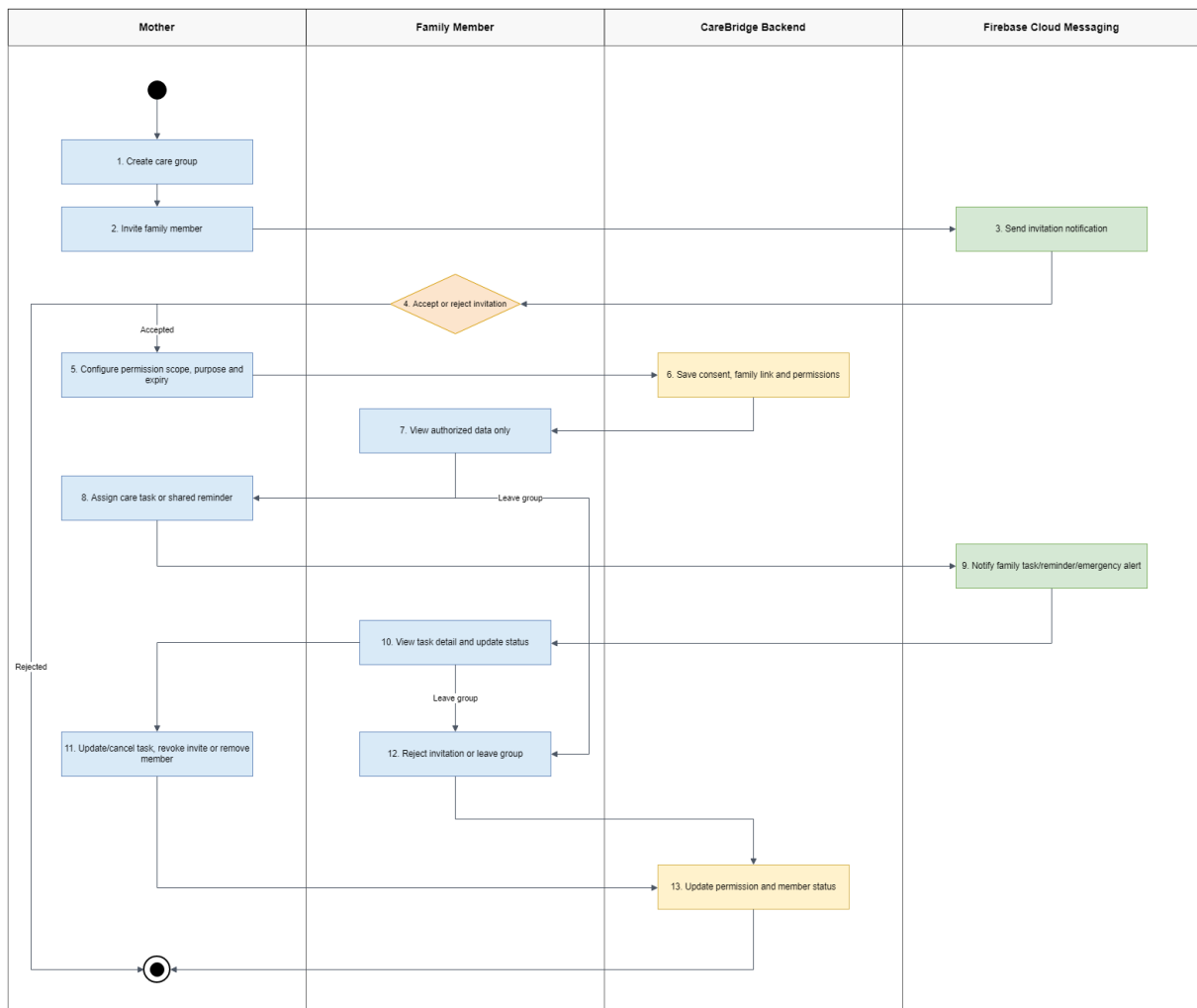


Figure 9: Family Sync, Permission & Care Task Workflow

Image Detail: [Family Sync, Permission & Care Task Workflow](#)

Description

#	Step Name	Detail Description	Role	Note
1	Create care group	Mother creates a care group.	Mother	-
2	Invite family member	Mother invites a family member.	Mother	-
3	Send invitation notification	Firebase Cloud Messaging sends invitation notifications.	Firebase Cloud Messaging	-
4	Accept or reject invitation	Family members accept or reject invitations.	Family Member	Accepted / rejected.
5	Configure permission scope, purpose and expiry	Mother sets permission scope, purpose and expiry.	Mother	-
6	Save consent, family link and permissions	Backend stores consent, family link and permission state.	CareBridge Backend	-
7	View authorized data only	Family members see only authorized data.	Family Member	-
8	Assign care task or shared reminder	Mother assigns a care task or shared reminder.	Mother	-
9	Notify family task/reminder/emergency alert	Firebase Cloud Messaging sends a task, reminder or emergency notification.	Firebase Cloud Messaging	-
10	View task detail and update status	Family members view task details and update status.	Family Member	-
11	Update/cancel task, revoke invite or remove member	Mother manages tasks, invitations and member lifecycle.	Mother	-
12	Reject invitation or leave group	Family member rejects the invitation or leaves the group.	Family Member	Leave the group branch.
13	Update permission and member status	Backend updates permission and member status after lifecycle changes.	CareBridge Backend	-

Table 10: Family Sync, Permission & Care Task Workflow Description

1.2.8 Smart Activity Monitoring & Safety Support Workflow

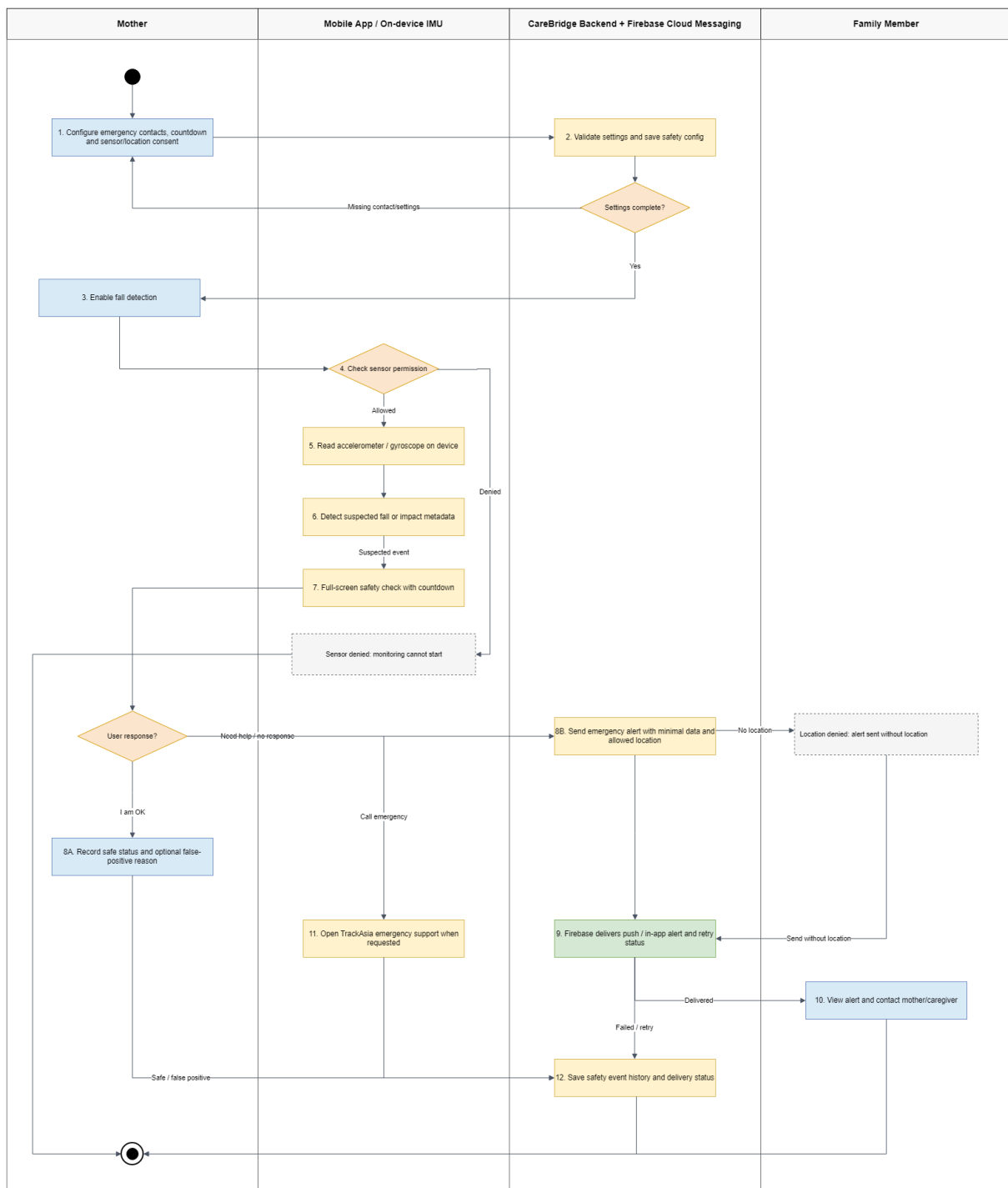


Figure 10: Smart Activity Monitoring & Safety Support Workflow

Image Detail: [Smart Activity Monitoring & Safety Support Workflow](#)

Description

#	Step Name	Detail Description	Role	Note
1	Configure emergency contacts, countdown and sensor/location consent	Mother configures emergency contacts, countdown and monitoring consent.	Mother	-
2	Validate settings and save safety config	Backend validates and saves safety configuration.	CareBridge Backend + Firebase Cloud Messaging	-
D-1	Settings complete?	Complete settings allow monitoring; missing contact/settings require correction.	CareBridge Backend + Firebase Cloud Messaging	Yes / missing contact/settings.
3	Enable fall detection	Mother enables fall detection.	Mother	-
4	Check sensor permission	App checks sensor permission before monitoring starts.	Mobile App / On-device IMU	Allowed / denied.
5	Read accelerometer / gyroscope on device	On-device IMU reads accelerometer and gyroscope signals.	Mobile App / On-device IMU	-
6	Detect suspected fall or impact metadata	On-device logic detects suspected fall/impact metadata.	Mobile App / On-device IMU	-
7	Full-screen safety check with countdown	App displays full-screen safety confirmation with countdown.	Mobile App / On-device IMU	-
D-2	User response?	I am OK records safe/false-positive result; need help or no response sends emergency alert.	Mother	I am OK / need help / no response.
8A	Record safe status and optional false-positive reason	Mother records safe status and optional false-positive reason.	Mother	Safe / false positive path.
8B	Send emergency alert with minimal data and allowed location	Backend sends emergency alert with minimal data and allowed location when permitted.	CareBridge Backend + Firebase Cloud Messaging	Location denied sends alert without location.
9	Firebase delivers push / in-app alert and retry status	Firebase delivers alert and records delivery/retry status.	CareBridge Backend + Firebase Cloud Messaging	Delivered / failed retry.
10	View alert and contact mother/caregiver	Family member receives alert and contacts mother/caregiver.	Family Member	-

#	Step Name	Detail Description	Role	Note
11	Open TrackAsia emergency support when requested	App opens TrackAsia emergency support when user requests emergency help.	Mobile App / On-device IMU	Call emergency branch.
12	Save safety event history and delivery status	Backend saves safety event history and notification delivery status.	CareBridge Backend + Firebase Cloud Messaging	-
EX-1	Sensor denied	Monitoring cannot start when sensor permission is denied.	Mobile App / On-device IMU	Exception path.
EX-2	Location denied	Alert is sent without location when location permission is denied.	Family Member	Exception path shown in diagram lane.

Table 11: Smart Activity Monitoring & Safety Support Workflow Description

1.2.9 Expert App, TrackAsia Nearby Support & Realtime Contact Workflow

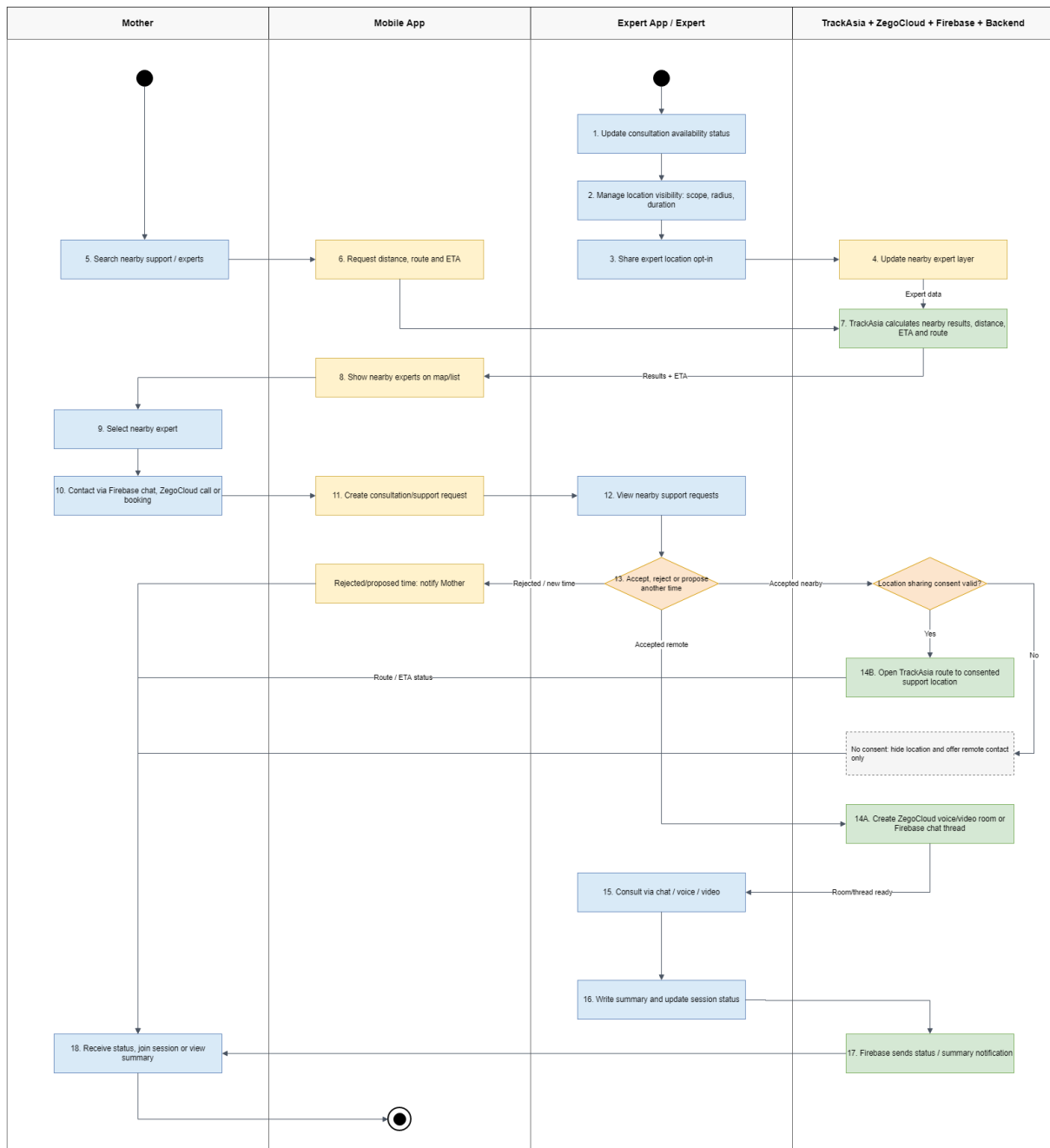


Figure 11: Expert App, TrackAsia Nearby Support & Realtime Contact Workflow

Image Detail: [Expert App, TrackAsia Nearby Support & Realtime Contact Workflow](#)

Description

#	Step Name	Detail Description	Role	Note
1	Update consultation availability status	Expert updates current consultation availability.	Expert App / Expert	-
2	Manage location visibility: scope, radius, duration	Expert controls location visibility scope, radius and duration.	Expert App / Expert	-
3	Share expert location opt-in	Expert opts in to sharing nearby support location.	Expert App / Expert	-
4	Update nearby expert layer	Backend/TrackAsia updates nearby expert map layer.	TrackAsia + ZegoCloud + Firebase + Backend	-
5	Search nearby support / experts	Mother searches nearby support or experts.	Mother	-
6	Request distance, route and ETA	Mobile app requests map distance, route and ETA.	Mobile App	-
7	TrackAsia calculates nearby results, distance, ETA and route	TrackAsia calculates nearby results and route metadata.	TrackAsia + ZegoCloud + Firebase + Backend	-
8	Show nearby experts on map/list	Mobile app displays nearby experts on map/list.	Mobile App	-
9	Select nearby expert	Mother selects a nearby expert.	Mother	-
10	Contact via Firebase chat, ZegoCloud call or booking	Mother chooses chat, voice/video call or booking.	Mother	-
11	Create consultation/support request	Mobile app creates consultation or support request.	Mobile App	-
12	View nearby support requests	Expert views incoming nearby support requests.	Expert App / Expert	-
13	Accept, reject or propose another time	Expert accepts, rejects or proposes another time.	Expert App / Expert	Accepted remote / rejected / new time.
D-1	Location sharing consent valid?	Valid consent opens TrackAsia route; no consent hides location and offers remote contact only.	TrackAsia + ZegoCloud + Firebase + Backend	Yes / No.
14A	Create ZegoCloud voice/video room or Firebase chat thread	Backend creates ZegoCloud room or Firebase chat thread for remote contact.	TrackAsia + ZegoCloud + Firebase + Backend	Accepted remote.

#	Step Name	Detail Description	Role	Note
14B	Open TrackAsia route to consented support location	TrackAsia opens route only to consented support location.	TrackAsia + ZegoCloud + Firebase + Backend	Accepted nearby with consent.
15	Consult via chat / voice / video	Expert and Mother conduct the remote consultation.	Expert App / Expert	-
16	Write summary and update session status	Expert writes summary and updates session status.	Expert App / Expert	-
17	Firebase sends status / summary notification	Firebase sends status or summary notification.	TrackAsia + ZegoCloud + Firebase + Backend	-
18	Receive status, join session or view summary	Mother receives status, joins session or views summary.	Mother	-
EX-1	Rejected/proposed time	Mother is notified when expert rejects or proposes another time.	Mobile App	-
EX-2	No consent	Location is hidden and remote contact is offered only.	TrackAsia + ZegoCloud + Firebase + Backend	-

Table 12: Expert App, TrackAsia Nearby Support & Realtime Contact Workflow Description

1.2.10 Pregnancy Exercise & MediaPipe Posture Workflow

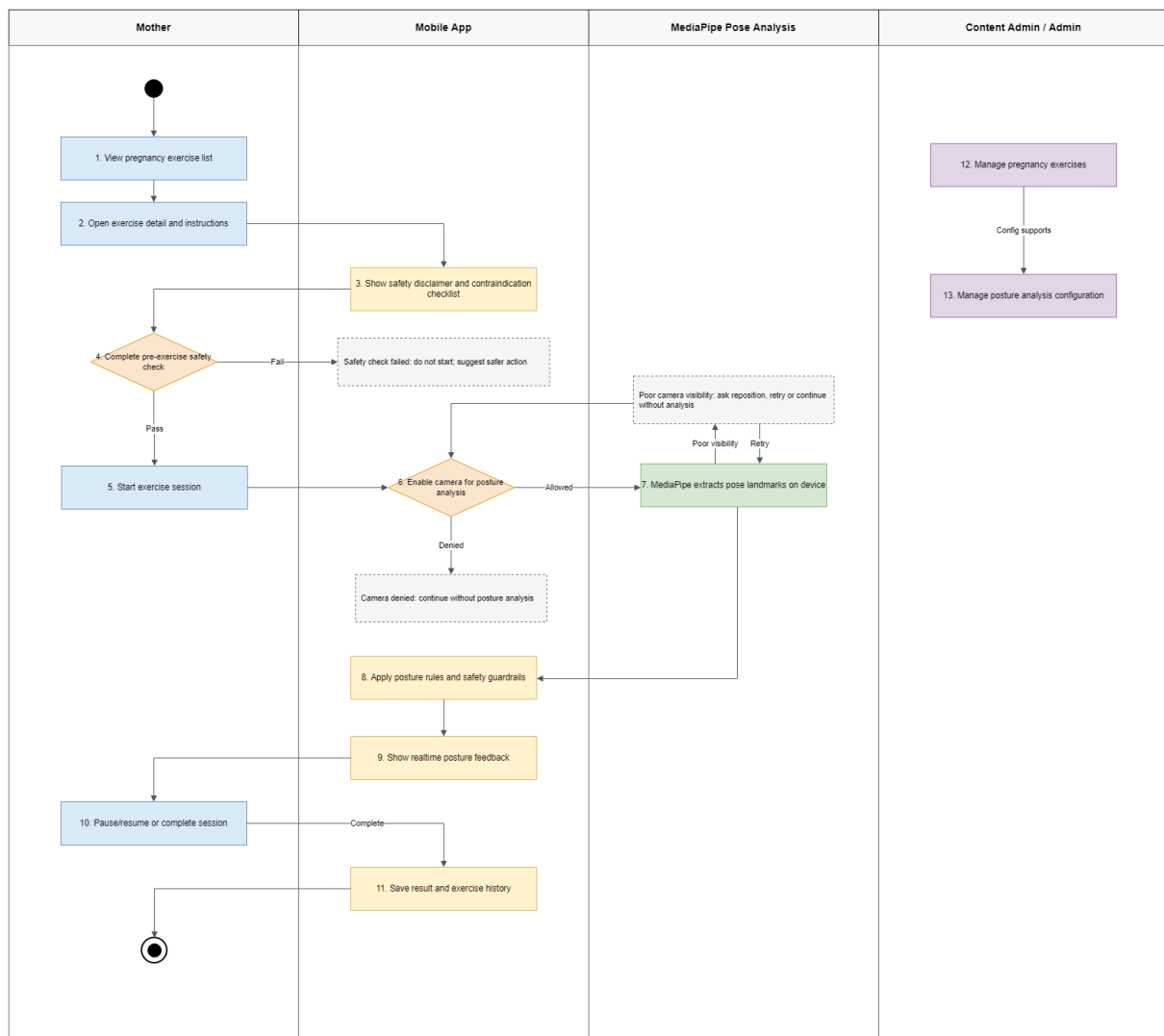


Figure 12: Pregnancy Exercise & MediaPipe Posture Workflow

Image Detail: [Pregnancy Exercise & MediaPipe Posture Workflow](#)

Description

#	Step Name	Detail Description	Role	Note
1	View pregnancy exercise list	Mother opens pregnancy exercise list.	Mother	-
2	Open exercise detail and instructions	Mother views exercise details and instructions.	Mother	-
3	Show safety disclaimer and contraindication checklist	Mobile app displays safety disclaimer and contraindication checklist.	Mobile App	-
4	Complete pre-exercise safety check	Mother completes pre-exercise safety check.	Mother	Pass / fail.

#	Step Name	Detail Description	Role	Note
5	Start exercise session	Mother starts exercise session after safety pass.	Mother	-
6	Enable camera for posture analysis	Mobile app requests/enables camera for posture analysis.	Mobile App	Allowed / denied.
7	MediaPipe extracts pose landmarks on device	MediaPipe extracts pose landmarks on device.	MediaPipe Pose Analysis	-
8	Apply posture rules and safety guardrails	Mobile app applies posture rules and safety guardrails.	Mobile App	-
9	Show realtime posture feedback	Mobile app shows realtime posture feedback.	Mobile App	-
10	Pause/resume or complete session	Mother pauses, resumes or completes exercise session.	Mother	-
11	Save result and exercise history	Mobile app saves result and exercise history.	Mobile App	-
12	Manage pregnancy exercises	Content Admin/Admin manages pregnancy exercise content.	Content Admin / Admin	-
13	Manage posture analysis configuration	Content Admin/Admin manages posture analysis configuration.	Content Admin / Admin	Config supports posture flow.
EX-1	Safety check failed	Session does not start; app suggests safer action.	Mobile App	Fail branch.
EX-2	Camera denied	User continues without posture analysis.	Mobile App	Denied branch.
EX-3	Poor camera visibility	App asks user to reposition, retry or continue without analysis.	MediaPipe Pose Analysis	Retry / continue without analysis.

Table 13: Pregnancy Exercise & MediaPipe Posture Workflow Description

1.2.11 Smartwatch Data Sync & Health Trend Workflow

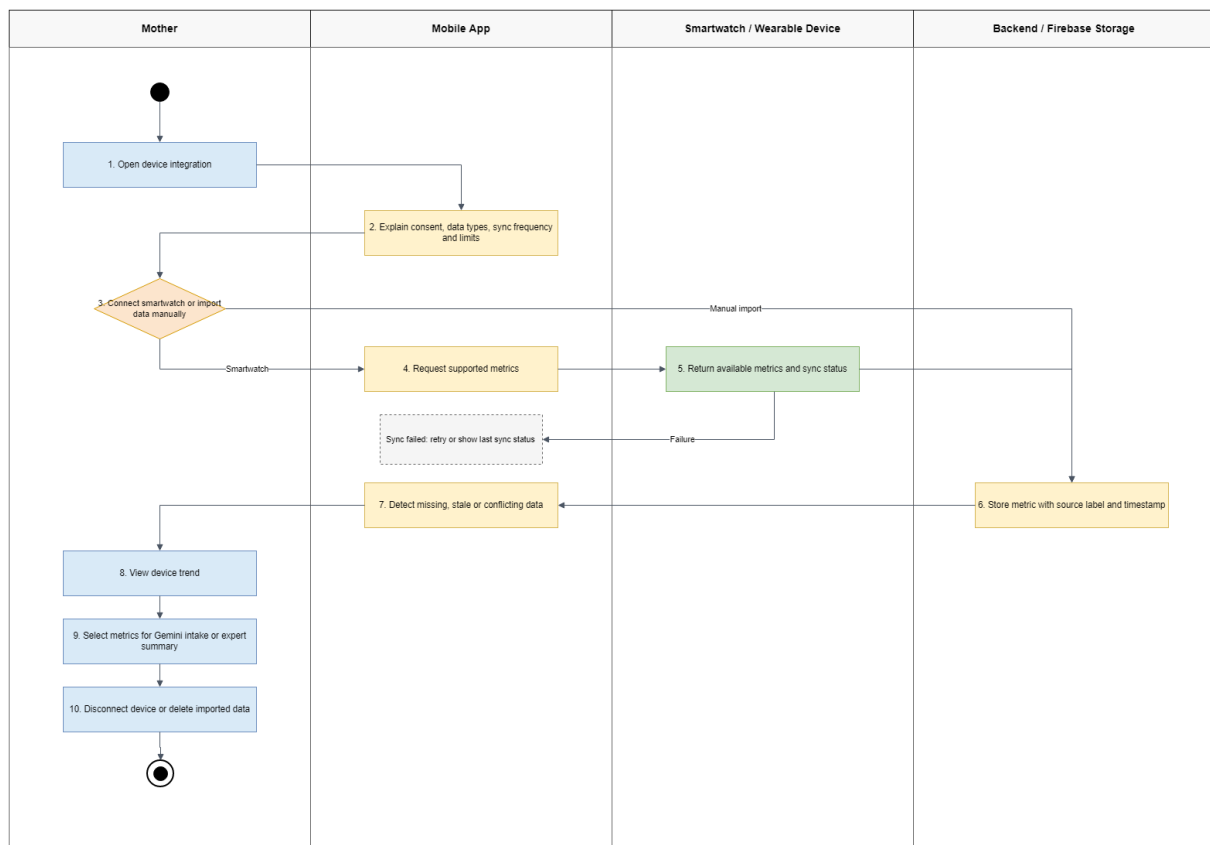


Figure 13: Smartwatch Data Sync & Health Trend Workflow

Image Detail: [Smartwatch Data Sync & Health Trend Workflow](#)

Description

#	Step Name	Detail Description	Role	Note
1	Open device integration	Mother opens device integration.	Mother	-
2	Explain consent, data types, sync frequency and limits	Mobile app explains consent, data types, sync frequency and limitations.	Mobile App	-
3	Connect smartwatch or import data manually	Mother connects smartwatch or imports data manually.	Mother	Smartwatch / manual import.
4	Request supported metrics	Mobile app requests supported metrics.	Mobile App	-
5	Return available metrics and sync status	Smartwatch/wearable returns metrics and sync status.	Smartwatch / Wearable Device	-
6	Store metric with source label and timestamp	Backend/Firebase Storage stores metric with source label and timestamp.	Backend / Firebase Storage	-

#	Step Name	Detail Description	Role	Note
7	Detect missing, stale or conflicting data	Mobile app detects missing, stale or conflicting data.	Mobile App	Sync failure path may retry/show last sync.
8	View device trend	Mother views device trend.	Mother	-
9	Select metrics for Gemini intake or expert summary	Mother selects metrics for Gemini intake or expert summary.	Mother	-
10	Disconnect device or delete imported data	Mother disconnects device or deletes imported data.	Mother	-
EX-1	Sync failed	Mobile app retries sync or shows last sync status.	Mobile App	Failure branch.

Table 14: Smartwatch Data Sync & Health Trend Workflow Description

1.2.12 Partner Clinic & Sponsored Content Governance Workflow

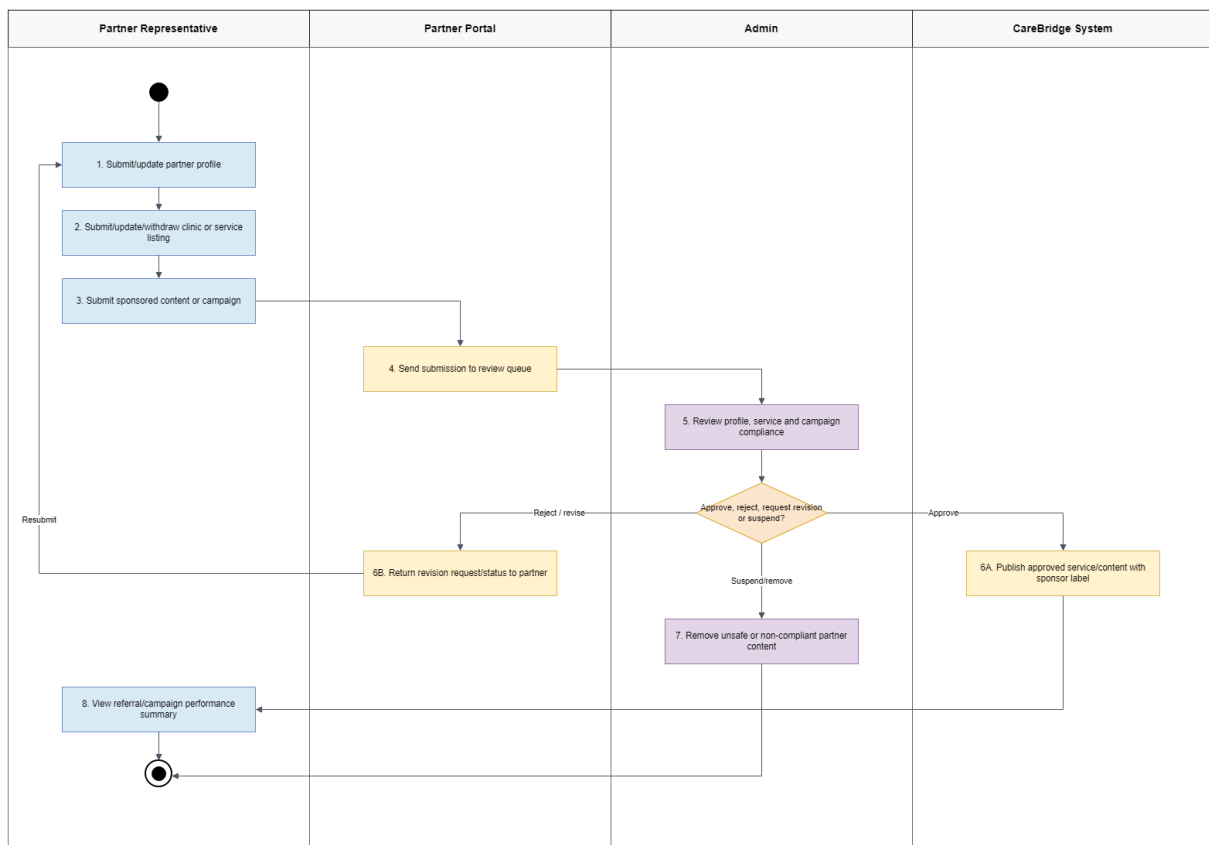


Figure 14: Partner Clinic & Sponsored Content Governance Workflow

Image Detail: [Partner Clinic & Sponsored Content Governance Workflow](#)

Description

#	Step Name	Detail Description	Role	Note
1	Submit/update partner profile	Partner representative submits or updates partner profile.	Partner Representative	-
2	Submit/update/withdraw clinic or service listing	Partner representative manages clinic/service listing submission.	Partner Representative	-
3	Submit sponsored content or campaign	Partner representative submits sponsored content or campaign.	Partner Representative	-
4	Send submission to review queue	Partner Portal sends submission to review queue.	Partner Portal	-
5	Review profile, service and campaign compliance	Admin reviews partner profile, service and campaign compliance.	Admin	-
D-1	Approve, reject, request revision or suspend?	Admin decision publishes, returns revision or removes/suspends content.	Admin	Approve / reject-revise / suspend-remove.
6A	Publish approved service/content with sponsor label	CareBridge System publishes approved service/content with sponsor label.	CareBridge System	Approve branch.
6B	Return revision request/status to partner	Partner Portal returns revision request or status to partner.	Partner Portal	Reject/revise branch; partner can resubmit.
7	Remove unsafe or non-compliant partner content	Admin removes unsafe or non-compliant partner content.	Admin	Suspend/remove branch.
8	View referral/campaign performance summary	Partner representative views referral/campaign performance summary.	Partner Representative	-

Table 15: Partner Clinic & Sponsored Content Governance Workflow Description

1.2.13 Verified Content, Checklist & Gemini RAG Knowledge Workflow

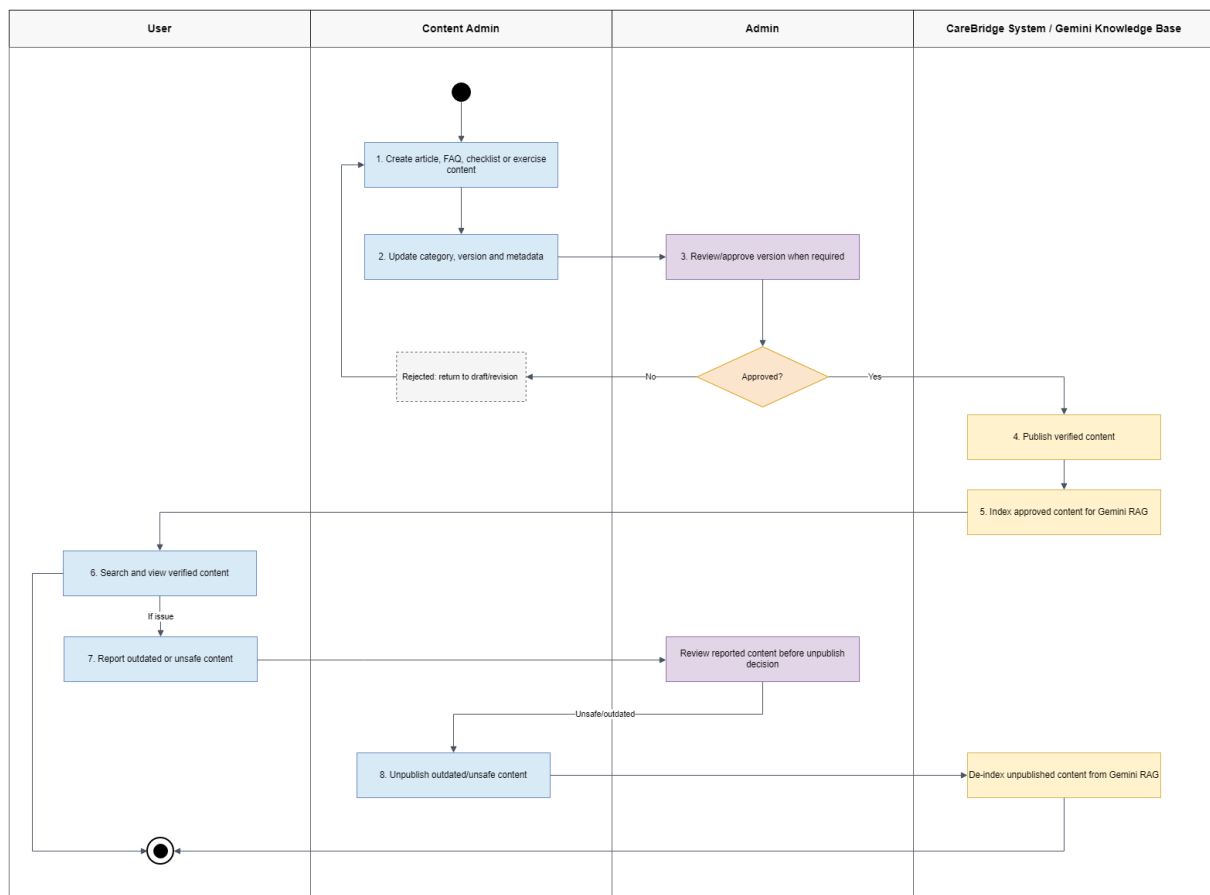


Figure 15: Verified Content, Checklist & Gemini RAG Knowledge Workflow

Image Detail: [Verified Content, Checklist & Gemini RAG Knowledge Workflow](#)

Description

#	Step Name	Detail Description	Role	Note
1	Create article, FAQ, checklist or exercise content	Content Admin creates verified content, FAQ, checklist or exercise content.	Content Admin	-
2	Update category, version and metadata	Content Admin updates category, version and metadata.	Content Admin	-
3	Review/approve version when required	Admin reviews and approves version when required.	Admin	-
D-1	Approved?	Approved content is published and indexed; rejected content returns to draft/revision.	Admin	Yes / No.

#	Step Name	Detail Description	Role	Note
4	Publish verified content	CareBridge System publishes verified content.	CareBridge System / Gemini Knowledge Base	-
5	Index approved content for Gemini RAG	CareBridge System indexes approved content for Gemini RAG.	CareBridge System / Gemini Knowledge Base	-
6	Search and view verified content	User searches and views verified content.	User	-
7	Report outdated or unsafe content	User reports outdated or unsafe content.	User	-
EX-1	Review reported content before unpublish decision	Admin reviews reported content before unpublish decision.	Admin	If issue.
8	Unpublish outdated/unsafe content	Content Admin unpublishes outdated or unsafe content.	Content Admin	Unsafe/outdated.
EX-2	Rejected: return to draft/revision	Rejected content returns to draft or revision state.	Content Admin	No branch.
EX-3	De-index unpublished content from Gemini RAG	CareBridge System removes unpublished content from Gemini RAG index.	CareBridge System / Gemini Knowledge Base	After unpublish.

Table 16: Verified Content, Checklist & Gemini RAG Knowledge Workflow Description

1.2.14 Audit, Security Incident & Access Review Workflow

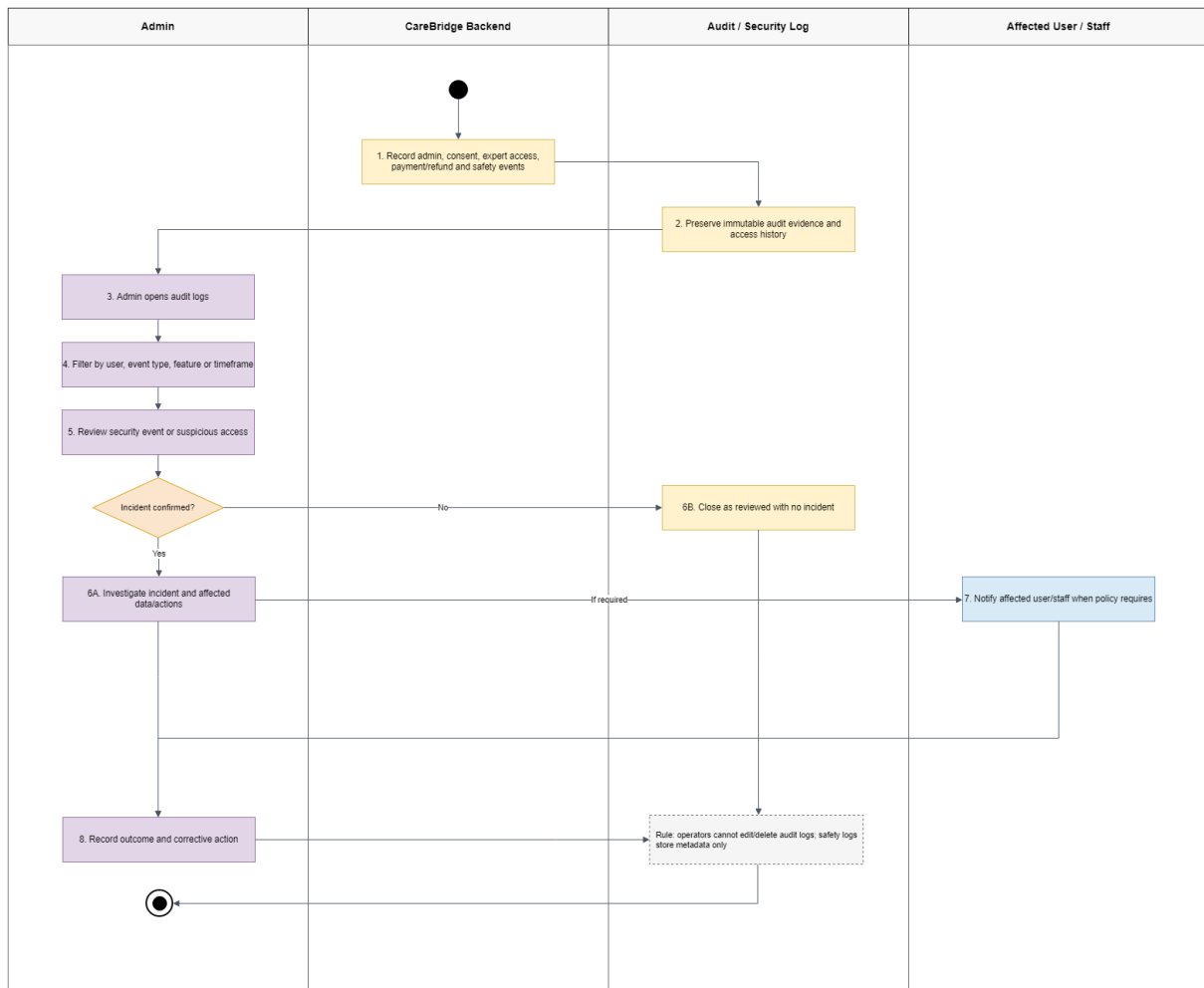


Figure 16: Audit, Security Incident & Access Review Workflow

Image Detail: [Audit, Security Incident & Access Review Workflow](#)

Description

#	Step Name	Detail Description	Role	Note
1	Record admin, consent, expert access, payment/refund and safety events	Backend records admin, consent, expert access, payment/refund and safety events.	CareBridge Backend	-
2	Preserve immutable audit evidence and access history	Audit/Security Log preserves immutable evidence and access history.	Audit / Security Log	-
3	Admin opens audit logs	Admin opens audit logs for review.	Admin	-
4	Filter by user, event type, feature or timeframe	Admin filters logs by user, event type, feature or timeframe.	Admin	-

#	Step Name	Detail Description	Role	Note
5	Review security event or suspicious access	Admin reviews security event or suspicious access.	Admin	-
D-1	Incident confirmed?	Confirmed incidents move to investigation and notification; non-incidents close as reviewed.	Admin	Yes / No.
6A	Investigate incident and affected data/actions	Admin investigates incident and affected data/actions.	Admin	-
6B	Close as reviewed with no incident	Audit/Security Log closes event as reviewed with no incident.	Audit / Security Log	No branch.
7	Notify affected user/staff when policy requires	Affected user or staff is notified when policy requires.	Affected User / Staff	If required.
8	Record outcome and corrective action	Admin records outcome and corrective action.	Admin	-
EX-1	Audit log immutability rule	Operators cannot edit/delete audit logs; safety logs store metadata only.	Audit / Security Log	Rule in diagram.

Table 17: Audit, Security Incident & Access Review Workflow Description

2. User Requirements

2.1 Actors

#	Actor	Description
1	Guest	A non-authenticated visitor who can access public content and start account registration. This actor has limited access and cannot use personalized functions until login or account verification is completed.
2	User	A general parent actor for all authenticated accounts in the CareBridge system. This actor shares common functions such as authentication, profile management, notification settings, search, reporting, and permission-based access.
3	Mother	The main end-user of the mobile application. This actor manages pregnancy, postpartum, mother and baby care journeys, asks community questions, stores personal health records, books consultations, receives reminders, and accesses emergency or safety support features.
4	Family Member	A family member or secondary caregiver invited to join a care group. This actor can view shared data, receive assigned tasks, get alerts, and support the care seeker according to granted permissions.
5	Verified Expert	A verified healthcare-related expert such as a doctor, nutritionist, lactation consultant, psychologist, or child development specialist.

#	Actor	Description
		This actor can manage an expert profile, answer questions, provide consultation support, and interact with users within professional and safety boundaries.
6	Moderator	A community moderation actor responsible for reviewing posts, comments, answers, reports, and community violations. This actor helps maintain safe, respectful, and reliable community interactions.
7	Content Admin	An administrative actor responsible for creating, editing, reviewing, and versioning educational content such as articles, FAQs, checklists, and verified guidance materials.
8	Admin	The system administrator who manages users, roles, expert verification, dashboards, partner accounts, safety rules, configuration, and overall platform operation.
9	Partner Representative	A representative of a partner clinic, healthcare facility, service provider, or sponsor. This actor submits partner profiles, manages service listings, and handles sponsored campaigns through the Partner Portal.
10	Payment Gateway	An external payment service that processes paid consultation transactions and returns payment status or transaction results to the CareBridge system.
11	TrackAsia Map Service	An external map service that supports geocoding, nearby healthcare place search, distance calculation, route guidance, and estimated travel time for emergency map and nearby care support.
12	Smartwatch / Wearable Device	An external wearable or smartwatch data source that can provide supported health metrics only when the user grants consent. Imported data must keep source labels and timestamps.
13	Gemini AI Service	The AI service used by CareBridge for intake structuring, normalization, summarization, FAQs, and RAG-supported responses; it does not make independent medical decisions.
14	Firebase Cloud Messaging	The notification service used to deliver push and in-app alerts for reminders, care tasks, consultation updates, moderation results, and emergency or safety events.
15	Phone Motion / IMU Sensor	The accelerometer and gyroscope on the phone provide motion data for safety monitoring. Raw data should preferably be processed on device and must not be considered medical-device data.
16	ZegoCloud + Firebase Realtime / Chat	Realtime communication services used for authenticated consultation sessions: ZegoCloud supports voice/video calls, while Firebase Realtime / Chat supports chat threads.
17	Firebase Storage	The file storage service used for uploaded health records, images, and supporting files. Access to stored health files is governed by CareBridge permissions and consent rules.
18	MediaPipe Pose Analysis	The on-device posture analysis capability used for pregnancy exercise posture support. It assists pose landmark analysis and must not be treated as medical diagnosis or clinical assessment.

Table 18: Actor Description

2.2 Use Cases

2.2.1 Diagram(s)

2.2.1.1 Use Cases for User

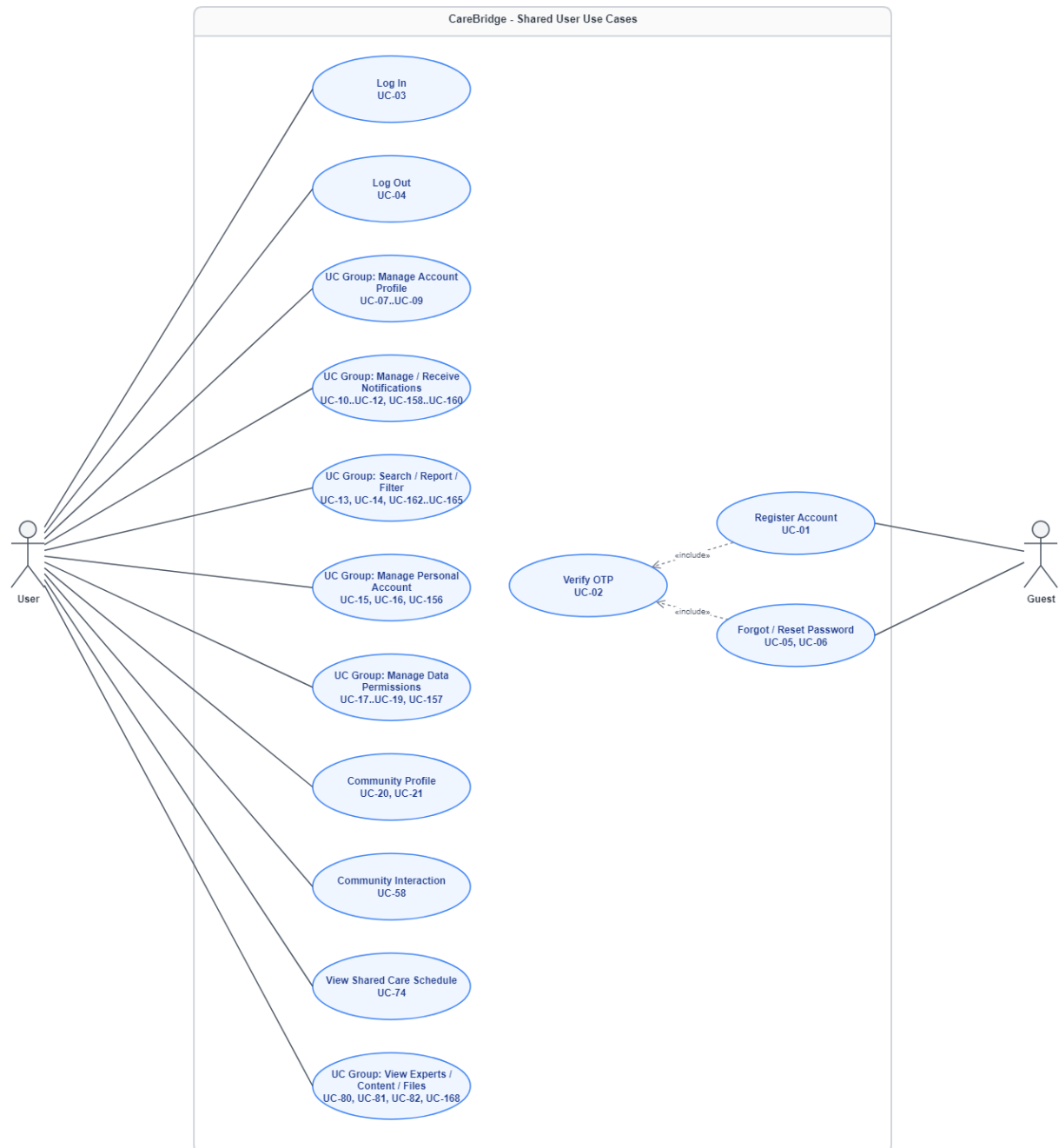


Figure 17: Shared User Use Cases

Image Detail: [Shared User Use Cases](#)

2.2.1.2 Use Cases for Mother

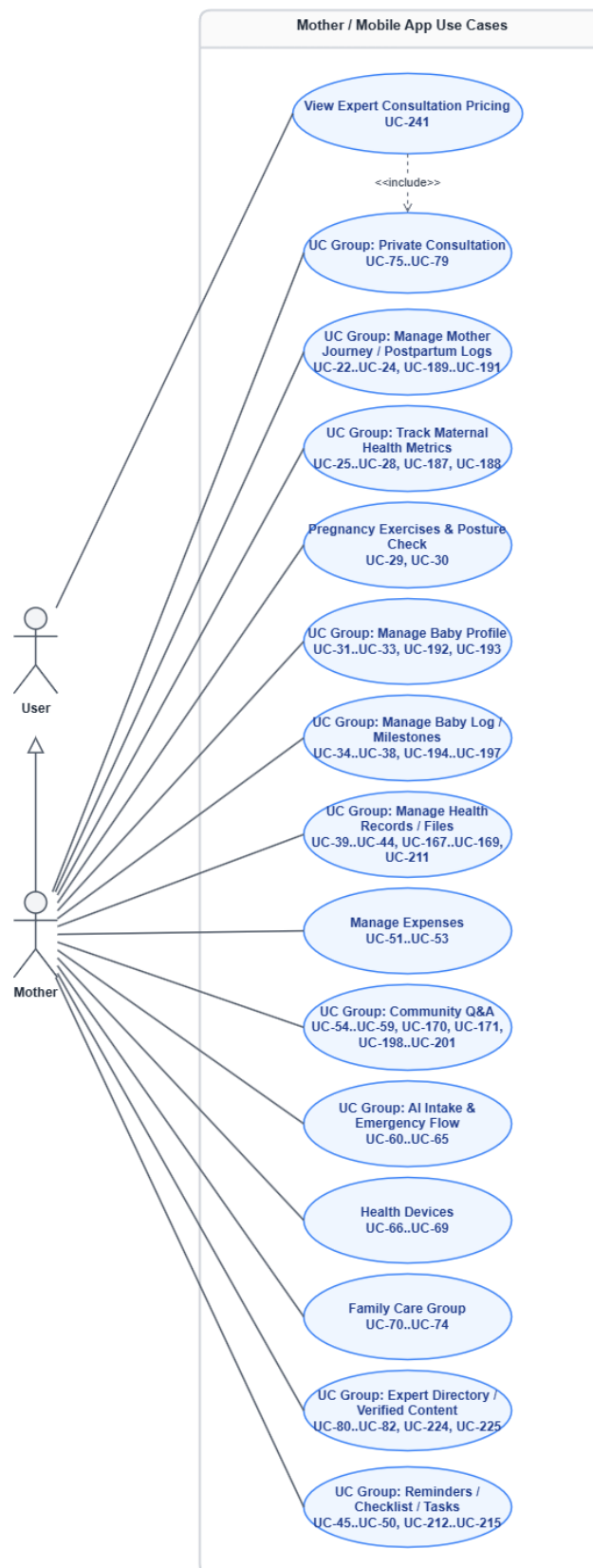


Figure 18: Mother Use Cases

Image Detail: [Mother Use Cases](#)

2.2.1.3 Use Cases for Family Member

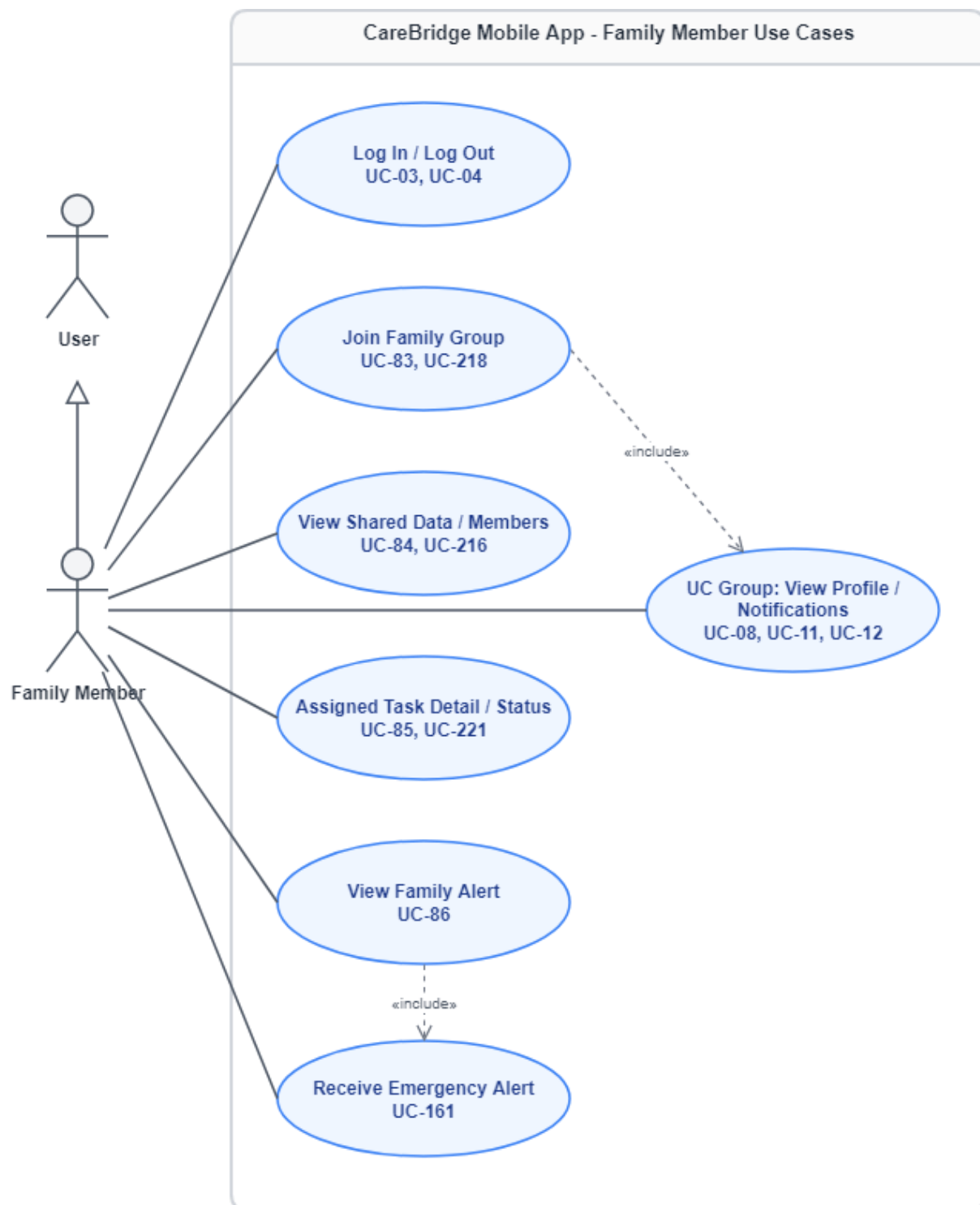


Figure 19: Family Member Use Cases

Image Detail: [Family Member Use Cases](#)

2.2.1.4 Use Cases for Expert

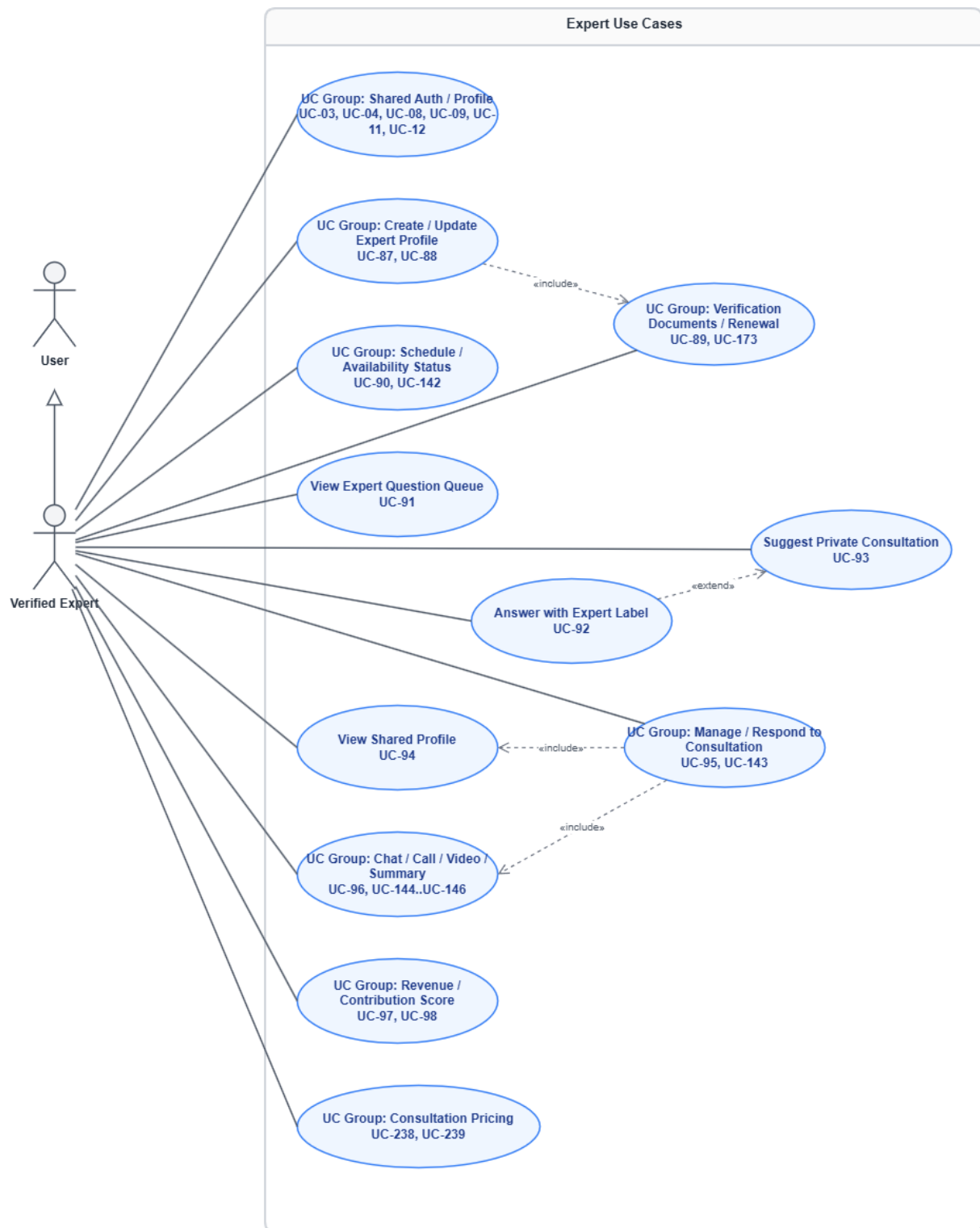


Figure 20: Expert Use Cases

Image Detail: [Expert Use Cases](#)

2.2.1.5 Use Cases for Content, Moderation & Operations

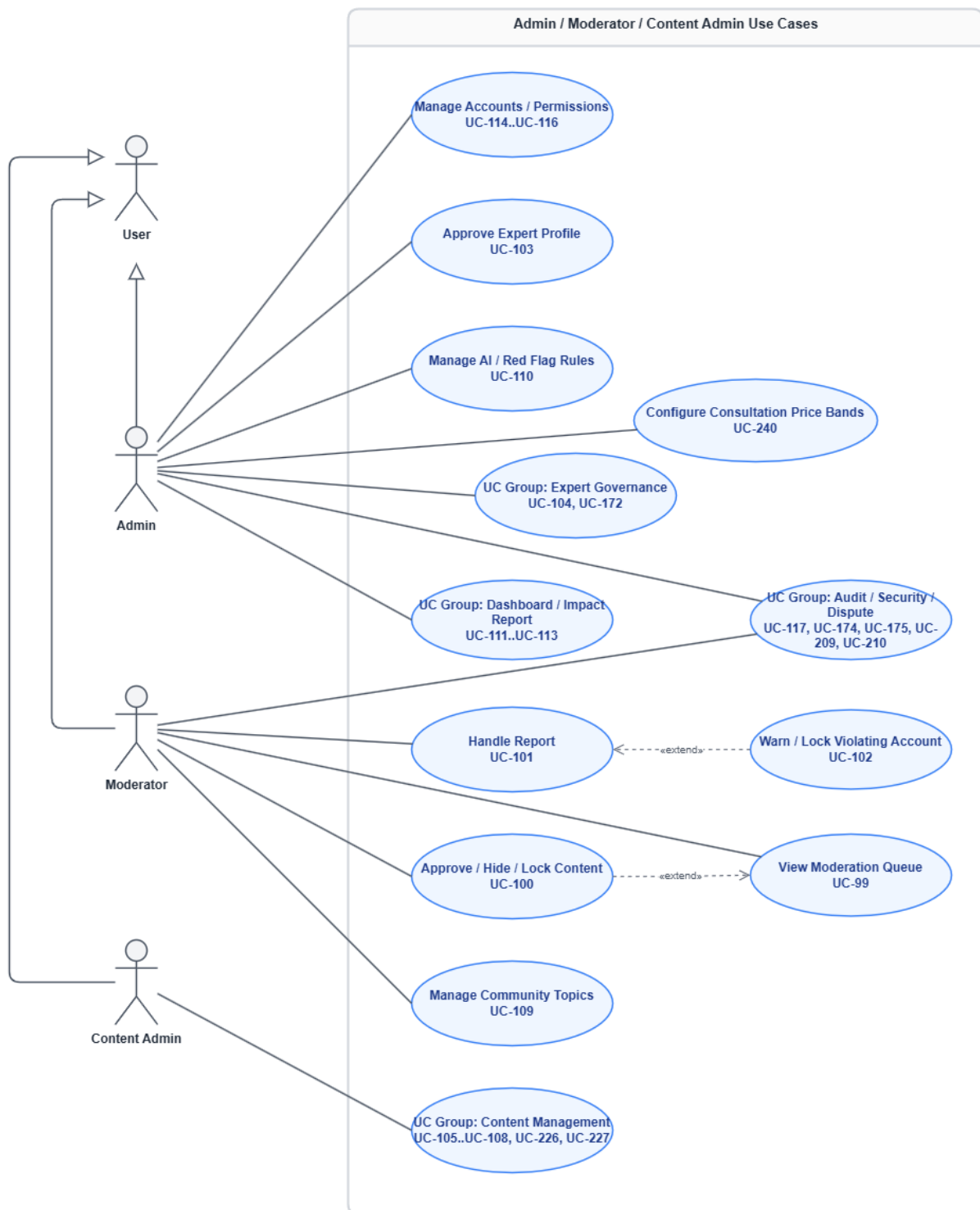


Figure 21: Content, Moderation & Operations Use Cases

Image Detail: [Content, Moderation & Operations Use Cases](#)

2.2.1.6 Use Cases for Partner

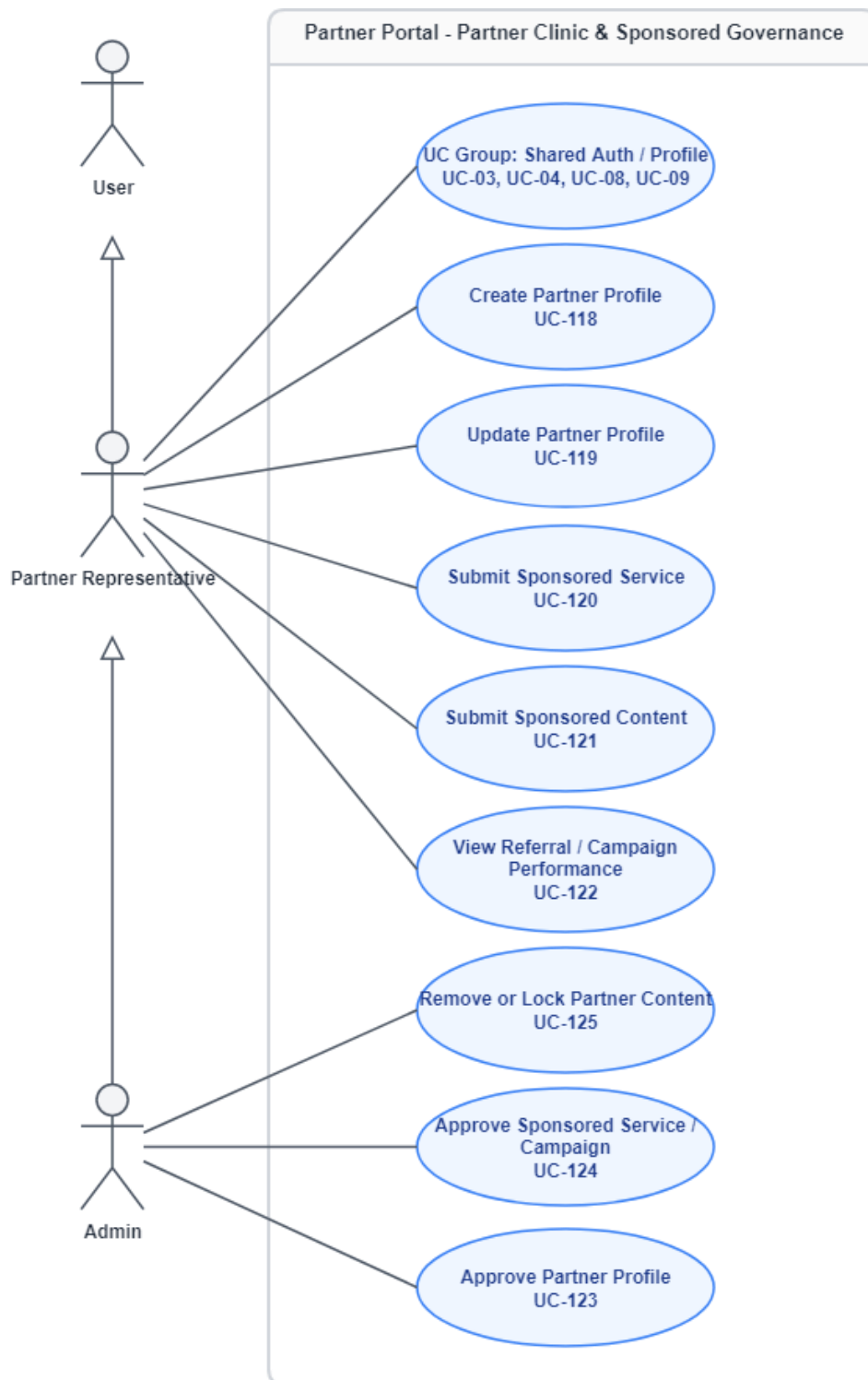


Figure 22: Partner Use Cases

Image Detail: [Partner Use Cases](#)

2.2.1.7 Use Cases for Supporting Services

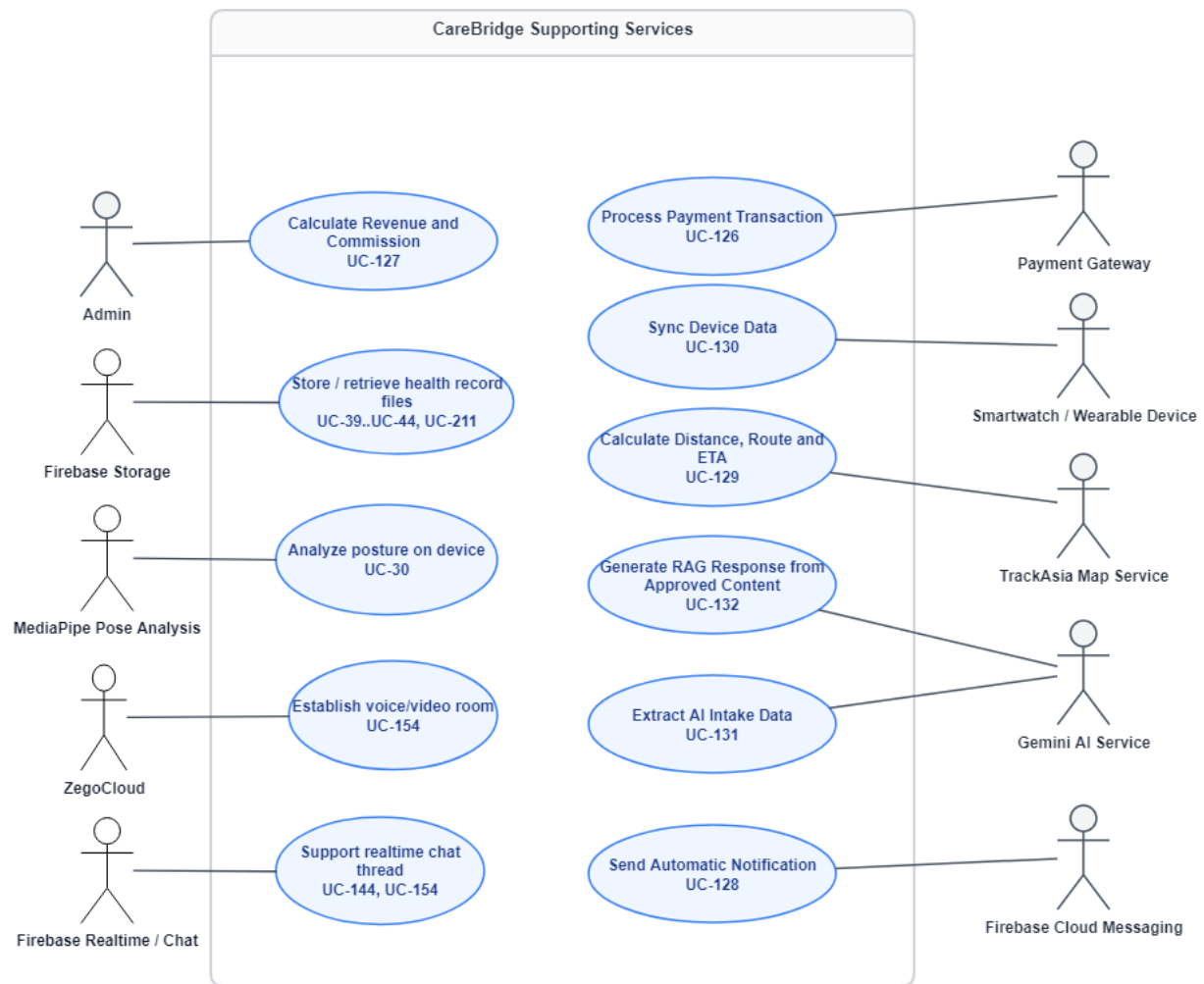


Figure 23: Supporting Services Use Cases

Image Detail: [Supporting Services Use Cases](#)

2.2.1.8 Use Cases for Smart Activity Monitoring & Safety Support

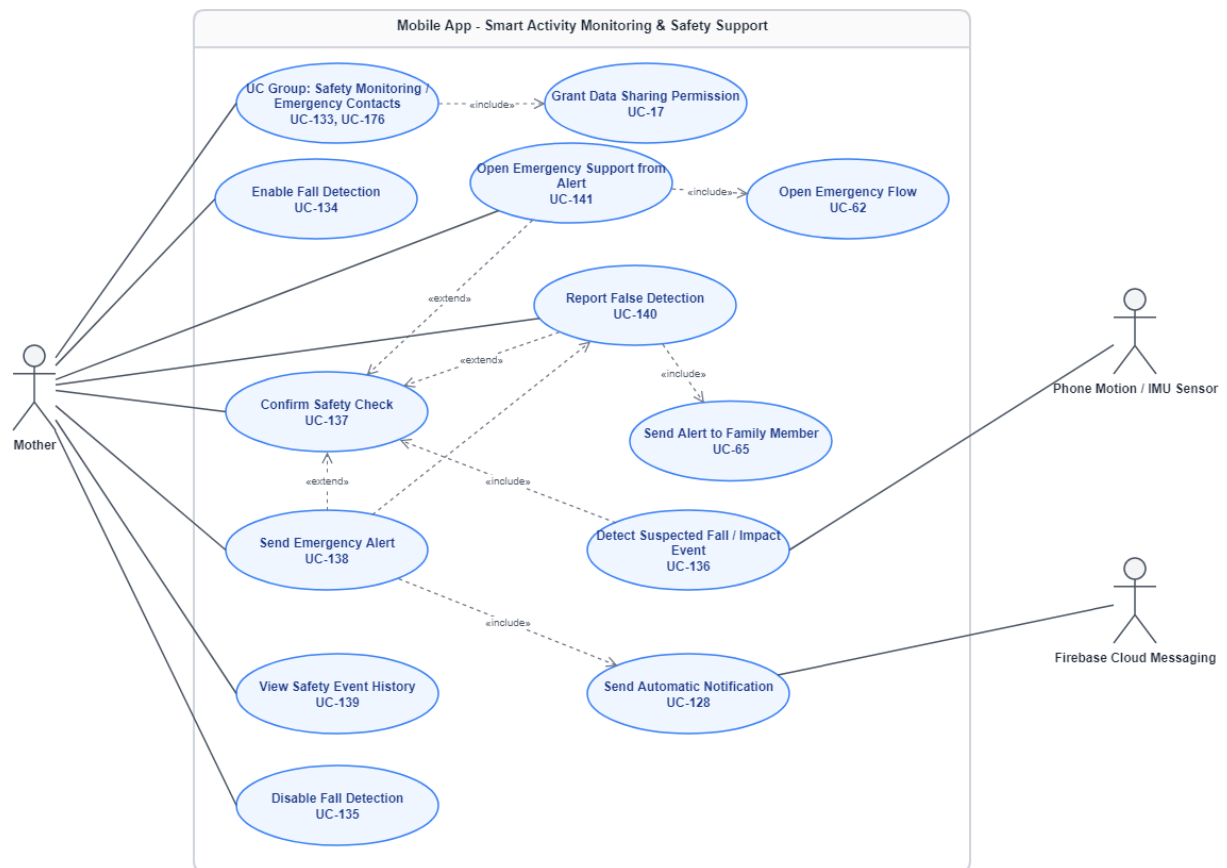


Figure 24: Smart Activity Monitoring & Safety Support Cases

Image Detail: [Smart Activity Monitoring & Safety Support Cases](#)

2.2.1.9 Use Cases for Community Q&A Extended Use Cases



Figure 25: Community Q&A Extended Use Cases

Image Detail: [Community Q&A Extended Use Cases](#)

2.2.1.10 Use Cases for Map, Expert App & Nearby Support

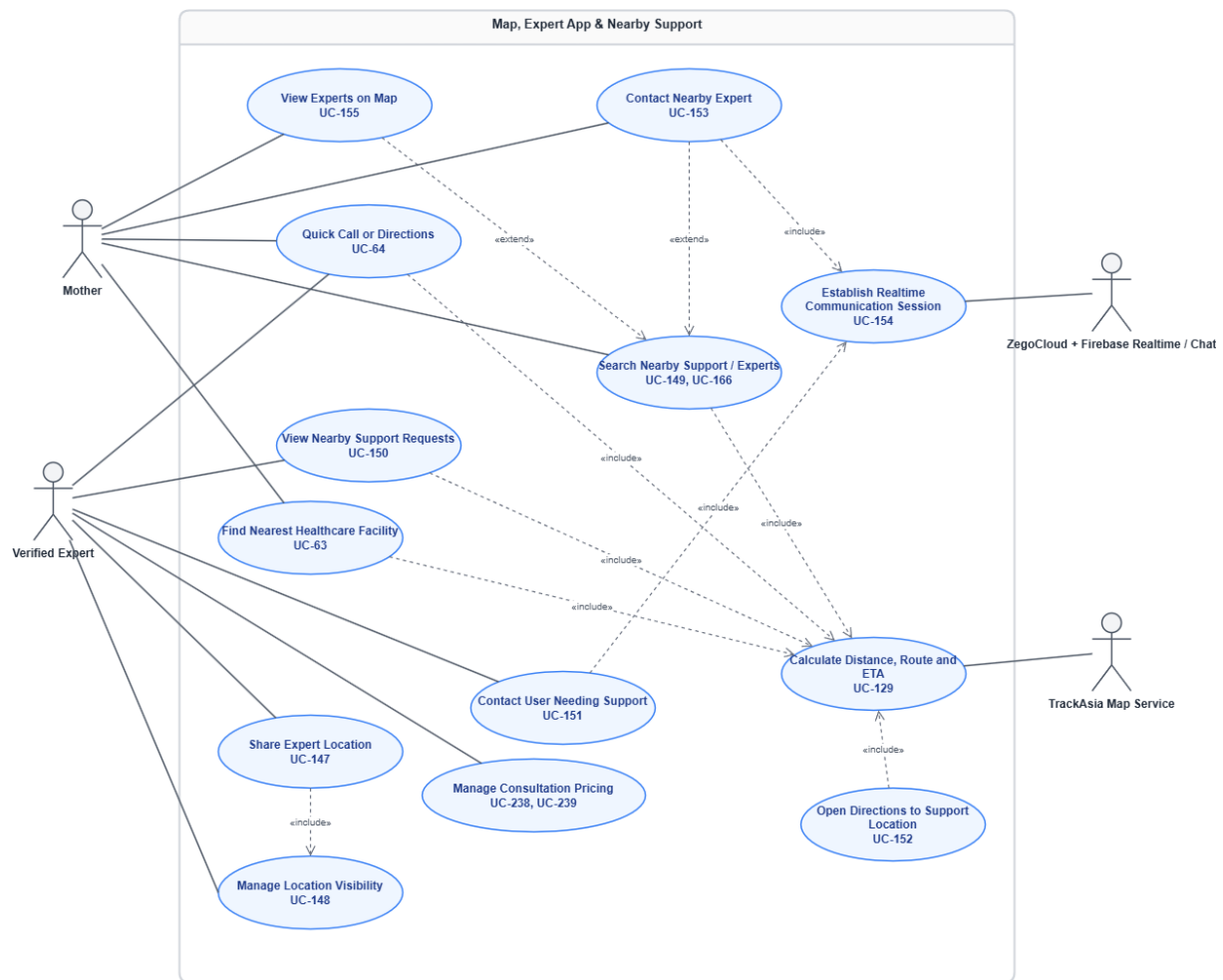


Figure 26: Map, Expert App & Nearby Support Use Cases

Image Detail: [Map, Expert App & Nearby Support Use Cases](#)

2.2.1.11 Use Cases for Pregnancy Exercise & MediaPipe Posture

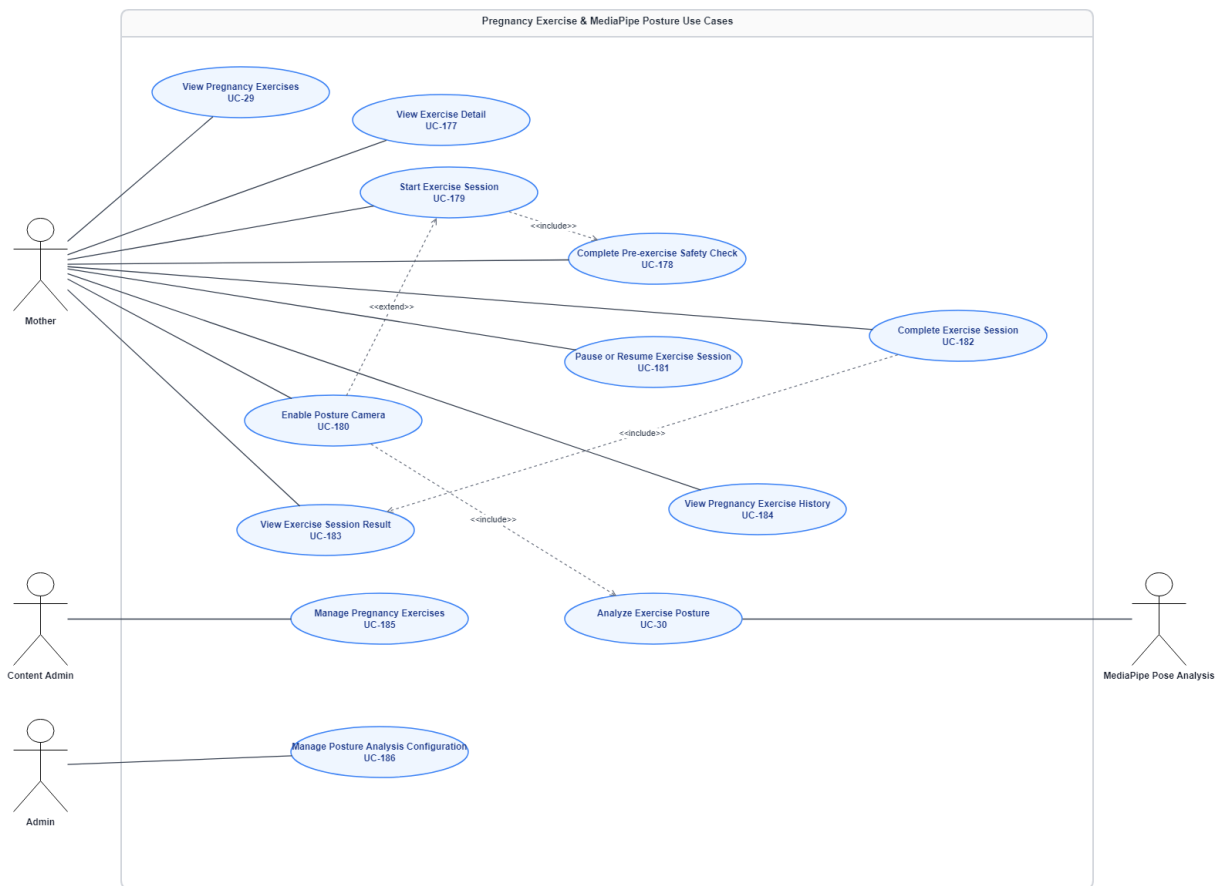


Figure 27: Pregnancy Exercise & MediaPipe Posture Use Cases

Image Detail: [Pregnancy Exercise & MediaPipe Posture Use Cases](#)

2.2.1.12 Use Cases for Consultation Lifecycle, Dispute & Refund

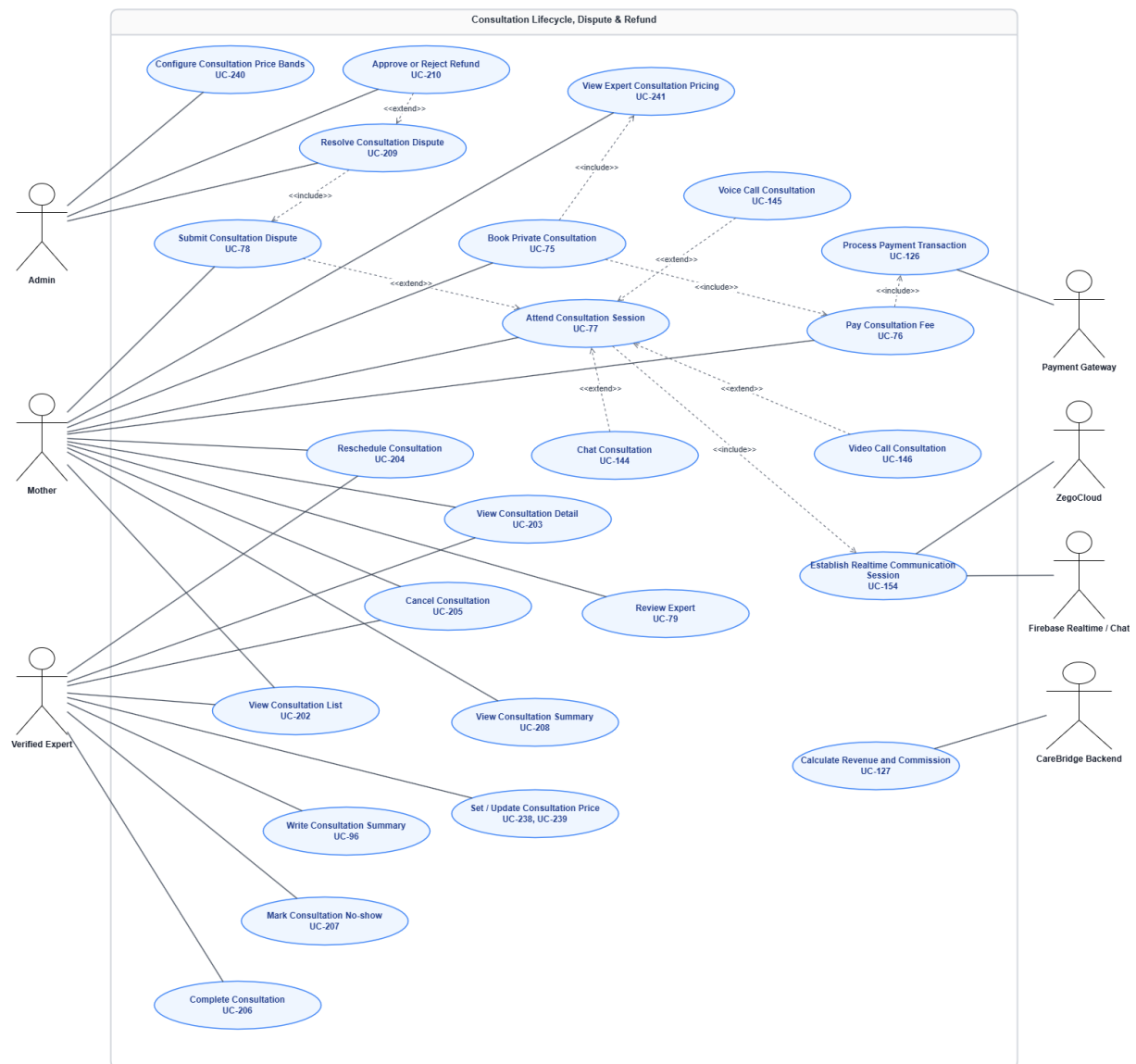


Figure 28: Consultation Lifecycle, Dispute & Refund Use Cases

Image Detail: [Consultation Lifecycle, Dispute & Refund Use Cases](#)

2.2.1.13 Use Cases for Family Sync & Care Task Management

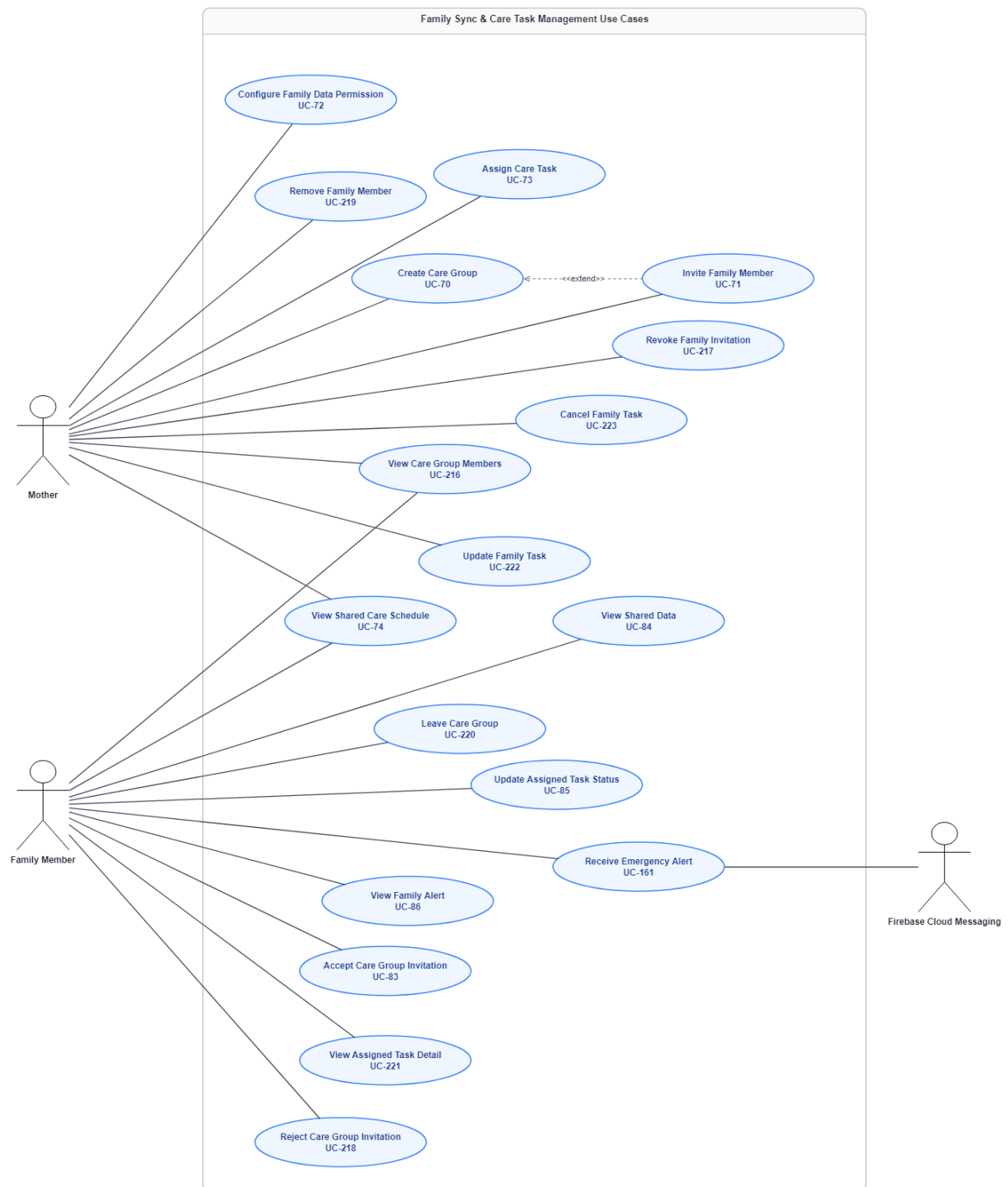


Figure 29: Family Sync & Care Task Management Use Cases

Image Detail: [Family Sync & Care Task Management Use Cases](#)

2.2.1.14 Use Cases for Verified Content, Checklist & RAG Knowledge

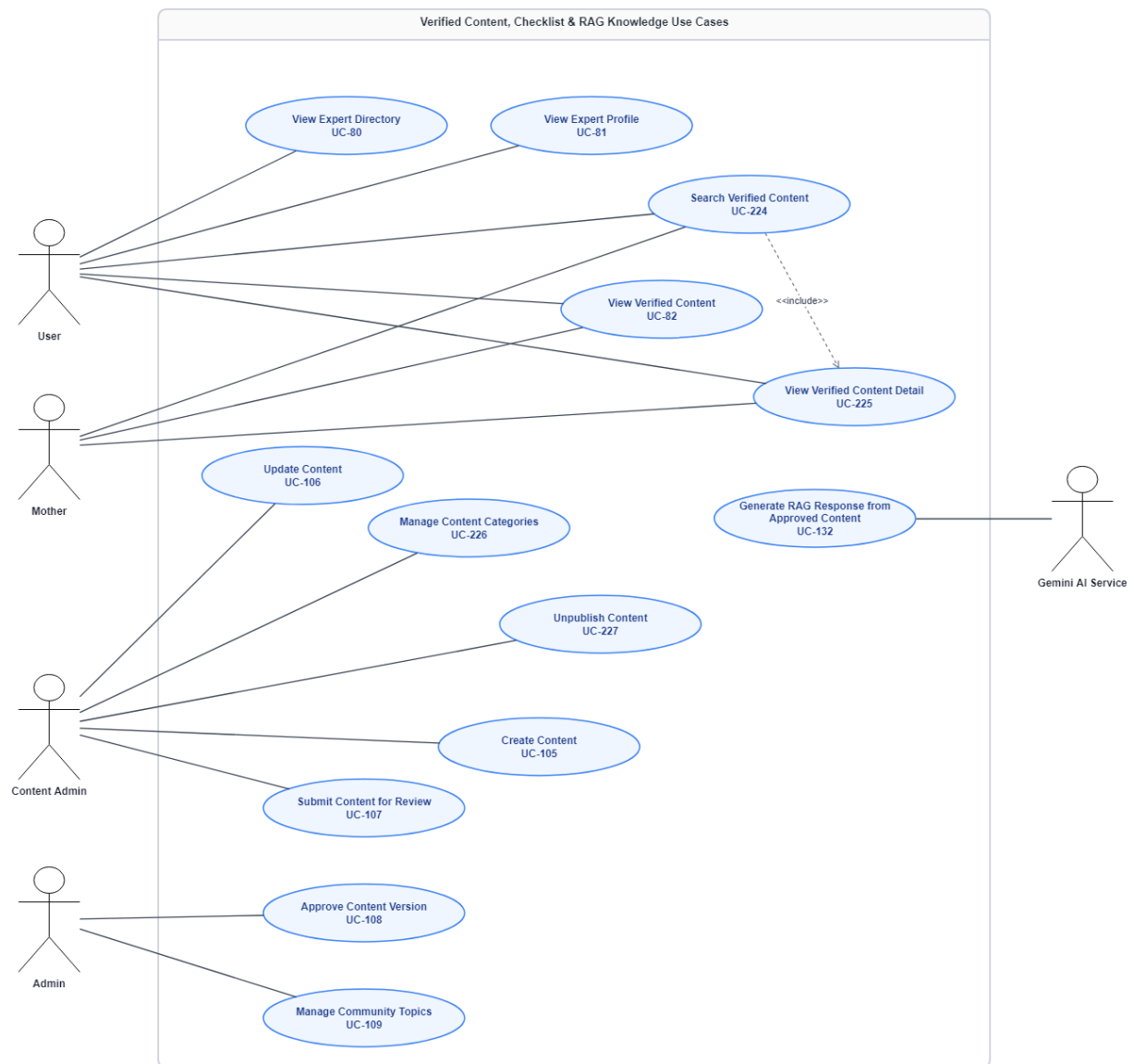


Figure 30: Verified Content, Checklist & RAG Knowledge Use Cases

Image Detail: [Verified Content, Checklist & RAG Knowledge Use Cases](#)

2.2.1.15 Use Cases for Vaccination & Growth Tracking

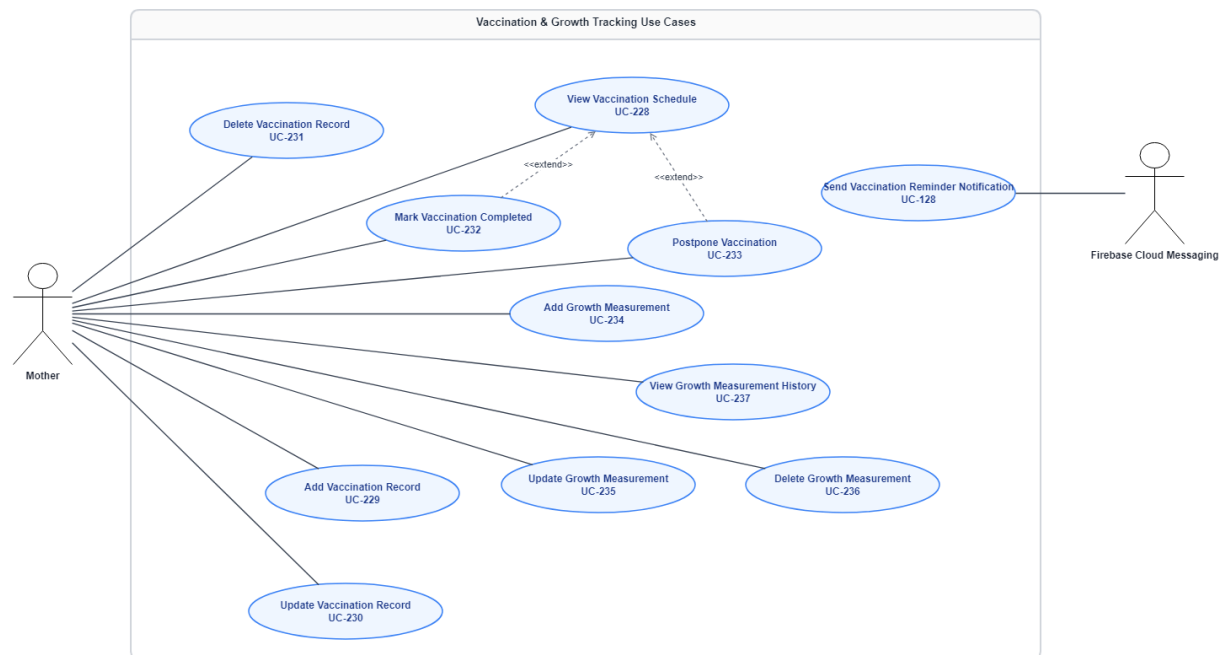


Figure 31: Vaccination & Growth Tracking Use Cases

Image Detail: [Vaccination & Growth Tracking Use Cases](#)

2.2.2 Descriptions

ID	Use Case	Actor	Description
UC-01	Register Account	Guest	Creates an account with a supported contact method and selected initial role. The new account remains unverified until OTP confirmation succeeds.
UC-02	Verify OTP	Guest	Validates the OTP sent for account activation or another sensitive account action.
UC-03	Log In	User	Authenticates a verified account and routes the user to the role-appropriate workspace.
UC-04	Log Out	User	Ends the current session on the active device without deleting the account.
UC-05	Request Password Reset	Guest	Starts password recovery by requesting a time-limited reset proof through the registered channel.
UC-06	Reset Password	Guest	Sets a new password only after a valid recovery proof is verified.

ID	Use Case	Actor	Description
UC-07	Change Password	User	Changes the signed-in user password after current-password validation.
UC-08	View Account Profile	User	Shows the user private account profile, role, status and permitted account controls.
UC-09	Update Account Profile	User	Updates non-sensitive private profile data used for account operation and contact preferences.
UC-10	Manage Community Identity	User	Creates or updates a public community identity that is separate from private maternal, baby and family information.
UC-11	Manage Notification Preferences	User	Lets the user choose supported delivery channels and categories for reminders, community replies, family alerts and account events.
UC-12	View and Mark Notifications	User	Shows authorized notifications and allows the user to mark individual or all eligible notifications as read.
UC-13	Manage Own Login Sessions	User	Lists active or recent own sessions and allows the user to revoke a selected device session.
UC-14	Deactivate or Delete Own Account	User	Requests account deactivation or deletion subject to retention, care-group and audit obligations.
UC-15	Grant Data Permission	User	Creates a purpose-specific, scope-based and time-limited permission for a family member or verified expert to access selected data.
UC-16	Review and Revoke Data Permission	User	Shows active and past sharing permissions and lets the owner revoke an active grant.
UC-17	Administer User Accounts and Role Access	System Admin	Allows authorized administrators to review account status, apply role/access changes and restrict misuse under separation-of-duties controls.
UC-18	Review Sensitive Access and Security Events	System Admin	Reviews sensitive-record access, abnormal login, permission change and file-access events to determine follow-up or investigation.

ID	Use Case	Actor	Description
UC-19	Initialize Mother Care Journey	Mother	Creates a mother journey for preconception, pregnancy or postpartum recovery using the minimum dates and stage context required for stage-based support.
UC-20	Update Mother Journey Stage and Dates	Mother	Updates stage, last menstrual period, expected due date, birth date or other permitted journey dates when circumstances change.
UC-21	View Mother Journey Dashboard	Mother	Shows current stage/week, relevant care priorities, reminders, checklists, reviewed content and safe shortcuts.
UC-22	Add Maternal Health Metric	Mother	Records a user-entered maternal indicator such as weight, blood pressure, glucose, fetal movement note or another supported observation.
UC-23	Update or Delete Maternal Health Metric	Mother	Corrects or removes a user-entered maternal metric when ownership and record state allow.
UC-24	View Maternal Health Trend	Mother	Displays time-based trends for recorded maternal indicators with source labels and clear non-diagnostic context.
UC-25	Add Postpartum Recovery Log	Mother	Records a postpartum recovery observation such as sleep, mood, milk-feeding context, pain note or recovery symptom note.
UC-26	Update or Delete Postpartum Recovery Log	Mother	Corrects or removes an owner-entered postpartum recovery log.
UC-27	Manage Preparation Checklist	Mother	Adds, edits and completes preparation items for pregnancy, hospital bag, postpartum recovery or early baby care.
UC-28	Browse Pregnancy Exercise Library	Mother	Shows reviewed pregnancy exercise content with embedded stage/difficulty/duration filters and exercise-specific safety notes.
UC-29	Complete Pre-exercise Safety Check	Mother	Collects required pre-exercise safety answers and stops the session entry when configured warning answers are present.

ID	Use Case	Actor	Description
UC-30	Conduct Pregnancy Exercise Session with Optional Posture Feedback	Mother	Runs an approved exercise session with start, pause, resume and complete actions and optional rule/ML-based posture feedback after camera consent.
UC-31	View Exercise History and Session Result	Mother	Shows completed exercise sessions, duration, completion state, aggregate posture feedback and safety notices.
UC-32	Create Baby Profile	Mother	Creates a baby profile with core identity and birth context needed for journals, growth and vaccination tracking.
UC-33	Update or Archive Baby Profile	Mother	Updates permitted baby profile fields or archives a profile that is no longer actively managed without destroying linked history.
UC-34	Switch Active Baby Profile	Mother	Selects the baby profile used by the current baby dashboard, journals, growth views and relevant reminders.
UC-35	View Baby Care Overview	Mother	Shows a baby overview with recent journals, growth, milestones, vaccination status and current care prompts.
UC-36	Add Baby Daily Log	Mother	Records feeding, sleep, diaper, symptom note or other supported daily care observation for a selected baby.
UC-37	Update or Delete Baby Daily Log	Mother	Corrects or removes an owner-entered baby journal entry.
UC-38	View Baby Log Summary	Mother	Shows recent 24-hour or 7-day feeding, sleep, diaper and logged observation patterns for the active baby.
UC-39	Record Development Milestone	Mother	Records a caregiver-observed development milestone with date and note.
UC-40	Update or Delete Development Milestone	Mother	Corrects or removes a caregiver-entered development milestone.
UC-41	Add Growth Measurement	Mother	Records supported growth measurements for a selected baby with measurement date and source.
UC-42	Update or Delete Growth Measurement	Mother	Corrects or removes a user-entered growth measurement.

ID	Use Case	Actor	Description
UC-43	View Growth Trend and Measurement History	Mother	Displays growth measurement history and charts with reference context and prompts to seek professional assessment when appropriate.
UC-44	Manage Vaccination Record	Mother	Adds, updates or removes a user-entered vaccination record for the selected baby.
UC-45	View Vaccination Reference Schedule and Reminder Status	Mother	Displays a reference vaccination schedule, completed/recorded status and linked reminder state for the selected baby.
UC-46	Browse Community Feed and Topics	User	Displays moderated questions and topic collections. Keyword search, topic chips, filters, sort and pagination are embedded controls of this browse outcome.
UC-47	View Community Question Detail	User	Shows a permitted question, its answer thread, source labels, moderation state and allowed interactions.
UC-48	Create Community Question	User	Creates a topic-based community question and optionally applies anonymous public display while preserving internal accountability.
UC-49	Edit or Delete Own Community Post	User	Edits or removes the user own community post while the post is not locked by moderation or investigation.
UC-50	Post Community Answer	User	Posts a personal-experience answer or a verified expert answer with the applicable source/role label.
UC-51	Edit or Delete Own Community Answer	User	Edits or removes the user own answer when no moderation lock or evidence hold applies.
UC-52	React to Community Answer	User	Adds or removes a permitted helpfulness reaction to a community answer.
UC-53	Bookmark Community Post	User	Saves or removes a permitted post from the user private bookmark list.
UC-54	Follow or Unfollow Community Topic	User	Subscribes or unsubscribes the user from a community topic for feed relevance and optional notifications.
UC-55	Report Unsafe Community Content or Account	User	Submits a report for misinformation, dangerous advice, hidden advertising,

ID	Use Case	Actor	Description
			harassment, privacy breach or suspicious account behavior.
UC-56	View Moderation Queue	Moderator	Shows reported, automatically flagged and pre-publication content awaiting authorized moderator review.
UC-57	Moderate Community Content	Moderator	Approves, hides, locks, labels or requests revision for a community post or answer based on policy and safety rules.
UC-58	Resolve Content or Account Report and Apply Enforcement	Moderator / System Admin	Resolves a report and, when authorized, applies warning, posting restriction, suspension escalation or no-action decision.
UC-59	Manage Community Topics and Visibility	Moderator / Content Admin	Creates, updates, reorders or disables community topics and their visibility rules.
UC-60	Submit Expert Profile	Expert Applicant	Creates an expert application profile with specialty, experience, support scope and public professional information.
UC-61	Update Expert Profile	Verified Expert	Updates approved professional profile fields and public service information subject to review rules.
UC-62	Submit or Replace Verification Documents	Expert Applicant / Verified Expert	Uploads required credentials, certificates and supporting evidence for verification or renewal.
UC-63	View Verification Status and Renew Submission	Expert Applicant / Verified Expert	Shows verification result, required corrections, expiry status and allows a renewal submission when eligible.
UC-64	Configure Expert Availability and Service Scope	Verified Expert	Configures availability, support methods, service area and public readiness state for the expert directory and nearby support eligibility.
UC-65	Browse Verified Expert Directory	User	Displays verified experts with embedded search and filters by specialty, availability, badge and consented service area.
UC-66	View Verified Expert Profile	User	Shows a verified expert public professional profile, specialties, supported scope, availability and contribution reputation.

ID	Use Case	Actor	Description
UC-67	View Expert Question Queue	Verified Expert	Shows community questions matched to the expert verified specialties and public support scope.
UC-68	Post Verified Expert Answer	Verified Expert	Posts an expert answer with verified badge, source/safety label and contribution attribution.
UC-69	View Contribution Points and Badges	Verified Expert	Shows contribution points, badges and qualifying community activities used to recognize responsible participation.
UC-70	Review Expert Verification Submission	System Admin	Approves, requests supplementation, rejects or renews an expert verification submission after credential review.
UC-71	Restrict, Suspend or Reinstatement Expert Trust Status	System Admin	Applies, updates or lifts an expert restriction/suspension when credentials expire or policy violations require it.
UC-72	Run AI Symptom Intake	Mother	Collects structured symptom information, timing, severity and permitted context through a guided intake flow.
UC-73	View Risk Triage Result	Mother	Shows non-diagnostic green, yellow or red orientation with safe next-step guidance and clearly labelled uncertainty.
UC-74	Open Emergency Support from a Red Risk Result	Mother	Transfers a user from a red-risk result to the Emergency Map and Nearby Care Support flow without treating the result as a diagnosis.
UC-75	Manage Approved AI Knowledge Sources	Content Admin / System Admin	Registers, reviews, versions or disables approved knowledge sources available to AI retrieval and safe answer composition.
UC-76	Configure AI Risk and Red-Flag Rules	System Admin	Maintains conservative risk thresholds, red-flag patterns, safe fallback actions and rule version history.
UC-77	Open Emergency Map	Mother	Opens the shared emergency map after informing the user about location permission and the non-dispatch scope.
UC-78	Find Nearby Care Facilities	Mother	Finds permitted nearby care facilities based on the user approved location or a selected area.

ID	Use Case	Actor	Description
UC-79	View Route, ETA and Quick Call or Navigate	Mother	Shows route/ETA and enables quick call or navigation to a selected care facility through the map/device capability.
UC-80	Share Time-Limited Location and Send Family Emergency Alert	Mother	Shares minimum necessary location/context for a time-limited emergency alert to selected family recipients.
UC-81	Create or Cancel Nearby Support Request	Mother	Creates or cancels a consented nearby-support request visible only to eligible verified experts; no expert arrival is guaranteed.
UC-82	Manage Expert Nearby Availability and Respond to Nearby Support Request	Verified Expert	Lets an eligible verified expert opt in to nearby availability, view minimal request context and accept, decline or stop responding to a request.
UC-83	Add Personal Health Record and Attachment	Mother	Creates a maternal or baby health record and attaches permitted documents or images with type, date and source metadata.
UC-84	Update Health Record Metadata	Mother	Updates title, category, date, source label, tags or note of an owner-entered health record without altering the protected file itself unless allowed.
UC-85	Archive or Delete User-entered Health Record	Mother	Archives or requests deletion of an owner-entered record when retention, sharing and safety-evidence obligations allow.
UC-86	View Health Record Timeline and Detail	Mother	Shows authorized maternal and baby records by time, category and source with embedded filtering inside the timeline.
UC-87	Generate Health Summary	Mother	Creates a selected-period summary from user-chosen maternal or baby records for personal review or permitted sharing.
UC-88	Share Health Summary or Selected Records Under Consent	Mother	Shares selected records or a summary only through an active MF-01 consent grant with scope and expiry.
UC-89	Create Appointment Reminder	Mother	Creates a reminder for a self-entered appointment, follow-up, examination, test or other care event.

ID	Use Case	Actor	Description
UC-90	Create Medicine or Vitamin Reminder	Mother	Creates a user-entered medicine or vitamin reminder based on prior professional advice or personal routine.
UC-91	Create Vaccination Reminder	Mother	Creates or confirms a reminder linked to the selected baby vaccination record or reference schedule.
UC-92	Update, Snooze, Complete, Skip or Delete Care Reminder	Mother	Manages the lifecycle of a user reminder, including time updates, snooze, complete, skip and deletion where permitted.
UC-93	View Today Tasks and Care Plan	Mother	Shows due reminders, checklist items and permitted family tasks in a prioritized today view.
UC-94	Create Care Group	Mother	Creates a care group for a mother, baby or household care context with the creator as initial owner.
UC-95	Invite or Revoke Family Member Invitation	Mother	Sends or revokes a care-group invitation using a permitted contact or invite link.
UC-96	Accept or Reject Care Group Invitation	Family Member	Lets a recipient accept or reject an active care-group invitation after authentication.
UC-97	Manage Care Group Membership	Mother / Family Member	Allows a group owner to remove a member and an individual member to leave a group where rules permit.
UC-98	Manage Family Permission Scope	Mother	Defines and updates family access to selected calendar, tasks, alerts, logs, summaries or records.
UC-99	Create, Update or Cancel Family Care Task	Mother	Creates, edits or cancels a care task assigned to an eligible care-group member.
UC-100	Update Assigned Task Status	Family Member	Updates the status of a task assigned to the signed-in family member, such as in progress, completed or needs help.
UC-101	View Shared Care Calendar, Data and Alerts	Family Member	Displays calendar items, shared care data and family alerts only within the member active permission scope.
UC-102	Browse Verified Content	User	Displays approved articles, FAQs and checklists with embedded keyword, stage and topic controls.

ID	Use Case	Actor	Description
UC-103	View Verified Content Detail	User	Shows the selected article, FAQ or checklist with source, review/version state, updated date and safety notes.
UC-104	Create Verified Content	Content Admin	Creates a draft article, FAQ or checklist with stage/topic mapping, sources and initial review state.
UC-105	Update Verified Content and Sources	Content Admin	Updates an existing content draft or creates a new version when published content needs changes to body, source, tag or safety notes.
UC-106	Review and Publish Content Version	System Admin / Content Admin	Reviews a content version and approves, rejects or requests revision before public publication.
UC-107	Unpublish or Archive Content	Content Admin / System Admin	Withdraws or archives content that is outdated, unsafe, incorrect or replaced by a newer version.
UC-108	Manage Content Categories and Stage/Topic Mapping	Content Admin	Maintains categories, stage mappings, topic tags and display order used by the content hub.
UC-109	Add Expense Entry	Mother	Records a household preparation expense with category, amount, date and optional experience note.
UC-110	Update or Delete Expense Entry	Mother	Corrects or deletes an owner-entered preparation expense.
UC-111	View Expense Summary	Mother	Shows simple expense totals by category, month or journey/baby period.
UC-112	Connect Health Device or Platform	Mother	Connects a supported device/platform after the user reviews permissions and selected indicator scope.
UC-113	Import or Synchronize Device Observations	Mother	Imports or synchronizes selected device observations such as heart rate, sleep, steps, SpO2, temperature or blood pressure where supported.
UC-114	View Device Data Trend and Quality	Mother	Displays imported data trends, source labels, data gaps and quality cautions for the selected indicator.
UC-115	Disconnect Device and Delete Imported Data	Mother	Stops future synchronization and allows the user to request deletion of imported data according to retention rules.

ID	Use Case	Actor	Description
UC-116	Manage Emergency Contacts	Mother	Adds, verifies, reprioritizes or removes emergency alert recipients used by safety monitoring and emergency flows.
UC-117	Configure Smart Activity Monitoring	Mother	Configures sensor consent, alert recipients, countdown, location-sharing option and monitoring conditions.
UC-118	Enable Smart Activity Monitoring	Mother	Starts IMU-based monitoring after configuration, permissions and consent are valid.
UC-119	Disable Smart Activity Monitoring	Mother	Stops IMU monitoring and prevents new candidate safety events from being created after the state change.
UC-120	Handle Suspected Fall or Impact and Safety Check	Mother / Phone IMU Sensor	Detects a candidate fall/impact, opens a safety countdown and records I am OK, Need Help or emergency-support handoff response.
UC-121	Review Safety Event History and Report False Positive	Mother	Shows prior safety events and lets the user mark a false detection with a reason to improve future threshold evaluation.

Table 19: Use Cases Description

3. Functional Specifications

3.1 Shared

3.1.1 Authentication & Profile

3.1.1.1 UC-01 Register Account

UC ID and Name	UC-01 Register Account		
Created By	PhuongNT	Date Created	01/07/2026
Primary Actor	Guest	Secondary Actors	OTP / Email Service
Trigger	Guest selects Register Account.		
Description	Creates an account with a supported contact method and selected initial role. The new account remains unverified until OTP confirmation succeeds.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.		
Postconditions	The requested account outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. Guest opens the Register Account function.		

	2. The system loads the required authorized context and input fields for the account. 3. Guest provides email or phone, password, initial role. 4. The system validates business rules, ownership/consent and data format. 5. The system creates the account, records required audit information and returns the resulting state.
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.
Priority	P0
Frequency of Use	Regular
Business Rules	BR-RBAC: access is limited by effective role and assigned permission. BR-PRIVACY: private maternal, baby and health data are separated from public community identity. BR-ACCOUNT-01: The contact identifier must be unique. BR-ACCOUNT-02: The account is not active until verification succeeds.
Other Information	Platform: Mobile App / Web Portal; Feature / Group: MF-01 - Shared - Account, Trust & Access Control. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.1.1.2 UC-02 Verify OTP

UC ID and Name	UC-02 Verify OTP		
Created By	PhuongNT	Date Created	01/07/2026
Primary Actor	Guest	Secondary Actors	OTP / Email Service
Trigger	Guest submits an OTP code.		
Description	Validates the OTP sent for account activation or another sensitive account action.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.		
Postconditions	The requested OTP proof outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. Guest starts Verify OTP. 2. The system verifies the current authorization, required context and applicable consent state. 3. Guest provides OTP code and target action. 4. The system applies the relevant domain and safety rules, including external-service checks when needed.		

	5. The system completes the OTP proof action, records required audit information and presents the safe result state.
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.
Priority	P0
Frequency of Use	Regular
Business Rules	BR-RBAC: access is limited by effective role and assigned permission. BR-PRIVACY: private maternal, baby and health data are separated from public community identity. BR-ACCOUNT-03: OTP attempts are rate-limited and expire. BR-ACCOUNT-04: Failed verification must not reveal whether unrelated accounts exist.
Other Information	Platform: Mobile App / Web Portal; Feature / Group: MF-01 - Shared - Account, Trust & Access Control. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.1.1.3 UC-03 Log In

UC ID and Name	UC-03 Log In		
Created By	PhuongNT	Date Created	01/07/2026
Primary Actor	User	Secondary Actors	None
Trigger	User submits credentials on the login screen.		
Description	Authenticates a verified account and routes the user to the role-appropriate workspace.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.		
Postconditions	The requested account session outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. User starts Log In. 2. The system verifies the current authorization, required context and applicable consent state. 3. User provides identifier and password. 4. The system applies the relevant domain and safety rules, including external-service checks when needed. 5. The system completes the account session action, records required audit information and presents the safe result state.		

Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.
Priority	P0
Frequency of Use	Regular
Business Rules	BR-RBAC: access is limited by effective role and assigned permission. BR-PRIVACY: private maternal, baby and health data are separated from public community identity. BR-RBAC-01: Routing is based on effective role and account status. BR-SECURITY-01: Suspended or deactivated accounts cannot create a session.
Other Information	Platform: Mobile App / Web Portal; Feature / Group: MF-01 - Shared - Account, Trust & Access Control. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.1.1.4 UC-04 Log Out

UC ID and Name	UC-04 Log Out		
Created By	PhuongNT	Date Created	01/07/2026
Primary Actor	User	Secondary Actors	None
Trigger	User selects Log Out.		
Description	Ends the current session on the active device without deleting the account.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.		
Postconditions	The requested current session outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. User starts Log Out. 2. The system verifies the current authorization, required context and applicable consent state. 3. User provides the required selection, confirmation or input. 4. The system applies the relevant domain and safety rules, including external-service checks when needed. 5. The system completes the current session action, records required audit information and presents the safe result state.		
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.		
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is		

	rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.
Priority	P0
Frequency of Use	Regular
Business Rules	BR-RBAC: access is limited by effective role and assigned permission. BR-PRIVACY: private maternal, baby and health data are separated from public community identity. BR-SECURITY-02: The current refresh token/session must be revoked.
Other Information	Platform: Mobile App / Web Portal; Feature / Group: MF-01 - Shared - Account, Trust & Access Control. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.1.1.5 UC-05 Request Password Reset

UC ID and Name	UC-05 Request Password Reset		
Created By	PhuongNT	Date Created	01/07/2026
Primary Actor	Guest	Secondary Actors	OTP / Email Service
Trigger	Guest selects Forgot Password.		
Description	Starts password recovery by requesting a time-limited reset proof through the registered channel.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.		
Postconditions	The requested password reset outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. Guest starts Request Password Reset. 2. The system verifies the current authorization, required context and applicable consent state. 3. Guest provides registered email or phone. 4. The system applies the relevant domain and safety rules, including external-service checks when needed. 5. The system completes the password reset action, records required audit information and presents the safe result state.		
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.		
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.		
Priority	P0		
Frequency of Use	Regular		
Business Rules	BR-RBAC: access is limited by effective role and assigned permission.		

	BR-PRIVACY: private maternal, baby and health data are separated from public community identity. BR-ACCOUNT-05: The system returns a neutral message whether or not the identifier exists. BR-SECURITY-03: Reset requests are rate-limited.
Other Information	Platform: Mobile App / Web Portal; Feature / Group: MF-01 - Shared - Account, Trust & Access Control. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.1.1.6 UC-06 Reset Password

UC ID and Name	UC-06 Reset Password		
Created By	PhuongNT	Date Created	01/07/2026
Primary Actor	Guest	Secondary Actors	None
Trigger	Guest opens a valid recovery link or submits a valid reset code.		
Description	Sets a new password only after a valid recovery proof is verified.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.		
Postconditions	The requested password outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. Guest opens the selected password. 2. The system verifies the actor authorization, ownership and current record state. 3. Guest selects an allowed action and provides reset proof and new password. 4. The system validates applicable policy, safety and consistency rules. 5. The system applies the change to the password, writes required audit data and refreshes related screens or notifications.		
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.		
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.		
Priority	P0		
Frequency of Use	Regular		
Business Rules	BR-RBAC: access is limited by effective role and assigned permission. BR-PRIVACY: private maternal, baby and health data are separated from public community identity. BR-SECURITY-04: Existing sessions are revoked after a successful password reset.		

Other Information	Platform: Mobile App / Web Portal; Feature / Group: MF-01 - Shared - Account, Trust & Access Control. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.1.1.7 UC-07 Change Password

UC ID and Name	UC-07 Change Password		
Created By	PhuongNT	Date Created	01/07/2026
Primary Actor	User	Secondary Actors	None
Trigger	User selects Change Password.		
Description	Changes the signed-in user password after current-password validation.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.		
Postconditions	The requested password outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. User opens the selected password. 2. The system verifies the actor authorization, ownership and current record state. 3. User selects an allowed action and provides current password and new password. 4. The system validates applicable policy, safety and consistency rules. 5. The system applies the change to the password, writes required audit data and refreshes related screens or notifications.		
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.		
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.		
Priority	P0		
Frequency of Use	Regular		
Business Rules	BR-RBAC: access is limited by effective role and assigned permission. BR-PRIVACY: private maternal, baby and health data are separated from public community identity. BR-SECURITY-05: The current password must be validated before the change.		
Other Information	Platform: Mobile App / Web Portal; Feature / Group: MF-01 - Shared - Account, Trust & Access Control. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.		
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.		

3.1.1.8 UC-08 View Account Profile

UC ID and Name	UC-08 View Account Profile		
Created By	PhuongNT	Date Created	01/07/2026
Primary Actor	User	Secondary Actors	None
Trigger	User opens Account Profile.		
Description	Shows the user private account profile, role, status and permitted account controls.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.		
Postconditions	The requested account profile outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. User opens the View Account Profile function. 2. The system validates role, ownership, active consent and visibility scope before retrieval. 3. The system retrieves the permitted account profile data using any embedded list controls where applicable. 4. User may select an authorized item, time range or permitted embedded control. 5. The system displays the authorized result with applicable source, safety or status labels.		
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.		
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.		
Priority	P0		
Frequency of Use	Frequent		
Business Rules	BR-RBAC: access is limited by effective role and assigned permission. BR-PRIVACY: private maternal, baby and health data are separated from public community identity. BR-PRIVACY-01: Private account data is never displayed in public community contexts.		
Other Information	Platform: Mobile App / Web Portal; Feature / Group: MF-01 - Shared - Account, Trust & Access Control. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.		
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.		

3.1.1.9 UC-09 Update Account Profile

UC ID and Name	UC-09 Update Account Profile		
Created By	PhuongNT	Date Created	01/07/2026

Primary Actor	User	Secondary Actors	None
Trigger	User saves profile changes.		
Description	Updates non-sensitive private profile data used for account operation and contact preferences.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.		
Postconditions	The requested account profile outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. User opens the selected account profile. 2. The system verifies the actor authorization, ownership and current record state. 3. User selects an allowed action and provides name, avatar, area, permitted contact fields. 4. The system validates applicable policy, safety and consistency rules. 5. The system applies the change to the account profile, writes required audit data and refreshes related screens or notifications.		
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.		
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.		
Priority	P0		
Frequency of Use	Regular		
Business Rules	BR-RBAC: access is limited by effective role and assigned permission. BR-PRIVACY: private maternal, baby and health data are separated from public community identity. BR-ACCOUNT-06: Sensitive identity changes may require re-verification.		
Other Information	Platform: Mobile App / Web Portal; Feature / Group: MF-01 - Shared - Account, Trust & Access Control. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.		
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.		

3.1.1.10 UC-10 Manage Community Identity

UC ID and Name	UC-10 Manage Community Identity		
Created By	PhuongNT	Date Created	01/07/2026
Primary Actor	User	Secondary Actors	None
Trigger	User opens Community Identity settings.		
Description	Creates or updates a public community identity that is separate from private maternal, baby and family information.		

Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.
Postconditions	The requested community identity outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.
Normal Flow	1. User opens the selected community identity. 2. The system verifies the actor authorization, ownership and current record state. 3. User selects an allowed action and provides display name, avatar, visibility settings. 4. The system validates applicable policy, safety and consistency rules. 5. The system applies the change to the community identity, writes required audit data and refreshes related screens or notifications.
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.
Priority	P0
Frequency of Use	Regular
Business Rules	BR-RBAC: access is limited by effective role and assigned permission. BR-PRIVACY: private maternal, baby and health data are separated from public community identity. BR-PRIVACY-02: Community identity must not expose private health records by default. BR-COMMUNITY-01: Display names and avatars remain subject to moderation.
Other Information	Platform: Mobile App / Web Portal; Feature / Group: MF-01 - Shared - Account, Trust & Access Control. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.1.1.11 UC-11 Manage Notification Preferences

UC ID and Name	UC-11 Manage Notification Preferences		
Created By	PhuongNT	Date Created	01/07/2026
Primary Actor	User	Secondary Actors	None
Trigger	User changes notification preferences.		
Description	Lets the user choose supported delivery channels and categories for reminders, community replies, family alerts and account events.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.		

Postconditions	The requested notification preference outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.
Normal Flow	1. User opens the selected notification preference. 2. The system verifies the actor authorization, ownership and current record state. 3. User selects an allowed action and provides channels and categories. 4. The system validates applicable policy, safety and consistency rules. 5. The system applies the change to the notification preference, writes required audit data and refreshes related screens or notifications.
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.
Priority	P0
Frequency of Use	Regular
Business Rules	BR-RBAC: access is limited by effective role and assigned permission. BR-PRIVACY: private maternal, baby and health data are separated from public community identity. BR-NOTIFY-01: Safety-critical alerts may remain enabled where required by active monitoring consent.
Other Information	Platform: Mobile App / Web Portal; Feature / Group: MF-01 - Shared - Account, Trust & Access Control. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.1.1.12 UC-12 View and Mark Notifications

UC ID and Name	UC-12 View and Mark Notifications		
Created By	PhuongNT	Date Created	01/07/2026
Primary Actor	User	Secondary Actors	Firebase Cloud Messaging
Trigger	User opens Notification Center.		
Description	Shows authorized notifications and allows the user to mark individual or all eligible notifications as read.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.		
Postconditions	The requested notifications outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. User opens View and Mark Notifications. 2. The system validates access and loads authorized context.		

	3. User provides read state. 4. The system applies business rules and records the result. 5. The system shows the outcome and refreshes related data.
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.
Priority	P0
Frequency of Use	Frequent
Business Rules	BR-RBAC: access is limited by effective role and assigned permission. BR-PRIVACY: private maternal, baby and health data are separated from public community identity. BR-NOTIFY-02: Notification deep links must re-check authorization before opening protected data.
Other Information	Platform: Mobile App / Web Portal; Feature / Group: MF-01 - Shared - Account, Trust & Access Control. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.1.1.13 UC-13 Manage Own Login Sessions

UC ID and Name	UC-13 Manage Own Login Sessions		
Created By	PhuongNT	Date Created	01/07/2026
Primary Actor	User	Secondary Actors	None
Trigger	User opens Login Sessions.		
Description	Lists active or recent own sessions and allows the user to revoke a selected device session.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.		
Postconditions	The requested login sessions outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. User opens the selected login sessions. 2. The system verifies the actor authorization, ownership and current record state. 3. User selects an allowed action and provides selected session. 4. The system validates applicable policy, safety and consistency rules. 5. The system applies the change to the login sessions, writes required audit data and refreshes related screens or notifications.		
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.		

Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.
Priority	P0
Frequency of Use	Regular
Business Rules	BR-RBAC: access is limited by effective role and assigned permission. BR-PRIVACY: private maternal, baby and health data are separated from public community identity. BR-SECURITY-06: A user may not revoke an unrelated account session.
Other Information	Platform: Mobile App / Web Portal; Feature / Group: MF-01 - Shared - Account, Trust & Access Control. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.1.1.14 UC-14 Deactivate or Delete Own Account

UC ID and Name	UC-14 Deactivate or Delete Own Account		
Created By	PhuongNT	Date Created	01/07/2026
Primary Actor	User	Secondary Actors	None
Trigger	User submits an account lifecycle request.		
Description	Requests account deactivation or deletion subject to retention, care-group and audit obligations.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.		
Postconditions	The requested own account outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. User opens the selected own account. 2. The system verifies the actor authorization, ownership and current record state. 3. User selects an allowed action and provides deactivation or deletion request. 4. The system validates applicable policy, safety and consistency rules. 5. The system applies the change to the own account, writes required audit data and refreshes related screens or notifications.		
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.		
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.		
Priority	P0		
Frequency of Use	Occasional		

Business Rules	BR-RBAC: access is limited by effective role and assigned permission. BR-PRIVACY: private maternal, baby and health data are separated from public community identity. BR-PRIVACY-03: Deletion follows retention and legal obligations. BR-FAMILY-01: The system explains ownership impact on active care groups.
Other Information	Platform: Mobile App / Web Portal; Feature / Group: MF-01 - Shared - Account, Trust & Access Control. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.1.2 Consent & Data Permission

3.1.2.1 UC-15 Grant Data Permission

UC ID and Name	UC-15 Grant Data Permission		
Created By	PhuongNT	Date Created	01/07/2026
Primary Actor	User	Secondary Actors	None
Trigger	User confirms a sharing grant.		
Description	Creates a purpose-specific, scope-based and time-limited permission for a family member or verified expert to access selected data.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.		
Postconditions	The requested data permission outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. User starts Grant Data Permission. 2. The system verifies the current authorization, required context and applicable consent state. 3. User provides recipient, data scope, purpose, expiry. 4. The system applies the relevant domain and safety rules, including external-service checks when needed. 5. The system completes the data permission action, records required audit information and presents the safe result state.		
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.		
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.		
Priority	P0		
Frequency of Use	Regular		
Business Rules	BR-RBAC: access is limited by effective role and assigned permission.		

	BR-PRIVACY: private maternal, baby and health data are separated from public community identity. BR-CONSENT-01: The grant must include recipient, scope, purpose and expiry. BR-CONSENT-02: Sharing is denied until the recipient identity is valid.
Other Information	Platform: Mobile App / Web Portal; Feature / Group: MF-01 - Shared - Account, Trust & Access Control. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.1.2.2 UC-16 Review and Revoke Data Permission

UC ID and Name	UC-16 Review and Revoke Data Permission		
Created By	PhuongNT	Date Created	01/07/2026
Primary Actor	User	Secondary Actors	None
Trigger	User opens Sharing History or selects Revoke.		
Description	Shows active and past sharing permissions and lets the owner revoke an active grant.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.		
Postconditions	The requested data permission outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. User starts Review and Revoke Data Permission. 2. The system verifies the current authorization, required context and applicable consent state. 3. User provides selected permission. 4. The system applies the relevant domain and safety rules, including external-service checks when needed. 5. The system completes the data permission action, records required audit information and presents the safe result state.		
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.		
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.		
Priority	P0		
Frequency of Use	Frequent		
Business Rules	BR-RBAC: access is limited by effective role and assigned permission. BR-PRIVACY: private maternal, baby and health data are separated from public community identity. BR-CONSENT-03: Revocation takes effect before the next protected read. BR-AUDIT-01: Grant and revocation events are auditable.		

Other Information	Platform: Mobile App / Web Portal; Feature / Group: MF-01 - Shared - Account, Trust & Access Control. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.2 WebApp

3.2.1 Account & Security Administration

3.2.1.1 UC-17 Administer User Accounts and Role Access

UC ID and Name	UC-17 Administer User Accounts and Role Access		
Created By	PhuongNT	Date Created	01/07/2026
Primary Actor	System Admin	Secondary Actors	None
Trigger	Admin opens user administration.		
Description	Allows authorized administrators to review account status, apply role/access changes and restrict misuse under separation-of-duties controls.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.		
Postconditions	The requested user account and role access outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. System Admin opens Administer User Accounts and Role Access. 2. The system validates access and loads authorized context. 3. System Admin provides target account, role/access decision, reason. 4. The system applies business rules and records the result. 5. The system shows the outcome and refreshes related data.		
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.		
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.		
Priority	P0		
Frequency of Use	Regular		
Business Rules	BR-RBAC: access is limited by effective role and assigned permission. BR-PRIVACY: private maternal, baby and health data are separated from public community identity. BR-RBAC-02: Administrators may only change roles within assigned authority. BR-AUDIT-02: Access changes require an auditable reason.		
Other Information	Platform: Admin Portal; Feature / Group: MF-01 - Shared - Account, Trust & Access Control. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.		

Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.
--------------------	--

3.2.1.2 UC-18 Review Sensitive Access and Security Events

UC ID and Name	UC-18 Review Sensitive Access and Security Events		
Created By	PhuongNT	Date Created	01/07/2026
Primary Actor	System Admin	Secondary Actors	None
Trigger	Admin opens audit and security review.		
Description	Reviews sensitive-record access, abnormal login, permission change and file-access events to determine follow-up or investigation.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.		
Postconditions	The requested audit/security event outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. System Admin opens the Review Sensitive Access and Security Events function. 2. The system validates role, ownership, active consent and visibility scope before retrieval. 3. The system retrieves the permitted audit/security event data using any embedded list controls where applicable. 4. System Admin may select an authorized item, time range or permitted embedded control. 5. The system displays the authorized result with applicable source, safety or status labels.		
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.		
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.		
Priority	P0		
Frequency of Use	Frequent		
Business Rules	BR-RBAC: access is limited by effective role and assigned permission. BR-PRIVACY: private maternal, baby and health data are separated from public community identity. BR-AUDIT-03: Audit records are append-only for normal operators. BR-PRIVACY-04: Audit viewers must not receive more health content than necessary.		
Other Information	Platform: Admin Portal; Feature / Group: MF-01 - Shared - Account, Trust & Access Control. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.		

Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.
--------------------	--

3.2.2 Community Moderation & Governance

3.2.2.1 UC-56 View Moderation Queue

UC ID and Name	UC-56 View Moderation Queue		
Created By	BachNQ	Date Created	04/07/2026
Primary Actor	Moderator	Secondary Actors	None
Trigger	Moderator opens Moderation Queue.		
Description	Shows reported, automatically flagged and pre-publication content awaiting authorized moderator review.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.		
Postconditions	The requested moderation queue outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. Moderator opens the View Moderation Queue function. 2. The system validates role, ownership, active consent and visibility scope before retrieval. 3. The system retrieves the permitted moderation queue data using any embedded list controls where applicable. 4. Moderator may select an authorized item, time range or permitted embedded control. 5. The system displays the authorized result with applicable source, safety or status labels.		
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.		
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.		
Priority	P0		
Frequency of Use	Frequent		
Business Rules	BR-MODERATION: visible community content is governed by moderation state, report workflow and safety policy. BR-SAFETY: community content may not be used to diagnose, prescribe or promote unsafe treatment. BR-MODERATION-02: Moderators see only evidence and private context needed to resolve the case.		
Other Information	Platform: Admin Portal; Feature / Group: MF-04 - Community & Admin - Community Q&A & Moderation. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.		

Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.
--------------------	--

3.2.2.2 UC-57 Moderate Community Content

UC ID and Name	UC-57 Moderate Community Content		
Created By	BachNQ	Date Created	04/07/2026
Primary Actor	Moderator	Secondary Actors	None
Trigger	Moderator selects a moderation item.		
Description	Approves, hides, locks, labels or requests revision for a community post or answer based on policy and safety rules.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.		
Postconditions	The requested community content outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. Moderator starts Moderate Community Content. 2. The system verifies the current authorization, required context and applicable consent state. 3. Moderator provides decision, reason, visibility state. 4. The system applies the relevant domain and safety rules, including external-service checks when needed. 5. The system completes the community content action, records required audit information and presents the safe result state.		
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.		
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.		
Priority	P0		
Frequency of Use	Occasional		
Business Rules	BR-MODERATION: visible community content is governed by moderation state, report workflow and safety policy. BR-SAFETY: community content may not be used to diagnose, prescribe or promote unsafe treatment. BR-MODERATION-03: Safety-sensitive content may be escalated rather than published. BR-AUDIT-04: Moderation decisions require a reason and audit trail.		
Other Information	Platform: Admin Portal; Feature / Group: MF-04 - Community & Admin - Community Q&A & Moderation. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.		
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.		

3.2.2.3 UC-58 Resolve Content or Account Report and Apply Enforcement

UC ID and Name	UC-58 Resolve Content or Account Report and Apply Enforcement		
Created By	BachNQ	Date Created	04/07/2026
Primary Actor	Moderator / System Admin	Secondary Actors	None
Trigger	Authorized reviewer resolves a report.		
Description	Resolves a report and, when authorized, applies warning, posting restriction, suspension escalation or no-action decision.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.		
Postconditions	The requested report and enforcement action outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. Moderator / System Admin starts Resolve Content or Account Report and Apply Enforcement. 2. The system verifies the current authorization, required context and applicable consent state. 3. Moderator / System Admin provides resolution, target state, reason. 4. The system applies the relevant domain and safety rules, including external-service checks when needed. 5. The system completes the report and enforcement action action, records required audit information and presents the safe result state.		
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.		
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.		
Priority	P0		
Frequency of Use	Regular		
Business Rules	BR-MODERATION: visible community content is governed by moderation state, report workflow and safety policy. BR-SAFETY: community content may not be used to diagnose, prescribe or promote unsafe treatment. BR-ENFORCEMENT-01: Enforcement severity follows configured policy and authority. BR-ENFORCEMENT-02: A restriction must identify scope and effective period.		
Other Information	Platform: Admin Portal; Feature / Group: MF-04 - Community & Admin - Community Q&A & Moderation. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.		
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.		

3.2.2.4 UC-59 Manage Community Topics and Visibility

UC ID and Name	UC-59 Manage Community Topics and Visibility		
Created By	BachNQ	Date Created	04/07/2026
Primary Actor	Moderator / Content Admin	Secondary Actors	None
Trigger	Authorized staff saves a topic configuration.		
Description	Creates, updates, reorders or disables community topics and their visibility rules.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.		
Postconditions	The requested community topic outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. Moderator / Content Admin opens the selected community topic. 2. The system verifies the actor authorization, ownership and current record state. 3. Moderator / Content Admin selects an allowed action and provides name, description, order, visibility. 4. The system validates applicable policy, safety and consistency rules. 5. The system applies the change to the community topic, writes required audit data and refreshes related screens or notifications.		
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.		
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.		
Priority	P0		
Frequency of Use	Regular		
Business Rules	BR-MODERATION: visible community content is governed by moderation state, report workflow and safety policy. BR-SAFETY: community content may not be used to diagnose, prescribe or promote unsafe treatment. BR-COMMUNITY-11: Topic changes do not alter the authorship of existing posts.		
Other Information	Platform: Admin Portal; Feature / Group: MF-04 - Community & Admin - Community Q&A & Moderation. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.		
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.		

3.2.3 Expert Network & Verification

3.2.3.1 UC-60 Submit Expert Profile

UC ID and Name	UC-60 Submit Expert Profile		
Created By	LamVH	Date Created	05/07/2026
Primary Actor	Expert Applicant	Secondary Actors	None
Trigger	Expert applicant starts expert onboarding.		
Description	Creates an expert application profile with specialty, experience, support scope and public professional information.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.		
Postconditions	The requested expert profile outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. Expert Applicant opens the Submit Expert Profile function. 2. The system loads the required authorized context and input fields for the expert profile. 3. Expert Applicant provides professional identity, specialty, experience, support scope. 4. The system validates business rules, ownership/consent and data format. 5. The system creates the expert profile, records required audit information and returns the resulting state.		
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.		
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.		
Priority	P0		
Frequency of Use	Regular		
Business Rules	BR-EXPERT: verified status is required for verified badge and verified expert response. BR-TRUST: contribution points and badges show participation, not clinical competence. BR-EXPERT-02: An applicant is not shown as verified until approval.		
Other Information	Platform: Expert Portal / Expert App; Feature / Group: MF-05 - Expert & Admin - Verified Expert Network & Contribution. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.		
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.		

3.2.3.2 UC-61 Update Expert Profile

UC ID and Name	UC-61 Update Expert Profile
----------------	-----------------------------

Created By	LamVH	Date Created	05/07/2026
Primary Actor	Verified Expert	Secondary Actors	None
Trigger	Expert saves profile changes.		
Description	Updates approved professional profile fields and public service information subject to review rules.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.		
Postconditions	The requested expert profile outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. Verified Expert opens the selected expert profile. 2. The system verifies the actor authorization, ownership and current record state. 3. Verified Expert selects an allowed action and provides profile fields and public display information. 4. The system validates applicable policy, safety and consistency rules. 5. The system applies the change to the expert profile, writes required audit data and refreshes related screens or notifications.		
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.		
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.		
Priority	P0		
Frequency of Use	Regular		
Business Rules	BR-EXPERT: verified status is required for verified badge and verified expert response. BR-TRUST: contribution points and badges show participation, not clinical competence. BR-EXPERT-03: Material credential/scope changes may return the profile to review.		
Other Information	Platform: Expert Portal / Expert App; Feature / Group: MF-05 - Expert & Admin - Verified Expert Network & Contribution. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.		
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.		

3.2.3.3 UC-62 Submit or Replace Verification Documents

UC ID and Name	UC-62 Submit or Replace Verification Documents		
Created By	LamVH	Date Created	05/07/2026
Primary Actor	Expert Applicant / Verified Expert	Secondary Actors	Protected Object Storage

Trigger	Expert submits verification documents.
Description	Uploads required credentials, certificates and supporting evidence for verification or renewal.
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.
Postconditions	The requested verification document set outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.
Normal Flow	1. Expert Applicant / Verified Expert opens the Submit or Replace Verification Documents function. 2. The system loads the required authorized context and input fields for the verification document set. 3. Expert Applicant / Verified Expert provides credential files, expiry dates, declarations. 4. The system validates business rules, ownership/consent and data format. 5. The system creates the verification document set, records required audit information and returns the resulting state.
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.
Priority	P0
Frequency of Use	Regular
Business Rules	BR-EXPERT: verified status is required for verified badge and verified expert response. BR-TRUST: contribution points and badges show participation, not clinical competence. BR-EXPERT-04: Documents are protected and visible only to authorized reviewers. BR-EXPERT-05: Expiry dates are retained where applicable.
Other Information	Platform: Expert Portal / Expert App; Feature / Group: MF-05 - Expert & Admin - Verified Expert Network & Contribution. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.2.3.4 UC-63 View Verification Status and Renew Submission

UC ID and Name	UC-63 View Verification Status and Renew Submission		
Created By	LamVH	Date Created	05/07/2026
Primary Actor	Expert Applicant / Verified Expert	Secondary Actors	None

Trigger	Expert opens Verification Status.
Description	Shows verification result, required corrections, expiry status and allows a renewal submission when eligible.
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.
Postconditions	The requested verification status outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.
Normal Flow	1. Expert Applicant / Verified Expert opens View Verification Status and Renew Submission. 2. The system validates access and loads authorized context. 3. Expert Applicant / Verified Expert provides selected renewal action. 4. The system applies business rules and records the result. 5. The system shows the outcome and refreshes related data.
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.
Priority	P0
Frequency of Use	Frequent
Business Rules	BR-EXPERT: verified status is required for verified badge and verified expert response. BR-TRUST: contribution points and badges show participation, not clinical competence. BR-EXPERT-06: A rejected or expired verification limits verified-only capabilities.
Other Information	Platform: Expert Portal / Expert App; Feature / Group: MF-05 - Expert & Admin - Verified Expert Network & Contribution. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.2.3.5 UC-64 Configure Expert Availability and Service Scope

UC ID and Name	UC-64 Configure Expert Availability and Service Scope		
Created By	LamVH	Date Created	05/07/2026
Primary Actor	Verified Expert	Secondary Actors	None
Trigger	Expert saves availability or scope settings.		
Description	Configures availability, support methods, service area and public readiness state for the expert directory and nearby support eligibility.		

Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.
Postconditions	The requested expert availability and service scope outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.
Normal Flow	1. Verified Expert opens the selected expert availability and service scope. 2. The system verifies the actor authorization, ownership and current record state. 3. Verified Expert selects an allowed action and provides availability slots, online state, support methods, service area. 4. The system validates applicable policy, safety and consistency rules. 5. The system applies the change to the expert availability and service scope, writes required audit data and refreshes related screens or notifications.
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.
Priority	P0
Frequency of Use	Regular
Business Rules	BR-EXPERT: verified status is required for verified badge and verified expert response. BR-TRUST: contribution points and badges show participation, not clinical competence. BR-EXPERT-07: Location-related settings remain consent-based and time-limited where precise sharing is used.
Other Information	Platform: Expert Portal / Expert App; Feature / Group: MF-05 - Expert & Admin - Verified Expert Network & Contribution. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.2.3.6 UC-67 View Expert Question Queue

UC ID and Name	UC-67 View Expert Question Queue		
Created By	LamVH	Date Created	05/07/2026
Primary Actor	Verified Expert	Secondary Actors	None
Trigger	Expert opens Expert Question Queue.		
Description	Shows community questions matched to the expert verified specialties and public support scope.		

Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.
Postconditions	The requested expert question queue outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.
Normal Flow	1. Verified Expert opens the View Expert Question Queue function. 2. The system validates role, ownership, active consent and visibility scope before retrieval. 3. The system retrieves the permitted expert question queue data using any embedded list controls where applicable. 4. Verified Expert may select an authorized item, time range or permitted embedded control. 5. The system displays the authorized result with applicable source, safety or status labels.
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.
Priority	P0
Frequency of Use	Frequent
Business Rules	BR-EXPERT: verified status is required for verified badge and verified expert response. BR-TRUST: contribution points and badges show participation, not clinical competence. BR-EXPERT-11: The queue cannot expose private reporter health data beyond publicly visible question context.
Other Information	Platform: Expert Portal / Expert App; Feature / Group: MF-05 - Expert & Admin - Verified Expert Network & Contribution. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.2.3.7 UC-68 Post Verified Expert Answer

UC ID and Name	UC-68 Post Verified Expert Answer		
Created By	LamVH	Date Created	05/07/2026
Primary Actor	Verified Expert	Secondary Actors	None
Trigger	Expert submits an answer.		
Description	Posts an expert answer with verified badge, source/safety label and contribution attribution.		

Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.
Postconditions	The requested verified expert answer outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.
Normal Flow	1. Verified Expert opens the Post Verified Expert Answer function. 2. The system loads the required authorized context and input fields for the verified expert answer. 3. Verified Expert provides answer, source label, safe scope note. 4. The system validates business rules, ownership/consent and data format. 5. The system creates the verified expert answer, records required audit information and returns the resulting state.
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.
Priority	P0
Frequency of Use	Regular
Business Rules	BR-EXPERT: verified status is required for verified badge and verified expert response. BR-TRUST: contribution points and badges show participation, not clinical competence. BR-EXPERT-12: The answer must remain non-diagnostic and non-prescriptive. BR-COMMUNITY-12: Promotional misuse is prohibited.
Other Information	Platform: Expert Portal / Expert App; Feature / Group: MF-05 - Expert & Admin - Verified Expert Network & Contribution. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.2.3.8 UC-69 View Contribution Points and Badges

UC ID and Name	UC-69 View Contribution Points and Badges		
Created By	LamVH	Date Created	05/07/2026
Primary Actor	Verified Expert	Secondary Actors	None
Trigger	Expert opens Contribution Activity.		
Description	Shows contribution points, badges and qualifying community activities used to recognize responsible participation.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.		

Postconditions	The requested contribution points and badges outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.
Normal Flow	1. Verified Expert opens the View Contribution Points and Badges function. 2. The system validates role, ownership, active consent and visibility scope before retrieval. 3. The system retrieves the permitted contribution points and badges data using any embedded list controls where applicable. 4. Verified Expert may select an authorized item, time range or permitted embedded control. 5. The system displays the authorized result with applicable source, safety or status labels.
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.
Priority	P0
Frequency of Use	Frequent
Business Rules	BR-EXPERT: verified status is required for verified badge and verified expert response. BR-TRUST: contribution points and badges show participation, not clinical competence. BR-EXPERT-13: Points and badges do not prove clinical competence or medical accuracy.
Other Information	Platform: Expert Portal / Expert App; Feature / Group: MF-05 - Expert & Admin - Verified Expert Network & Contribution. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.2.3.9 UC-70 Review Expert Verification Submission

UC ID and Name	UC-70 Review Expert Verification Submission		
Created By	LamVH	Date Created	05/07/2026
Primary Actor	System Admin	Secondary Actors	None
Trigger	Admin opens a verification submission.		
Description	Approves, requests supplementation, rejects or renews an expert verification submission after credential review.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.		

Postconditions	The requested expert verification submission outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.
Normal Flow	1. System Admin opens the Review Expert Verification Submission function. 2. The system validates role, ownership, active consent and visibility scope before retrieval. 3. The system retrieves the permitted expert verification submission data using any embedded list controls where applicable. 4. System Admin may select an authorized item, time range or permitted embedded control. 5. The system displays the authorized result with applicable source, safety or status labels.
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.
Priority	P0
Frequency of Use	Frequent
Business Rules	BR-EXPERT: verified status is required for verified badge and verified expert response. BR-TRUST: contribution points and badges show participation, not clinical competence. BR-EXPERT-14: Verification decisions require documented evidence and rationale. BR-AUDIT-05: Verification decisions are auditable.
Other Information	Platform: Admin Portal; Feature / Group: MF-05 - Expert & Admin - Verified Expert Network & Contribution. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.2.3.10 UC-71 Restrict, Suspend or Reinstate Expert Trust Status

UC ID and Name	UC-71 Restrict, Suspend or Reinstate Expert Trust Status		
Created By	LamVH	Date Created	05/07/2026
Primary Actor	System Admin	Secondary Actors	None
Trigger	Admin submits an expert trust-status decision.		
Description	Applies, updates or lifts an expert restriction/suspension when credentials expire or policy violations require it.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.		

Postconditions	The requested expert trust status outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.
Normal Flow	1. System Admin opens the selected expert trust status. 2. The system verifies the actor authorization, ownership and current record state. 3. System Admin selects an allowed action and provides restriction scope, effective period, reason. 4. The system validates applicable policy, safety and consistency rules. 5. The system applies the change to the expert trust status, writes required audit data and refreshes related screens or notifications.
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.
Priority	P0
Frequency of Use	Occasional
Business Rules	BR-EXPERT: verified status is required for verified badge and verified expert response. BR-TRUST: contribution points and badges show participation, not clinical competence. BR-EXPERT-15: Public badge, question response and nearby availability must reflect the active trust state.
Other Information	Platform: Admin Portal; Feature / Group: MF-05 - Expert & Admin - Verified Expert Network & Contribution. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.2.4 AI Knowledge & Risk Configuration

3.2.4.1 UC-75 Manage Approved AI Knowledge Sources

UC ID and Name	UC-75 Manage Approved AI Knowledge Sources		
Created By	PhuongNT	Date Created	06/07/2026
Primary Actor	Content Admin / System Admin	Secondary Actors	None
Trigger	Authorized staff saves a knowledge source decision.		
Description	Registers, reviews, versions or disables approved knowledge sources available to AI retrieval and safe answer composition.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.		

Postconditions	The requested approved AI knowledge source outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.
Normal Flow	1. Content Admin / System Admin opens the selected approved AI knowledge source. 2. The system verifies the actor authorization, ownership and current record state. 3. Content Admin / System Admin selects an allowed action and provides source, reviewer, version, validity, status. 4. The system validates applicable policy, safety and consistency rules. 5. The system applies the change to the approved AI knowledge source, writes required audit data and refreshes related screens or notifications.
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.
Priority	P0
Frequency of Use	Regular
Business Rules	BR-AI: outputs are non-diagnostic and non-prescriptive. BR-SAFETY: red-flag routing must provide conservative safe next actions. BR-AI-06: AI retrieval must use approved, traceable source material only.
Other Information	Platform: Admin Portal; Feature / Group: MF-06 - Mother & Admin - AI Nurse Assistant & Risk Triage. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.2.4.2 UC-76 Configure AI Risk and Red-Flag Rules

UC ID and Name	UC-76 Configure AI Risk and Red-Flag Rules		
Created By	PhuongNT	Date Created	06/07/2026
Primary Actor	System Admin	Secondary Actors	None
Trigger	Admin publishes a risk-rule configuration.		
Description	Maintains conservative risk thresholds, red-flag patterns, safe fallback actions and rule version history.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.		
Postconditions	The requested AI risk/red-flag rule outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. System Admin opens the selected AI risk/red-flag rule. 2. The system verifies the actor authorization, ownership and current record state.		

	3. System Admin selects an allowed action and provides rule version, threshold, action, effective status. 4. The system validates applicable policy, safety and consistency rules. 5. The system applies the change to the AI risk/red-flag rule, writes required audit data and refreshes related screens or notifications.
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.
Priority	P0
Frequency of Use	Regular
Business Rules	BR-AI: outputs are non-diagnostic and non-prescriptive. BR-SAFETY: red-flag routing must provide conservative safe next actions. BR-AI-07: Rule changes require versioning and audit history. BR-AI-08: A safer fallback takes precedence when confidence is insufficient.
Other Information	Platform: Admin Portal; Feature / Group: MF-06 - Mother & Admin - AI Nurse Assistant & Risk Triage. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.2.5 Nearby Expert Support

3.2.5.1 UC-82 Manage Expert Nearby Availability and Respond to Nearby Support Request

UC ID and Name	UC-82 Manage Expert Nearby Availability and Respond to Nearby Support Request		
Created By	ChuongBD	Date Created	07/07/2026
Primary Actor	Verified Expert	Secondary Actors	TrackAsia Map Service
Trigger	Expert manages nearby availability or opens an eligible request.		
Description	Lets an eligible verified expert opt in to nearby availability, view minimal request context and accept, decline or stop responding to a request.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.		
Postconditions	The requested nearby support request outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. Verified Expert starts Manage Expert Nearby Availability and Respond to Nearby Support Request. 2. The system verifies the current authorization, required context and applicable consent state.		

	3. Verified Expert provides availability, response decision, location consent. 4. The system applies the relevant domain and safety rules, including external-service checks when needed. 5. The system completes the nearby support request action, records required audit information and presents the safe result state.
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.
Priority	P0
Frequency of Use	Regular
Business Rules	BR-EMERGENCY: the system supports faster action but does not dispatch emergency services or guarantee arrival. BR-LOCATION: location is consented, minimum necessary and time-limited. BR-EMERGENCY-05: Expert location sharing is opt-in, configurable and time-limited. BR-EMERGENCY-06: Acceptance does not create a medical emergency dispatch obligation.
Other Information	Platform: Expert App / Expert Portal; Feature / Group: MF-07 - Mother/Expert Mobile - Emergency Map & Nearby Care Support. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.2.6 Verified Content Administration

3.2.6.1 UC-104 Create Verified Content

UC ID and Name	UC-104 Create Verified Content		
Created By	LamVH	Date Created	11/07/2026
Primary Actor	Content Admin	Secondary Actors	None
Trigger	Content Admin saves a new content draft.		
Description	Creates a draft article, FAQ or checklist with stage/topic mapping, sources and initial review state.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; the relevant P0 identity/consent context exists; and required reference data exists.		
Postconditions	The requested verified content draft outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. Content Admin opens the Create Verified Content function.		

	2. The system loads the required authorized context and input fields for the verified content draft. 3. Content Admin provides content type, title, body, sources, stage/topic. 4. The system validates business rules, ownership/consent and data format. 5. The system creates the verified content draft, records required audit information and returns the resulting state.
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.
Priority	P1
Frequency of Use	Regular
Business Rules	BR-CONTENT: public content must have review state, source labeling and version traceability. BR-CONTENT-04: Drafts are not publicly visible before approval.
Other Information	Platform: Admin Portal; Feature / Group: MF-11 - Mobile & Admin - Verified Content & Checklist Hub. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.2.6.2 UC-105 Update Verified Content and Sources

UC ID and Name	UC-105 Update Verified Content and Sources		
Created By	LamVH	Date Created	11/07/2026
Primary Actor	Content Admin	Secondary Actors	None
Trigger	Content Admin saves content changes.		
Description	Updates an existing content draft or creates a new version when published content needs changes to body, source, tag or safety notes.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; the relevant P0 identity/consent context exists; and required reference data exists.		
Postconditions	The requested verified content and sources outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. Content Admin opens the selected verified content and sources. 2. The system verifies the actor authorization, ownership and current record state. 3. Content Admin selects an allowed action and provides content changes, source changes, version note. 4. The system validates applicable policy, safety and consistency rules. 5. The system applies the change to the verified content and sources, writes required audit data and refreshes related screens or notifications.		

Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.
Priority	P1
Frequency of Use	Regular
Business Rules	BR-CONTENT: public content must have review state, source labeling and version traceability. BR-CONTENT-05: Published content changes must preserve version traceability.
Other Information	Platform: Admin Portal; Feature / Group: MF-11 - Mobile & Admin - Verified Content & Checklist Hub. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.2.6.3 UC-106 Review and Publish Content Version

UC ID and Name	UC-106 Review and Publish Content Version		
Created By	LamVH	Date Created	11/07/2026
Primary Actor	System Admin / Content Admin	Secondary Actors	None
Trigger	Authorized reviewer submits a review decision.		
Description	Reviews a content version and approves, rejects or requests revision before public publication.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; the relevant P0 identity/consent context exists; and required reference data exists.		
Postconditions	The requested content version outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. System Admin / Content Admin opens the Review and Publish Content Version function. 2. The system validates role, ownership, active consent and visibility scope before retrieval. 3. The system retrieves the permitted content version data using any embedded list controls where applicable. 4. System Admin / Content Admin may select an authorized item, time range or permitted embedded control. 5. The system displays the authorized result with applicable source, safety or status labels.		
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.		

Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.
Priority	P1
Frequency of Use	Frequent
Business Rules	BR-CONTENT: public content must have review state, source labeling and version traceability. BR-CONTENT-06: A published version must identify responsible review and effective date.
Other Information	Platform: Admin Portal; Feature / Group: MF-11 - Mobile & Admin - Verified Content & Checklist Hub. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.2.6.4 UC-107 Unpublish or Archive Content

UC ID and Name	UC-107 Unpublish or Archive Content		
Created By	LamVH	Date Created	11/07/2026
Primary Actor	Content Admin / System Admin	Secondary Actors	None
Trigger	Authorized staff unpublishes or archives content.		
Description	Withdraws or archives content that is outdated, unsafe, incorrect or replaced by a newer version.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; the relevant P0 identity/consent context exists; and required reference data exists.		
Postconditions	The requested verified content outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. Content Admin / System Admin opens the selected verified content. 2. The system verifies the actor authorization, ownership and current record state. 3. Content Admin / System Admin selects an allowed action and provides selected item, reason, state. 4. The system validates applicable policy, safety and consistency rules. 5. The system applies the change to the verified content, writes required audit data and refreshes related screens or notifications.		
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.		
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.		
Priority	P1		

Frequency of Use	Occasional
Business Rules	BR-CONTENT: public content must have review state, source labeling and version traceability. BR-CONTENT-07: Withdrawal preserves a historical record for audit while removing public visibility.
Other Information	Platform: Admin Portal; Feature / Group: MF-11 - Mobile & Admin - Verified Content & Checklist Hub. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.2.6.5 UC-108 Manage Content Categories and Stage/Topic Mapping

UC ID and Name	UC-108 Manage Content Categories and Stage/Topic Mapping		
Created By	LamVH	Date Created	11/07/2026
Primary Actor	Content Admin	Secondary Actors	None
Trigger	Content Admin saves category configuration.		
Description	Maintains categories, stage mappings, topic tags and display order used by the content hub.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; the relevant P0 identity/consent context exists; and required reference data exists.		
Postconditions	The requested content category mapping outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. Content Admin opens the selected content category mapping. 2. The system verifies the actor authorization, ownership and current record state. 3. Content Admin selects an allowed action and provides category, stage, topic, display order, status. 4. The system validates applicable policy, safety and consistency rules. 5. The system applies the change to the content category mapping, writes required audit data and refreshes related screens or notifications.		
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.		
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.		
Priority	P1		
Frequency of Use	Regular		
Business Rules	BR-CONTENT: public content must have review state, source labeling and version traceability. BR-CONTENT-08: Category changes do not silently delete content history.		

Other Information	Platform: Admin Portal; Feature / Group: MF-11 - Mobile & Admin - Verified Content & Checklist Hub. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.3 MobileApp

3.3.1 Mother Care Journey

3.3.1.1 UC-19 Initialize Mother Care Journey

UC ID and Name	UC-19 Initialize Mother Care Journey		
Created By	ChuongBD	Date Created	02/07/2026
Primary Actor	Mother	Secondary Actors	None
Trigger	Mother starts Journey Setup.		
Description	Creates a mother journey for preconception, pregnancy or postpartum recovery using the minimum dates and stage context required for stage-based support.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.		
Postconditions	The requested mother care journey outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. Mother opens the Initialize Mother Care Journey function. 2. The system loads the required authorized context and input fields for the mother care journey. 3. Mother provides journey stage and key dates. 4. The system validates business rules, ownership/consent and data format. 5. The system creates the mother care journey, records required audit information and returns the resulting state.		
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.		
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.		
Priority	P0		
Frequency of Use	Regular		
Business Rules	BR-HEALTH-BOUNDARY: journey and exercise support does not diagnose, prescribe or replace professional assessment. BR-OWNERSHIP: mother journey data is private by default. BR-JOURNEY-01: Estimated dates are supportive planning values, not medical diagnosis.		

Other Information	Platform: Mother Mobile App; Feature / Group: MF-02 - Mother Mobile - Mother Care Journey. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.3.1.2 UC-20 Update Mother Journey Stage and Dates

UC ID and Name	UC-20 Update Mother Journey Stage and Dates		
Created By	ChuongBD	Date Created	02/07/2026
Primary Actor	Mother	Secondary Actors	None
Trigger	Mother saves journey changes.		
Description	Updates stage, last menstrual period, expected due date, birth date or other permitted journey dates when circumstances change.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.		
Postconditions	The requested mother care journey outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. Mother opens the selected mother care journey. 2. The system verifies the actor authorization, ownership and current record state. 3. Mother selects an allowed action and provides stage and date fields. 4. The system validates applicable policy, safety and consistency rules. 5. The system applies the change to the mother care journey, writes required audit data and refreshes related screens or notifications.		
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.		
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.		
Priority	P0		
Frequency of Use	Regular		
Business Rules	BR-HEALTH-BOUNDARY: journey and exercise support does not diagnose, prescribe or replace professional assessment. BR-OWNERSHIP: mother journey data is private by default. BR-JOURNEY-02: The system recalculates stage-dependent suggestions after valid changes.		
Other Information	Platform: Mother Mobile App; Feature / Group: MF-02 - Mother Mobile - Mother Care Journey. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.		
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.		

3.3.1.3 UC-21 View Mother Journey Dashboard

UC ID and Name	UC-21 View Mother Journey Dashboard		
Created By	ChuongBD	Date Created	02/07/2026
Primary Actor	Mother	Secondary Actors	None
Trigger	Mother opens Mother Journey.		
Description	Shows current stage/week, relevant care priorities, reminders, checklists, reviewed content and safe shortcuts.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.		
Postconditions	The requested mother journey dashboard outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. Mother opens the View Mother Journey Dashboard function. 2. The system validates role, ownership, active consent and visibility scope before retrieval. 3. The system retrieves the permitted mother journey dashboard data using any embedded list controls where applicable. 4. Mother may select an authorized item, time range or permitted embedded control. 5. The system displays the authorized result with applicable source, safety or status labels.		
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.		
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.		
Priority	P0		
Frequency of Use	Frequent		
Business Rules	BR-HEALTH-BOUNDARY: journey and exercise support does not diagnose, prescribe or replace professional assessment. BR-OWNERSHIP: mother journey data is private by default. BR-JOURNEY-03: The dashboard is informational and does not make clinical conclusions.		
Other Information	Platform: Mother Mobile App; Feature / Group: MF-02 - Mother Mobile - Mother Care Journey. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.		
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.		

3.3.1.4 UC-22 Add Maternal Health Metric

UC ID and Name	UC-22 Add Maternal Health Metric		
Created By	ChuongBD	Date Created	02/07/2026

Primary Actor	Mother	Secondary Actors	None
Trigger	Mother saves a new maternal metric.		
Description	Records a user-entered maternal indicator such as weight, blood pressure, glucose, fetal movement note or another supported observation.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.		
Postconditions	The requested maternal health metric outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. Mother opens the Add Maternal Health Metric function. 2. The system loads the required authorized context and input fields for the maternal health metric. 3. Mother provides metric type, value, time, source, note. 4. The system validates business rules, ownership/consent and data format. 5. The system creates the maternal health metric, records required audit information and returns the resulting state.		
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.		
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.		
Priority	P0		
Frequency of Use	Regular		
Business Rules	BR-HEALTH-BOUNDARY: journey and exercise support does not diagnose, prescribe or replace professional assessment. BR-OWNERSHIP: mother journey data is private by default. BR-METRIC-01: The source and timestamp are required for trend interpretation. BR-HEALTH-BOUNDARY-01: Values are not used to diagnose disease.		
Other Information	Platform: Mother Mobile App; Feature / Group: MF-02 - Mother Mobile - Mother Care Journey. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.		
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.		

3.3.1.5 UC-23 Update or Delete Maternal Health Metric

UC ID and Name	UC-23 Update or Delete Maternal Health Metric		
Created By	ChuongBD	Date Created	02/07/2026
Primary Actor	Mother	Secondary Actors	None
Trigger	Mother edits or removes an existing metric.		
Description	Corrects or removes a user-entered maternal metric when ownership and record state allow.		

Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.
Postconditions	The requested maternal health metric outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.
Normal Flow	1. Mother opens the selected maternal health metric. 2. The system verifies the actor authorization, ownership and current record state. 3. Mother selects an allowed action and provides selected metric and changes. 4. The system validates applicable policy, safety and consistency rules. 5. The system applies the change to the maternal health metric, writes required audit data and refreshes related screens or notifications.
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.
Priority	P0
Frequency of Use	Occasional
Business Rules	BR-HEALTH-BOUNDARY: journey and exercise support does not diagnose, prescribe or replace professional assessment. BR-OWNERSHIP: mother journey data is private by default. BR-METRIC-02: Delete is limited to owner-entered data and retains audit metadata when required.
Other Information	Platform: Mother Mobile App; Feature / Group: MF-02 - Mother Mobile - Mother Care Journey. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.3.1.6 UC-24 View Maternal Health Trend

UC ID and Name	UC-24 View Maternal Health Trend		
Created By	ChuongBD	Date Created	02/07/2026
Primary Actor	Mother	Secondary Actors	None
Trigger	Mother selects a metric trend.		
Description	Displays time-based trends for recorded maternal indicators with source labels and clear non-diagnostic context.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.		

Postconditions	The requested maternal health trend outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.
Normal Flow	1. Mother opens the View Maternal Health Trend function. 2. The system validates role, ownership, active consent and visibility scope before retrieval. 3. The system retrieves the permitted maternal health trend data using any embedded list controls where applicable. 4. Mother may select an authorized item, time range or permitted embedded control. 5. The system displays the authorized result with applicable source, safety or status labels.
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.
Priority	P0
Frequency of Use	Frequent
Business Rules	BR-HEALTH-BOUNDARY: journey and exercise support does not diagnose, prescribe or replace professional assessment. BR-OWNERSHIP: mother journey data is private by default. BR-HEALTH-BOUNDARY-02: Charts must not state or imply a medical diagnosis.
Other Information	Platform: Mother Mobile App; Feature / Group: MF-02 - Mother Mobile - Mother Care Journey. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.3.1.7 UC-25 Add Postpartum Recovery Log

UC ID and Name	UC-25 Add Postpartum Recovery Log		
Created By	ChuongBD	Date Created	02/07/2026
Primary Actor	Mother	Secondary Actors	None
Trigger	Mother saves a postpartum log.		
Description	Records a postpartum recovery observation such as sleep, mood, milk-feeding context, pain note or recovery symptom note.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.		
Postconditions	The requested postpartum recovery log outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. Mother opens the Add Postpartum Recovery Log function.		

	2. The system loads the required authorized context and input fields for the postpartum recovery log. 3. Mother provides log type, time, note, optional severity. 4. The system validates business rules, ownership/consent and data format. 5. The system creates the postpartum recovery log, records required audit information and returns the resulting state.
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.
Priority	P0
Frequency of Use	Regular
Business Rules	BR-HEALTH-BOUNDARY: journey and exercise support does not diagnose, prescribe or replace professional assessment. BR-OWNERSHIP: mother journey data is private by default. BR-POSTPARTUM-01: Sensitive entries remain private unless actively shared.
Other Information	Platform: Mother Mobile App; Feature / Group: MF-02 - Mother Mobile - Mother Care Journey. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.3.1.8 UC-26 Update or Delete Postpartum Recovery Log

UC ID and Name	UC-26 Update or Delete Postpartum Recovery Log		
Created By	ChuongBD	Date Created	02/07/2026
Primary Actor	Mother	Secondary Actors	None
Trigger	Mother edits or removes a postpartum log.		
Description	Corrects or removes an owner-entered postpartum recovery log.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.		
Postconditions	The requested postpartum recovery log outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. Mother opens the selected postpartum recovery log. 2. The system verifies the actor authorization, ownership and current record state. 3. Mother selects an allowed action and provides selected log and changes. 4. The system validates applicable policy, safety and consistency rules. 5. The system applies the change to the postpartum recovery log, writes required audit data and refreshes related screens or notifications.		

Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.
Priority	P0
Frequency of Use	Occasional
Business Rules	BR-HEALTH-BOUNDARY: journey and exercise support does not diagnose, prescribe or replace professional assessment. BR-OWNERSHIP: mother journey data is private by default. BR-POSTPARTUM-02: Deleted entries follow retention/audit policy.
Other Information	Platform: Mother Mobile App; Feature / Group: MF-02 - Mother Mobile - Mother Care Journey. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.3.1.9 UC-27 Manage Preparation Checklist

UC ID and Name	UC-27 Manage Preparation Checklist		
Created By	ChuongBD	Date Created	02/07/2026
Primary Actor	Mother	Secondary Actors	None
Trigger	Mother opens a preparation checklist.		
Description	Adds, edits and completes preparation items for pregnancy, hospital bag, postpartum recovery or early baby care.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.		
Postconditions	The requested preparation checklist outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. Mother opens the selected preparation checklist. 2. The system verifies the actor authorization, ownership and current record state. 3. Mother selects an allowed action and provides item, status, due context, note. 4. The system validates applicable policy, safety and consistency rules. 5. The system applies the change to the preparation checklist, writes required audit data and refreshes related screens or notifications.		
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.		
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is		

	rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.
Priority	P0
Frequency of Use	Regular
Business Rules	BR-HEALTH-BOUNDARY: journey and exercise support does not diagnose, prescribe or replace professional assessment. BR-OWNERSHIP: mother journey data is private by default. BR-CHECKLIST-01: Checklist suggestions are organizational prompts, not clinical orders.
Other Information	Platform: Mother Mobile App; Feature / Group: MF-02 - Mother Mobile - Mother Care Journey. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.3.1.10 UC-28 Browse Pregnancy Exercise Library

UC ID and Name	UC-28 Browse Pregnancy Exercise Library		
Created By	ChuongBD	Date Created	02/07/2026
Primary Actor	Mother	Secondary Actors	None
Trigger	Mother opens Exercise Library.		
Description	Shows reviewed pregnancy exercise content with embedded stage/difficulty/duration filters and exercise-specific safety notes.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.		
Postconditions	The requested pregnancy exercise library outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. Mother opens the Browse Pregnancy Exercise Library function. 2. The system validates role, ownership, active consent and visibility scope before retrieval. 3. The system retrieves the permitted pregnancy exercise library data using any embedded list controls where applicable. 4. Mother may select an authorized item, time range or permitted embedded control. 5. The system displays the authorized result with applicable source, safety or status labels.		
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.		
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.		
Priority	P0		
Frequency of Use	Frequent		

Business Rules	BR-HEALTH-BOUNDARY: journey and exercise support does not diagnose, prescribe or replace professional assessment. BR-OWNERSHIP: mother journey data is private by default. BR-EXERCISE-01: Only reviewed and active exercises are visible. BR-EXERCISE-02: Raw camera media is not retained.
Other Information	Platform: Mother Mobile App; Feature / Group: MF-02 - Mother Mobile - Mother Care Journey. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.3.1.11 UC-29 Complete Pre-exercise Safety Check

UC ID and Name	UC-29 Complete Pre-exercise Safety Check		
Created By	ChuongBD	Date Created	02/07/2026
Primary Actor	Mother	Secondary Actors	None
Trigger	Mother selects Start Exercise.		
Description	Collects required pre-exercise safety answers and stops the session entry when configured warning answers are present.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.		
Postconditions	The requested pre-exercise safety check outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. Mother starts Complete Pre-exercise Safety Check. 2. The system verifies the current authorization, required context and applicable consent state. 3. Mother provides safety responses. 4. The system applies the relevant domain and safety rules, including external-service checks when needed. 5. The system completes the pre-exercise safety check action, records required audit information and presents the safe result state.		
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.		
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.		
Priority	P0		
Frequency of Use	Regular		
Business Rules	BR-HEALTH-BOUNDARY: journey and exercise support does not diagnose, prescribe or replace professional assessment. BR-OWNERSHIP: mother journey data is private by default. BR-EXERCISE-03: A configured warning result blocks exercise session start and gives safe next-step guidance.		

	BR-EXERCISE-04: Passing the check does not certify medical fitness.
Other Information	Platform: Mother Mobile App; Feature / Group: MF-02 - Mother Mobile - Mother Care Journey. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.3.1.12 UC-30 Conduct Pregnancy Exercise Session with Optional Posture Feedback

UC ID and Name	UC-30 Conduct Pregnancy Exercise Session with Optional Posture Feedback		
Created By	ChuongBD	Date Created	02/07/2026
Primary Actor	Mother	Secondary Actors	MediaPipe / On-device Pose Service
Trigger	Mother starts an eligible exercise session.		
Description	Runs an approved exercise session with start, pause, resume and complete actions and optional rule/ML-based posture feedback after camera consent.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.		
Postconditions	The requested pregnancy exercise session outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. Mother starts Conduct Pregnancy Exercise Session with Optional Posture Feedback. 2. The system verifies the current authorization, required context and applicable consent state. 3. Mother provides exercise selection, optional camera consent, session state. 4. The system applies the relevant domain and safety rules, including external-service checks when needed. 5. The system completes the pregnancy exercise session action, records required audit information and presents the safe result state.		
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.		
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.		
Priority	P0		
Frequency of Use	Regular		
Business Rules	BR-HEALTH-BOUNDARY: journey and exercise support does not diagnose, prescribe or replace professional assessment. BR-OWNERSHIP: mother journey data is private by default. BR-EXERCISE-05: Posture feedback is supportive and not professional exercise supervision.		

	BR-EXERCISE-06: No raw image/video is retained.
Other Information	Platform: Mother Mobile App; Feature / Group: MF-02 - Mother Mobile - Mother Care Journey. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.3.1.13 UC-31 View Exercise History and Session Result

UC ID and Name	UC-31 View Exercise History and Session Result		
Created By	ChuongBD	Date Created	02/07/2026
Primary Actor	Mother	Secondary Actors	None
Trigger	Mother opens exercise history or a completed session.		
Description	Shows completed exercise sessions, duration, completion state, aggregate posture feedback and safety notices.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.		
Postconditions	The requested exercise session history outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. Mother opens the View Exercise History and Session Result function. 2. The system validates role, ownership, active consent and visibility scope before retrieval. 3. The system retrieves the permitted exercise session history data using any embedded list controls where applicable. 4. Mother may select an authorized item, time range or permitted embedded control. 5. The system displays the authorized result with applicable source, safety or status labels.		
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.		
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.		
Priority	P0		
Frequency of Use	Frequent		
Business Rules	BR-HEALTH-BOUNDARY: journey and exercise support does not diagnose, prescribe or replace professional assessment. BR-OWNERSHIP: mother journey data is private by default. BR-EXERCISE-07: Results describe session feedback only and must not assess health fitness.		
Other Information	Platform: Mother Mobile App; Feature / Group: MF-02 - Mother Mobile - Mother Care Journey. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.		

Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.
--------------------	--

3.3.2 Baby Care Journey, Growth & Vaccination

3.3.2.1 UC-32 Create Baby Profile

UC ID and Name	UC-32 Create Baby Profile		
Created By	HuyND	Date Created	03/07/2026
Primary Actor	Mother	Secondary Actors	None
Trigger	Mother selects Add Baby.		
Description	Creates a baby profile with core identity and birth context needed for journals, growth and vaccination tracking.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.		
Postconditions	The requested baby profile outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. Mother opens the Create Baby Profile function. 2. The system loads the required authorized context and input fields for the baby profile. 3. Mother provides nickname, date of birth, sex, birth measurements. 4. The system validates business rules, ownership/consent and data format. 5. The system creates the baby profile, records required audit information and returns the resulting state.		
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.		
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.		
Priority	P0		
Frequency of Use	Regular		
Business Rules	BR-HEALTH-BOUNDARY: baby records, trends and milestones support observation rather than diagnosis. BR-OWNERSHIP: baby data is private by default. BR-BABY-01: Only the caregiver owner may create a profile in their care space.		
Other Information	Platform: Mother Mobile App; Feature / Group: MF-03 - Mother Mobile - Baby Care Journey, Growth & Vaccination. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.		
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.		

3.3.2.2 UC-33 Update or Archive Baby Profile

UC ID and Name	UC-33 Update or Archive Baby Profile		
Created By	HuyND	Date Created	03/07/2026
Primary Actor	Mother	Secondary Actors	None
Trigger	Mother saves changes or selects Archive.		
Description	Updates permitted baby profile fields or archives a profile that is no longer actively managed without destroying linked history.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.		
Postconditions	The requested baby profile outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. Mother opens the selected baby profile. 2. The system verifies the actor authorization, ownership and current record state. 3. Mother selects an allowed action and provides selected profile and changes. 4. The system validates applicable policy, safety and consistency rules. 5. The system applies the change to the baby profile, writes required audit data and refreshes related screens or notifications.		
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.		
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.		
Priority	P0		
Frequency of Use	Occasional		
Business Rules	BR-HEALTH-BOUNDARY: baby records, trends and milestones support observation rather than diagnosis. BR-OWNERSHIP: baby data is private by default. BR-BABY-02: Archive preserves linked records and historical audit references.		
Other Information	Platform: Mother Mobile App; Feature / Group: MF-03 - Mother Mobile - Baby Care Journey, Growth & Vaccination. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.		
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.		

3.3.2.3 UC-34 Switch Active Baby Profile

UC ID and Name	UC-34 Switch Active Baby Profile		
Created By	HuyND	Date Created	03/07/2026
Primary Actor	Mother	Secondary Actors	None

Trigger	Mother selects a different baby profile.
Description	Selects the baby profile used by the current baby dashboard, journals, growth views and relevant reminders.
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.
Postconditions	The requested active baby profile outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.
Normal Flow	1. Mother starts Switch Active Baby Profile. 2. The system verifies the current authorization, required context and applicable consent state. 3. Mother provides selected baby. 4. The system applies the relevant domain and safety rules, including external-service checks when needed. 5. The system completes the active baby profile action, records required audit information and presents the safe result state.
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.
Priority	P0
Frequency of Use	Regular
Business Rules	BR-HEALTH-BOUNDARY: baby records, trends and milestones support observation rather than diagnosis. BR-OWNERSHIP: baby data is private by default. BR-BABY-03: The active selection changes display context only, not ownership or permission.
Other Information	Platform: Mother Mobile App; Feature / Group: MF-03 - Mother Mobile - Baby Care Journey, Growth & Vaccination. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.3.2.4 UC-35 View Baby Care Overview

UC ID and Name	UC-35 View Baby Care Overview		
Created By	HuyND	Date Created	03/07/2026
Primary Actor	Mother	Secondary Actors	None
Trigger	Mother opens Baby Care.		
Description	Shows a baby overview with recent journals, growth, milestones, vaccination status and current care prompts.		

Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.
Postconditions	The requested baby care overview outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.
Normal Flow	1. Mother opens the View Baby Care Overview function. 2. The system validates role, ownership, active consent and visibility scope before retrieval. 3. The system retrieves the permitted baby care overview data using any embedded list controls where applicable. 4. Mother may select an authorized item, time range or permitted embedded control. 5. The system displays the authorized result with applicable source, safety or status labels.
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.
Priority	P0
Frequency of Use	Frequent
Business Rules	BR-HEALTH-BOUNDARY: baby records, trends and milestones support observation rather than diagnosis. BR-OWNERSHIP: baby data is private by default. BR-BABY-04: The overview is descriptive and does not diagnose illness or developmental disorder.
Other Information	Platform: Mother Mobile App; Feature / Group: MF-03 - Mother Mobile - Baby Care Journey, Growth & Vaccination. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.3.2.5 UC-36 Add Baby Daily Log

UC ID and Name	UC-36 Add Baby Daily Log		
Created By	HuyND	Date Created	03/07/2026
Primary Actor	Mother	Secondary Actors	None
Trigger	Mother saves a baby journal entry.		
Description	Records feeding, sleep, diaper, symptom note or other supported daily care observation for a selected baby.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.		

Postconditions	The requested baby daily log outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.
Normal Flow	1. Mother opens the Add Baby Daily Log function. 2. The system loads the required authorized context and input fields for the baby daily log. 3. Mother provides log type, time, value, note, source. 4. The system validates business rules, ownership/consent and data format. 5. The system creates the baby daily log, records required audit information and returns the resulting state.
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.
Priority	P0
Frequency of Use	Regular
Business Rules	BR-HEALTH-BOUNDARY: baby records, trends and milestones support observation rather than diagnosis. BR-OWNERSHIP: baby data is private by default. BR-BABY-05: Medicine information may be recorded only as user-entered context; the system does not prescribe or set dosage.
Other Information	Platform: Mother Mobile App; Feature / Group: MF-03 - Mother Mobile - Baby Care Journey, Growth & Vaccination. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.3.2.6 UC-37 Update or Delete Baby Daily Log

UC ID and Name	UC-37 Update or Delete Baby Daily Log		
Created By	HuyND	Date Created	03/07/2026
Primary Actor	Mother	Secondary Actors	None
Trigger	Mother edits or deletes a baby journal entry.		
Description	Corrects or removes an owner-entered baby journal entry.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.		
Postconditions	The requested baby daily log outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. Mother opens the selected baby daily log. 2. The system verifies the actor authorization, ownership and current record state. 3. Mother selects an allowed action and provides selected log and changes.		

	4. The system validates applicable policy, safety and consistency rules. 5. The system applies the change to the baby daily log, writes required audit data and refreshes related screens or notifications.
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.
Priority	P0
Frequency of Use	Occasional
Business Rules	BR-HEALTH-BOUNDARY: baby records, trends and milestones support observation rather than diagnosis. BR-OWNERSHIP: baby data is private by default. BR-BABY-06: Deletion is limited to owner-entered data and follows retention rules.
Other Information	Platform: Mother Mobile App; Feature / Group: MF-03 - Mother Mobile - Baby Care Journey, Growth & Vaccination. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.3.2.7 UC-38 View Baby Log Summary

UC ID and Name	UC-38 View Baby Log Summary		
Created By	HuyND	Date Created	03/07/2026
Primary Actor	Mother	Secondary Actors	None
Trigger	Mother opens Baby Log Summary.		
Description	Shows recent 24-hour or 7-day feeding, sleep, diaper and logged observation patterns for the active baby.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.		
Postconditions	The requested baby log summary outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. Mother opens the View Baby Log Summary function. 2. The system validates role, ownership, active consent and visibility scope before retrieval. 3. The system retrieves the permitted baby log summary data using any embedded list controls where applicable. 4. Mother may select an authorized item, time range or permitted embedded control. 5. The system displays the authorized result with applicable source, safety or status labels.		

Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.
Priority	P0
Frequency of Use	Frequent
Business Rules	BR-HEALTH-BOUNDARY: baby records, trends and milestones support observation rather than diagnosis. BR-OWNERSHIP: baby data is private by default. BR-BABY-07: Summary patterns are observational and not diagnostic.
Other Information	Platform: Mother Mobile App; Feature / Group: MF-03 - Mother Mobile - Baby Care Journey, Growth & Vaccination. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.3.2.8 UC-39 Record Development Milestone

UC ID and Name	UC-39 Record Development Milestone		
Created By	HuyND	Date Created	03/07/2026
Primary Actor	Mother	Secondary Actors	None
Trigger	Mother saves a milestone.		
Description	Records a caregiver-observed development milestone with date and note.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.		
Postconditions	The requested development milestone outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. Mother opens the Record Development Milestone function. 2. The system loads the required authorized context and input fields for the development milestone. 3. Mother provides milestone type, observed date, note. 4. The system validates business rules, ownership/consent and data format. 5. The system creates the development milestone, records required audit information and returns the resulting state.		
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.		
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.		
Priority	P0		

Frequency of Use	Regular
Business Rules	BR-HEALTH-BOUNDARY: baby records, trends and milestones support observation rather than diagnosis. BR-OWNERSHIP: baby data is private by default. BR-DEVELOPMENT-01: Milestone records support observation only and do not determine developmental disorders.
Other Information	Platform: Mother Mobile App; Feature / Group: MF-03 - Mother Mobile - Baby Care Journey, Growth & Vaccination. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.3.2.9 UC-40 Update or Delete Development Milestone

UC ID and Name	UC-40 Update or Delete Development Milestone		
Created By	HuyND	Date Created	03/07/2026
Primary Actor	Mother	Secondary Actors	None
Trigger	Mother edits or deletes a milestone.		
Description	Corrects or removes a caregiver-entered development milestone.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.		
Postconditions	The requested development milestone outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. Mother opens the selected development milestone. 2. The system verifies the actor authorization, ownership and current record state. 3. Mother selects an allowed action and provides selected milestone and changes. 4. The system validates applicable policy, safety and consistency rules. 5. The system applies the change to the development milestone, writes required audit data and refreshes related screens or notifications.		
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.		
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.		
Priority	P0		
Frequency of Use	Occasional		
Business Rules	BR-HEALTH-BOUNDARY: baby records, trends and milestones support observation rather than diagnosis. BR-OWNERSHIP: baby data is private by default. BR-DEVELOPMENT-02: The record owner may edit only their own entry.		

Other Information	Platform: Mother Mobile App; Feature / Group: MF-03 - Mother Mobile - Baby Care Journey, Growth & Vaccination. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.3.2.10 UC-41 Add Growth Measurement

UC ID and Name	UC-41 Add Growth Measurement		
Created By	HuyND	Date Created	03/07/2026
Primary Actor	Mother	Secondary Actors	None
Trigger	Mother saves a growth measurement.		
Description	Records supported growth measurements for a selected baby with measurement date and source.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.		
Postconditions	The requested growth measurement outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. Mother opens the Add Growth Measurement function. 2. The system loads the required authorized context and input fields for the growth measurement. 3. Mother provides weight, height/length, head circumference, date, source. 4. The system validates business rules, ownership/consent and data format. 5. The system creates the growth measurement, records required audit information and returns the resulting state.		
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.		
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.		
Priority	P0		
Frequency of Use	Regular		
Business Rules	BR-HEALTH-BOUNDARY: baby records, trends and milestones support observation rather than diagnosis. BR-OWNERSHIP: baby data is private by default. BR-GROWTH-01: Measurement source and time are required to interpret trends.		
Other Information	Platform: Mother Mobile App; Feature / Group: MF-03 - Mother Mobile - Baby Care Journey, Growth & Vaccination. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.		
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.		

3.3.2.11 UC-42 Update or Delete Growth Measurement

UC ID and Name	UC-42 Update or Delete Growth Measurement		
Created By	HuyND	Date Created	03/07/2026
Primary Actor	Mother	Secondary Actors	None
Trigger	Mother edits or deletes a growth measurement.		
Description	Corrects or removes a user-entered growth measurement.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.		
Postconditions	The requested growth measurement outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. Mother opens the selected growth measurement. 2. The system verifies the actor authorization, ownership and current record state. 3. Mother selects an allowed action and provides selected measurement and changes. 4. The system validates applicable policy, safety and consistency rules. 5. The system applies the change to the growth measurement, writes required audit data and refreshes related screens or notifications.		
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.		
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.		
Priority	P0		
Frequency of Use	Occasional		
Business Rules	BR-HEALTH-BOUNDARY: baby records, trends and milestones support observation rather than diagnosis. BR-OWNERSHIP: baby data is private by default. BR-GROWTH-02: The system retains audit information where a measurement affected a shared summary.		
Other Information	Platform: Mother Mobile App; Feature / Group: MF-03 - Mother Mobile - Baby Care Journey, Growth & Vaccination. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.		
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.		

3.3.2.12 UC-43 View Growth Trend and Measurement History

UC ID and Name	UC-43 View Growth Trend and Measurement History		
Created By	HuyND	Date Created	03/07/2026
Primary Actor	Mother	Secondary Actors	None
Trigger	Mother opens Growth Tracking.		

Description	Displays growth measurement history and charts with reference context and prompts to seek professional assessment when appropriate.
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.
Postconditions	The requested growth trend outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.
Normal Flow	1. Mother opens the View Growth Trend and Measurement History function. 2. The system validates role, ownership, active consent and visibility scope before retrieval. 3. The system retrieves the permitted growth trend data using any embedded list controls where applicable. 4. Mother may select an authorized item, time range or permitted embedded control. 5. The system displays the authorized result with applicable source, safety or status labels.
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.
Priority	P0
Frequency of Use	Frequent
Business Rules	BR-HEALTH-BOUNDARY: baby records, trends and milestones support observation rather than diagnosis. BR-OWNERSHIP: baby data is private by default. BR-GROWTH-03: Charts do not diagnose disease or determine a clinical growth status.
Other Information	Platform: Mother Mobile App; Feature / Group: MF-03 - Mother Mobile - Baby Care Journey, Growth & Vaccination. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.3.2.13 UC-44 Manage Vaccination Record

UC ID and Name	UC-44 Manage Vaccination Record		
Created By	HuyND	Date Created	03/07/2026
Primary Actor	Mother	Secondary Actors	None
Trigger	Mother saves or modifies a vaccination record.		
Description	Adds, updates or removes a user-entered vaccination record for the selected baby.		

Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.
Postconditions	The requested vaccination record outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.
Normal Flow	1. Mother opens the selected vaccination record. 2. The system verifies the actor authorization, ownership and current record state. 3. Mother selects an allowed action and provides vaccine, dose, date, provider note, attachment. 4. The system validates applicable policy, safety and consistency rules. 5. The system applies the change to the vaccination record, writes required audit data and refreshes related screens or notifications.
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.
Priority	P0
Frequency of Use	Regular
Business Rules	BR-HEALTH-BOUNDARY: baby records, trends and milestones support observation rather than diagnosis. BR-OWNERSHIP: baby data is private by default. BR-VACCINE-01: The digital record is a personal support record and does not replace official vaccination documentation.
Other Information	Platform: Mother Mobile App; Feature / Group: MF-03 - Mother Mobile - Baby Care Journey, Growth & Vaccination. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.3.2.14 UC-45 View Vaccination Reference Schedule and Reminder Status

UC ID and Name	UC-45 View Vaccination Reference Schedule and Reminder Status		
Created By	HuyND	Date Created	03/07/2026
Primary Actor	Mother	Secondary Actors	None
Trigger	Mother opens Vaccination Tracking.		
Description	Displays a reference vaccination schedule, completed/recorded status and linked reminder state for the selected baby.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.		

Postconditions	The requested vaccination reference schedule outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.
Normal Flow	1. Mother opens the View Vaccination Reference Schedule and Reminder Status function. 2. The system validates role, ownership, active consent and visibility scope before retrieval. 3. The system retrieves the permitted vaccination reference schedule data using any embedded list controls where applicable. 4. Mother may select an authorized item, time range or permitted embedded control. 5. The system displays the authorized result with applicable source, safety or status labels.
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.
Priority	P0
Frequency of Use	Frequent
Business Rules	BR-HEALTH-BOUNDARY: baby records, trends and milestones support observation rather than diagnosis. BR-OWNERSHIP: baby data is private by default. BR-VACCINE-02: Reference schedules are informational and users should confirm with qualified providers.
Other Information	Platform: Mother Mobile App; Feature / Group: MF-03 - Mother Mobile - Baby Care Journey, Growth & Vaccination. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.3.3 Community Q&A

3.3.3.1 UC-46 Browse Community Feed and Topics

UC ID and Name	UC-46 Browse Community Feed and Topics		
Created By	BachNQ	Date Created	04/07/2026
Primary Actor	User	Secondary Actors	None
Trigger	User opens Community.		
Description	Displays moderated questions and topic collections. Keyword search, topic chips, filters, sort and pagination are embedded controls of this browse outcome.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.		

Postconditions	The requested community feed and topics outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.
Normal Flow	1. User opens the Browse Community Feed and Topics function. 2. The system validates role, ownership, active consent and visibility scope before retrieval. 3. The system retrieves the permitted community feed and topics data using any embedded list controls where applicable. 4. User may select an authorized item, time range or permitted embedded control. 5. The system displays the authorized result with applicable source, safety or status labels.
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.
Priority	P0
Frequency of Use	Frequent
Business Rules	BR-MODERATION: visible community content is governed by moderation state, report workflow and safety policy. BR-SAFETY: community content may not be used to diagnose, prescribe or promote unsafe treatment. BR-COMMUNITY-02: Only content visible under moderation state is returned. BR-COMMUNITY-03: Search/filter is not a separate use case.
Other Information	Platform: Mobile App; Feature / Group: MF-04 - Community & Admin - Community Q&A & Moderation. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.3.3.2 UC-47 View Community Question Detail

UC ID and Name	UC-47 View Community Question Detail		
Created By	BachNQ	Date Created	04/07/2026
Primary Actor	User	Secondary Actors	None
Trigger	User opens a question from the community feed.		
Description	Shows a permitted question, its answer thread, source labels, moderation state and allowed interactions.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.		

Postconditions	The requested community question detail outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.
Normal Flow	1. User opens the View Community Question Detail function. 2. The system validates role, ownership, active consent and visibility scope before retrieval. 3. The system retrieves the permitted community question detail data using any embedded list controls where applicable. 4. User may select an authorized item, time range or permitted embedded control. 5. The system displays the authorized result with applicable source, safety or status labels.
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.
Priority	P0
Frequency of Use	Frequent
Business Rules	BR-MODERATION: visible community content is governed by moderation state, report workflow and safety policy. BR-SAFETY: community content may not be used to diagnose, prescribe or promote unsafe treatment. BR-COMMUNITY-04: Hidden or restricted content remains inaccessible to ordinary users.
Other Information	Platform: Mobile App; Feature / Group: MF-04 - Community & Admin - Community Q&A & Moderation. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.3.3.3 UC-48 Create Community Question

UC ID and Name	UC-48 Create Community Question		
Created By	BachNQ	Date Created	04/07/2026
Primary Actor	User	Secondary Actors	None
Trigger	User submits a new community question.		
Description	Creates a topic-based community question and optionally applies anonymous public display while preserving internal accountability.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.		
Postconditions	The requested community question outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		

Normal Flow	1. User opens the Create Community Question function. 2. The system loads the required authorized context and input fields for the community question. 3. User provides topic, content, optional anonymous display, stage context. 4. The system validates business rules, ownership/consent and data format. 5. The system creates the community question, records required audit information and returns the resulting state.
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.
Priority	P0
Frequency of Use	Regular
Business Rules	BR-MODERATION: visible community content is governed by moderation state, report workflow and safety policy. BR-SAFETY: community content may not be used to diagnose, prescribe or promote unsafe treatment. BR-COMMUNITY-05: Anonymous display does not remove internal traceability. BR-SAFETY-01: High-risk wording may be routed for review or safe escalation.
Other Information	Platform: Mobile App; Feature / Group: MF-04 - Community & Admin - Community Q&A & Moderation. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.3.3.4 UC-49 Edit or Delete Own Community Post

UC ID and Name	UC-49 Edit or Delete Own Community Post		
Created By	BachNQ	Date Created	04/07/2026
Primary Actor	User	Secondary Actors	None
Trigger	User chooses Edit or Delete on own post.		
Description	Edits or removes the user own community post while the post is not locked by moderation or investigation.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.		
Postconditions	The requested community post outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. User opens the selected community post. 2. The system verifies the actor authorization, ownership and current record state.		

	3. User selects an allowed action and provides selected post and changes. 4. The system validates applicable policy, safety and consistency rules. 5. The system applies the change to the community post, writes required audit data and refreshes related screens or notifications.
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.
Priority	P0
Frequency of Use	Occasional
Business Rules	BR-MODERATION: visible community content is governed by moderation state, report workflow and safety policy. BR-SAFETY: community content may not be used to diagnose, prescribe or promote unsafe treatment. BR-COMMUNITY-06: Locked, escalated or evidence-preserved posts cannot be silently altered.
Other Information	Platform: Mobile App; Feature / Group: MF-04 - Community & Admin - Community Q&A & Moderation. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.3.3.5 UC-50 Post Community Answer

UC ID and Name	UC-50 Post Community Answer		
Created By	BachNQ	Date Created	04/07/2026
Primary Actor	User	Secondary Actors	None
Trigger	User submits an answer to a permitted question.		
Description	Posts a personal-experience answer or a verified expert answer with the applicable source/role label.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.		
Postconditions	The requested community answer outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. User opens the Post Community Answer function. 2. The system loads the required authorized context and input fields for the community answer. 3. User provides answer content and optional source context. 4. The system validates business rules, ownership/consent and data format. 5. The system creates the community answer, records required audit information and returns the resulting state.		

Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.
Priority	P0
Frequency of Use	Regular
Business Rules	BR-MODERATION: visible community content is governed by moderation state, report workflow and safety policy. BR-SAFETY: community content may not be used to diagnose, prescribe or promote unsafe treatment. BR-COMMUNITY-07: Personal answers must not be presented as diagnosis or prescription. BR-EXPERT-01: Verified labels are applied only to active verified experts.
Other Information	Platform: Mobile App / Expert Portal; Feature / Group: MF-04 - Community & Admin - Community Q&A & Moderation. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.3.3.6 UC-51 Edit or Delete Own Community Answer

UC ID and Name	UC-51 Edit or Delete Own Community Answer		
Created By	BachNQ	Date Created	04/07/2026
Primary Actor	User	Secondary Actors	None
Trigger	User chooses Edit or Delete on own answer.		
Description	Edits or removes the user own answer when no moderation lock or evidence hold applies.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.		
Postconditions	The requested community answer outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. User opens the selected community answer. 2. The system verifies the actor authorization, ownership and current record state. 3. User selects an allowed action and provides selected answer and changes. 4. The system validates applicable policy, safety and consistency rules. 5. The system applies the change to the community answer, writes required audit data and refreshes related screens or notifications.		
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.		

Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.
Priority	P0
Frequency of Use	Occasional
Business Rules	BR-MODERATION: visible community content is governed by moderation state, report workflow and safety policy. BR-SAFETY: community content may not be used to diagnose, prescribe or promote unsafe treatment. BR-COMMUNITY-08: Changes to expert answers remain attributable and auditable.
Other Information	Platform: Mobile App / Expert Portal; Feature / Group: MF-04 - Community & Admin - Community Q&A & Moderation. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.3.3.7 UC-52 React to Community Answer

UC ID and Name	UC-52 React to Community Answer		
Created By	BachNQ	Date Created	04/07/2026
Primary Actor	User	Secondary Actors	None
Trigger	User taps a reaction control.		
Description	Adds or removes a permitted helpfulness reaction to a community answer.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.		
Postconditions	The requested community answer outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. User starts React to Community Answer. 2. The system verifies the current authorization, required context and applicable consent state. 3. User provides selected answer and reaction state. 4. The system applies the relevant domain and safety rules, including external-service checks when needed. 5. The system completes the community answer action, records required audit information and presents the safe result state.		
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.		
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.		

Priority	P0
Frequency of Use	Regular
Business Rules	BR-MODERATION: visible community content is governed by moderation state, report workflow and safety policy. BR-SAFETY: community content may not be used to diagnose, prescribe or promote unsafe treatment. BR-COMMUNITY-09: Reactions are participation signals, not proof of medical accuracy.
Other Information	Platform: Mobile App; Feature / Group: MF-04 - Community & Admin - Community Q&A & Moderation. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.3.3.8 UC-53 Bookmark Community Post

UC ID and Name	UC-53 Bookmark Community Post		
Created By	BachNQ	Date Created	04/07/2026
Primary Actor	User	Secondary Actors	None
Trigger	User toggles a bookmark.		
Description	Saves or removes a permitted post from the user private bookmark list.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.		
Postconditions	The requested community post outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. User starts Bookmark Community Post. 2. The system verifies the current authorization, required context and applicable consent state. 3. User provides selected post. 4. The system applies the relevant domain and safety rules, including external-service checks when needed. 5. The system completes the community post action, records required audit information and presents the safe result state.		
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.		
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.		
Priority	P0		
Frequency of Use	Regular		
Business Rules	BR-MODERATION: visible community content is governed by moderation state, report workflow and safety policy.		

	BR-SAFETY: community content may not be used to diagnose, prescribe or promote unsafe treatment. BR-PRIVACY-05: Bookmarks are private to the account unless the owner explicitly shares them.
Other Information	Platform: Mobile App; Feature / Group: MF-04 - Community & Admin - Community Q&A & Moderation. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.3.3.9 UC-54 Follow or Unfollow Community Topic

UC ID and Name	UC-54 Follow or Unfollow Community Topic		
Created By	BachNQ	Date Created	04/07/2026
Primary Actor	User	Secondary Actors	None
Trigger	User toggles topic follow state.		
Description	Subscribes or unsubscribes the user from a community topic for feed relevance and optional notifications.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.		
Postconditions	The requested community topic outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. User starts Follow or Unfollow Community Topic. 2. The system verifies the current authorization, required context and applicable consent state. 3. User provides selected topic. 4. The system applies the relevant domain and safety rules, including external-service checks when needed. 5. The system completes the community topic action, records required audit information and presents the safe result state.		
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.		
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.		
Priority	P0		
Frequency of Use	Regular		
Business Rules	BR-MODERATION: visible community content is governed by moderation state, report workflow and safety policy. BR-SAFETY: community content may not be used to diagnose, prescribe or promote unsafe treatment. BR-COMMUNITY-10: Topic notifications remain subject to notification preferences.		

Other Information	Platform: Mobile App; Feature / Group: MF-04 - Community & Admin - Community Q&A & Moderation. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.3.3.10 UC-55 Report Unsafe Community Content or Account

UC ID and Name	UC-55 Report Unsafe Community Content or Account		
Created By	BachNQ	Date Created	04/07/2026
Primary Actor	User	Secondary Actors	None
Trigger	User selects Report from a permitted target.		
Description	Submits a report for misinformation, dangerous advice, hidden advertising, harassment, privacy breach or suspicious account behavior.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.		
Postconditions	The requested content or account outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. User starts Report Unsafe Community Content or Account. 2. The system verifies the current authorization, required context and applicable consent state. 3. User provides reason, optional evidence/note. 4. The system applies the relevant domain and safety rules, including external-service checks when needed. 5. The system completes the content or account action, records required audit information and presents the safe result state.		
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.		
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.		
Priority	P0		
Frequency of Use	Regular		
Business Rules	BR-MODERATION: visible community content is governed by moderation state, report workflow and safety policy. BR-SAFETY: community content may not be used to diagnose, prescribe or promote unsafe treatment. BR-MODERATION-01: Reports must preserve the target evidence reference and reporter context. BR-PRIVACY-06: Report visibility is restricted to authorized reviewers.		
Other Information	Platform: Mobile App / Web Portal; Feature / Group: MF-04 - Community & Admin - Community Q&A & Moderation. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.		

Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.
--------------------	--

3.3.4 Verified Expert Discovery

3.3.4.1 UC-65 Browse Verified Expert Directory

UC ID and Name	UC-65 Browse Verified Expert Directory		
Created By	LamVH	Date Created	05/07/2026
Primary Actor	User	Secondary Actors	None
Trigger	User opens Expert Directory.		
Description	Displays verified experts with embedded search and filters by specialty, availability, badge and consented service area.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.		
Postconditions	The requested verified expert directory outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. User opens the Browse Verified Expert Directory function. 2. The system validates role, ownership, active consent and visibility scope before retrieval. 3. The system retrieves the permitted verified expert directory data using any embedded list controls where applicable. 4. User may select an authorized item, time range or permitted embedded control. 5. The system displays the authorized result with applicable source, safety or status labels.		
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.		
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.		
Priority	P0		
Frequency of Use	Frequent		
Business Rules	BR-EXPERT: verified status is required for verified badge and verified expert response. BR-TRUST: contribution points and badges show participation, not clinical competence. BR-EXPERT-08: Only currently verified and permitted public profiles are listed. BR-EXPERT-09: Search/filter is embedded in the directory, not a standalone use case.		
Other Information	Platform: Mobile App; Feature / Group: MF-05 - Expert & Admin - Verified Expert Network & Contribution. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.		

Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.
--------------------	--

3.3.4.2 UC-66 View Verified Expert Profile

UC ID and Name	UC-66 View Verified Expert Profile		
Created By	LamVH	Date Created	05/07/2026
Primary Actor	User	Secondary Actors	None
Trigger	User opens an expert profile.		
Description	Shows a verified expert public professional profile, specialties, supported scope, availability and contribution reputation.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.		
Postconditions	The requested verified expert profile outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. User opens the View Verified Expert Profile function. 2. The system validates role, ownership, active consent and visibility scope before retrieval. 3. The system retrieves the permitted verified expert profile data using any embedded list controls where applicable. 4. User may select an authorized item, time range or permitted embedded control. 5. The system displays the authorized result with applicable source, safety or status labels.		
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.		
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.		
Priority	P0		
Frequency of Use	Frequent		
Business Rules	BR-EXPERT: verified status is required for verified badge and verified expert response. BR-TRUST: contribution points and badges show participation, not clinical competence. BR-EXPERT-10: Contribution/review signals must not be represented as medical competence guarantees.		
Other Information	Platform: Mobile App; Feature / Group: MF-05 - Expert & Admin - Verified Expert Network & Contribution. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.		
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.		

3.3.5 AI Nurse Assistant & Risk Triage

3.3.5.1 UC-72 Run AI Symptom Intake

UC ID and Name	UC-72 Run AI Symptom Intake		
Created By	PhuongNT	Date Created	06/07/2026
Primary Actor	Mother	Secondary Actors	Gemini / Rule Engine
Trigger	Mother starts AI Symptom Intake.		
Description	Collects structured symptom information, timing, severity and permitted context through a guided intake flow.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.		
Postconditions	The requested symptom intake session outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. Mother starts Run AI Symptom Intake. 2. The system verifies the current authorization, required context and applicable consent state. 3. Mother provides symptoms, duration, severity, optional context. 4. The system applies the relevant domain and safety rules, including external-service checks when needed. 5. The system completes the symptom intake session action, records required audit information and presents the safe result state.		
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.		
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.		
Priority	P0		
Frequency of Use	Regular		
Business Rules	BR-AI: outputs are non-diagnostic and non-prescriptive. BR-SAFETY: red-flag routing must provide conservative safe next actions. BR-AI-01: Intake must state that the assistant is not diagnosing or prescribing. BR-AI-02: Only necessary consented context may be used.		
Other Information	Platform: Mother Mobile App; Feature / Group: MF-06 - Mother & Admin - AI Nurse Assistant & Risk Triage. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.		
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.		

3.3.5.2 UC-73 View Risk Triage Result

UC ID and Name	UC-73 View Risk Triage Result		
Created By	PhuongNT	Date Created	06/07/2026

Primary Actor	Mother	Secondary Actors	Gemini / Rule Engine
Trigger	The system completes a symptom intake evaluation.		
Description	Shows non-diagnostic green, yellow or red orientation with safe next-step guidance and clearly labelled uncertainty.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.		
Postconditions	The requested risk triage result outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. Mother opens the View Risk Triage Result function. 2. The system validates role, ownership, active consent and visibility scope before retrieval. 3. The system retrieves the permitted risk triage result data using any embedded list controls where applicable. 4. Mother may select an authorized item, time range or permitted embedded control. 5. The system displays the authorized result with applicable source, safety or status labels.		
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.		
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.		
Priority	P0		
Frequency of Use	Frequent		
Business Rules	BR-AI: outputs are non-diagnostic and non-prescriptive. BR-SAFETY: red-flag routing must provide conservative safe next actions. BR-AI-03: Results must not name a disease, prescribe medication or set dosage. BR-AI-04: Red results must prioritize emergency/offline care guidance.		
Other Information	Platform: Mother Mobile App; Feature / Group: MF-06 - Mother & Admin - AI Nurse Assistant & Risk Triage. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.		
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.		

3.3.5.3 UC-74 Open Emergency Support from a Red Risk Result

UC ID and Name	UC-74 Open Emergency Support from a Red Risk Result		
Created By	PhuongNT	Date Created	06/07/2026
Primary Actor	Mother	Secondary Actors	None
Trigger	Mother selects emergency support from a red-risk result.		
Description	Transfers a user from a red-risk result to the Emergency Map and Nearby Care Support flow without treating the result as a diagnosis.		

Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.
Postconditions	The requested red risk result outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.
Normal Flow	1. Mother starts Open Emergency Support from a Red Risk Result. 2. The system verifies the current authorization, required context and applicable consent state. 3. Mother provides current result context. 4. The system applies the relevant domain and safety rules, including external-service checks when needed. 5. The system completes the red risk result action, records required audit information and presents the safe result state.
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.
Priority	P0
Frequency of Use	Occasional
Business Rules	BR-AI: outputs are non-diagnostic and non-prescriptive. BR-SAFETY: red-flag routing must provide conservative safe next actions. BR-AI-05: The emergency handoff does not dispatch emergency services.
Other Information	Platform: Mother Mobile App; Feature / Group: MF-06 - Mother & Admin - AI Nurse Assistant & Risk Triage. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.3.6 Emergency Map & Nearby Care Support

3.3.6.1 UC-77 Open Emergency Map

UC ID and Name	UC-77 Open Emergency Map		
Created By	ChuongBD	Date Created	07/07/2026
Primary Actor	Mother	Secondary Actors	TrackAsia Map Service
Trigger	Mother selects Emergency Support.		
Description	Opens the shared emergency map after informing the user about location permission and the non-dispatch scope.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.		

Postconditions	The requested emergency map outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.
Normal Flow	1. Mother starts Open Emergency Map. 2. The system verifies the current authorization, required context and applicable consent state. 3. Mother provides location permission choice. 4. The system applies the relevant domain and safety rules, including external-service checks when needed. 5. The system completes the emergency map action, records required audit information and presents the safe result state.
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.
Priority	P0
Frequency of Use	Occasional
Business Rules	BR-EMERGENCY: the system supports faster action but does not dispatch emergency services or guarantee arrival. BR-LOCATION: location is consented, minimum necessary and time-limited. BR-MAP-01: Location permission is optional and purpose-specific. BR-EMERGENCY-01: CareBridge does not dispatch ambulances.
Other Information	Platform: Mother Mobile App; Feature / Group: MF-07 - Mother/Expert Mobile - Emergency Map & Nearby Care Support. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.3.6.2 UC-78 Find Nearby Care Facilities

UC ID and Name	UC-78 Find Nearby Care Facilities		
Created By	ChuongBD	Date Created	07/07/2026
Primary Actor	Mother	Secondary Actors	TrackAsia Map Service
Trigger	Mother selects Nearby Facilities.		
Description	Finds permitted nearby care facilities based on the user approved location or a selected area.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.		
Postconditions	The requested nearby care facilities outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		

Normal Flow	1. Mother opens the Find Nearby Care Facilities function. 2. The system validates role, ownership, active consent and visibility scope before retrieval. 3. The system retrieves the permitted nearby care facilities data using any embedded list controls where applicable. 4. Mother may select an authorized item, time range or permitted embedded control. 5. The system displays the authorized result with applicable source, safety or status labels.
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.
Priority	P0
Frequency of Use	Regular
Business Rules	BR-EMERGENCY: the system supports faster action but does not dispatch emergency services or guarantee arrival. BR-LOCATION: location is consented, minimum necessary and time-limited. BR-MAP-02: Results are informational and facility availability is not guaranteed.
Other Information	Platform: Mother Mobile App; Feature / Group: MF-07 - Mother/Expert Mobile - Emergency Map & Nearby Care Support. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.3.6.3 UC-79 View Route, ETA and Quick Call or Navigate

UC ID and Name	UC-79 View Route, ETA and Quick Call or Navigate		
Created By	ChuongBD	Date Created	07/07/2026
Primary Actor	Mother	Secondary Actors	TrackAsia Map Service
Trigger	Mother selects a facility action.		
Description	Shows route/ETA and enables quick call or navigation to a selected care facility through the map/device capability.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.		
Postconditions	The requested selected care facility outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. Mother starts View Route, ETA and Quick Call or Navigate.		

	2. The system verifies the current authorization, required context and applicable consent state. 3. Mother provides selected facility and destination mode. 4. The system applies the relevant domain and safety rules, including external-service checks when needed. 5. The system completes the selected care facility action, records required audit information and presents the safe result state.
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.
Priority	P0
Frequency of Use	Frequent
Business Rules	BR-EMERGENCY: the system supports faster action but does not dispatch emergency services or guarantee arrival. BR-LOCATION: location is consented, minimum necessary and time-limited. BR-MAP-03: Route/ETA is an estimate and does not guarantee arrival time.
Other Information	Platform: Mother Mobile App; Feature / Group: MF-07 - Mother/Expert Mobile - Emergency Map & Nearby Care Support. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.3.6.4 UC-80 Share Time-Limited Location and Send Family Emergency Alert

UC ID and Name	UC-80 Share Time-Limited Location and Send Family Emergency Alert		
Created By	ChuongBD	Date Created	07/07/2026
Primary Actor	Mother	Secondary Actors	Firebase Cloud Messaging
Trigger	Mother confirms Send Family Alert.		
Description	Shares minimum necessary location/context for a time-limited emergency alert to selected family recipients.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.		
Postconditions	The requested family emergency alert outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. Mother starts Share Time-Limited Location and Send Family Emergency Alert.		

	2. The system verifies the current authorization, required context and applicable consent state. 3. Mother provides selected recipients, location precision, expiry. 4. The system applies the relevant domain and safety rules, including external-service checks when needed. 5. The system completes the family emergency alert action, records required audit information and presents the safe result state.
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.
Priority	P0
Frequency of Use	Occasional
Business Rules	BR-EMERGENCY: the system supports faster action but does not dispatch emergency services or guarantee arrival. BR-LOCATION: location is consented, minimum necessary and time-limited. BR-EMERGENCY-02: Alerts contain minimum necessary context. BR-CONSENT-04: Location sharing expires automatically or on user stop.
Other Information	Platform: Mother Mobile App; Feature / Group: MF-07 - Mother/Expert Mobile - Emergency Map & Nearby Care Support. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.3.6.5 UC-81 Create or Cancel Nearby Support Request

UC ID and Name	UC-81 Create or Cancel Nearby Support Request		
Created By	ChuongBD	Date Created	07/07/2026
Primary Actor	Mother	Secondary Actors	None
Trigger	Mother submits or cancels a nearby-support request.		
Description	Creates or cancels a consented nearby-support request visible only to eligible verified experts; no expert arrival is guaranteed.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; required reference data exists; and role, ownership or consent conditions are met.		
Postconditions	The requested nearby support request outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. Mother starts Create or Cancel Nearby Support Request. 2. The system verifies the current authorization, required context and applicable consent state. 3. Mother provides support reason, location scope, expiry.		

	4. The system applies the relevant domain and safety rules, including external-service checks when needed. 5. The system completes the nearby support request action, records required audit information and presents the safe result state.
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.
Priority	P0
Frequency of Use	Regular
Business Rules	BR-EMERGENCY: the system supports faster action but does not dispatch emergency services or guarantee arrival. BR-LOCATION: location is consented, minimum necessary and time-limited. BR-EMERGENCY-03: The request is not an ambulance dispatch or guaranteed service booking. BR-EMERGENCY-04: Request visibility uses minimum necessary context.
Other Information	Platform: Mother Mobile App; Feature / Group: MF-07 - Mother/Expert Mobile - Emergency Map & Nearby Care Support. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.3.7 Personal Health Records & Source Labeling

3.3.7.1 UC-83 Add Personal Health Record and Attachment

UC ID and Name	UC-83 Add Personal Health Record and Attachment		
Created By	HuyND	Date Created	08/07/2026
Primary Actor	Mother	Secondary Actors	Firebase Storage
Trigger	Mother saves a health record.		
Description	Creates a maternal or baby health record and attaches permitted documents or images with type, date and source metadata.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; the relevant P0 identity/consent context exists; and required reference data exists.		
Postconditions	The requested health record outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. Mother opens the Add Personal Health Record and Attachment function. 2. The system loads the required authorized context and input fields for the health record. 3. Mother provides owner, category, date, source, note, attachment.		

	4. The system validates business rules, ownership/consent and data format. 5. The system creates the health record, records required audit information and returns the resulting state.
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.
Priority	P1
Frequency of Use	Regular
Business Rules	BR-RECORD: records and attachments carry source labels and are protected by ownership/consent. BR-PRIVACY: sharing is explicit, scoped and time-limited. BR-RECORD-01: Source label and owner context are required. BR-RECORD-02: File type/size and access policy are validated.
Other Information	Platform: Mother Mobile App; Feature / Group: MF-08 - Mother Mobile - Personal Health Records & Source Labeling. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.3.7.2 UC-84 Update Health Record Metadata

UC ID and Name	UC-84 Update Health Record Metadata		
Created By	HuyND	Date Created	08/07/2026
Primary Actor	Mother	Secondary Actors	None
Trigger	Mother edits a health record.		
Description	Updates title, category, date, source label, tags or note of an owner-entered health record without altering the protected file itself unless allowed.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; the relevant P0 identity/consent context exists; and required reference data exists.		
Postconditions	The requested health record metadata outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. Mother opens the selected health record metadata. 2. The system verifies the actor authorization, ownership and current record state. 3. Mother selects an allowed action and provides selected record and metadata changes. 4. The system validates applicable policy, safety and consistency rules. 5. The system applies the change to the health record metadata, writes required audit data and refreshes related screens or notifications.		

Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.
Priority	P1
Frequency of Use	Regular
Business Rules	BR-RECORD: records and attachments carry source labels and are protected by ownership/consent. BR-PRIVACY: sharing is explicit, scoped and time-limited. BR-RECORD-03: Edits are limited to records the user owns or is authorized to manage.
Other Information	Platform: Mother Mobile App; Feature / Group: MF-08 - Mother Mobile - Personal Health Records & Source Labeling. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.3.7.3 UC-85 Archive or Delete User-entered Health Record

UC ID and Name	UC-85 Archive or Delete User-entered Health Record		
Created By	HuyND	Date Created	08/07/2026
Primary Actor	Mother	Secondary Actors	None
Trigger	Mother selects Archive or Delete.		
Description	Archives or requests deletion of an owner-entered record when retention, sharing and safety-evidence obligations allow.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; the relevant P0 identity/consent context exists; and required reference data exists.		
Postconditions	The requested health record outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. Mother opens the selected health record. 2. The system verifies the actor authorization, ownership and current record state. 3. Mother selects an allowed action and provides selected record and lifecycle action. 4. The system validates applicable policy, safety and consistency rules. 5. The system applies the change to the health record, writes required audit data and refreshes related screens or notifications.		
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.		
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is		

	rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.
Priority	P1
Frequency of Use	Occasional
Business Rules	BR-RECORD: records and attachments carry source labels and are protected by ownership/consent. BR-PRIVACY: sharing is explicit, scoped and time-limited. BR-RECORD-04: Evidence/retention locks prevent silent deletion.
Other Information	Platform: Mother Mobile App; Feature / Group: MF-08 - Mother Mobile - Personal Health Records & Source Labeling. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.3.7.4 UC-86 View Health Record Timeline and Detail

UC ID and Name	UC-86 View Health Record Timeline and Detail		
Created By	HuyND	Date Created	08/07/2026
Primary Actor	Mother	Secondary Actors	None
Trigger	Mother opens Health Records.		
Description	Shows authorized maternal and baby records by time, category and source with embedded filtering inside the timeline.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; the relevant P0 identity/consent context exists; and required reference data exists.		
Postconditions	The requested health record timeline outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. Mother opens the View Health Record Timeline and Detail function. 2. The system validates role, ownership, active consent and visibility scope before retrieval. 3. The system retrieves the permitted health record timeline data using any embedded list controls where applicable. 4. Mother may select an authorized item, time range or permitted embedded control. 5. The system displays the authorized result with applicable source, safety or status labels.		
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.		
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.		
Priority	P1		
Frequency of Use	Frequent		

Business Rules	BR-RECORD: records and attachments carry source labels and are protected by ownership/consent. BR-PRIVACY: sharing is explicit, scoped and time-limited. BR-RECORD-05: Timeline queries enforce ownership and active consent at read time.
Other Information	Platform: Mother Mobile App; Feature / Group: MF-08 - Mother Mobile - Personal Health Records & Source Labeling. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.3.7.5 UC-87 Generate Health Summary

UC ID and Name	UC-87 Generate Health Summary		
Created By	HuyND	Date Created	08/07/2026
Primary Actor	Mother	Secondary Actors	None
Trigger	Mother selects Generate Summary.		
Description	Creates a selected-period summary from user-chosen maternal or baby records for personal review or permitted sharing.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; the relevant P0 identity/consent context exists; and required reference data exists.		
Postconditions	The requested health summary outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. Mother starts Generate Health Summary. 2. The system verifies the current authorization, required context and applicable consent state. 3. Mother provides record set and time range. 4. The system applies the relevant domain and safety rules, including external-service checks when needed. 5. The system completes the health summary action, records required audit information and presents the safe result state.		
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.		
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.		
Priority	P1		
Frequency of Use	Regular		
Business Rules	BR-RECORD: records and attachments carry source labels and are protected by ownership/consent. BR-PRIVACY: sharing is explicit, scoped and time-limited. BR-RECORD-06: Summaries preserve source labels and do not create diagnoses.		

Other Information	Platform: Mother Mobile App; Feature / Group: MF-08 - Mother Mobile - Personal Health Records & Source Labeling. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.3.7.6 UC-88 Share Health Summary or Selected Records Under Consent

UC ID and Name	UC-88 Share Health Summary or Selected Records Under Consent		
Created By	HuyND	Date Created	08/07/2026
Primary Actor	Mother	Secondary Actors	None
Trigger	Mother selects Share under an active or new grant.		
Description	Shares selected records or a summary only through an active MF-01 consent grant with scope and expiry.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; the relevant P0 identity/consent context exists; and required reference data exists.		
Postconditions	The requested health summary/record set outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. Mother starts Share Health Summary or Selected Records Under Consent. 2. The system verifies the current authorization, required context and applicable consent state. 3. Mother provides recipient, selected data, purpose, expiry. 4. The system applies the relevant domain and safety rules, including external-service checks when needed. 5. The system completes the health summary/record set action, records required audit information and presents the safe result state.		
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.		
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.		
Priority	P1		
Frequency of Use	Regular		
Business Rules	BR-RECORD: records and attachments carry source labels and are protected by ownership/consent. BR-PRIVACY: sharing is explicit, scoped and time-limited. BR-CONSENT-05: Sharing requires a valid consent scope and re-checks it for every access.		
Other Information	Platform: Mother Mobile App; Feature / Group: MF-08 - Mother Mobile - Personal Health Records & Source Labeling. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.		

Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.
--------------------	--

3.3.8 Reminders, Tasks & Care Plan

3.3.8.1 UC-89 Create Appointment Reminder

UC ID and Name	UC-89 Create Appointment Reminder		
Created By	BachNQ	Date Created	09/07/2026
Primary Actor	Mother	Secondary Actors	Firebase Cloud Messaging
Trigger	Mother saves an appointment reminder.		
Description	Creates a reminder for a self-entered appointment, follow-up, examination, test or other care event.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; the relevant P0 identity/consent context exists; and required reference data exists.		
Postconditions	The requested appointment reminder outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. Mother opens the Create Appointment Reminder function. 2. The system loads the required authorized context and input fields for the appointment reminder. 3. Mother provides title, date/time, recurrence, note. 4. The system validates business rules, ownership/consent and data format. 5. The system creates the appointment reminder, records required audit information and returns the resulting state.		
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.		
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.		
Priority	P1		
Frequency of Use	Regular		
Business Rules	BR-REMINDER: reminders organize user-entered/approved schedules and do not create prescriptions or formal treatment plans. BR-REMINDER-01: Reminder content is organizational and not a care order.		
Other Information	Platform: Mother Mobile App; Feature / Group: MF-09 - Mother Mobile - Reminders, Tasks & Care Plan. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.		
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.		

3.3.8.2 UC-90 Create Medicine or Vitamin Reminder

UC ID and Name	UC-90 Create Medicine or Vitamin Reminder		
Created By	BachNQ	Date Created	09/07/2026
Primary Actor	Mother	Secondary Actors	Firestore Cloud Messaging
Trigger	Mother saves a medicine/vitamin reminder.		
Description	Creates a user-entered medicine or vitamin reminder based on prior professional advice or personal routine.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; the relevant P0 identity/consent context exists; and required reference data exists.		
Postconditions	The requested medicine/vitamin reminder outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. Mother opens the Create Medicine or Vitamin Reminder function. 2. The system loads the required authorized context and input fields for the medicine/vitamin reminder. 3. Mother provides item name, time, recurrence, note. 4. The system validates business rules, ownership/consent and data format. 5. The system creates the medicine/vitamin reminder, records required audit information and returns the resulting state.		
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.		
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.		
Priority	P1		
Frequency of Use	Regular		
Business Rules	BR-REMINDER: reminders organize user-entered/approved schedules and do not create prescriptions or formal treatment plans. BR-REMINDER-02: The system does not prescribe medicine, determine dosage or alter existing advice.		
Other Information	Platform: Mother Mobile App; Feature / Group: MF-09 - Mother Mobile - Reminders, Tasks & Care Plan. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.		
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.		

3.3.8.3 UC-91 Create Vaccination Reminder

UC ID and Name	UC-91 Create Vaccination Reminder		
Created By	BachNQ	Date Created	09/07/2026
Primary Actor	Mother	Secondary Actors	Firestore Cloud Messaging

Trigger	Mother saves a vaccination reminder.
Description	Creates or confirms a reminder linked to the selected baby vaccination record or reference schedule.
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; the relevant P0 identity/consent context exists; and required reference data exists.
Postconditions	The requested vaccination reminder outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.
Normal Flow	1. Mother opens the Create Vaccination Reminder function. 2. The system loads the required authorized context and input fields for the vaccination reminder. 3. Mother provides baby, vaccine/dose label, date/time, note. 4. The system validates business rules, ownership/consent and data format. 5. The system creates the vaccination reminder, records required audit information and returns the resulting state.
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.
Priority	P1
Frequency of Use	Regular
Business Rules	BR-REMINDER: reminders organize user-entered/approved schedules and do not create prescriptions or formal treatment plans. BR-REMINDER-03: The reminder does not replace official vaccination advice or documentation.
Other Information	Platform: Mother Mobile App; Feature / Group: MF-09 - Mother Mobile - Reminders, Tasks & Care Plan. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.3.8.4 UC-92 Update, Snooze, Complete, Skip or Delete Care Reminder

UC ID and Name	UC-92 Update, Snooze, Complete, Skip or Delete Care Reminder		
Created By	BachNQ	Date Created	09/07/2026
Primary Actor	Mother	Secondary Actors	None
Trigger	Mother changes a reminder state.		
Description	Manages the lifecycle of a user reminder, including time updates, snooze, complete, skip and deletion where permitted.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; the relevant P0 identity/consent context exists; and required reference data exists.		

Postconditions	The requested care reminder outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.
Normal Flow	1. Mother opens the selected care reminder. 2. The system verifies the actor authorization, ownership and current record state. 3. Mother selects an allowed action and provides selected reminder and lifecycle state. 4. The system validates applicable policy, safety and consistency rules. 5. The system applies the change to the care reminder, writes required audit data and refreshes related screens or notifications.
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.
Priority	P1
Frequency of Use	Occasional
Business Rules	BR-REMINDER: reminders organize user-entered/approved schedules and do not create prescriptions or formal treatment plans. BR-REMINDER-04: Completion or skip records an outcome but does not imply clinical adherence.
Other Information	Platform: Mother Mobile App; Feature / Group: MF-09 - Mother Mobile - Reminders, Tasks & Care Plan. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.3.8.5 UC-93 View Today Tasks and Care Plan

UC ID and Name	UC-93 View Today Tasks and Care Plan		
Created By	BachNQ	Date Created	09/07/2026
Primary Actor	Mother	Secondary Actors	None
Trigger	Mother opens Today Tasks.		
Description	Shows due reminders, checklist items and permitted family tasks in a prioritized today view.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; the relevant P0 identity/consent context exists; and required reference data exists.		
Postconditions	The requested today tasks and care plan outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. Mother opens the View Today Tasks and Care Plan function. 2. The system validates role, ownership, active consent and visibility scope before retrieval.		

	3. The system retrieves the permitted today tasks and care plan data using any embedded list controls where applicable. 4. Mother may select an authorized item, time range or permitted embedded control. 5. The system displays the authorized result with applicable source, safety or status labels.
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.
Priority	P1
Frequency of Use	Frequent
Business Rules	BR-REMINDER: reminders organize user-entered/approved schedules and do not create prescriptions or formal treatment plans. BR-REMINDER-05: The view aggregates organizational tasks and does not create a formal treatment plan.
Other Information	Platform: Mother Mobile App; Feature / Group: MF-09 - Mother Mobile - Reminders, Tasks & Care Plan. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.3.9 Family Sync & Cooperative Care

3.3.9.1 UC-94 Create Care Group

UC ID and Name	UC-94 Create Care Group		
Created By	LamVH	Date Created	10/07/2026
Primary Actor	Mother	Secondary Actors	None
Trigger	Mother selects Create Care Group.		
Description	Creates a care group for a mother, baby or household care context with the creator as initial owner.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; the relevant P0 identity/consent context exists; and required reference data exists.		
Postconditions	The requested care group outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. Mother opens the Create Care Group function. 2. The system loads the required authorized context and input fields for the care group. 3. Mother provides group name and care context. 4. The system validates business rules, ownership/consent and data format.		

	5. The system creates the care group, records required audit information and returns the resulting state.
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.
Priority	P1
Frequency of Use	Regular
Business Rules	BR-FAMILY: family access is limited by active membership and permission scope. BR-PRIVACY: shared views show minimum necessary context. BR-FAMILY-02: The group owner controls membership and default permission scope.
Other Information	Platform: Mother Mobile App; Feature / Group: MF-10 - Mother/Family Mobile - Family Sync & Cooperative Care. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.3.9.2 UC-95 Invite or Revoke Family Member Invitation

UC ID and Name	UC-95 Invite or Revoke Family Member Invitation		
Created By	LamVH	Date Created	10/07/2026
Primary Actor	Mother	Secondary Actors	Firebase Cloud Messaging
Trigger	Mother sends or revokes an invitation.		
Description	Sends or revokes a care-group invitation using a permitted contact or invite link.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; the relevant P0 identity/consent context exists; and required reference data exists.		
Postconditions	The requested care group invitation outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. Mother opens the Invite or Revoke Family Member Invitation function. 2. The system loads the required authorized context and input fields for the care group invitation. 3. Mother provides recipient contact/link and group. 4. The system validates business rules, ownership/consent and data format. 5. The system creates the care group invitation, records required audit information and returns the resulting state.		
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.		

Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.
Priority	P1
Frequency of Use	Regular
Business Rules	BR-FAMILY: family access is limited by active membership and permission scope. BR-PRIVACY: shared views show minimum necessary context. BR-FAMILY-03: An invitation has an expiry and does not grant access before acceptance.
Other Information	Platform: Mother Mobile App; Feature / Group: MF-10 - Mother/Family Mobile - Family Sync & Cooperative Care. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.3.9.3 UC-96 Accept or Reject Care Group Invitation

UC ID and Name	UC-96 Accept or Reject Care Group Invitation		
Created By	LamVH	Date Created	10/07/2026
Primary Actor	Family Member	Secondary Actors	None
Trigger	Family member opens an invitation.		
Description	Lets a recipient accept or reject an active care-group invitation after authentication.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; the relevant P0 identity/consent context exists; and required reference data exists.		
Postconditions	The requested care group invitation outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. Family Member starts Accept or Reject Care Group Invitation. 2. The system verifies the current authorization, required context and applicable consent state. 3. Family Member provides accept/reject decision. 4. The system applies the relevant domain and safety rules, including external-service checks when needed. 5. The system completes the care group invitation action, records required audit information and presents the safe result state.		
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.		
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.		
Priority	P1		

Frequency of Use	Regular
Business Rules	BR-FAMILY: family access is limited by active membership and permission scope. BR-PRIVACY: shared views show minimum necessary context. BR-FAMILY-04: Acceptance creates membership only after group and invitation validation.
Other Information	Platform: Family Mobile App; Feature / Group: MF-10 - Mother/Family Mobile - Family Sync & Cooperative Care. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.3.9.4 UC-97 Manage Care Group Membership

UC ID and Name	UC-97 Manage Care Group Membership		
Created By	LamVH	Date Created	10/07/2026
Primary Actor	Mother / Family Member	Secondary Actors	None
Trigger	A member or owner selects a membership action.		
Description	Allows a group owner to remove a member and an individual member to leave a group where rules permit.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; the relevant P0 identity/consent context exists; and required reference data exists.		
Postconditions	The requested care group membership outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. Mother / Family Member opens the selected care group membership. 2. The system verifies the actor authorization, ownership and current record state. 3. Mother / Family Member selects an allowed action and provides selected member and action. 4. The system validates applicable policy, safety and consistency rules. 5. The system applies the change to the care group membership, writes required audit data and refreshes related screens or notifications.		
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.		
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.		
Priority	P1		
Frequency of Use	Regular		
Business Rules	BR-FAMILY: family access is limited by active membership and permission scope. BR-PRIVACY: shared views show minimum necessary context.		

	BR-FAMILY-05: Removing a member immediately ends access derived solely from group membership.
Other Information	Platform: Mother Mobile App / Family Mobile App; Feature / Group: MF-10 - Mother/Family Mobile - Family Sync & Cooperative Care. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.3.9.5 UC-98 Manage Family Permission Scope

UC ID and Name	UC-98 Manage Family Permission Scope		
Created By	LamVH	Date Created	10/07/2026
Primary Actor	Mother	Secondary Actors	None
Trigger	Mother saves a permission scope.		
Description	Defines and updates family access to selected calendar, tasks, alerts, logs, summaries or records.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; the relevant P0 identity/consent context exists; and required reference data exists.		
Postconditions	The requested family permission scope outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. Mother opens the selected family permission scope. 2. The system verifies the actor authorization, ownership and current record state. 3. Mother selects an allowed action and provides member, data scope, purpose, expiry. 4. The system validates applicable policy, safety and consistency rules. 5. The system applies the change to the family permission scope, writes required audit data and refreshes related screens or notifications.		
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.		
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.		
Priority	P1		
Frequency of Use	Regular		
Business Rules	BR-FAMILY: family access is limited by active membership and permission scope. BR-PRIVACY: shared views show minimum necessary context. BR-FAMILY-06: Family access follows the same minimum-necessary, scope and expiry principles as consent grants.		

Other Information	Platform: Mother Mobile App; Feature / Group: MF-10 - Mother/Family Mobile - Family Sync & Cooperative Care. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.3.9.6 UC-99 Create, Update or Cancel Family Care Task

UC ID and Name	UC-99 Create, Update or Cancel Family Care Task		
Created By	LamVH	Date Created	10/07/2026
Primary Actor	Mother	Secondary Actors	Firebase Cloud Messaging
Trigger	Mother saves a family care task.		
Description	Creates, edits or cancels a care task assigned to an eligible care-group member.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; the relevant P0 identity/consent context exists; and required reference data exists.		
Postconditions	The requested family care task outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. Mother opens the selected family care task. 2. The system verifies the actor authorization, ownership and current record state. 3. Mother selects an allowed action and provides assignee, due date, task details, status. 4. The system validates applicable policy, safety and consistency rules. 5. The system applies the change to the family care task, writes required audit data and refreshes related screens or notifications.		
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.		
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.		
Priority	P1		
Frequency of Use	Regular		
Business Rules	BR-FAMILY: family access is limited by active membership and permission scope. BR-PRIVACY: shared views show minimum necessary context. BR-FAMILY-07: A task reveals only the context allowed by the assignee permission scope.		
Other Information	Platform: Mother Mobile App; Feature / Group: MF-10 - Mother/Family Mobile - Family Sync & Cooperative Care. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.		

Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.
--------------------	--

3.3.9.7 UC-100 Update Assigned Task Status

UC ID and Name	UC-100 Update Assigned Task Status		
Created By	LamVH	Date Created	10/07/2026
Primary Actor	Family Member	Secondary Actors	None
Trigger	Family member changes assigned task status.		
Description	Updates the status of a task assigned to the signed-in family member, such as in progress, completed or needs help.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; the relevant P0 identity/consent context exists; and required reference data exists.		
Postconditions	The requested assigned task status outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. Family Member opens the selected assigned task status. 2. The system verifies the actor authorization, ownership and current record state. 3. Family Member selects an allowed action and provides selected task and status/note. 4. The system validates applicable policy, safety and consistency rules. 5. The system applies the change to the assigned task status, writes required audit data and refreshes related screens or notifications.		
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.		
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.		
Priority	P1		
Frequency of Use	Regular		
Business Rules	BR-FAMILY: family access is limited by active membership and permission scope. BR-PRIVACY: shared views show minimum necessary context. BR-FAMILY-08: Members may update only tasks assigned to them unless the group role permits otherwise.		
Other Information	Platform: Family Mobile App; Feature / Group: MF-10 - Mother/Family Mobile - Family Sync & Cooperative Care. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.		
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.		

3.3.9.8 UC-101 View Shared Care Calendar, Data and Alerts

UC ID and Name	UC-101 View Shared Care Calendar, Data and Alerts		
Created By	LamVH	Date Created	10/07/2026
Primary Actor	Family Member	Secondary Actors	None
Trigger	Family member opens Shared Care.		
Description	Displays calendar items, shared care data and family alerts only within the member active permission scope.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; the relevant P0 identity/consent context exists; and required reference data exists.		
Postconditions	The requested shared care context outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. Family Member opens the View Shared Care Calendar, Data and Alerts function. 2. The system validates role, ownership, active consent and visibility scope before retrieval. 3. The system retrieves the permitted shared care context data using any embedded list controls where applicable. 4. Family Member may select an authorized item, time range or permitted embedded control. 5. The system displays the authorized result with applicable source, safety or status labels.		
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.		
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.		
Priority	P1		
Frequency of Use	Frequent		
Business Rules	BR-FAMILY: family access is limited by active membership and permission scope. BR-PRIVACY: shared views show minimum necessary context. BR-FAMILY-09: Every protected read re-checks active membership and scope.		
Other Information	Platform: Family Mobile App; Feature / Group: MF-10 - Mother/Family Mobile - Family Sync & Cooperative Care. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.		
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.		

3.3.10 Verified Content & Checklist Hub

3.3.10.1 UC-102 Browse Verified Content

UC ID and Name	UC-102 Browse Verified Content		
Created By	LamVH	Date Created	11/07/2026
Primary Actor	User	Secondary Actors	None
Trigger	User opens Verified Content Hub.		
Description	Displays approved articles, FAQs and checklists with embedded keyword, stage and topic controls.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; the relevant P0 identity/consent context exists; and required reference data exists.		
Postconditions	The requested verified content hub outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. User opens the Browse Verified Content function. 2. The system validates role, ownership, active consent and visibility scope before retrieval. 3. The system retrieves the permitted verified content hub data using any embedded list controls where applicable. 4. User may select an authorized item, time range or permitted embedded control. 5. The system displays the authorized result with applicable source, safety or status labels.		
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.		
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.		
Priority	P1		
Frequency of Use	Frequent		
Business Rules	BR-CONTENT: public content must have review state, source labeling and version traceability. BR-CONTENT-01: Only approved, active content is visible to end users. BR-CONTENT-02: Search/filter remains embedded in this browse use case.		
Other Information	Platform: Mobile App; Feature / Group: MF-11 - Mobile & Admin - Verified Content & Checklist Hub. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.		
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.		

3.3.10.2 UC-103 View Verified Content Detail

UC ID and Name	UC-103 View Verified Content Detail
-----------------------	--

Created By	LamVH	Date Created	11/07/2026
Primary Actor	User	Secondary Actors	None
Trigger	User opens a verified content item.		
Description	Shows the selected article, FAQ or checklist with source, review/version state, updated date and safety notes.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; the relevant P0 identity/consent context exists; and required reference data exists.		
Postconditions	The requested verified content detail outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. User opens the View Verified Content Detail function. 2. The system validates role, ownership, active consent and visibility scope before retrieval. 3. The system retrieves the permitted verified content detail data using any embedded list controls where applicable. 4. User may select an authorized item, time range or permitted embedded control. 5. The system displays the authorized result with applicable source, safety or status labels.		
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.		
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.		
Priority	P1		
Frequency of Use	Frequent		
Business Rules	BR-CONTENT: public content must have review state, source labeling and version traceability. BR-CONTENT-03: The content view must display applicable source/review information.		
Other Information	Platform: Mobile App; Feature / Group: MF-11 - Mobile & Admin - Verified Content & Checklist Hub. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.		
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.		

3.3.11 Expense & Preparation Planner

3.3.11.1 UC-109 Add Expense Entry

UC ID and Name	UC-109 Add Expense Entry		
Created By	ChuongBD	Date Created	12/07/2026
Primary Actor	Mother	Secondary Actors	None
Trigger	Mother saves an expense entry.		

Description	Records a household preparation expense with category, amount, date and optional experience note.
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; the relevant P0 identity/consent context exists; and required reference data exists.
Postconditions	The requested expense entry outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.
Normal Flow	1. Mother opens the Add Expense Entry function. 2. The system loads the required authorized context and input fields for the expense entry. 3. Mother provides category, amount, date, note. 4. The system validates business rules, ownership/consent and data format. 5. The system creates the expense entry, records required audit information and returns the resulting state.
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.
Priority	P1
Frequency of Use	Regular
Business Rules	BR-EXPENSE: expense tracking is household preparation only, not formal accounting or billing. BR-EXPENSE-01: Expense data is for household planning only and is not a bill, claim or accounting record.
Other Information	Platform: Mother Mobile App; Feature / Group: MF-12 - Mother Mobile - Expense & Preparation Planner. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.3.11.2 UC-110 Update or Delete Expense Entry

UC ID and Name	UC-110 Update or Delete Expense Entry		
Created By	ChuongBD	Date Created	12/07/2026
Primary Actor	Mother	Secondary Actors	None
Trigger	Mother edits or deletes an expense entry.		
Description	Corrects or deletes an owner-entered preparation expense.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; the relevant P0 identity/consent context exists; and required reference data exists.		
Postconditions	The requested expense entry outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		

Normal Flow	1. Mother opens the selected expense entry. 2. The system verifies the actor authorization, ownership and current record state. 3. Mother selects an allowed action and provides selected entry and changes. 4. The system validates applicable policy, safety and consistency rules. 5. The system applies the change to the expense entry, writes required audit data and refreshes related screens or notifications.
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.
Priority	P1
Frequency of Use	Occasional
Business Rules	BR-EXPENSE: expense tracking is household preparation only, not formal accounting or billing. BR-EXPENSE-02: Only the owner can modify personal expense entries.
Other Information	Platform: Mother Mobile App; Feature / Group: MF-12 - Mother Mobile - Expense & Preparation Planner. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.3.11.3 UC-111 View Expense Summary

UC ID and Name	UC-111 View Expense Summary		
Created By	ChuongBD	Date Created	12/07/2026
Primary Actor	Mother	Secondary Actors	None
Trigger	Mother opens Expense Summary.		
Description	Shows simple expense totals by category, month or journey/baby period.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; the relevant P0 identity/consent context exists; and required reference data exists.		
Postconditions	The requested expense summary outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. Mother opens the View Expense Summary function. 2. The system validates role, ownership, active consent and visibility scope before retrieval. 3. The system retrieves the permitted expense summary data using any embedded list controls where applicable. 4. Mother may select an authorized item, time range or permitted embedded control.		

	5. The system displays the authorized result with applicable source, safety or status labels.
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.
Priority	P1
Frequency of Use	Frequent
Business Rules	BR-EXPENSE: expense tracking is household preparation only, not formal accounting or billing. BR-EXPENSE-03: Summaries are informational and not financial advice.
Other Information	Platform: Mother Mobile App; Feature / Group: MF-12 - Mother Mobile - Expense & Preparation Planner. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.3.12 Connected Device & Health Data Integration

3.3.12.1 UC-112 Connect Health Device or Platform

UC ID and Name	UC-112 Connect Health Device or Platform		
Created By	HuyND	Date Created	13/07/2026
Primary Actor	Mother	Secondary Actors	Health Connect / Gadgetbridge / Controlled Import
Trigger	Mother selects Connect Device.		
Description	Connects a supported device/platform after the user reviews permissions and selected indicator scope.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; the relevant P0 identity/consent context exists; and required reference data exists.		
Postconditions	The requested device/platform connection outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. Mother starts Connect Health Device or Platform. 2. The system verifies the current authorization, required context and applicable consent state. 3. Mother provides provider, selected permissions, indicator scope. 4. The system applies the relevant domain and safety rules, including external-service checks when needed. 5. The system completes the device/platform connection action, records required audit information and presents the safe result state.		

Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.
Priority	P1
Frequency of Use	Regular
Business Rules	BR-DEVICE: device observations are source-labeled, quality-qualified and non-clinical. BR-PRIVACY: users can disconnect and request deletion of imported data. BR-DEVICE-01: Connection requires explicit user consent. BR-DEVICE-02: Only supported indicators are requested.
Other Information	Platform: Mother Mobile App; Feature / Group: MF-13 - Mother Mobile - Connected Device & Health Data Integration. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.3.12.2 UC-113 Import or Synchronize Device Observations

UC ID and Name	UC-113 Import or Synchronize Device Observations		
Created By	HuyND	Date Created	13/07/2026
Primary Actor	Mother	Secondary Actors	Health Connect / Gadgetbridge / Controlled Import
Trigger	Mother initiates sync/import or authorized background sync runs.		
Description	Imports or synchronizes selected device observations such as heart rate, sleep, steps, SpO2, temperature or blood pressure where supported.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; the relevant P0 identity/consent context exists; and required reference data exists.		
Postconditions	The requested device observations outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. Mother starts Import or Synchronize Device Observations. 2. The system verifies the current authorization, required context and applicable consent state. 3. Mother provides selected indicators and source timestamps. 4. The system applies the relevant domain and safety rules, including external-service checks when needed. 5. The system completes the device observations action, records required audit information and presents the safe result state.		
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.		

Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.
Priority	P1
Frequency of Use	Regular
Business Rules	BR-DEVICE: device observations are source-labeled, quality-qualified and non-clinical. BR-PRIVACY: users can disconnect and request deletion of imported data. BR-DEVICE-03: Every observation retains source, time and quality metadata. BR-DEVICE-04: Device data is not treated as a clinical measurement.
Other Information	Platform: Mother Mobile App; Feature / Group: MF-13 - Mother Mobile - Connected Device & Health Data Integration. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.3.12.3 UC-114 View Device Data Trend and Quality

UC ID and Name	UC-114 View Device Data Trend and Quality		
Created By	HuyND	Date Created	13/07/2026
Primary Actor	Mother	Secondary Actors	None
Trigger	Mother opens Device Trends.		
Description	Displays imported data trends, source labels, data gaps and quality cautions for the selected indicator.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; the relevant P0 identity/consent context exists; and required reference data exists.		
Postconditions	The requested device data trend and quality outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. Mother opens the View Device Data Trend and Quality function. 2. The system validates role, ownership, active consent and visibility scope before retrieval. 3. The system retrieves the permitted device data trend and quality data using any embedded list controls where applicable. 4. Mother may select an authorized item, time range or permitted embedded control. 5. The system displays the authorized result with applicable source, safety or status labels.		
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.		
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is		

	rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.
Priority	P1
Frequency of Use	Frequent
Business Rules	BR-DEVICE: device observations are source-labeled, quality-qualified and non-clinical. BR-PRIVACY: users can disconnect and request deletion of imported data. BR-DEVICE-05: Trend views must not diagnose illness or claim medical accuracy.
Other Information	Platform: Mother Mobile App; Feature / Group: MF-13 - Mother Mobile - Connected Device & Health Data Integration. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.3.12.4 UC-115 Disconnect Device and Delete Imported Data

UC ID and Name	UC-115 Disconnect Device and Delete Imported Data		
Created By	HuyND	Date Created	13/07/2026
Primary Actor	Mother	Secondary Actors	None
Trigger	Mother selects Disconnect or Delete Imported Data.		
Description	Stops future synchronization and allows the user to request deletion of imported data according to retention rules.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; the relevant P0 identity/consent context exists; and required reference data exists.		
Postconditions	The requested device connection/imported data outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. Mother starts Disconnect Device and Delete Imported Data. 2. The system verifies the current authorization, required context and applicable consent state. 3. Mother provides connection and deletion scope. 4. The system applies the relevant domain and safety rules, including external-service checks when needed. 5. The system completes the device connection/imported data action, records required audit information and presents the safe result state.		
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.		
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.		
Priority	P1		
Frequency of Use	Occasional		

Business Rules	BR-DEVICE: device observations are source-labeled, quality-qualified and non-clinical. BR-PRIVACY: users can disconnect and request deletion of imported data. BR-DEVICE-06: Disconnect stops future import immediately. BR-DEVICE-07: Deletion respects audit/retention obligations.
Other Information	Platform: Mother Mobile App; Feature / Group: MF-13 - Mother Mobile - Connected Device & Health Data Integration. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.3.13 Smart Activity Monitoring & Safety Support

3.3.13.1 UC-116 Manage Emergency Contacts

UC ID and Name	UC-116 Manage Emergency Contacts		
Created By	BachNQ	Date Created	14/07/2026
Primary Actor	Mother	Secondary Actors	Firestore Cloud Messaging
Trigger	Mother saves emergency-contact settings.		
Description	Adds, verifies, reprioritizes or removes emergency alert recipients used by safety monitoring and emergency flows.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; the relevant P0 identity/consent context exists; and required reference data exists.		
Postconditions	The requested emergency contacts outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. Mother opens the selected emergency contacts. 2. The system verifies the actor authorization, ownership and current record state. 3. Mother selects an allowed action and provides contact, priority, verification state. 4. The system validates applicable policy, safety and consistency rules. 5. The system applies the change to the emergency contacts, writes required audit data and refreshes related screens or notifications.		
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.		
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.		
Priority	P1		
Frequency of Use	Occasional		
Business Rules	BR-SAFETY-MONITORING: monitoring is optional, consent-based and cannot guarantee fall detection or emergency response.		

	BR-SAFETY-01: Contacts must be explicitly selected for safety alert delivery.
Other Information	Platform: Mother Mobile App; Feature / Group: MF-14 - Mother Mobile - Smart Activity Monitoring & Safety Support. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.

3.3.13.2 UC-117 Configure Smart Activity Monitoring

UC ID and Name	UC-117 Configure Smart Activity Monitoring		
Created By	BachNQ	Date Created	14/07/2026
Primary Actor	Mother	Secondary Actors	None
Trigger	Mother opens Smart Safety settings.		
Description	Configures sensor consent, alert recipients, countdown, location-sharing option and monitoring conditions.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; the relevant P0 identity/consent context exists; and required reference data exists.		
Postconditions	The requested smart activity monitoring outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. Mother opens the selected smart activity monitoring. 2. The system verifies the actor authorization, ownership and current record state. 3. Mother selects an allowed action and provides consent, contacts, countdown, location option. 4. The system validates applicable policy, safety and consistency rules. 5. The system applies the change to the smart activity monitoring, writes required audit data and refreshes related screens or notifications.		
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.		
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.		
Priority	P1		
Frequency of Use	Regular		
Business Rules	BR-SAFETY-MONITORING: monitoring is optional, consent-based and cannot guarantee fall detection or emergency response. BR-SAFETY-02: Monitoring is optional and can be configured or disabled by the user. BR-SAFETY-03: The system explains device/OS limitations.		
Other Information	Platform: Mother Mobile App; Feature / Group: MF-14 - Mother Mobile - Smart Activity Monitoring & Safety Support. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.		

Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.
--------------------	--

3.3.13.3 UC-118 Enable Smart Activity Monitoring

UC ID and Name	UC-118 Enable Smart Activity Monitoring		
Created By	BachNQ	Date Created	14/07/2026
Primary Actor	Mother	Secondary Actors	Phone IMU Sensor
Trigger	Mother selects Enable Monitoring.		
Description	Starts IMU-based monitoring after configuration, permissions and consent are valid.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; the relevant P0 identity/consent context exists; and required reference data exists.		
Postconditions	The requested smart activity monitoring outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. Mother starts Enable Smart Activity Monitoring. 2. The system verifies the current authorization, required context and applicable consent state. 3. Mother provides active monitoring state. 4. The system applies the relevant domain and safety rules, including external-service checks when needed. 5. The system completes the smart activity monitoring action, records required audit information and presents the safe result state.		
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.		
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.		
Priority	P1		
Frequency of Use	Regular		
Business Rules	BR-SAFETY-MONITORING: monitoring is optional, consent-based and cannot guarantee fall detection or emergency response. BR-SAFETY-04: Sensor access is requested only while needed and permitted. BR-SAFETY-05: The system cannot guarantee continuous background monitoring.		
Other Information	Platform: Mother Mobile App; Feature / Group: MF-14 - Mother Mobile - Smart Activity Monitoring & Safety Support. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.		
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.		

3.3.13.4 UC-119 Disable Smart Activity Monitoring

UC ID and Name	UC-119 Disable Smart Activity Monitoring		
Created By	BachNQ	Date Created	14/07/2026
Primary Actor	Mother	Secondary Actors	Phone IMU Sensor
Trigger	Mother selects Disable Monitoring.		
Description	Stops IMU monitoring and prevents new candidate safety events from being created after the state change.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; the relevant P0 identity/consent context exists; and required reference data exists.		
Postconditions	The requested smart activity monitoring outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. Mother starts Disable Smart Activity Monitoring. 2. The system verifies the current authorization, required context and applicable consent state. 3. Mother provides inactive monitoring state. 4. The system applies the relevant domain and safety rules, including external-service checks when needed. 5. The system completes the smart activity monitoring action, records required audit information and presents the safe result state.		
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.		
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.		
Priority	P1		
Frequency of Use	Regular		
Business Rules	BR-SAFETY-MONITORING: monitoring is optional, consent-based and cannot guarantee fall detection or emergency response. BR-SAFETY-06: Disabling monitoring must be immediately reflected in future sensor-event processing.		
Other Information	Platform: Mother Mobile App; Feature / Group: MF-14 - Mother Mobile - Smart Activity Monitoring & Safety Support. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.		
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.		

3.3.13.5 UC-120 Handle Suspected Fall or Impact and Safety Check

UC ID and Name	UC-120 Handle Suspected Fall or Impact and Safety Check		
Created By	BachNQ	Date Created	14/07/2026

Primary Actor	Mother / Phone IMU Sensor	Secondary Actors	Phone IMU Sensor / Firebase Cloud Messaging
Trigger	A configured sensor pattern creates a candidate event.		
Description	Detects a candidate fall/impact, opens a safety countdown and records I am OK, Need Help or emergency-support handoff response.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; the relevant P0 identity/consent context exists; and required reference data exists.		
Postconditions	The requested suspected fall/impact event outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. Mother / Phone IMU Sensor starts Handle Suspected Fall or Impact and Safety Check. 2. The system verifies the current authorization, required context and applicable consent state. 3. Mother / Phone IMU Sensor provides candidate event, safety response, optional location. 4. The system applies the relevant domain and safety rules, including external-service checks when needed. 5. The system completes the suspected fall/impact event action, records required audit information and presents the safe result state.		
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.		
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.		
Priority	P1		
Frequency of Use	Regular		
Business Rules	BR-SAFETY-MONITORING: monitoring is optional, consent-based and cannot guarantee fall detection or emergency response. BR-SAFETY-07: A candidate event is not proof of a fall or medical emergency. BR-SAFETY-08: Automatic alerts use configured recipients and minimum necessary context.		
Other Information	Platform: Mother Mobile App; Feature / Group: MF-14 - Mother Mobile - Smart Activity Monitoring & Safety Support. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.		
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.		

3.3.13.6 UC-121 Review Safety Event History and Report False Positive

UC ID and Name	UC-121 Review Safety Event History and Report False Positive		
Created By	BachNQ	Date Created	14/07/2026

Primary Actor	Mother	Secondary Actors	None
Trigger	Mother opens Safety Event History or marks a false positive.		
Description	Shows prior safety events and lets the user mark a false detection with a reason to improve future threshold evaluation.		
Preconditions	The user is on the relevant authorized screen; the CareBridge backend is available; the relevant P0 identity/consent context exists; and required reference data exists.		
Postconditions	The requested safety event history outcome is completed or a clear safe result state is shown. Related records, visibility, permissions and notifications are updated where applicable. Sensitive actions are auditable.		
Normal Flow	1. Mother opens the Review Safety Event History and Report False Positive function. 2. The system validates role, ownership, active consent and visibility scope before retrieval. 3. The system retrieves the permitted safety event history data using any embedded list controls where applicable. 4. Mother may select an authorized item, time range or permitted embedded control. 5. The system displays the authorized result with applicable source, safety or status labels.		
Alternative Flows	The actor may cancel before final submission. If no eligible record or result exists, the system shows an empty state and the next allowed action without exposing unauthorized data.		
Exceptions	Access is denied for an unauthenticated, unauthorized, expired-consent or out-of-scope request. Invalid, missing, expired or conflicting data is rejected with a clear message. External-service failure must preserve a safe state and provide retry guidance without duplicating the action.		
Priority	P1		
Frequency of Use	Frequent		
Business Rules	BR-SAFETY-MONITORING: monitoring is optional, consent-based and cannot guarantee fall detection or emergency response. BR-SAFETY-09: Feedback improves system quality but does not alter retained audit evidence improperly.		
Other Information	Platform: Mother Mobile App; Feature / Group: MF-14 - Mother Mobile - Smart Activity Monitoring & Safety Support. Embedded search/filter/list controls, if needed, are part of this use case and not separately modelled.		
Assumptions	The actor has a legitimate purpose for using this function; required role, ownership and consent conditions remain valid; and CareBridge retains data according to privacy, safety and audit policies.		

4. Functional Requirements

4.1 System Functional Overview

4.1.1 Screens Flow

Explain:



Figure 32: Screens Flow

4.1.1.1 Mother Mobile App Screens Flow

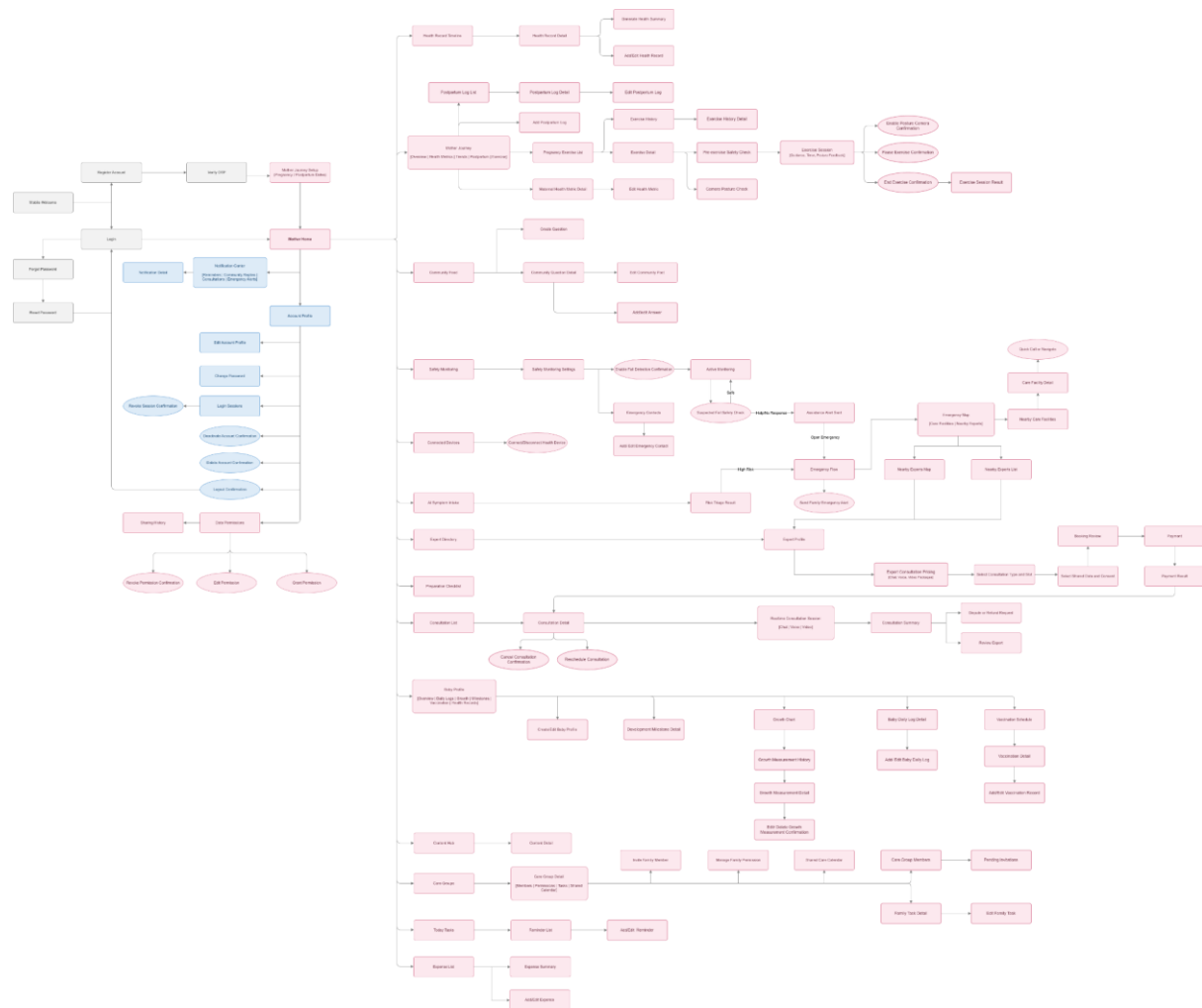


Figure 33: Mother Mobile App Screens Flow

Mother Mobile App Screens Flow

4.1.1.2 Family Member Mobile App Screens Flow

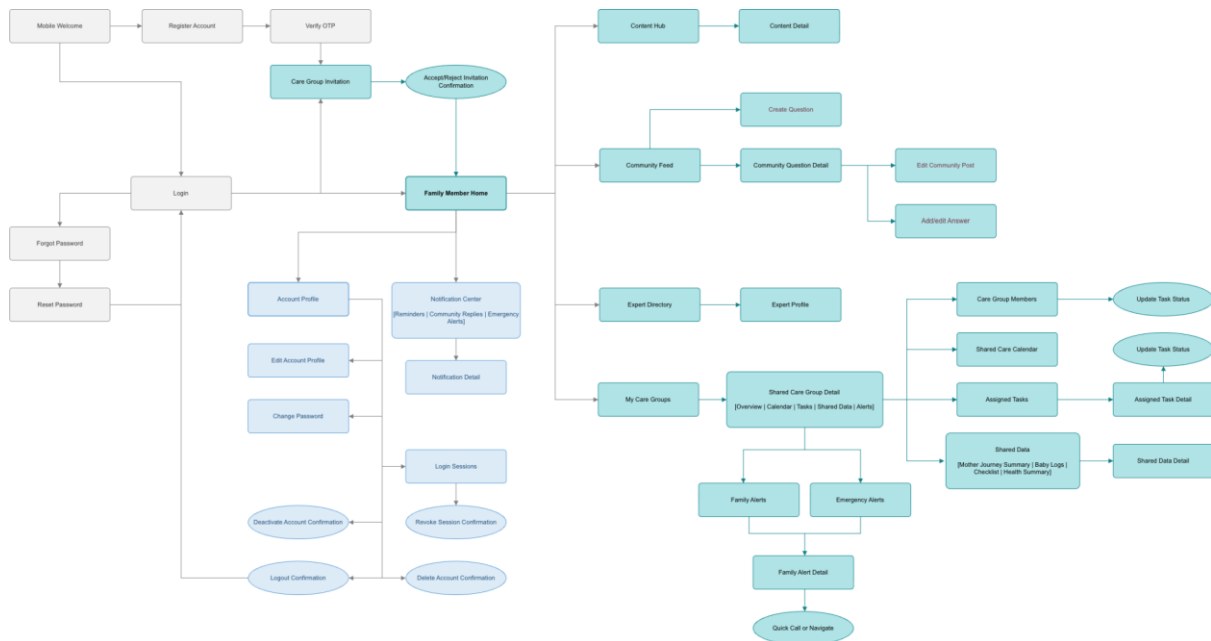


Figure 34: Family Member Mobile App Screens Flow

Family Member Mobile App Screens Flow

4.1.1.3 Verified Expert App Screens Flow

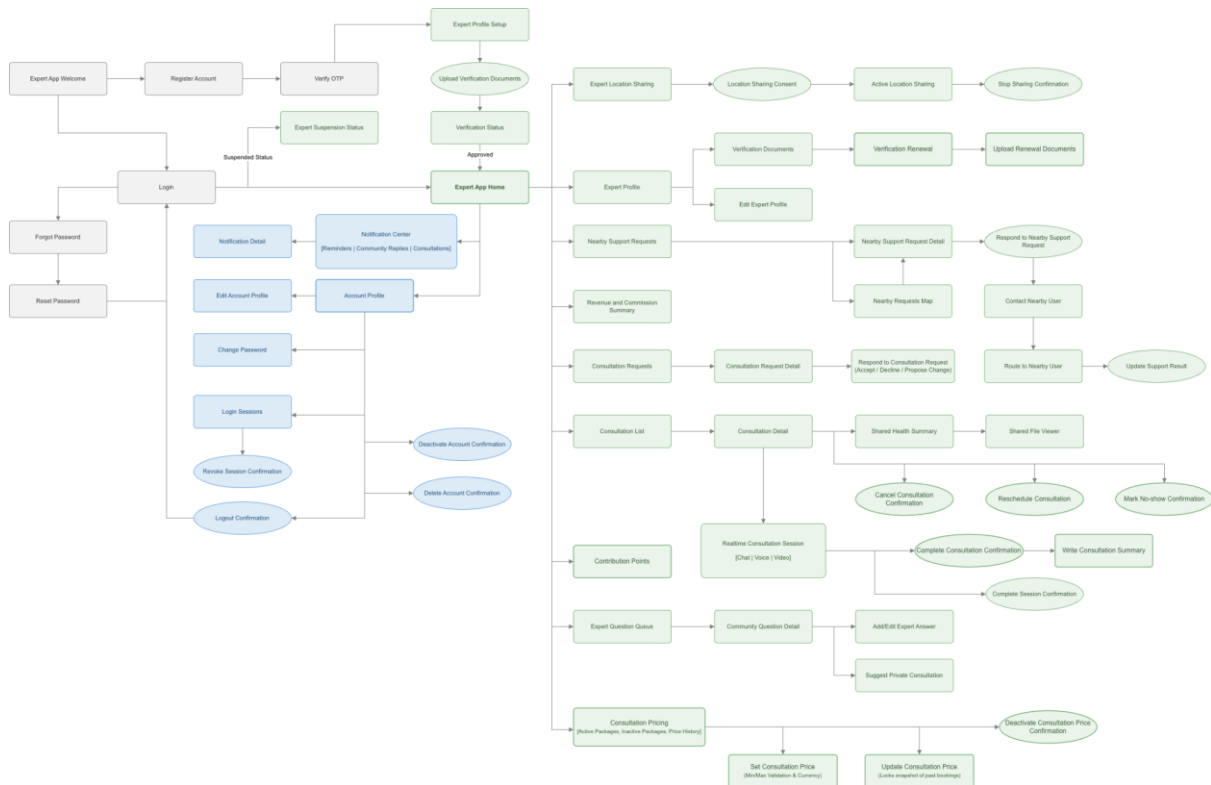


Figure 35: Verified Expert App Screens Flow

4.1.1.4 Verified Expert Web Portal Screens Flow

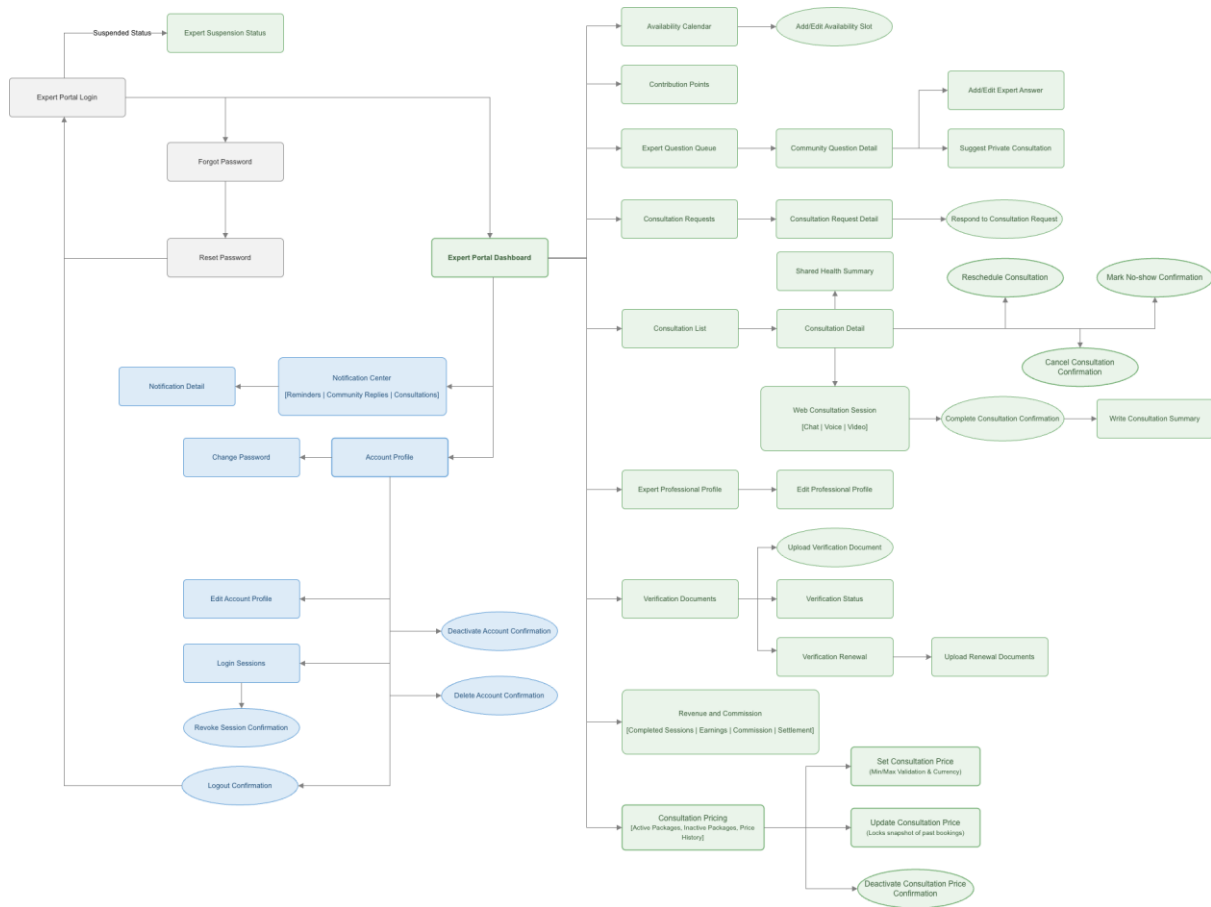


Figure 36: Verified Expert Web Portal Screens Flow

Verified Expert Web Portal Screens Flow

4.1.1.5 System Admin Web Portal Screens Flow

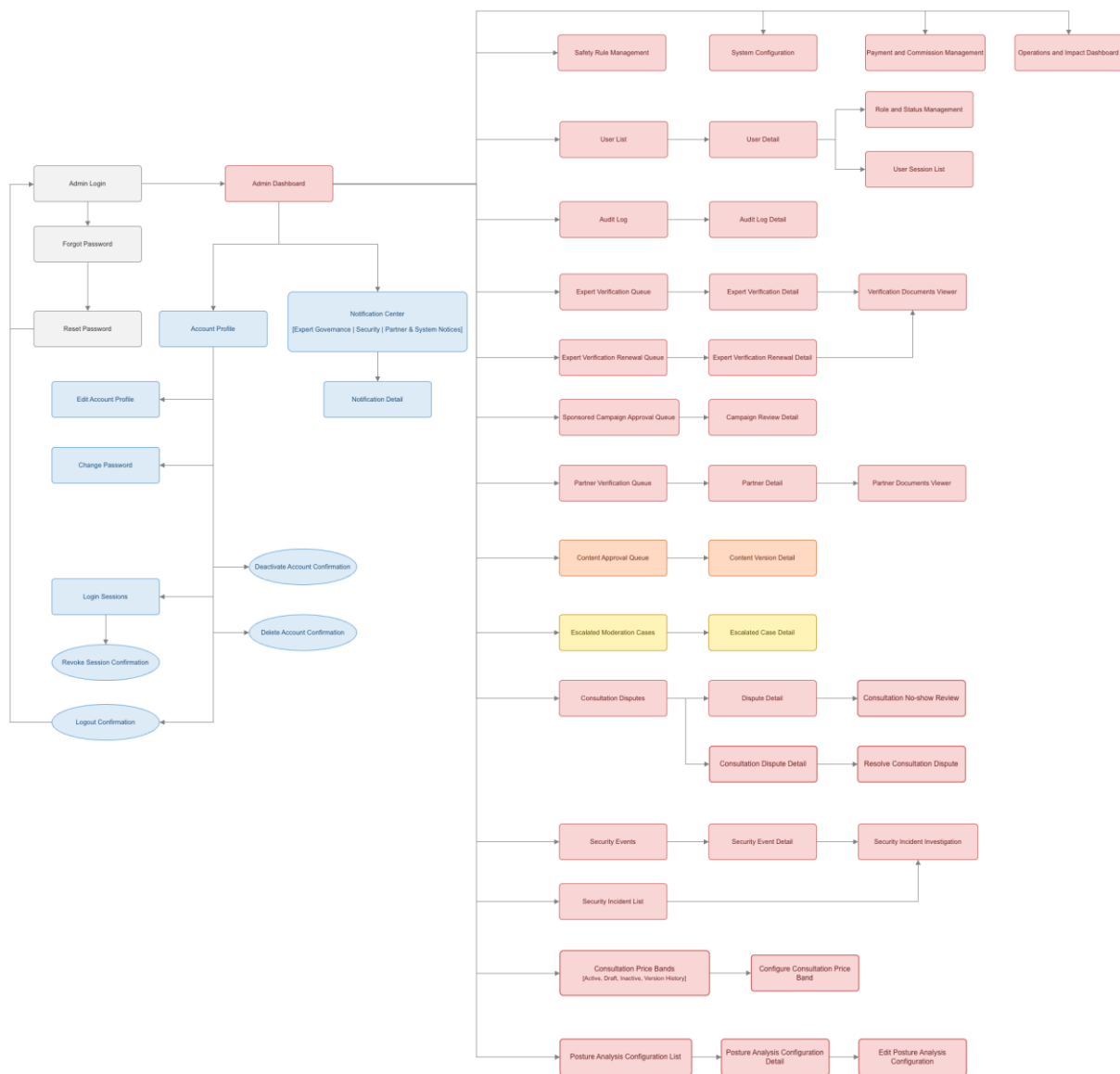


Figure 37: System Admin Web Portal Screens Flow

System Admin Web Portal Screens Flow

4.1.1.6 Content Admin Web Portal Screens Flow

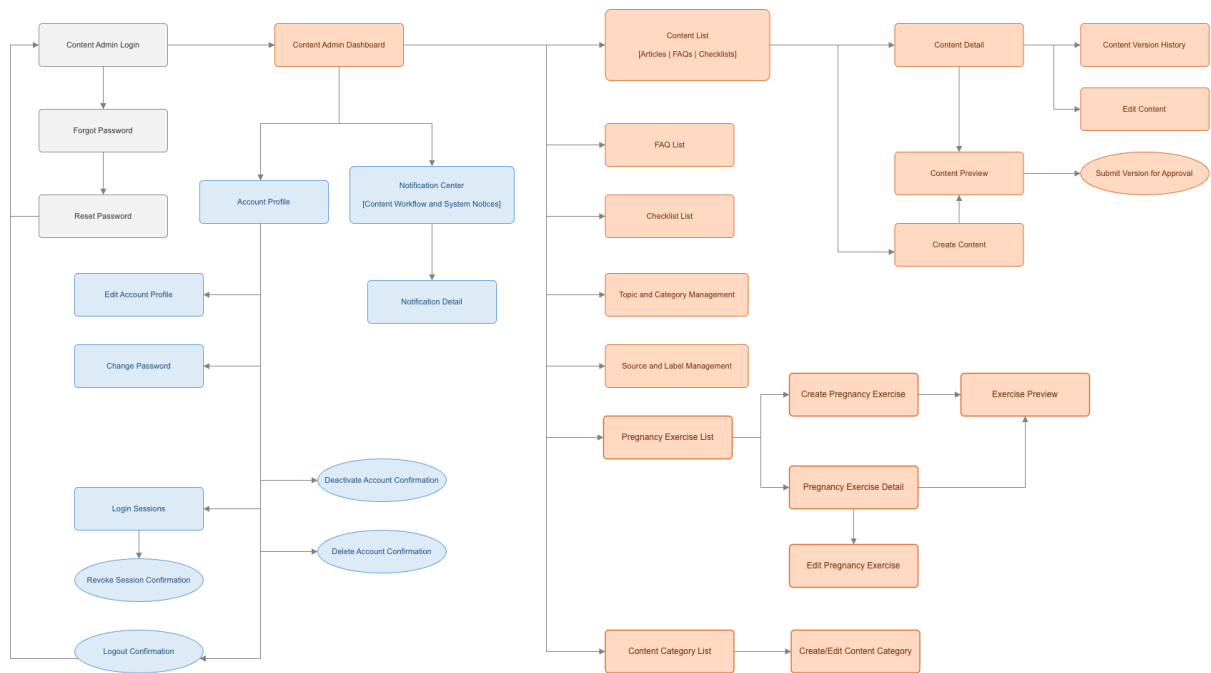


Figure 38: Content Admin Web Portal Screens Flow

Content Admin Web Portal Screens Flow

4.1.1.7 Moderator Web Portal Screens Flow

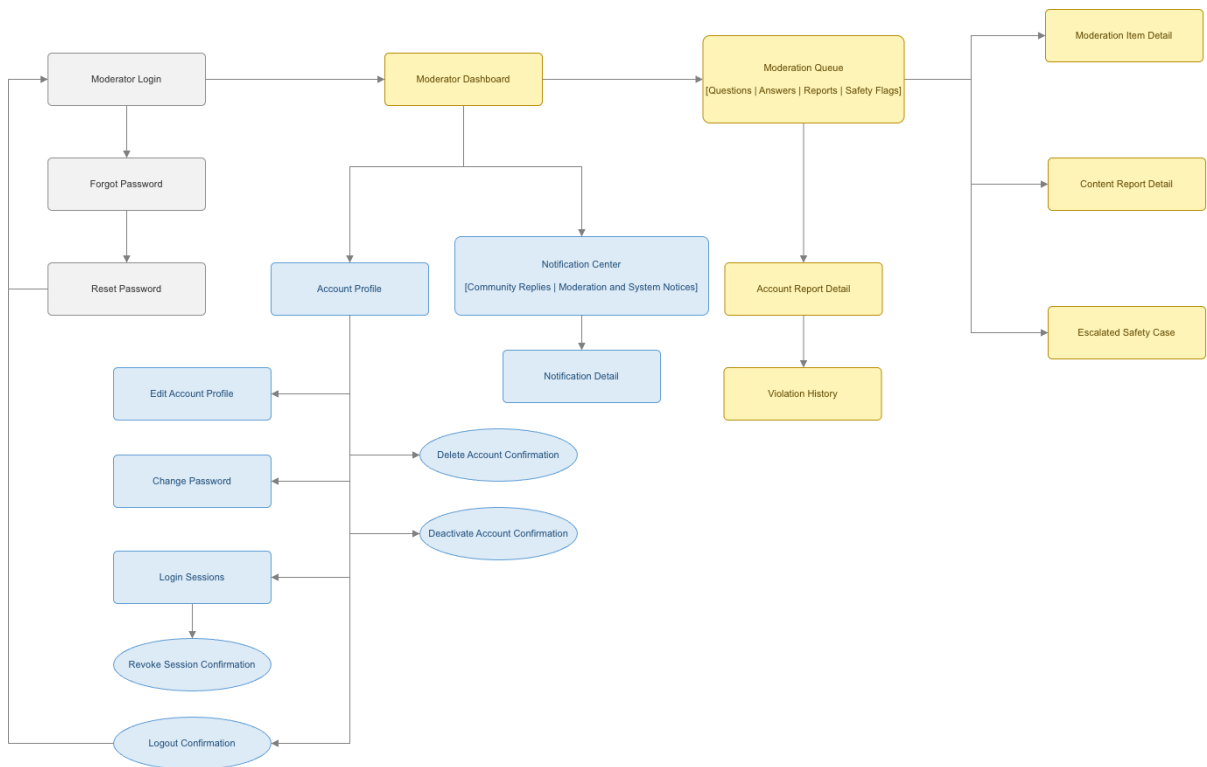


Figure 39: Moderator Web Portal Screens Flow

4.1.1.8 Partner Representative Web Portal Screens Flow

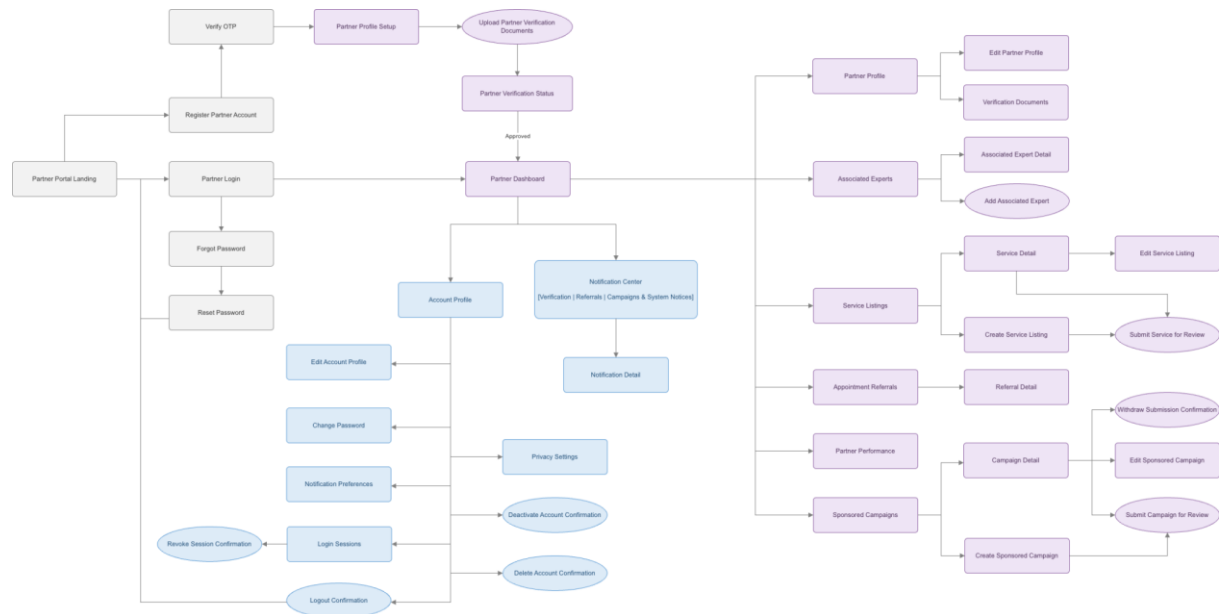


Figure 40: Partner Representative Web Portal Screens Flow

Partner Representative Web Portal Screens Flow

4.1.2 Screen Descriptions

STT	Platform	Feature	Screen	Description
1	Shared Mobile Apps	Authentication & setup	Mobile Welcome	Shared entry screen for Mother, Family Member, and Verified Expert mobile apps; role-specific routing is handled after registration or login.
2	Shared Mobile Apps	Authentication & setup	Register Account	Collects common account credentials and registration details for mobile users, with role-specific fields shown only when required.
3	Shared Mobile Apps	Authentication & setup	Verify OTP	Validates a one-time password for account activation, recovery, or sensitive account actions across mobile roles.
4	Shared Mobile Apps	Authentication & setup	Login	Authenticates mobile users and routes each role to the correct home screen.

5	Shared Mobile Apps	Authentication & setup	Forgot Password	Starts password recovery for mobile users by collecting the registered email or phone number.
6	Shared Mobile Apps	Authentication & setup	Reset Password	Allows mobile users to set a new password after a valid recovery code or link is confirmed.
7	Mother Mobile App	Authentication & setup	Mother Journey Setup	Collects the mother's current stage and key dates to initialize a personalized care journey.
8	Mother Mobile App	Core experience	Mother Home	Shows the mother's personalized overview, priorities, reminders, alerts, and shortcuts.
9	Mother Mobile App	Core experience	Mother Journey	Displays and manages pre-pregnancy, pregnancy, and postpartum journey information.
10	Mother Mobile App	Core experience	Baby Profiles	Lists the baby profiles managed by the mother and provides access to create or open a profile.
11	Mother Mobile App	Core experience	Baby Profile	Shows one baby's overview, daily logs, growth, milestones, vaccination, and health records.
12	Mother Mobile App	Core experience	Health Record Timeline	Displays maternal and baby health records chronologically with source and category filters.
13	Mother Mobile App	Core experience	Today Tasks	Shows reminders, appointments, checklist items, and family tasks due today.
14	Mother Mobile App	Core experience	Community Feed	Displays moderated community questions, answers, topics, and saved content.
15	Mother Mobile App	Core experience	AI Symptom Intake	Collects symptoms and relevant context

				through a guided, structured intake flow.
16	Mother Mobile App	Core experience	Risk Triage Result	Shows the non-diagnostic risk level and the recommended next safe action.
17	Mother Mobile App	Core experience	Emergency Map	Shows nearby care facilities, available experts, routes, ETA, quick call, and emergency actions.
18	Mother Mobile App	Core experience	Expert Directory	Lists and filters verified experts by specialty, availability, location, and rating.
19	Mother Mobile App	Core experience	Consultation Detail	Shows consultation booking information, status, shared data scope, payment, and session actions.
20	Mother Mobile App	Core experience	Realtime Consultation Session	Provides authenticated chat, voice call, or video call for an active consultation.
21	Mother Mobile App	Core experience	Care Groups	Lists and manages family care groups, members, permissions, invitations, and assigned tasks.
22	Mother Mobile App	Core experience	Connected Devices	Manages connected health platforms or devices and displays imported data status.
23	Mother Mobile App	Core experience	Safety Monitoring	Controls IMU-based activity monitoring and shows safety status, detected events, and alerts.
24	Family Member Mobile App	Authentication & invitation	Care Group Invitation	Displays an invitation and allows the family member to accept or decline joining a care group.
25	Family Member Mobile App	Authentication & invitation	Family Member Home	Shows shared care groups, assigned tasks, calendar items, and family alerts.
26	Family Member Mobile App	Shared care	My Care Groups	Lists the care groups that the family member has joined.

27	Family Member Mobile App	Shared care	Shared Care Group Detail	Shows the selected group, members, granted permissions, shared data, and tasks.
28	Family Member Mobile App	Shared care	Shared Care Calendar	Displays appointments, reminders, and care tasks shared with the family member.
29	Family Member Mobile App	Shared care	Assigned Tasks	Lists assigned care tasks and allows the family member to update task status.
30	Family Member Mobile App	Shared care	Shared Data	Displays only the maternal or baby data included in the family member's active permission scope.
31	Family Member Mobile App	Shared care	Family Alerts	Lists safety, emergency, and important care alerts shared with the family member.
32	Family Member Mobile App	Shared care	Family Alert Detail	Shows the selected alert, permitted location or context, time, and response actions.
33	Expert App	Setup & mobile operations	Expert Profile Setup	Collects the expert's specialty, experience, service scope, and public profile information.
34	Expert App	Setup & mobile operations	Upload Verification Documents	Allows the expert to upload credentials and supporting documents for verification.
35	Expert App	Setup & mobile operations	Verification Status	Shows the expert verification result, missing information, rejection reason, or approval state.
36	Expert App	Setup & mobile operations	Expert App Home	Shows mobile consultation activity, question queue, availability, notifications, and nearby support requests.
37	Expert App	Setup & mobile operations	Availability Status	Allows the expert to set online status, support methods,

				location-sharing status, and availability duration.
38	Expert App	Setup & mobile operations	Consultation Requests	Lists new mobile consultation requests and provides quick accept or decline actions.
39	Expert App	Setup & mobile operations	Consultation Detail	Shows the request, user consent scope, schedule, status, and actions to accept, reject, or join the consultation.
40	Expert App	Setup & mobile operations	Shared Health Summary	Displays the consented health summary needed for the selected mobile consultation.
41	Expert App	Setup & mobile operations	Realtime Consultation Session	Provides authenticated chat, voice call, or video call for an active consultation.
42	Expert App	Setup & mobile operations	Expert Question Queue	Lists community questions matched to the expert's verified specialties.
43	Expert App	Setup & mobile operations	Expert Location Sharing	Allows the expert to start, pause, update, or stop controlled location sharing.
44	Expert App	Setup & mobile operations	Nearby Support Requests	Lists active nearby support requests that match the expert's specialty and availability.
45	Expert App	Setup & mobile operations	Nearby Requests Map	Displays eligible nearby support requests on the shared MF-19 map.
46	Expert App	Setup & mobile operations	Nearby Support Request Detail	Shows the request summary, distance, consent scope, urgency label, and response options.
47	Expert App	Setup & mobile operations	Contact Nearby User	Provides approved chat or call actions after the expert accepts a nearby support request.

48	Expert App	Setup & mobile operations	Route to Nearby User	Displays the route and ETA to the consented user location using the map service.
49	Shared Mobile Apps	Account & privacy	Notification Center	Lists alerts and account notifications across mobile roles, filtered by the user's permissions and care-group context.
50	Shared Mobile Apps	Account & privacy	Privacy Settings	Lets mobile users manage consent, data sharing, notification privacy, and account visibility settings.
51	Shared Mobile Apps	Account & privacy	Delete Account Confirmation	Confirms permanent account deletion and explains role-specific data retention or care-group impact.
52	Shared Mobile Apps	Account & privacy	Deactivate Account Confirmation	Confirms temporary account deactivation where supported and explains access limitations while inactive.
53	Shared Web Portals	Authentication & account access	Web Login	Authenticates portal users and routes them to the correct role dashboard.
54	Expert Web Portal	Professional operations	Expert Portal Dashboard	Summarizes consultations, availability, questions, documents, revenue, and contribution activity.
55	Expert Web Portal	Professional operations	Expert Professional Profile	Displays and edits the expert's professional profile and approved service scope.
56	Expert Web Portal	Professional operations	Verification Documents	Lists submitted credentials, verification status, expiry dates, and document update actions.
57	Expert Web Portal	Professional operations	Availability Calendar	Manages recurring availability, consultation slots, exceptions, and support methods.

58	Expert Web Portal	Professional operations	Consultation Requests	Lists new consultation requests that the expert can review, accept, or decline.
59	Expert Web Portal	Professional operations	Consultation List	Lists scheduled, active, completed, cancelled, and disputed consultations.
60	Expert Web Portal	Professional operations	Consultation Detail	Shows the consultation schedule, participant, consented records, payment state, session access, and summary actions.
61	Expert Web Portal	Professional operations	Shared Health Summary	Displays consented maternal or baby health summaries and records for the selected consultation.
62	Expert Web Portal	Professional operations	Web Consultation Session	Provides the web interface for authenticated chat, voice, or video consultation.
63	Expert Web Portal	Professional operations	Expert Question Queue	Lists matched community questions and supports detailed expert response drafting.
64	Expert Web Portal	Professional operations	Revenue and Commission	Shows completed paid consultations, gross revenue, platform commission, and settlement status.
65	Expert Web Portal	Professional operations	Contribution Points	Shows contribution points, badges, and qualifying community activities.
66	Admin Web Portal	Moderation	Moderator Dashboard	Summarizes pending reports, moderation queues, violations, and escalated safety cases.
67	Admin Web Portal	Moderation	Moderation Queue	Lists reported or automatically flagged content awaiting moderator review.
68	Admin Web Portal	Moderation	Moderation Item Detail	Shows the flagged content, context, evidence, history, and available moderation actions.

69	Admin Web Portal	Moderation	Content Report Detail	Shows a report against a post or answer and supports resolution actions.
70	Admin Web Portal	Moderation	Account Report Detail	Shows a report against an account, related evidence, and enforcement options.
71	Admin Web Portal	Moderation	Violation History	Displays prior warnings, restrictions, suspensions, and resolved violations for an account.
72	Admin Web Portal	Moderation	Escalated Safety Case	Shows a high-risk content or safety case requiring urgent review and escalation.
73	Admin Web Portal	Content management	Content Admin Dashboard	Summarizes drafts, pending updates, published content, FAQs, checklists, and categories.
74	Admin Web Portal	Content management	Content List	Lists articles, FAQs, and checklists with filters for type, status, topic, and version.
75	Admin Web Portal	Content management	Content Detail	Shows the full content, metadata, source labels, publication status, and version information.
76	Admin Web Portal	Content management	Create Content	Provides the form for creating a new article, FAQ, or checklist.
77	Admin Web Portal	Content management	Edit Content	Allows authorized fields, metadata, sources, tags, and status of existing content to be updated.
78	Admin Web Portal	Content management	Content Preview	Shows how content will appear to end users before submission or publication.
79	Admin Web Portal	Content management	Content Version History	Lists previous content versions, editors, timestamps, and change summaries.
80	Admin Web Portal	Content management	FAQ List	Lists and manages verified frequently asked questions and answers.

81	Admin Web Portal	Content management	Checklist List	Lists and manages verified preparation and care checklists.
82	Admin Web Portal	Content management	Topic and Category Management	Creates and maintains content topics, categories, tags, and display order.
83	Admin Web Portal	System administration	Admin Dashboard	Summarizes users, expert and partner verification, operations, disputes, safety, and payments.
84	Admin Web Portal	System administration	User List	Lists system accounts with filters for role, status, verification, and risk indicators.
85	Admin Web Portal	System administration	User Detail	Shows account information, roles, status, sessions, reports, and permitted administration actions.
86	Admin Web Portal	System administration	Expert Verification Queue	Lists expert applications and credential updates awaiting administrative review.
87	Admin Web Portal	System administration	Content Approval Queue	Lists content versions awaiting approval, rejection, or a request for revision.
88	Admin Web Portal	System administration	Escalated Moderation Cases	Lists moderation cases escalated to the system administrator for final action.
89	Admin Web Portal	System administration	Partner Verification Queue	Lists partner applications and submitted evidence awaiting verification.
90	Admin Web Portal	System administration	Operations and Impact Dashboard	Displays operational KPIs and aggregated social-impact metrics without exposing personal health data.
91	Admin Web Portal	System administration	Safety Rule Management	Manages controlled safety, triage, escalation, and abuse-prevention rules.
92	Admin Web Portal	System administration	System Configuration	Manages system-wide settings, limits, reference data, and integration configuration.

93	Admin Web Portal	System administration	Consultation Disputes	Lists consultation complaints and supports investigation, resolution, and refund decisions.
94	Admin Web Portal	System administration	Payment and Commission Management	Manages payment records, platform commission, refunds, and settlement status.
95	Admin Web Portal	System administration	Audit Log	Displays immutable records of sensitive access, permission, moderation, and administration events.
96	Partner Web Portal	Partner operations	Partner Portal Landing	Introduces the partner portal and provides registration or login entry points.
97	Partner Web Portal	Partner operations	Register Partner Account	Collects organization and representative information to create a partner application.
98	Partner Web Portal	Partner operations	Verify OTP	Validates the OTP used to verify the partner representative's contact information.
99	Partner Web Portal	Partner operations	Partner Profile Setup	Collects organization identity, contacts, facilities, services, and verification evidence.
100	Partner Web Portal	Partner operations	Partner Verification Status	Shows the partner application status and any requested corrections or documents.
101	Partner Web Portal	Partner operations	Partner Dashboard	Summarizes profile status, associated experts, services, referrals, campaigns, and performance.
102	Partner Web Portal	Partner operations	Partner Profile	Displays and updates the approved organization profile and public partner information.
103	Partner Web Portal	Partner operations	Associated Experts	Lists experts linked to the partner and supports submission or removal requests.

104	Partner Web Portal	Partner operations	Service Listings	Lists and manages partner service or appointment-referral listings submitted for review.
105	Partner Web Portal	Partner operations	Appointment Referrals	Shows appointment referrals, their status, and permitted operational details.
106	Partner Web Portal	Partner operations	Sponsored Campaigns	Lists and manages sponsored education or support campaigns submitted for approval.
107	Partner Web Portal	Partner operations	Partner Performance	Displays aggregated referral, service, and campaign performance metrics.
108	Each actor platform	Shared account	Notifications	Lists role-specific notifications with read and unread status.
109	Each actor platform	Shared account	Notification Detail	Shows the full notification and opens the related authorized screen or action.
110	Each actor platform	Shared account	Account Profile	Displays the signed-in user's account information and account-setting shortcuts.
111	Each actor platform	Shared account	Edit Account Profile	Updates the signed-in user's own non-sensitive account information.
112	Each actor platform	Shared account	Change Password	Changes the current account password after validating the existing password.
113	Each actor platform	Shared account	Notification Preferences	Manages the notification channels and categories available to the current role.
114	Each actor platform	Shared account	Login Sessions	Lists active and recent login sessions for the current account.
115	Each actor platform	Shared account	Revoke Session Confirmation	Confirms signing out a selected device by revoking its stored session or refresh token.

116	Each actor platform	Shared account	Logout Confirmation	Confirms ending the current login session and returning to the login screen.
117	Mother Mobile App	Community Q&A	Community Search	Searches community questions by keyword, stage, topic, answer status and expert label.
118	Mother Mobile App	Community Q&A	Topic Directory	Lists and searches community topics for pregnancy, postpartum, child care, nutrition, psychology and safety.
119	Mother Mobile App	Community Q&A	Topic Detail	Shows a topic feed and allows the user to follow or unfollow that topic.
120	Mother Mobile App	Community Q&A	Edit Community Post	Edits the user's own community post while it remains editable and is not locked by moderation.
121	Mother Mobile App	Community Q&A	Delete Community Post Confirmation	Confirms deletion or archival of the user's own post when no moderation or investigation lock applies.
122	Mother Mobile App	File management	File Manager	Lists uploaded ultrasound images, medical records, vaccination files and child photos with ownership and access status.
123	Mother Mobile App	File management	File Viewer	Previews or downloads an authorized file through a protected access link.
124	Mother Mobile App	File management	Upload File	Collects file, category, owner, date and metadata before secure upload validation.
125	Mother Mobile App	File management	Delete File Confirmation	Confirms soft deletion of a user-owned file when retention or record-link rules do not block deletion.

126	Mother Mobile App	Expert discovery	Expert Search and Filters	Searches and filters verified experts by name, specialty, channel, availability, fee, rating, online state and consented distance.
127	Mother Mobile App	Safety monitoring	Emergency Contacts	Lists verified emergency contacts and their priority for MF-21 and MF-19 alert delivery.
128	Mother Mobile App	Safety monitoring	Edit Emergency Contact	Adds, verifies, reprioritizes or removes an emergency contact.
129	Mother Mobile App	Safety monitoring	Enable Fall Detection Confirmation	Confirms consent, sensor permission and monitoring conditions before fall detection starts.
130	Mother Mobile App	Safety monitoring	Disable Fall Detection Confirmation	Confirms stopping sensor monitoring and optionally retaining the configuration for later use.
131	Family Member Mobile App	Community Q&A	Community Search	Searches community questions and topics available to the family member.
132	Family Member Mobile App	Community Q&A	Topic Detail	Shows posts within a topic and allows follow or unfollow actions.
133	Family Member Mobile App	Emergency alerts	Emergency Alert Detail	Shows the minimum consented alert context, location and response actions for an emergency notification.
134	Expert App	Expert governance	Verification Renewal	Starts renewal of expert verification before credential expiry.
135	Expert App	Expert governance	Renewal Status	Shows submitted renewal documents, review status, expiry date and required follow-up.
136	Expert App	Expert governance	Expert Suspension Status	Explains a suspension, restricted capabilities, effective period and

				permitted appeal or support action.
137	Expert Web Portal	Expert governance	Verification Renewal	Submits updated credentials and tracks renewal before the current verification expires.
138	Expert Web Portal	Expert governance	Expert Suspension Status	Displays suspension reason, restricted capabilities, evidence and available appeal or support actions.
139	Admin Web Portal	Expert governance	Expert Verification Renewal Queue	Lists expert renewal submissions ordered by expiry risk, completeness and review status.
140	Admin Web Portal	Expert governance	Expert Verification Renewal Detail	Shows renewal credentials, previous verification, expiry history and approve, supplement or reject actions.
141	Admin Web Portal	Expert governance	Suspend Expert Confirmation	Confirms suspension scope, reason, effective period and impact on public listing and consultation access.
142	Admin Web Portal	Audit & security	Security Events	Lists anomalous logins, permission changes, unusual file access and sensitive-record access alerts.
143	Admin Web Portal	Audit & security	Security Event Detail	Shows event evidence, affected account or resource, risk indicators and related audit records.
144	Admin Web Portal	Audit & security	Security Incident Investigation	Manages an investigation case, evidence, timeline, assigned reviewer and containment actions.
145	Admin Web Portal	Audit & security	Security Incident Resolution	Confirms the incident outcome, corrective actions, notifications and case closure.
146	Expert App	Shared files	Shared File Viewer	Previews a health-record or consultation file only while the

				expert's consent scope and access period remain valid.
147	Expert App	Community Q&A	Search Community Questions	Searches community questions by keyword, topic, status and specialty relevance.
148	Expert App	Community Q&A	Search Community Topics	Searches and browses community topics available to the expert.
149	Expert App	Community Q&A	Topic Detail	Shows posts in a selected topic and supports following or unfollowing it.
150	Expert Web Portal	Shared files	Shared File Viewer	Displays an authorized shared file within the expert's valid consent scope.
151	Admin Web Portal	Audit & security	Security Incident List	Lists opened and resolved security investigations by severity, status, affected resource and assigned reviewer.
152	Mother Mobile App	Pregnancy Exercise & Posture Support	Pre-exercise Safety Check	Collects mandatory safety answers and blocks starting the exercise when a configured warning condition is present.
153	Mother Mobile App	Pregnancy Exercise & Posture Support	Exercise Session	Runs the selected exercise with instructions, timer, pause/resume controls and optional live posture feedback.
154	Mother Mobile App	Pregnancy Exercise & Posture Support	Enable Posture Camera Confirmation	Requests explicit consent before enabling camera-based posture analysis for the current exercise session.
155	Mother Mobile App	Pregnancy Exercise & Posture Support	Exercise Session Result	Shows duration, completion status, aggregate posture score, common posture issues and safety warnings from the completed session.
156	Mother Mobile App	Pregnancy Exercise & Posture Support	Exercise History	Lists stored exercise sessions and opens the

				result of a selected session.
157	Mother Mobile App	Mother Care Journey	Maternal Health Metric Detail	Shows one maternal metric with value, time, source and note, and provides authorized edit or delete actions.
158	Mother Mobile App	Mother Care Journey	Postpartum Log List	Lists postpartum recovery logs chronologically with filters by symptom or log type.
159	Mother Mobile App	Mother Care Journey	Postpartum Log Detail	Shows one postpartum log and provides edit or soft-delete actions for the owner.
160	Mother Mobile App	Baby Care Journey	Switch Active Baby Selector	Allows a mother who manages multiple babies to select the profile used by dashboards, logs and reminders.
161	Mother Mobile App	Baby Care Journey	Baby Daily Log Detail	Shows one feeding, sleep, diaper, symptom or care entry and provides authorized update or delete actions.
162	Mother Mobile App	Baby Care Journey	Development Milestone Detail	Shows a recorded milestone, date, note and status and provides update or soft-delete actions.
163	Mother Mobile App	Community Q&A	My Answer Detail	Shows the mother's own answer and provides edit or soft-delete actions while moderation rules allow changes.
164	Mother Mobile App	Direct Consultation	Reschedule Consultation	Collects a proposed replacement slot and confirms the schedule change under the consultation policy.
165	Mother Mobile App	Direct Consultation	Cancel Consultation Confirmation	Confirms cancellation, reason and any fee or timing consequence before changing consultation status.

166	Mother Mobile App	Direct Consultation	Consultation Summary	Shows the expert's post-session summary and safe follow-up steps for the completed consultation.
167	Mother Mobile App	Health Records	Reminder Detail	Shows reminder type, recurrence, next time, status and notes with complete, skip, edit and delete actions.
168	Mother Mobile App	Family Sync	Care Group Members	Lists members, roles, invitation status and permission summary for a selected care group.
169	Mother Mobile App	Family Sync	Pending Invitations	Lists unaccepted family invitations and allows the mother to revoke an invitation.
170	Mother Mobile App	Family Sync	Family Task Detail	Shows assignee, due date, status and notes of a family task and supports update or cancellation.
171	Mother Mobile App	Verified Content	Verified Content Search	Searches approved articles, FAQs and checklists by keyword, care stage and topic.
172	Mother Mobile App	Verified Content	Verified Content Detail	Shows approved content, source, version, update date and applicable safety notes.
173	Mother Mobile App	Vaccination & Growth	Vaccination Detail	Shows one scheduled or completed vaccination item with date, status, facility, notes and supporting file.
174	Mother Mobile App	Vaccination & Growth	Add Vaccination Record	Collects vaccination date, vaccine, facility, notes and optional proof file.
175	Mother Mobile App	Vaccination & Growth	Growth Measurement History	Lists weight, height and head-circumference measurements used by the growth chart.

176	Mother Mobile App	Vaccination & Growth	Growth Measurement Detail	Shows one growth measurement and provides authorized edit or soft-delete actions.
177	Family Member Mobile App	Family Sync	Reject Invitation Confirmation	Confirms declining a care-group invitation without creating membership or data access.
178	Family Member Mobile App	Family Sync	Leave Care Group Confirmation	Confirms leaving a care group and ending future access to its shared data.
179	Family Member Mobile App	Community Q&A	My Answer Detail	Shows the family member's own answer and allows editing or soft deletion when permitted.
180	Family Member Mobile App	Verified Content	Verified Content Search	Searches approved articles, FAQs and checklists available to the family member.
181	Family Member Mobile App	Verified Content	Verified Content Detail	Shows approved content with source, version and safety information.
182	Expert App	Community Q&A	My Expert Answer Detail	Shows the expert's own community answer and supports authorized edit or soft-delete actions.
183	Expert App	Direct Consultation	Reschedule Consultation	Allows the expert to propose or confirm a new consultation time before the policy deadline.
184	Expert App	Direct Consultation	Mark No-show Confirmation	Confirms recording a participant no-show after the required waiting period and session evidence.
185	Expert App	Direct Consultation	Complete Consultation Confirmation	Confirms ending the active consultation before writing the post-session summary.
186	Expert Web Portal	Community Q&A	My Expert Answer Detail	Shows an expert-authored answer and provides edit or soft-delete actions when permitted.

187	Expert Web Portal	Direct Consultation	Reschedule Consultation	Allows the expert to propose or confirm a replacement consultation slot.
188	Expert Web Portal	Direct Consultation	Mark No-show Confirmation	Confirms a no-show decision after checking timing and technical session evidence.
189	Content Admin Web Portal	Pregnancy Exercise Management	Pregnancy Exercise List	Lists pregnancy exercises with trimester, difficulty, status and posture-analysis filters.
190	Content Admin Web Portal	Pregnancy Exercise Management	Pregnancy Exercise Detail	Shows exercise instructions, duration, safety notes, eligibility metadata and publication state.
191	Content Admin Web Portal	Pregnancy Exercise Management	Create Pregnancy Exercise	Creates a pregnancy exercise with instructions, trimester, difficulty, duration and safety guidance.
192	Content Admin Web Portal	Pregnancy Exercise Management	Edit Pregnancy Exercise	Updates an existing exercise without adding medical auto-prescription or video-storage behavior.
193	Content Admin Web Portal	Pregnancy Exercise Management	Exercise Preview	Shows the mobile presentation of an exercise before activation.
194	Content Admin Web Portal	Verified Content	Content Category List	Lists, orders and manages active or hidden verified-content categories.
195	Content Admin Web Portal	Verified Content	Unpublish Content Confirmation	Confirms removal of an approved content version from user visibility while retaining version history.
196	Admin Web Portal	Pregnancy Exercise Management	Posture Analysis Configuration List	Lists rule/model configurations, versions, status and linked exercises.
197	Admin Web Portal	Pregnancy Exercise Management	Posture Analysis Configuration Detail	Shows analysis mode, version, thresholds, feedback severity and linked exercise scope.

198	Admin Web Portal	Pregnancy Exercise Management	Edit Posture Analysis Configuration	Updates posture-analysis rules, model version and confidence thresholds under controlled versioning.
199	Admin Web Portal	Direct Consultation	Consultation No-show Review	Shows consultation timing and technical evidence used to review or confirm a no-show.
200	Admin Web Portal	Direct Consultation	Resolve Consultation Dispute	Records the administrative dispute outcome and routes the approved financial decision.
201	Admin Web Portal	Direct Consultation	Approve Refund Confirmation	Confirms the approved refund amount and transaction action.
202	Admin Web Portal	Direct Consultation	Reject Refund Confirmation	Confirms rejection of the refund request with a recorded reason.
203	Mother Mobile App	Expert Consultation & Pricing	Expert Profile	Shows the expert's verified identity, specialties, service channels, availability, reviews and entry point to the effective consultation pricing.
204	Mother Mobile App	Expert Consultation & Pricing	Expert Consultation Pricing	Lists the expert's active chat, voice and video packages with duration, effective price, total payable estimate and cancellation/refund policy before booking.
205	Mother Mobile App	Expert Consultation & Pricing	Booking Review	Shows the selected expert, slot, shared-data consent and immutable booking price snapshot before final confirmation.
206	Mother Mobile App	Expert Consultation & Pricing	Payment	Collects the payment method and charges the amount stored in the booking price snapshot rather than the expert's current public price.

207	Mother Mobile App	Expert Consultation & Pricing	Payment Result	Shows transaction status, booking identifier, locked amount and safe retry or return actions.
208	Expert App	Consultation Pricing	Consultation Pricing	Lists the expert's active and inactive consultation packages and read-only price history by channel and duration.
209	Expert App	Consultation Pricing	Set Consultation Price	Creates the initial price for a supported channel and duration after validating it against the active CareBridge price band.
210	Expert App	Consultation Pricing	Update Consultation Price	Creates a new future-effective price version while preserving price history and existing booking snapshots.
211	Expert App	Consultation Pricing	Deactivate Consultation Price Confirmation	Confirms stopping new bookings for a selected package without changing existing confirmed or paid bookings.
212	Expert Web Portal	Consultation Pricing	Consultation Pricing	Provides detailed package management, status filters and price-version history for the signed-in expert.
213	Expert Web Portal	Consultation Pricing	Set Consultation Price	Creates a package price within the active minimum and maximum for the selected channel and duration.
214	Expert Web Portal	Consultation Pricing	Update Consultation Price	Schedules a new effective price for future bookings and retains the previous version for audit and reconciliation.
215	Expert Web Portal	Consultation Pricing	Deactivate Consultation Price Confirmation	Confirms deactivation of a package for future bookings while preserving existing booking obligations.

216	Admin Web Portal	Consultation Pricing & Commission	Consultation Price Bands	Lists active, draft, inactive and historical price-band versions by channel, duration, limits, commission rate and effective period.
217	Admin Web Portal	Consultation Pricing & Commission	Configure Consultation Price Band	Creates or updates a versioned minimum, maximum, commission rate and effective period for a consultation channel and duration.
218	Admin Web Portal	Consultation Pricing & Commission	Deactivate Price Band Confirmation	Confirms deactivation of a price band for new price changes without retroactively changing locked bookings.
219	Shared Web Portals	Account & privacy	Notification Center	Lists operational, moderation, account, and partner notifications according to the signed-in portal role.
220	Shared Web Portals	Account & privacy	Privacy Settings	Lets portal users manage account privacy and security preferences available to their role.
221	Shared Web Portals	Account & privacy	Delete Account Confirmation	Confirms permanent account deletion for portal roles where deletion is allowed.
222	Shared Web Portals	Account & privacy	Deactivate Account Confirmation	Confirms temporary account deactivation for portal roles where deactivation is allowed.

Table 261: Screen Descriptions

4.1.3 Screen Authorization

Screen	Guest	Mother	Family Member	Verified Expert	Moderator	Content Admin	Admin	Partner
Mobile Welcome (Shared Mobile Apps)	X	X	X	X				

Register Account (Shared Mobile Apps)	X	X	X	X				
Verify OTP (Shared Mobile Apps)	X	X	X	X				
Login (Shared Mobile Apps)	X	X	X	X				
Forgot Password (Shared Mobile Apps)	X	X	X	X				
Reset Password (Shared Mobile Apps)	X	X	X	X				
Mother Journey Setup (Mother Mobile App)		X						
Mother Home (Mother Mobile App)		X						
Mother Journey (Mother Mobile App)		X						
Baby Profiles (Mother Mobile App)		X						
Baby Profile (Mother Mobile App)		X						
Health Record Timeline (Mother Mobile App)		X						
Today Tasks (Mother Mobile App)		X						
Community Feed (Mother Mobile App)		X						
AI Symptom Intake (Mother Mobile App)		X						
Risk Triage Result (Mother Mobile App)		X						
Emergency Map (Mother Mobile App)		X						
Expert Directory (Mother Mobile App)		X						

Consultation Detail (Mother Mobile App)		X						
Realtime Consultation Session (Mother Mobile App)		X						
Care Groups (Mother Mobile App)		X						
Connected Devices (Mother Mobile App)		X						
Safety Monitoring (Mother Mobile App)		X						
Care Group Invitation (Family Member Mobile App)	X		X					
Family Member Home (Family Member Mobile App)			X					
My Care Groups (Family Member Mobile App)			X					
Shared Care Group Detail (Family Member Mobile App)			X					
Shared Care Calendar (Family Member Mobile App)			X					
Assigned Tasks (Family Member Mobile App)			X					
Shared Data (Family Member Mobile App)			X					
Family Alerts (Family Member Mobile App)			X					
Family Alert Detail (Family Member Mobile App)			X					
Expert Profile Setup (Expert App)				X				

Upload Verification Documents (Expert App)				X				
Verification Status (Expert App)				X				
Expert App Home (Expert App)				X				
Availability Status (Expert App)				X				
Consultation Requests (Expert App)				X				
Consultation Detail (Expert App)				X				
Shared Health Summary (Expert App)				X				
Realtime Consultation Session (Expert App)				X				
Expert Question Queue (Expert App)				X				
Expert Location Sharing (Expert App)				X				
Nearby Support Requests (Expert App)				X				
Nearby Requests Map (Expert App)				X				
Nearby Support Request Detail (Expert App)				X				
Contact Nearby User (Expert App)				X				
Route to Nearby User (Expert App)				X				
Notification Center (Shared Mobile Apps)		X	X	X				
Privacy Settings (Shared Mobile Apps)		X	X	X				
Delete Account Confirmation		X	X	X				

(Shared Mobile Apps)								
Deactivate Account Confirmation (Shared Mobile Apps)				X				
Web Login (Shared Web Portals)	X			X	X	X	X	X
Expert Portal Dashboard (Expert Web Portal)				X				
Expert Professional Profile (Expert Web Portal)				X				
Verification Documents (Expert Web Portal)				X				
Availability Calendar (Expert Web Portal)				X				
Consultation Requests (Expert Web Portal)				X				
Consultation List (Expert Web Portal)				X				
Consultation Detail (Expert Web Portal)				X				
Shared Health Summary (Expert Web Portal)				X				
Web Consultation Session (Expert Web Portal)				X				
Expert Question Queue (Expert Web Portal)				X				
Revenue and Commission (Expert Web Portal)				X				
Contribution Points (Expert Web Portal)				X				

Moderator Dashboard (Admin Web Portal)					X		X	
Moderation Queue (Admin Web Portal)					X		X	
Moderation Item Detail (Admin Web Portal)					X		X	
Content Report Detail (Admin Web Portal)					X		X	
Account Report Detail (Admin Web Portal)					X		X	
Violation History (Admin Web Portal)					X		X	
Escalated Safety Case (Admin Web Portal)					X		X	
Content Admin Dashboard (Admin Web Portal)						X	X	
Content List (Admin Web Portal)						X	X	
Content Detail (Admin Web Portal)						X	X	
Create Content (Admin Web Portal)						X	X	
Edit Content (Admin Web Portal)						X	X	
Content Preview (Admin Web Portal)						X	X	
Content Version History (Admin Web Portal)						X	X	
FAQ List (Admin Web Portal)						X	X	
Checklist List (Admin Web Portal)						X	X	
Topic and Category						X	X	

Management (Admin Web Portal)								
Admin Dashboard (Admin Web Portal)							X	
User List (Admin Web Portal)							X	
User Detail (Admin Web Portal)							X	
Expert Verification Queue (Admin Web Portal)							X	
Content Approval Queue (Admin Web Portal)							X	
Escalated Moderation Cases (Admin Web Portal)							X	
Partner Verification Queue (Admin Web Portal)							X	
Operations and Impact Dashboard (Admin Web Portal)							X	
Safety Rule Management (Admin Web Portal)							X	
System Configuration (Admin Web Portal)							X	
Consultation Disputes (Admin Web Portal)							X	
Payment and Commission Management (Admin Web Portal)							X	
Audit Log (Admin Web Portal)							X	

Partner Portal Landing (Partner Web Portal)	X							X
Register Partner Account (Partner Web Portal)	X							X
Verify OTP (Partner Web Portal)	X							X
Partner Profile Setup (Partner Web Portal)	X							X
Partner Verification Status (Partner Web Portal)	X							X
Partner Dashboard (Partner Web Portal)								X
Partner Profile (Partner Web Portal)								X
Associated Experts (Partner Web Portal)								X
Service Listings (Partner Web Portal)								X
Appointment Referrals (Partner Web Portal)								X
Sponsored Campaigns (Partner Web Portal)								X
Partner Performance (Partner Web Portal)								X
Notifications (Each actor platform)		X	X	X	X	X	X	X
Notification Detail (Each actor platform)		X	X	X	X	X	X	X
Account Profile (Each actor platform)		X	X	X	X	X	X	X
Edit Account Profile (Each actor platform)		X	X	X	X	X	X	X

Change Password (Each actor platform)		X	X	X	X	X	X	X
Notification Preferences (Each actor platform)		X	X	X	X	X	X	X
Login Sessions (Each actor platform)		X	X	X	X	X	X	X
Revoke Session Confirmation (Each actor platform)		X	X	X	X	X	X	X
Logout Confirmation (Each actor platform)		X	X	X	X	X	X	X
Community Search (Mother Mobile App)		X						
Topic Directory (Mother Mobile App)		X						
Topic Detail (Mother Mobile App)		X						
Edit Community Post (Mother Mobile App)		X						
Delete Community Post Confirmation (Mother Mobile App)		X						
File Manager (Mother Mobile App)		X						
File Viewer (Mother Mobile App)		X						
Upload File (Mother Mobile App)		X						
Delete File Confirmation (Mother Mobile App)		X						
Expert Search and Filters (Mother Mobile App)		X						

Emergency Contacts (Mother Mobile App)		X						
Edit Emergency Contact (Mother Mobile App)		X						
Enable Fall Detection Confirmation (Mother Mobile App)		X						
Disable Fall Detection Confirmation (Mother Mobile App)		X						
Community Search (Family Member Mobile App)			X					
Topic Detail (Family Member Mobile App)			X					
Emergency Alert Detail (Family Member Mobile App)			X					
Verification Renewal (Expert App)				X				
Renewal Status (Expert App)				X				
Expert Suspension Status (Expert App)				X				
Verification Renewal (Expert Web Portal)				X				
Expert Suspension Status (Expert Web Portal)				X				
Expert Verification Renewal Queue (Admin Web Portal)							X	
Expert Verification Renewal Detail (Admin Web Portal)							X	

Suspend Expert Confirmation (Admin Web Portal)							X	
Security Events (Admin Web Portal)							X	
Security Event Detail (Admin Web Portal)							X	
Security Incident Investigation (Admin Web Portal)							X	
Security Incident Resolution (Admin Web Portal)							X	
Shared File Viewer (Expert App)				X				
Search Community Questions (Expert App)				X				
Search Community Topics (Expert App)				X				
Topic Detail (Expert App)				X				
Shared File Viewer (Expert Web Portal)				X				
Security Incident List (Admin Web Portal)							X	
Pre-exercise Safety Check (Mother Mobile App)		X						
Exercise Session (Mother Mobile App)		X						
Enable Posture Camera Confirmation (Mother Mobile App)		X						

Exercise Session Result (Mother Mobile App)		X						
Exercise History (Mother Mobile App)		X						
Maternal Health Metric Detail (Mother Mobile App)		X						
Postpartum Log List (Mother Mobile App)		X						
Postpartum Log Detail (Mother Mobile App)		X						
Switch Active Baby Selector (Mother Mobile App)		X						
Baby Daily Log Detail (Mother Mobile App)		X						
Development Milestone Detail (Mother Mobile App)		X						
My Answer Detail (Mother Mobile App)		X						
Reschedule Consultation (Mother Mobile App)		X						
Cancel Consultation Confirmation (Mother Mobile App)		X						
Consultation Summary (Mother Mobile App)		X						
Reminder Detail (Mother Mobile App)		X						
Care Group Members (Mother Mobile App)		X						
Pending Invitations		X						

(Mother Mobile App)								
Family Task Detail (Mother Mobile App)		X						
Verified Content Search (Mother Mobile App)		X						
Verified Content Detail (Mother Mobile App)		X						
Vaccination Detail (Mother Mobile App)		X						
Add Vaccination Record (Mother Mobile App)		X						
Growth Measurement History (Mother Mobile App)		X						
Growth Measurement Detail (Mother Mobile App)		X						
Reject Invitation Confirmation (Family Member Mobile App)			X					
Leave Care Group Confirmation (Family Member Mobile App)			X					
My Answer Detail (Family Member Mobile App)			X					
Verified Content Search (Family Member Mobile App)			X					
Verified Content Detail (Family Member Mobile App)			X					
My Expert Answer Detail (Expert App)				X				
Reschedule Consultation (Expert App)				X				

Mark No-show Confirmation (Expert App)				X				
Complete Consultation Confirmation (Expert App)				X				
My Expert Answer Detail (Expert Web Portal)				X				
Reschedule Consultation (Expert Web Portal)				X				
Mark No-show Confirmation (Expert Web Portal)				X				
Pregnancy Exercise List (Content Admin Web Portal)						X	X	
Pregnancy Exercise Detail (Content Admin Web Portal)						X	X	
Create Pregnancy Exercise (Content Admin Web Portal)						X	X	
Edit Pregnancy Exercise (Content Admin Web Portal)						X	X	
Exercise Preview (Content Admin Web Portal)						X	X	
Content Category List (Content Admin Web Portal)						X	X	
Unpublish Content Confirmation (Content Admin Web Portal)						X	X	
Posture Analysis Configuration List (Admin Web Portal)							X	

Posture Analysis Configuration Detail (Admin Web Portal)							X	
Edit Posture Analysis Configuration (Admin Web Portal)							X	
Consultation No-show Review (Admin Web Portal)							X	
Resolve Consultation Dispute (Admin Web Portal)							X	
Approve Refund Confirmation (Admin Web Portal)							X	
Reject Refund Confirmation (Admin Web Portal)							X	
Expert Profile (Mother Mobile App)		X						
Expert Consultation Pricing (Mother Mobile App)		X						
Booking Review (Mother Mobile App)		X						
Payment (Mother Mobile App)		X						
Payment Result (Mother Mobile App)		X						
Consultation Pricing (Expert App)				X				
Set Consultation Price (Expert App)				X				
Update Consultation Price (Expert App)				X				
Deactivate Consultation Price				X				

Confirmation (Expert App)								
Consultation Pricing (Expert Web Portal)				X				
Set Consultation Price (Expert Web Portal)				X				
Update Consultation Price (Expert Web Portal)				X				
Deactivate Consultation Price Confirmation (Expert Web Portal)				X				
Consultation Price Bands (Admin Web Portal)							X	
Configure Consultation Price Band (Admin Web Portal)							X	
Deactivate Price Band Confirmation (Admin Web Portal)							X	
Notification Center (Shared Web Portals)				X	X	X	X	X
Privacy Settings (Shared Web Portals)				X	X	X	X	X
Delete Account Confirmation (Shared Web Portals)				X	X	X		X
Deactivate Account Confirmation (Shared Web Portals)				X	X	X	X	X

Table 262: Screen Authorization

4.1.4 Non-Screen Functions

STT	Feature	System Function	Description
-----	---------	-----------------	-------------

1	MF-01 Account, Trust & Access Control	OTP Issuance and Verification Service	Generates expiring OTP challenges, enforces resend and attempt limits, verifies the submitted code and records security-relevant outcomes.
2	MF-01 Account, Trust & Access Control	Authentication Token and Session Service	Issues short-lived access tokens and device-specific refresh sessions, rotates tokens, detects revoked or expired sessions and supports sign-out from selected devices.
3	MF-01 Account, Trust & Access Control	RBAC Authorization Service	Evaluates role, permission, ownership, account state and platform context before protected data or operations are allowed.
4	MF-01 Consent & Privacy	Consent Scope and Expiry Service	Creates, validates, edits, expires and revokes permissions by recipient, purpose, data scope and time limit; access is denied when any condition is no longer valid.
5	MF-01 Consent & Privacy	Privacy Preference Enforcement Service	Applies profile visibility, online-status, location and data-use settings consistently to search, community, expert discovery and notification processing.
6	MF-01 Account Lifecycle	Account Deactivation Service	Temporarily disables sign-in and nonessential notifications, revokes active sessions and preserves data for possible reactivation under policy.
7	MF-01 Account Lifecycle	Account Deletion Workflow Job	Creates a deletion request, checks waiting period and unresolved obligations, anonymizes or deletes eligible data and preserves only legally or operationally required evidence.
8	MF-01 Expert Governance	Expert Credential Expiry Monitor	Checks credential expiry dates, notifies experts before expiry and changes verification state when evidence is no longer valid.
9	MF-01 Expert Governance	Expert Verification Renewal Service	Validates renewal submissions, links new evidence to prior verification, routes review to Admin and applies approved expiry dates or required-supplement status.
10	MF-01 Expert Governance	Expert Suspension Enforcement Service	Disables expert badge, public discovery, new consultation acceptance and restricted expert actions while preserving audit, appeal and settlement access as configured.
11	MF-02 Mother Care Journey	Journey Stage Calculation Service	Calculates pre-pregnancy, pregnancy week or postpartum stage from user-provided dates and refreshes stage-dependent dashboard data.
12	MF-03 Baby Care Journey & Growth	Baby Log Aggregation Service	Aggregates feeding, sleep, diaper, symptom and care logs into 24-hour and seven-day summaries without producing a medical diagnosis.
13	MF-03 / MF-13 Growth & Vaccination	Growth and Vaccination Schedule Service	Calculates age-based growth points and reference vaccination dates, marks sources and creates due or overdue states.
14	MF-04 Community Q&A	Community Safety Screening Service	Screens posts and answers for spam, abuse, prohibited advertising and high-risk medical wording, then labels or routes suspicious content for review.
15	MF-04 Community Q&A	Community Post Lifecycle Service	Enforces ownership and moderation locks for create, edit, soft-delete and restore states while retaining evidence required for reports or investigations.

16	MF-04 Community Q&A	Community Engagement Service	Creates and removes bookmarks, answer likes and topic follows idempotently, updates counters and excludes hidden or deleted content.
17	MF-04 / MF-17 Moderation	Moderation Queue Builder	Prioritizes cases from user reports, automated safety flags, repeated violations and escalated health-risk content.
18	MF-04 / MF-17 Moderation	Moderation Action Propagation Service	Applies approved hide, lock, warn or suspend actions to related content and accounts, sends notices and records the result.
19	MF-05 Verified Expert Network	Expert Availability and Presence Service	Maintains online, available and consultation-channel status and automatically expires stale presence.
20	MF-05 Verified Expert Network	Expert Search and Filter Service	Searches verified experts and applies specialty, channel, availability, fee, rating, online-state and consented-distance filters.
21	MF-05 / MF-19 Expert Discovery	Expert Matching and Nearby Search Service	Returns only verified, eligible and opted-in nearby experts using privacy-safe location precision and current availability.
22	MF-06 Free Expert Contribution	Contribution Point Calculation Job	Calculates expert contribution points and badges from eligible answers, likes and moderation outcomes while excluding removed activity.
23	MF-07 Direct Consultation	Consultation Booking and Slot Reservation Service	Validates availability and the selected active pricing package, reserves a slot, creates an immutable booking price snapshot, prevents duplicate booking and releases unpaid or expired reservations.
24	MF-07 Direct Consultation	Realtime Communication Session Service	Creates protected chat, voice or video rooms and short-lived participant tokens for a valid booking or accepted support request.
25	MF-07 Direct Consultation	Consultation Status and Timeout Job	Updates scheduled, accepted, active, completed, cancelled and no-show states according to events and configured time windows.
26	MF-07 Direct Consultation	Payment Callback Processing Service	Validates signed and idempotent payment callbacks and updates the transaction against the booking price snapshot; it never recalculates the payable amount from the expert's current price.
27	MF-07 Direct Consultation	Payment Reconciliation and Commission Job	Reconciles eligible completed consultations and stores gross amount, effective commission rate, platform commission, payment-gateway fee, refunded amount, expert net amount and settlement status.
28	MF-07 / MF-17 Disputes	Refund and Dispute Workflow Service	Collects consultation and technical evidence, routes cases to Admin and applies the approved refund to the booking/payment breakdown without changing the expert's public pricing or historical price versions.
29	MF-08 AI Nurse Assistant & Risk Triage	Structured Symptom Intake Service	Transforms user answers and selected context into a structured intake record for controlled safety processing.
30	MF-08 AI Nurse Assistant & Risk Triage	Risk Triage and Red-Flag Rule Engine	Applies controlled green, yellow and red rules; red flags prioritize emergency or in-person care and never diagnose or prescribe.

31	MF-08 / MF-12 AI and Verified Content	Verified RAG Answer Service	Retrieves approved content, FAQ and checklist sources and returns bounded answers with citations and safe fallback.
32	MF-09 Personal Health Records	Health Summary Generation Service	Builds a time-bounded summary from records explicitly selected by the owner and applies source labels and permission scope.
33	MF-09 File Management	Secure File Upload Processing Service	Validates ownership, file type and size, scans content, stores metadata and protects storage references.
34	MF-09 File Management	Protected File Access Service	Authorizes preview or download by ownership, consent scope and expiry and issues a short-lived protected file URL.
35	MF-09 File Management	File Deletion and Retention Service	Soft-deletes eligible owner files, blocks deletion when retention or record links apply and schedules physical purge when allowed.
36	MF-10 Reminders, Tasks & Care Plan	Reminder Scheduling and Dispatch Job	Finds due appointment, medication, vaccination, checklist and family-task reminders and creates reminder notifications.
37	MF-10 Reminders, Tasks & Care Plan	Reminder Recurrence and Snooze Service	Calculates the next occurrence, snooze time, completion and missed state while preventing duplicate delivery.
38	MF-11 Family Sync & Cooperative Care	Care Group Invitation Service	Creates expiring invitations, validates acceptance, prevents duplicate membership and initializes granted permissions.
39	MF-11 Family Sync & Cooperative Care	Family Alert Distribution Service	Sends consent-limited alerts and task updates only to active members whose scope includes the event type.
40	MF-12 Verified Content & Checklist Hub	Content Versioning and Publication Service	Stores immutable versions, manages draft and approval state and exposes only the latest approved version.
41	MF-12 Verified Content & Checklist Hub	Scheduled Content Publication Job	Publishes or unpublishes approved content at configured times and invalidates affected search indexes.
42	MF-16 Operation, Impact & Partner Dashboard	Operational Metrics Aggregation Job	Aggregates non-identifying usage, moderation, consultation and social-impact indicators for authorized dashboards.
43	MF-17 Safety, Compliance & Abuse Prevention	Audit Logging Service	Records authentication, sensitive access, permission changes, file activity, expert views, moderation, safety and financial actions in append-only logs.
44	MF-17 Audit & Security	Security Event Detection Service	Detects anomalous login, permission, file-access and sensitive-record patterns and creates reviewable security events with supporting evidence.
45	MF-17 Audit & Security	Security Incident Investigation Service	Links security events, audit records and evidence into an investigation case, tracks containment and resolution and records the final decision.
46	MF-17 Safety, Compliance & Abuse Prevention	Account and Content Risk Scoring Service	Combines report frequency, prohibited patterns and enforcement history to prioritize human review without making medical decisions.
47	MF-18 Partner Management	Partner Verification and Listing Publication Service	Controls partner review states and publishes only approved profiles, services and sponsored content with required labels and validity periods.

48	MF-18 Partner Management	Referral and Campaign Metric Aggregation Job	Aggregates listing views, referral events and campaign results without exposing personal health records.
49	MF-19 Map & Nearby Care Support	Map Geocoding, Nearby Search and ETA Service	Geocodes approved locations, searches nearby care facilities or opted-in experts and calculates distance, route and ETA.
50	MF-19 Map & Nearby Care Support	Location Privacy and Expiry Service	Limits location precision and visibility by role and purpose and automatically ends sharing when its TTL expires.
51	MF-19 Map & Nearby Care Support	Nearby Support Request Routing Service	Matches consented support requests with eligible experts and shares only minimum necessary information until expiry or cancellation.
52	MF-20 Connected Device Integration	Health Platform Synchronization Service	Imports supported observations, normalizes units, attaches source and timestamp labels and prevents duplicate records.
53	MF-20 Connected Device Integration	Device Synchronization Retry Job	Retries recoverable failures with backoff and marks stale or permanently failed imports without fabricating values.
54	MF-21 Smart Activity Monitoring	Fall Detection Configuration Service	Stores enable or disable state, sensor permission, monitoring window and device capability while keeping detection off until explicit consent.
55	MF-21 Smart Activity Monitoring	On-Device IMU Event Detection Service	Processes accelerometer and gyroscope signals on-device to identify a suspected fall or strong impact and starts a safety check.
56	MF-21 Smart Activity Monitoring	Safety Confirmation Countdown Service	Tracks the response window, cancels escalation when safety is confirmed and continues when help is requested or no response occurs.
57	MF-21 Smart Activity Monitoring	Emergency Contact Verification Service	Adds, verifies, prioritizes and removes emergency contacts and checks that at least one active recipient is available before alerting.
58	MF-21 Smart Activity Monitoring	Safety Alert Escalation Service	Sends the minimum necessary emergency alert and consented location, records delivery and offers MF-19 actions while prioritizing emergency services for red flags.
59	Cross-cutting Notifications	Notification Event Classification Service	Classifies events into reminder, community reply, consultation or emergency categories and applies category-specific urgency, deep link and audience rules.
60	Cross-cutting Notifications	Notification Queue and Delivery Service	Queues in-app, push and email messages, applies preferences and mandatory-safety rules, retries failures and stores delivery status.
61	Cross-cutting Notifications	Emergency Notification Delivery Service	Bypasses noncritical muting where policy permits, sends to verified emergency contacts, tracks acknowledgements and escalates failed delivery safely.
62	Cross-cutting Search	Community Question Search Index Job	Indexes eligible community questions and answers after create, edit, moderation or deletion and supports keyword and structured filters.
63	Cross-cutting Search	Community Topic Search and Follow Service	Indexes active topics, returns topic discovery results and maintains follow subscriptions used for feed personalization and notifications.

64	Cross-cutting Search	Expert Search Index Synchronization Job	Updates searchable expert fields, badge, availability and suspension state so ineligible experts are removed promptly.
65	Cross-cutting Data Protection	Data Retention and Soft-Deletion Job	Applies retention, anonymization and purge rules to deleted or deactivated accounts and expired data while preserving required evidence.
66	MF-02 Pregnancy Exercise & Posture Support	Exercise Catalog Eligibility Service	Returns only active exercises matching the selected trimester, difficulty and supported duration; it does not automatically prescribe an exercise as medical treatment.
67	MF-02 Pregnancy Exercise & Posture Support	Pre-exercise Safety Rule Service	Evaluates required safety answers against configured blocking and caution rules before an exercise session may start.
68	MF-02 Pregnancy Exercise & Posture Support	Exercise Session State Service	Creates and transitions exercise sessions through started, paused, resumed, completed and abandoned states and stores duration and completion metadata.
69	MF-02 Pregnancy Exercise & Posture Support	Posture Analysis Runtime Service	Processes live camera landmarks using the configured rule-based or ML-based mode and returns bounded posture feedback without saving raw camera video.
70	MF-02 Pregnancy Exercise & Posture Support	Exercise Result Aggregation Service	Aggregates session duration, completion level, posture score, recurring issues and safety warnings into the stored session result.
71	MF-02 Pregnancy Exercise & Posture Support	Exercise Content Publication Service	Validates exercise metadata, safety notes and status before an exercise is activated for the mobile catalog.
72	MF-02 Pregnancy Exercise & Posture Support	Posture Configuration Version Service	Versions analysis mode, rule/model identifier, confidence thresholds and feedback severity and applies only published configurations.
73	MF-02 Mother Care Journey	Maternal Record Soft-Deletion Service	Soft-deletes owner-entered maternal metrics or postpartum logs while retaining audit evidence and excluding deleted records from trends.
74	MF-03 Baby Care Journey & Growth	Active Baby Context Service	Stores and validates the currently selected baby profile so logs, reminders and dashboards use the intended child.
75	MF-03 Baby Care Journey & Growth	Baby Record Lifecycle Service	Applies ownership, edit and soft-delete rules to baby daily logs and development milestones.
76	MF-04 Community Q&A	Community Answer Lifecycle Service	Enforces author ownership and moderation locks for editing and soft-deleting the user's own answers.
77	MF-07 Direct Consultation	Consultation Rescheduling Service	Validates the policy deadline, slot availability and both-party confirmation before replacing a consultation schedule.
78	MF-07 Direct Consultation	Consultation Cancellation Service	Calculates allowed cancellation transitions, reason, payment consequence and notifications for both participants.
79	MF-07 Direct Consultation	Consultation No-show Evaluation Service	Checks waiting time and technical session evidence before allowing an expert or administrator to record a no-show.
80	MF-07 Direct Consultation	Consultation Completion Service	Transitions an active session to completed, unlocks summary and review steps and queues settlement processing.

81	MF-09 Personal Health Records	Health Record Detail Authorization Service	Returns metadata and protected attachments only when ownership or a valid consent scope permits access.
82	MF-10 Reminders, Tasks & Care Plan	Reminder Lifecycle Service	Applies complete, skip and delete actions, updates recurrence and prevents a skipped occurrence from deleting the recurring configuration.
83	MF-11 Family Sync & Cooperative Care	Care Group Membership Lifecycle Service	Revokes pending invitations, removes members, handles voluntary leave and immediately terminates associated data permissions.
84	MF-11 Family Sync & Cooperative Care	Family Task Lifecycle Service	Validates task ownership and status transitions for detail view, update, cancellation and assignee notification.
85	MF-12 Verified Content & Checklist Hub	Verified Content Search Service	Indexes and searches only approved, currently published articles, FAQs and checklists by keyword, stage and category.
86	MF-12 Verified Content & Checklist Hub	Content Category Management Service	Creates, reorders, updates and hides content categories while preserving references from historical content versions.
87	MF-12 Verified Content & Checklist Hub	Content Unpublication Service	Removes an approved version from user visibility without deleting its version history or audit trail.
88	MF-13 Vaccination & Growth Tracking	Vaccination Record Lifecycle Service	Creates, updates, soft-deletes, completes or postpones user-entered vaccination records and recalculates schedule status.
89	MF-13 Vaccination & Growth Tracking	Growth Measurement Lifecycle Service	Creates, updates and soft-deletes growth measurements and rebuilds chart data from the remaining valid history.
90	MF-07 Paid Direct Consultation & Commission	Consultation Price Band Versioning Service	Creates draft or effective price-band versions by consultation channel and duration, validates minimum/maximum and commission-rate rules, rejects overlapping effective periods and preserves historical versions.
91	MF-07 Paid Direct Consultation & Commission	Expert Consultation Price Management Service	Creates, updates and deactivates expert package prices only within the currently effective band, records effective time and price history and rejects out-of-band prices automatically without manual Admin approval.
92	MF-07 Paid Direct Consultation & Commission	Expert Pricing Read Service	Returns only active expert packages with channel, duration, effective price, estimated total payable amount and cancellation/refund information, and revalidates the selected version before booking.
93	MF-07 Paid Direct Consultation & Commission	Booking Price Snapshot Service	Copies the selected package, price-band version, gross amount, commission rate and applicable fee policy into an immutable booking snapshot used by payment, refund and settlement processing.

Table 263: Non-Screen Functions

4.1.5 Entity Relationship Diagram

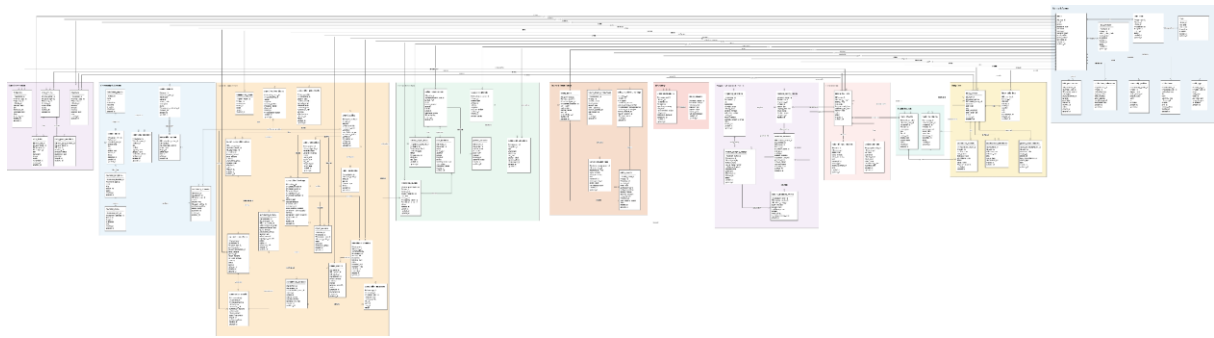


Figure 41: Entity Relationship Diagram

Entity Relationship Diagram

Entities Description

#	Entity	Description
1	roles	Catalog of system roles used for RBAC, including USER, MOTHER, FAMILY_MEMBER, EXPERT, MODERATOR, CONTENT_ADMIN, ADMIN, and PARTNER.
2	users	Root account entity for all CareBridge users; stores identity and account-status data without mixing health records or community-profile data.
3	user_roles	Junction entity that resolves the many-to-many relationship between users and roles.
4	user_sessions	User login sessions and refresh tokens used to support viewing and revoking sessions on other devices.
5	community_profiles	Public community display profile separated from identity records and health profiles.
6	notification_preferences	Per-user settings for notification channels and notification categories.
7	notifications	Personal notifications generated by the system from reminders, community activity, consultations, safety events, or operations.
8	data_permissions	Data-sharing permissions granted by users to family members or experts by scope, purpose, and expiry time.
9	audit_logs	Immutable log of important actions such as login, profile viewing, permission changes, moderation, consultation, payment, and safety events.
10	mother_journeys	Mother care journeys by pre-pregnancy, pregnancy, or postpartum state; one user may have multiple journeys over time.
11	maternal_health_metrics	Maternal health metrics entered manually or synchronized from another source, with source labeling.

#	Entity	Description
12	postpartum_logs	Postpartum recovery logs covering discharge, pain, breast milk, sleep, mood, and symptoms.
13	baby_profiles	Baby profiles created and managed by the mother or authorized caregiver.
14	baby_daily_logs	Daily baby logs for feeding, sleep, diaper, temperature, reflux/vomiting, prescribed medicine, and notes.
15	development_milestones	Baby development milestones such as rolling, crawling, walking, speaking, teething, or weaning.
16	growth_measurements	Weight, height, and head-circumference measurements used to build growth charts.
17	vaccination_records	Baby vaccination schedule and status data; used for personal tracking and not as a replacement for the official vaccination book.
18	health_records	Personal health-record repository for mothers and babies, including ultrasound images, lab results, prescriptions, vaccination slips, and clinical results, with source labels.
19	health_summaries	Time-bounded health-data summaries generated by users for viewing or temporary sharing.
20	reminders	Reminders for appointments, medicine or vitamins, vaccination, checklists, and care tasks.
21	care_groups	Care groups created by the mother or manager for a journey or baby.
22	care_group_members	Junction entity between care_groups and users; stores invitation status and member permissions.
23	care_tasks	Care tasks assigned to group members, with due dates and completion status.
24	expenses	Expenses related to pregnancy, postpartum care, or baby care.
25	community_topics	Community topic catalog used to classify questions and content.
26	community_questions	Questions posted by users; may appear anonymous publicly while remaining internally traceable.
27	community_answers	Answers from users or experts; expert answers carry badge and professional-scope information.
28	content_reports	Reports submitted by users for violating content or accounts.
29	moderation_actions	Moderation actions for reports or content, such as approve, hide, lock, warn, or suspend.
30	content_items	Articles, FAQs, or educational content created and moderated by content administrators.
31	checklist_templates	Moderated checklist templates for preparation by care stage.
32	checklist_items	Child items belonging to a checklist template.
33	expert_profiles	Professional profiles for registered doctors or experts; verification badge is granted only after admin approval.
34	expert_credentials	Uploaded expert documents and certificates used for verification.
35	expert_availability	Expert availability slots for consultations, including support mode and slot status; prices are versioned separately in expert_consultation_prices.
36	expert_location_shares	Expert-shared location or service area, controlled by consent and expiry, to support nearby search.

#	Entity	Description
37	consultation_bookings	Consultation booking request between a user and expert; locks the price snapshot, commission rate, and cancellation policy at confirmation or payment time.
38	consultation_sessions	Chat, voice-call, or video-call sessions created from valid bookings.
39	consultation_messages	Messages exchanged during consultation sessions, protected so only authorized parties can access them.
40	payment_transactions	Payment transactions for bookings based on the price snapshot, including gross price, gateway fee, refund amount, net amount, and transaction status.
41	commission_records	Commission calculation result after an eligible completed session, including original price, commission rate or amount, gateway fee, refund, expert net amount, and reconciliation status.
42	expert_reviews	Expert reviews created only after a valid completed booking.
43	contribution_points	Contribution point and badge history for users or experts based on community activity.
44	consultation_price_bands	Versioned CareBridge price bands by consultation mode, duration, and professional scope, including minimum price, maximum price, commission rate, and effective period.
45	expert_consultation_prices	Versioned expert-specific consultation prices within allowed bands; new prices apply only to future bookings and never retroactively change locked bookings.
46	consultation_disputes	Consultation dispute records with reason, evidence, handling status, and operations decision.
47	refund_records	Refund records linked to transactions and disputes, storing refund amount, processing status, and payment-gateway refund code.
48	settlement_records	Expert settlement records by period or commission item, including gross value, commission, gateway fee, refund, net payout, and payment status.
49	triage_assessments	Structured symptom-intake sessions with green, yellow, or red risk classification; not a medical diagnosis.
50	triage_answers	Step-by-step answers submitted during a symptom-intake session.
51	partner_organizations	Profiles for registered or approved clinics, healthcare facilities, social organizations, or sponsors.
52	partner_expert_links	Junction entity between partner_organizations and expert_profiles; links are effective only after approval.
53	partner_services	Reference clinic or service listings submitted by partners and approved by admins for public display.
54	sponsored_campaigns	Partner-submitted sponsored campaigns or content that require clear labeling and moderation.
55	care_facilities	Healthcare-facility catalog used for emergency maps and nearby search; records may come from partners or map sources.
56	emergency_events	Emergency or support-search flow opened by the user or system from triage, emergency button, or IMU event.
57	location_snapshots	Minimal time-limited location snapshots used for emergency, family alerts, or nearby search.

#	Entity	Description
58	health_device_connections	Optional smartwatch or health-platform connections storing secure token references and consent status.
59	device_measurements	Device or health-platform measurements with source and confidence labels; used for trend viewing and not as a substitute for clinical measurement.
60	safety_monitoring_settings	Configuration for IMU monitoring, countdown, alert recipients, and location consent.
61	safety_events	Suspected fall or strong-impact events detected from IMU data; stores minimal metadata and does not assert injury.
62	safety_alerts	Minimal alerts sent to family members after the user requests help or does not respond.
63	pregnancy_exercises	Moderated pregnancy-exercise catalog with trimester scope, difficulty, duration, instructions, safety warnings, and version status.
64	exercise_safety_checks	Mother's pre-exercise safety-check results for a specific exercise and journey, including structured answers, red flags, and allow-or-block decision.
65	exercise_sessions	Mother's pregnancy-exercise sessions, storing time, pause/completion status, completion level, aggregate posture score, and alerts.
66	posture_analysis_configs	Rule-based or ML-based posture-analysis configuration for each exercise, including rule or model version, confidence threshold, and effective period.
67	posture_feedback_events	Posture feedback events generated during an exercise session; stores only necessary result/keypoint summaries and no raw camera video.

Table 264: Entities Description

Data Dictionary

1. roles

No	Field	Name	Type	Constraint	Description
1	role_id	Role ID	UUID	PK, NOT NULL	Primary key that uniquely identifies one role record.
2	role_code	Role Code	VARCHAR(50)	UNIQUE	Display or reference value for role code.
3	role_name	Role Name	VARCHAR(100)	NULL allowed / business validation	Display or reference value for role name.
4	description	Description	VARCHAR(500)	NULL allowed / business validation	Descriptive text for the role record.

5	is_active	Is Active	BOOLEAN	NULL allowed / business validation	Boolean flag indicating is active.
6	created_at	Created At	TIMESTAMP	NOT NULL	Timestamp when the record was created.
7	updated_at	Updated At	TIMESTAMP	NULL allowed / business validation	Timestamp when the record was last updated.

Table 265: roles Data Dictionary

2. users

No	Field	Name	Type	Constraint	Description
1	user_id	User ID	UUID	PK, NOT NULL	Primary key that uniquely identifies one user record.
2	email	Email	VARCHAR(255)	UNIQUE	Email address associated with this record.
3	phone	Phone	VARCHAR(30)	NULL allowed / business validation	Phone number associated with this record.
4	password_hash	Password Hash	VARCHAR(255)	NULL allowed / business validation	Secure password-related value stored for authentication.
5	full_name	Full Name	VARCHAR(150)	NULL allowed / business validation	Display or reference value for full name.
6	avatar_url	Avatar URL	VARCHAR(500)	NULL allowed / business validation	URL or URI value for avatar URL.
7	account_status	Account Status	VARCHAR(30)	NULL allowed / business validation	Current status value for account status.
8	email_verified	Email Verified	BOOLEAN	NULL allowed / business validation	Boolean flag indicating email verified.
9	phone_verified	Phone Verified	BOOLEAN	NULL allowed / business validation	Boolean flag indicating phone verified.
10	last_login_at	Last Login At	TIMESTAMP	NULL allowed / business validation	Business timestamp or

					date value for last login at.
11	created_at	Created At	TIMESTAMP	NOT NULL	Timestamp when the record was created.
12	updated_at	Updated At	TIMESTAMP	NULL allowed / business validation	Timestamp when the record was last updated.

Table 266: users Data Dictionary

3. user_roles

No	Field	Name	Type	Constraint	Description
1	user_role_id	User Role ID	UUID	PK, NOT NULL	Primary key that uniquely identifies one user_role record.
2	user_id	User ID	UUID	FK	Foreign key that links this record to the related entity through user_id.
3	role_id	Role ID	UUID	FK	Foreign key that links this record to the related entity through role_id.
4	assigned_by	Assigned By	UUID	FK	Foreign key that links this record to the related entity through assigned_by.
5	assigned_at	Assigned At	TIMESTAMP	NULL allowed / business validation	Business timestamp or date value for assigned at.
6	expires_at	Expires At	TIMESTAMP	NULL allowed / business validation	Business timestamp or date value for expires at.
7	status	Status	VARCHAR(20)	NOT NULL	Current status value for status.

8	created_at	Created At	TIMESTAMP	NOT NULL	Timestamp when the record was created.
9	updated_at	Updated At	TIMESTAMP	NULL allowed / business validation	Timestamp when the record was last updated.

Table 267: user_roles Data Dictionary

4. user_sessions

No	Field	Name	Type	Constraint	Description
1	session_id	Session ID	UUID	PK, NOT NULL	Primary key that uniquely identifies one session record.
2	user_id	User ID	UUID	FK	Foreign key that links this record to the related entity through user_id.
3	refresh_token_hash	Refresh Token Hash	VARCHAR(255)	NULL allowed / business validation	Secure token-related value used by the workflow.
4	device_name	Device Name	VARCHAR(150)	NULL allowed / business validation	Display or reference value for device name.
5	ip_address	IP Address	VARCHAR(64)	NULL allowed / business validation	Stores IP address for this session record.
6	expires_at	Expires At	TIMESTAMP	NULL allowed / business validation	Business timestamp or date value for expires at.
7	revoked_at	Revoked At	TIMESTAMP	NULL allowed / business validation	Business timestamp or date value for revoked at.
8	status	Status	VARCHAR(20)	NOT NULL	Current status value for status.
9	created_at	Created At	TIMESTAMP	NOT NULL	Timestamp when the record was created.

10	updated_at	Updated At	TIMESTAMP	NULL allowed / business validation	Timestamp when the record was last updated.
----	------------	------------	-----------	------------------------------------	---

Table 268: user_sessions Data Dictionary

5. community_profiles

No	Field	Name	Type	Constraint	Description
1	community_profile_id	Community Profile ID	UUID	PK, NOT NULL	Primary key that uniquely identifies one community_profile record.
2	user_id	User ID	UUID	FK	Foreign key that links this record to the related entity through user_id.
3	display_name	Display Name	VARCHAR(100)	NULL allowed / business validation	Display or reference value for display name.
4	public_avatar_url	Public Avatar URL	VARCHAR(500)	NULL allowed / business validation	URL or URI value for public avatar URL.
5	interest_stage	Interest Stage	VARCHAR(30)	NULL allowed / business validation	Stores interest stage for this community_profile record.
6	region	Region	VARCHAR(120)	NULL allowed / business validation	Stores region for this community_profile record.
7	bio	Bio	VARCHAR(500)	NULL allowed / business validation	Stores bio for this community_profile record.
8	is_visible	Is Visible	BOOLEAN	NULL allowed / business validation	Boolean flag indicating is visible.
9	created_at	Created At	TIMESTAMP	NOT NULL	Timestamp when the record was created.
10	updated_at	Updated At	TIMESTAMP	NULL allowed / business validation	Timestamp when the record

					was last updated.
--	--	--	--	--	-------------------

Table 269: community_profiles Data Dictionary

6. notification_preferences

No	Field	Name	Type	Constraint	Description
1	preference_id	Preference ID	UUID	PK, NOT NULL	Primary key that uniquely identifies one preference record.
2	user_id	User ID	UUID	FK	Foreign key that links this record to the related entity through user_id.
3	notification_type	Notification Type	VARCHAR(50)	NULL allowed / business validation	Classification value for notification type.
4	in_app_enabled	In App Enabled	BOOLEAN	NULL allowed / business validation	Boolean flag indicating in app enabled.
5	push_enabled	Push Enabled	BOOLEAN	NULL allowed / business validation	Boolean flag indicating push enabled.
6	email_enabled	Email Enabled	BOOLEAN	NULL allowed / business validation	Boolean flag indicating email enabled.
7	quiet_hours_start	Quiet Hours Start	TIME	NULL allowed / business validation	Stores quiet hours start for this preference record.
8	quiet_hours_end	Quiet Hours End	TIME	NULL allowed / business validation	Stores quiet hours end for this preference record.
9	created_at	Created At	TIMESTAMP	NOT NULL	Timestamp when the record was created.
10	updated_at	Updated At	TIMESTAMP	NULL allowed / business validation	Timestamp when the record was last updated.

Table 270: notification_preferences Data Dictionary

7. notifications

No	Field	Name	Type	Constraint	Description
1	notification_id	Notification ID	UUID	PK, NOT NULL	Primary key that uniquely identifies one notification record.
2	recipient_user_id	Recipient User ID	UUID	FK	Foreign key that links this record to the related entity through recipient_user_id.
3	notification_type	Notification Type	VARCHAR(50)	NULL allowed / business validation	Classification value for notification type.
4	title	Title	VARCHAR(200)	NULL allowed / business validation	Display or reference value for title.
5	body	Body	TEXT	NULL allowed / business validation	Stores body for this notification record.
6	reference_type	Reference Type	VARCHAR(50)	NULL allowed / business validation	Classification value for reference type.
7	reference_id	Reference ID	UUID	NULL allowed / business validation	Stores reference ID for this notification record.
8	is_read	Is Read	BOOLEAN	NULL allowed / business validation	Boolean flag indicating is read.
9	sent_at	Sent At	TIMESTAMP	NULL allowed / business validation	Business timestamp or date value for sent at.
10	delivery_status	Delivery Status	VARCHAR(20)	NULL allowed / business validation	Current status value for delivery status.
11	created_at	Created At	TIMESTAMP	NOT NULL	Timestamp when the record was created.
12	updated_at	Updated At	TIMESTAMP	NULL allowed / business validation	Timestamp when the record

					was last updated.
--	--	--	--	--	-------------------

Table 271: notifications Data Dictionary

8. data_permissions

No	Field	Name	Type	Constraint	Description
1	permission_id	Permission ID	UUID	PK, NOT NULL	Primary key that uniquely identifies one permission record.
2	owner_user_id	Owner User ID	UUID	FK	Foreign key that links this record to the related entity through owner_user_id.
3	grantee_user_id	Grantee User ID	UUID	FK	Foreign key that links this record to the related entity through grantee_user_id.
4	scope_type	Scope Type	VARCHAR(50)	NULL allowed / business validation	Classification value for scope type.
5	scope_reference_id	Scope Reference ID	UUID	NULL allowed / business validation	Stores scope reference ID for this permission record.
6	purpose	Purpose	VARCHAR(255)	NULL allowed / business validation	Stores purpose for this permission record.
7	granted_at	Granted At	TIMESTAMP	NULL allowed / business validation	Business timestamp or date value for granted at.
8	expires_at	Expires At	TIMESTAMP	NULL allowed / business validation	Business timestamp or date value for expires at.
9	revoked_at	Revoked At	TIMESTAMP	NULL allowed / business validation	Business timestamp or date value for revoked at.

10	status	Status	VARCHAR(20)	NOT NULL	Current status value for status.
11	created_at	Created At	TIMESTAMP	NOT NULL	Timestamp when the record was created.
12	updated_at	Updated At	TIMESTAMP	NULL allowed / business validation	Timestamp when the record was last updated.

Table 272: data_permissions Data Dictionary

9. audit_logs

No	Field	Name	Type	Constraint	Description
1	audit_log_id	Audit Log ID	UUID	PK, NOT NULL	Primary key that uniquely identifies one audit_log record.
2	actor_user_id	Actor User ID	UUID	FK	Foreign key that links this record to the related entity through actor_user_id.
3	action	Action	VARCHAR(100)	NULL allowed / business validation	Stores action for this audit_log record.
4	entity_type	Entity Type	VARCHAR(80)	NULL allowed / business validation	Classification value for entity type.
5	entity_id	Entity ID	UUID	NULL allowed / business validation	Stores entity ID for this audit_log record.
6	old_value_json	Old Value Json	JSON	NULL allowed / business validation	Stores old value json for this audit_log record.
7	new_value_json	New Value Json	JSON	NULL allowed / business validation	Stores new value json for this audit_log record.
8	ip_address	IP Address	VARCHAR(64)	NULL allowed / business validation	Stores IP address for this audit_log record.

9	created_at	Created At	TIMESTAMP	NOT NULL	Timestamp when the record was created.
---	------------	------------	-----------	----------	--

Table 273: audit_logs Data Dictionary

10. mother_journeys

No	Field	Name	Type	Constraint	Description
1	journey_id	Journey ID	UUID	PK, NOT NULL	Primary key that uniquely identifies one journey record.
2	owner_user_id	Owner User ID	UUID	FK	Foreign key that links this record to the related entity through owner_user_id.
3	journey_type	Journey Type	VARCHAR(30)	NULL allowed / business validation	Classification value for journey type.
4	start_date	Start Date	DATE	NULL allowed / business validation	Business timestamp or date value for start date.
5	last_menstrual_date	Last Menstrual Date	DATE	NULL allowed / business validation	Business timestamp or date value for last menstrual date.
6	estimated_due_date	Estimated Due Date	DATE	NULL allowed / business validation	Business timestamp or date value for estimated due date.
7	delivery_date	Delivery Date	DATE	NULL allowed / business validation	Business timestamp or date value for delivery date.
8	status	Status	VARCHAR(20)	NOT NULL	Current status value for status.
9	notes	Notes	TEXT	NULL allowed / business validation	Text content recorded for notes.
10	created_at	Created At	TIMESTAMP	NOT NULL	Timestamp when the record was created.

11	updated_at	Updated At	TIMESTAMP	NULL allowed / business validation	Timestamp when the record was last updated.
----	------------	------------	-----------	------------------------------------	---

Table 274: mother_journeys Data Dictionary

11. maternal_health_metrics

No	Field	Name	Type	Constraint	Description
1	metric_id	Metric ID	UUID	PK, NOT NULL	Primary key that uniquely identifies one metric record.
2	journey_id	Journey ID	UUID	FK	Foreign key that links this record to the related entity through journey_id.
3	metric_type	Metric Type	VARCHAR(40)	NULL allowed / business validation	Classification value for metric type.
4	value_numeric	Value Numeric	DECIMAL(12, 3)	NULL allowed / business validation	Stores value numeric for this metric record.
5	value_secondary	Value Secondary	DECIMAL(12, 3)	NULL allowed / business validation	Stores value secondary for this metric record.
6	unit	Unit	VARCHAR(30)	NULL allowed / business validation	Stores unit for this metric record.
7	measured_at	Measured At	TIMESTAMP	NULL allowed / business validation	Business timestamp or date value for measured at.
8	source_type	Source Type	VARCHAR(30)	NULL allowed / business validation	Classification value for source type.
9	source_reference_id	Source Reference ID	UUID	NULL allowed / business validation	Stores source reference ID for this metric record.
10	note	Note	VARCHAR(500)	NULL allowed / business validation	Text content recorded for note.

11	created_at	Created At	TIMESTAMP	NOT NULL	Timestamp when the record was created.
12	updated_at	Updated At	TIMESTAMP	NULL allowed / business validation	Timestamp when the record was last updated.

Table 275: maternal_health_metrics Data Dictionary

12. postpartum_logs

No	Field	Name	Type	Constraint	Description
1	postpartum_log_id	Postpartum Log ID	UUID	PK, NOT NULL	Primary key that uniquely identifies one postpartum_log record.
2	journey_id	Journey ID	UUID	FK	Foreign key that links this record to the related entity through journey_id.
3	log_date	Log Date	DATE	NULL allowed / business validation	Business timestamp or date value for log date.
4	pain_level	Pain Level	SMALLINT	CHECK	Classification value for pain level.
5	bleeding_level	Bleeding Level	VARCHAR(20)	NULL allowed / business validation	Classification value for bleeding level.
6	mood_level	Mood Level	SMALLINT	CHECK	Classification value for mood level.
7	sleep_hours	Sleep Hours	DECIMAL(4,1)	NULL allowed / business validation	Stores sleep hours for this postpartum_log record.
8	breastfeeding_note	Breastfeeding Note	TEXT	NULL allowed / business validation	Monetary value used for breastfeeding note.
9	symptom_note	Symptom Note	TEXT	NULL allowed / business validation	Text content recorded for symptom note.

10	created_at	Created At	TIMESTAMP	NOT NULL	Timestamp when the record was created.
11	updated_at	Updated At	TIMESTAMP	NULL allowed / business validation	Timestamp when the record was last updated.

Table 276: postpartum_logs Data Dictionary

13. baby_profiles

No	Field	Name	Type	Constraint	Description
1	baby_id	Baby ID	UUID	PK, NOT NULL	Primary key that uniquely identifies one baby record.
2	owner_user_id	Owner User ID	UUID	FK	Foreign key that links this record to the related entity through owner_user_id.
3	related_journey_id	Related Journey ID	UUID	FK	Foreign key that links this record to the related entity through related_journey_id.
4	nickname	Nickname	VARCHAR(100)	NULL allowed / business validation	Display or reference value for nickname.
5	birth_date	Birth Date	DATE	NULL allowed / business validation	Business timestamp or date value for birth date.
6	sex	Sex	VARCHAR(20)	NULL allowed / business validation	Stores sex for this baby record.
7	birth_weight_kg	Birth Weight Kg	DECIMAL(5,2)	NULL allowed / business validation	Stores birth weight kg for this baby record.
8	birth_length_cm	Birth Length Cm	DECIMAL(5,2)	NULL allowed / business validation	Stores birth length cm for this baby record.

9	status	Status	VARCHAR(20)	NOT NULL	Current status value for status.
10	created_at	Created At	TIMESTAMP	NOT NULL	Timestamp when the record was created.
11	updated_at	Updated At	TIMESTAMP	NULL allowed / business validation	Timestamp when the record was last updated.

Table 277: baby_profiles Data Dictionary

14. baby_daily_logs

No	Field	Name	Type	Constraint	Description
1	baby_log_id	Baby Log ID	UUID	PK, NOT NULL	Primary key that uniquely identifies one baby_log record.
2	baby_id	Baby ID	UUID	FK	Foreign key that links this record to the related entity through baby_id.
3	log_type	Log Type	VARCHAR(30)	NULL allowed / business validation	Classification value for log type.
4	started_at	Started At	TIMESTAMP	NULL allowed / business validation	Business timestamp or date value for started at.
5	ended_at	Ended At	TIMESTAMP	NULL allowed / business validation	Business timestamp or date value for ended at.
6	quantity	Quantity	DECIMAL(10, 2)	NULL allowed / business validation	Stores quantity for this baby_log record.
7	unit	Unit	VARCHAR(20)	NULL allowed / business validation	Stores unit for this baby_log record.
8	note	Note	TEXT	NULL allowed / business validation	Text content recorded for note.
9	recorded_by	Recorded By	UUID	FK	Foreign key that links this

					record to the related entity through recorded_by.
10	created_at	Created At	TIMESTAMP	NOT NULL	Timestamp when the record was created.
11	updated_at	Updated At	TIMESTAMP	NULL allowed / business validation	Timestamp when the record was last updated.

Table 278: baby_daily_logs Data Dictionary

15. development_milestones

No	Field	Name	Type	Constraint	Description
1	milestone_id	Milestone ID	UUID	PK, NOT NULL	Primary key that uniquely identifies one milestone record.
2	baby_id	Baby ID	UUID	FK	Foreign key that links this record to the related entity through baby_id.
3	milestone_type	Milestone Type	VARCHAR(60)	NULL allowed / business validation	Classification value for milestone type.
4	achieved_date	Achieved Date	DATE	NULL allowed / business validation	Business timestamp or date value for achieved date.
5	note	Note	TEXT	NULL allowed / business validation	Text content recorded for note.
6	source_type	Source Type	VARCHAR(30)	NULL allowed / business validation	Classification value for source type.
7	recorded_by	Recorded By	UUID	FK	Foreign key that links this record to the related entity through recorded_by.

8	created_at	Created At	TIMESTAMP	NOT NULL	Timestamp when the record was created.
9	updated_at	Updated At	TIMESTAMP	NULL allowed / business validation	Timestamp when the record was last updated.

Table 279: development_milestones Data Dictionary

16. growth_measurements

No	Field	Name	Type	Constraint	Description
1	growth_measurement_id	Growth Measurement ID	UUID	PK, NOT NULL	Primary key that uniquely identifies one growth_measurement record.
2	baby_id	Baby ID	UUID	FK	Foreign key that links this record to the related entity through baby_id.
3	measured_date	Measured Date	DATE	NULL allowed / business validation	Business timestamp or date value for measured date.
4	weight_kg	Weight Kg	DECIMAL(5,2)	NULL allowed / business validation	Stores weight kg for this growth_measurement record.
5	height_cm	Height Cm	DECIMAL(5,2)	NULL allowed / business validation	Stores height cm for this growth_measurement record.
6	head_circumference_cm	Head Circumference Cm	DECIMAL(5,2)	NULL allowed / business validation	Stores head circumference cm for this growth_measurement record.
7	source_type	Source Type	VARCHAR(30)	NULL allowed / business validation	Classification value for source type.
8	note	Note	VARCHAR(500)	NULL allowed / business validation	Text content recorded for note.

9	created_at	Created At	TIMESTAMP	NOT NULL	Timestamp when the record was created.
10	updated_at	Updated At	TIMESTAMP	NULL allowed / business validation	Timestamp when the record was last updated.

Table 280: growth_measurements Data Dictionary

17. vaccination_records

No	Field	Name	Type	Constraint	Description
1	vaccination_record_id	Vaccination Record ID	UUID	PK, NOT NULL	Primary key that uniquely identifies one vaccination_record record.
2	baby_id	Baby ID	UUID	FK	Foreign key that links this record to the related entity through baby_id.
3	vaccine_name	Vaccine Name	VARCHAR(200)	NULL allowed / business validation	Display or reference value for vaccine name.
4	dose_number	Dose Number	VARCHAR(30)	NULL allowed / business validation	Count or numeric reference used for dose number.
5	scheduled_date	Scheduled Date	DATE	NULL allowed / business validation	Business timestamp or date value for scheduled date.
6	administered_date	Administered Date	DATE	NULL allowed / business validation	Business timestamp or date value for administered date.
7	status	Status	VARCHAR(20)	NOT NULL	Current status value for status.
8	facility_name	Facility Name	VARCHAR(255)	NULL allowed / business validation	Display or reference value for facility name.

9	proof_record_id	Proof Record ID	UUID	FK	Foreign key that links this record to the related entity through proof_record_id.
10	created_at	Created At	TIMESTAMP	NOT NULL	Timestamp when the record was created.
11	updated_at	Updated At	TIMESTAMP	NULL allowed / business validation	Timestamp when the record was last updated.

Table 281: vaccination_records Data Dictionary

18. health_records

No	Field	Name	Type	Constraint	Description
1	health_record_id	Health Record ID	UUID	PK, NOT NULL	Primary key that uniquely identifies one health_record record.
2	owner_user_id	Owner User ID	UUID	FK	Foreign key that links this record to the related entity through owner_user_id.
3	journey_id	Journey ID	UUID	FK	Foreign key that links this record to the related entity through journey_id.
4	baby_id	Baby ID	UUID	FK	Foreign key that links this record to the related entity through baby_id.
5	record_type	Record Type	VARCHAR(50)	NULL allowed / business validation	Classification value for record type.
6	title	Title	VARCHAR(200)	NULL allowed / business validation	Display or reference value for title.

7	file_url	File URL	VARCHAR(500)	NULL allowed / business validation	URL or URI value for file URL.
8	record_date	Record Date	DATE	NULL allowed / business validation	Business timestamp or date value for record date.
9	source_type	Source Type	VARCHAR(30)	NULL allowed / business validation	Classification value for source type.
10	source_name	Source Name	VARCHAR(255)	NULL allowed / business validation	Display or reference value for source name.
11	status	Status	VARCHAR(20)	NOT NULL	Current status value for status.
12	created_at	Created At	TIMESTAMP	NOT NULL	Timestamp when the record was created.
13	updated_at	Updated At	TIMESTAMP	NULL allowed / business validation	Timestamp when the record was last updated.

Table 282: health_records Data Dictionary

19. health_summaries

No	Field	Name	Type	Constraint	Description
1	summary_id	Summary ID	UUID	PK, NOT NULL	Primary key that uniquely identifies one summary record.
2	owner_user_id	Owner User ID	UUID	FK	Foreign key that links this record to the related entity through owner_user_id.
3	journey_id	Journey ID	UUID	FK	Foreign key that links this record to the related entity through journey_id.
4	baby_id	Baby ID	UUID	FK	Foreign key that links this record to the

					related entity through baby_id.
5	summary_period	Summary Period	VARCHAR(30)	NULL allowed / business validation	Stores summary period for this summary record.
6	period_start	Period Start	TIMESTAMP	NULL allowed / business validation	Business timestamp or date value for period start.
7	period_end	Period End	TIMESTAMP	NULL allowed / business validation	Business timestamp or date value for period end.
8	summary_json	Summary Json	JSON	NULL allowed / business validation	Stores summary json for this summary record.
9	generated_by	Generated By	VARCHAR(20)	NULL allowed / business validation	Numeric value used for generated by.
10	status	Status	VARCHAR(20)	NOT NULL	Current status value for status.
11	created_at	Created At	TIMESTAMP	NOT NULL	Timestamp when the record was created.
12	updated_at	Updated At	TIMESTAMP	NULL allowed / business validation	Timestamp when the record was last updated.

Table 283: health_summaries Data Dictionary

20. reminders

No	Field	Name	Type	Constraint	Description
1	reminder_id	Reminder ID	UUID	PK, NOT NULL	Primary key that uniquely identifies one reminder record.
2	owner_user_id	Owner User ID	UUID	FK	Foreign key that links this record to the related entity through owner_user_id.

3	journey_id	Journey ID	UUID	FK	Foreign key that links this record to the related entity through journey_id.
4	baby_id	Baby ID	UUID	FK	Foreign key that links this record to the related entity through baby_id.
5	reminder_type	Reminder Type	VARCHAR(40)	NULL allowed / business validation	Classification value for reminder type.
6	title	Title	VARCHAR(200)	NULL allowed / business validation	Display or reference value for title.
7	scheduled_at	Scheduled At	TIMESTAMP	NULL allowed / business validation	Business timestamp or date value for scheduled at.
8	recurrence_rule	Recurrence Rule	VARCHAR(255)	NULL allowed / business validation	Stores recurrence rule for this reminder record.
9	status	Status	VARCHAR(20)	NOT NULL	Current status value for status.
10	snoozed_until	Snoozed Until	TIMESTAMP	NULL allowed / business validation	Business timestamp or date value for snoozed until.
11	created_at	Created At	TIMESTAMP	NOT NULL	Timestamp when the record was created.
12	updated_at	Updated At	TIMESTAMP	NULL allowed / business validation	Timestamp when the record was last updated.

Table 284: reminders Data Dictionary

21. care_groups

No	Field	Name	Type	Constraint	Description
1	care_group_id	Care Group ID	UUID	PK, NOT NULL	Primary key that uniquely identifies one

					care_group record.
2	owner_user_id	Owner User ID	UUID	FK	Foreign key that links this record to the related entity through owner_user_id.
3	journey_id	Journey ID	UUID	FK	Foreign key that links this record to the related entity through journey_id.
4	baby_id	Baby ID	UUID	FK	Foreign key that links this record to the related entity through baby_id.
5	group_name	Group Name	VARCHAR(150)	NULL allowed / business validation	Display or reference value for group name.
6	status	Status	VARCHAR(20)	NOT NULL	Current status value for status.
7	created_at	Created At	TIMESTAMP	NOT NULL	Timestamp when the record was created.
8	updated_at	Updated At	TIMESTAMP	NULL allowed / business validation	Timestamp when the record was last updated.

Table 285: care_groups Data Dictionary

22. care_group_members

No	Field	Name	Type	Constraint	Description
1	care_group_member_id	Care Group Member ID	UUID	PK, NOT NULL	Primary key that uniquely identifies one care_group_member record.
2	care_group_id	Care Group ID	UUID	FK	Foreign key that links this record to the related entity through care_group_id.

3	user_id	User ID	UUID	FK	Foreign key that links this record to the related entity through user_id.
4	member_role	Member Role	VARCHAR(30)	NULL allowed / business validation	Stores member role for this care_group_member record.
5	invitation_status	Invitation Status	VARCHAR(20)	NULL allowed / business validation	Current status value for invitation status.
6	permission_json	Permission Json	JSON	NULL allowed / business validation	Stores permission json for this care_group_member record.
7	joined_at	Joined At	TIMESTAMP	NULL allowed / business validation	Business timestamp or date value for joined at.
8	created_at	Created At	TIMESTAMP	NOT NULL	Timestamp when the record was created.
9	updated_at	Updated At	TIMESTAMP	NULL allowed / business validation	Timestamp when the record was last updated.

Table 286: care_group_members Data Dictionary

23. care_tasks

No	Field	Name	Type	Constraint	Description
1	care_task_id	Care Task ID	UUID	PK, NOT NULL	Primary key that uniquely identifies one care_task record.
2	care_group_id	Care Group ID	UUID	FK	Foreign key that links this record to the related entity through care_group_id.
3	assigned_by	Assigned By	UUID	FK	Foreign key that links this

					record to the related entity through assigned_by.
4	assigned_to	Assigned To	UUID	FK	Foreign key that links this record to the related entity through assigned_to.
5	title	Title	VARCHAR(200)	NULL allowed / business validation	Display or reference value for title.
6	description	Description	TEXT	NULL allowed / business validation	Descriptive text for the care_task record.
7	due_at	Due At	TIMESTAMP	NULL allowed / business validation	Business timestamp or date value for due at.
8	status	Status	VARCHAR(20)	NOT NULL	Current status value for status.
9	completed_at	Completed At	TIMESTAMP	NULL allowed / business validation	Business timestamp or date value for completed at.
10	created_at	Created At	TIMESTAMP	NOT NULL	Timestamp when the record was created.
11	updated_at	Updated At	TIMESTAMP	NULL allowed / business validation	Timestamp when the record was last updated.

Table 287: care_tasks Data Dictionary

24. expenses

No	Field	Name	Type	Constraint	Description
1	expense_id	Expense ID	UUID	PK, NOT NULL	Primary key that uniquely identifies one expense record.
2	owner_user_id	Owner User ID	UUID	FK	Foreign key that links this record to the related entity

					through owner_user_id.
3	journey_id	Journey ID	UUID	FK	Foreign key that links this record to the related entity through journey_id.
4	baby_id	Baby ID	UUID	FK	Foreign key that links this record to the related entity through baby_id.
5	category	Category	VARCHAR(40)	NULL allowed / business validation	Classification value for category.
6	amount	Amount	DECIMAL(15, 2)	NULL allowed / business validation	Monetary value used for amount.
7	currency	Currency	CHAR(3)	NULL allowed / business validation	Stores currency for this expense record.
8	expense_date	Expense Date	DATE	NULL allowed / business validation	Business timestamp or date value for expense date.
9	note	Note	VARCHAR(500)	NULL allowed / business validation	Text content recorded for note.
10	created_at	Created At	TIMESTAMP	NOT NULL	Timestamp when the record was created.
11	updated_at	Updated At	TIMESTAMP	NULL allowed / business validation	Timestamp when the record was last updated.

Table 288: expenses Data Dictionary

25. community_topics

No	Field	Name	Type	Constraint	Description
1	topic_id	Topic ID	UUID	PK, NOT NULL	Primary key that uniquely identifies one topic record.

2	name	Name	VARCHAR(150)	NULL allowed / business validation	Display or reference value for name.
3	slug	Slug	VARCHAR(160)	NULL allowed / business validation	Stores slug for this topic record.
4	description	Description	VARCHAR(500)	NULL allowed / business validation	Descriptive text for the topic record.
5	risk_level	Risk Level	VARCHAR(20)	NULL allowed / business validation	Classification value for risk level.
6	is_active	Is Active	BOOLEAN	NULL allowed / business validation	Boolean flag indicating is active.
7	created_at	Created At	TIMESTAMP	NOT NULL	Timestamp when the record was created.
8	updated_at	Updated At	TIMESTAMP	NULL allowed / business validation	Timestamp when the record was last updated.

Table 289: community_topics Data Dictionary

26. community_questions

No	Field	Name	Type	Constraint	Description
1	question_id	Question ID	UUID	PK, NOT NULL	Primary key that uniquely identifies one question record.
2	author_user_id	Author User ID	UUID	FK	Foreign key that links this record to the related entity through author_user_id.
3	topic_id	Topic ID	UUID	FK	Foreign key that links this record to the related entity through topic_id.
4	title	Title	VARCHAR(250)	NULL allowed / business validation	Display or reference value for title.

5	content	Content	TEXT	NULL allowed / business validation	Text content recorded for content.
6	is_anonymous	Is Anonymous	BOOLEAN	NULL allowed / business validation	Boolean flag indicating is anonymous.
7	urgency_level	Urgency Level	VARCHAR(20)	NULL allowed / business validation	Classification value for urgency level.
8	moderation_status	Moderation Status	VARCHAR(20)	NULL allowed / business validation	Current status value for moderation status.
9	published_at	Published At	TIMESTAMP	NULL allowed / business validation	Business timestamp or date value for published at.
10	created_at	Created At	TIMESTAMP	NOT NULL	Timestamp when the record was created.
11	updated_at	Updated At	TIMESTAMP	NULL allowed / business validation	Timestamp when the record was last updated.

Table 290: community_questions Data Dictionary

27. community_answers

No	Field	Name	Type	Constraint	Description
1	answer_id	Answer ID	UUID	PK, NOT NULL	Primary key that uniquely identifies one answer record.
2	question_id	Question ID	UUID	FK	Foreign key that links this record to the related entity through question_id.
3	author_user_id	Author User ID	UUID	FK	Foreign key that links this record to the related entity through author_user_id.
4	content	Content	TEXT	NULL allowed / business validation	Text content recorded for content.

5	answer_type	Answer Type	VARCHAR(20)	NULL allowed / business validation	Classification value for answer type.
6	expert_profile_id	Expert Profile ID	UUID	FK	Foreign key that links this record to the related entity through expert_profile_id.
7	moderation_status	Moderation Status	VARCHAR(20)	NULL allowed / business validation	Current status value for moderation status.
8	helpful_count	Helpful Count	INT	NULL allowed / business validation	Count or numeric reference used for helpful count.
9	published_at	Published At	TIMESTAMP	NULL allowed / business validation	Business timestamp or date value for published at.
10	created_at	Created At	TIMESTAMP	NOT NULL	Timestamp when the record was created.
11	updated_at	Updated At	TIMESTAMP	NULL allowed / business validation	Timestamp when the record was last updated.

Table 291: community_answers Data Dictionary

28. content_reports

No	Field	Name	Type	Constraint	Description
1	report_id	Report ID	UUID	PK, NOT NULL	Primary key that uniquely identifies one report record.
2	reporter_user_id	Reporter User ID	UUID	FK	Foreign key that links this record to the related entity through reporter_user_id.

3	target_type	Target Type	VARCHAR(30)	NULL allowed / business validation	Classification value for target type.
4	target_id	Target ID	UUID	NULL allowed / business validation	Stores target ID for this report record.
5	reason_code	Reason Code	VARCHAR(40)	NULL allowed / business validation	Text content recorded for reason code.
6	description	Description	TEXT	NULL allowed / business validation	Descriptive text for the report record.
7	status	Status	VARCHAR(20)	NOT NULL	Current status value for status.
8	assigned_moderator_id	Assigned Moderator ID	UUID	FK	Foreign key that links this record to the related entity through assigned_moderator_id.
9	resolved_at	Resolved At	TIMESTAMP	NULL allowed / business validation	Business timestamp or date value for resolved at.
10	created_at	Created At	TIMESTAMP	NOT NULL	Timestamp when the record was created.
11	updated_at	Updated At	TIMESTAMP	NULL allowed / business validation	Timestamp when the record was last updated.

Table 292: content_reports Data Dictionary

29. moderation_actions

No	Field	Name	Type	Constraint	Description
1	moderation_action_id	Moderation Action ID	UUID	PK, NOT NULL	Primary key that uniquely identifies one moderation_action record.
2	report_id	Report ID	UUID	FK	Foreign key that links this record to the related entity through report_id.

3	moderator_user_id	Moderator User ID	UUID	FK	Foreign key that links this record to the related entity through moderator_user_id.
4	target_type	Target Type	VARCHAR(30)	NULL allowed / business validation	Classification value for target type.
5	target_id	Target ID	UUID	NULL allowed / business validation	Stores target ID for this moderation_action record.
6	action_type	Action Type	VARCHAR(30)	NULL allowed / business validation	Classification value for action type.
7	reason	Reason	TEXT	NULL allowed / business validation	Text content recorded for reason.
8	action_at	Action At	TIMESTAMP	NULL allowed / business validation	Business timestamp or date value for action at.
9	expires_at	Expires At	TIMESTAMP	NULL allowed / business validation	Business timestamp or date value for expires at.

Table 293: moderation_actions Data Dictionary

30. content_items

No	Field	Name	Type	Constraint	Description
1	content_item_id	Content Item ID	UUID	PK, NOT NULL	Primary key that uniquely identifies one content_item record.
2	topic_id	Topic ID	UUID	FK	Foreign key that links this record to the related entity through topic_id.
3	author_user_id	Author User ID	UUID	FK	Foreign key that links this record to the related entity

					through author_user_id.
4	content_type	Content Type	VARCHAR(30)	NULL allowed / business validation	Text content recorded for content type.
5	title	Title	VARCHAR(250)	NULL allowed / business validation	Display or reference value for title.
6	body	Body	TEXT	NULL allowed / business validation	Stores body for this content_item record.
7	source_label	Source Label	VARCHAR(255)	NULL allowed / business validation	Stores source label for this content_item record.
8	version_no	Version No	INT	NULL allowed / business validation	Stores version no for this content_item record.
9	status	Status	VARCHAR(20)	NOT NULL	Current status value for status.
10	published_at	Published At	TIMESTAMP	NULL allowed / business validation	Business timestamp or date value for published at.
11	created_at	Created At	TIMESTAMP	NOT NULL	Timestamp when the record was created.
12	updated_at	Updated At	TIMESTAMP	NULL allowed / business validation	Timestamp when the record was last updated.

Table 294: content_items Data Dictionary

31. checklist_templates

No	Field	Name	Type	Constraint	Description
1	checklist_template_id	Checklist Template ID	UUID	PK, NOT NULL	Primary key that uniquely identifies one checklist_template record.
2	content_item_id	Content Item ID	UUID	FK	Foreign key that links this record to the related entity through

					content_item_id.
3	name	Name	VARCHAR(200)	NULL allowed / business validation	Display or reference value for name.
4	stage	Stage	VARCHAR(30)	NULL allowed / business validation	Stores stage for this checklist_template record.
5	version_no	Version No	INT	NULL allowed / business validation	Stores version no for this checklist_template record.
6	status	Status	VARCHAR(20)	NOT NULL	Current status value for status.
7	created_at	Created At	TIMESTAMP	NOT NULL	Timestamp when the record was created.
8	updated_at	Updated At	TIMESTAMP	NULL allowed / business validation	Timestamp when the record was last updated.

Table 295: checklist_templates Data Dictionary

32. checklist_items

No	Field	Name	Type	Constraint	Description
1	checklist_item_id	Checklist Item ID	UUID	PK, NOT NULL	Primary key that uniquely identifies one checklist_item record.
2	checklist_template_id	Checklist Template ID	UUID	FK	Foreign key that links this record to the related entity through checklist_template_id.
3	item_order	Item Order	INT	NULL allowed / business validation	Stores item order for this checklist_item record.
4	item_text	Item Text	VARCHAR(500)	NULL allowed / business validation	Stores item text for this checklist_item record.

5	is_required	Is Required	BOOLEAN	NULL allowed / business validation	Boolean flag indicating is required.
6	note	Note	VARCHAR(500)	NULL allowed / business validation	Text content recorded for note.
7	created_at	Created At	TIMESTAMP	NOT NULL	Timestamp when the record was created.
8	updated_at	Updated At	TIMESTAMP	NULL allowed / business validation	Timestamp when the record was last updated.

Table 296: checklist_items Data Dictionary

33. expert_profiles

No	Field	Name	Type	Constraint	Description
1	expert_profile_id	Expert Profile ID	UUID	PK, NOT NULL	Primary key that uniquely identifies one expert_profile record.
2	user_id	User ID	UUID	FK	Foreign key that links this record to the related entity through user_id.
3	specialty	Specialty	VARCHAR(150)	NULL allowed / business validation	Stores specialty for this expert_profile record.
4	professional_title	Professional Title	VARCHAR(150)	NULL allowed / business validation	Display or reference value for professional title.
5	experience_years	Experience Years	INT	NULL allowed / business validation	Stores experience years for this expert_profile record.
6	workplace	Workplace	VARCHAR(255)	NULL allowed / business validation	Stores workplace for this expert_profile record.

7	consultation_scope	Consultation Scope	TEXT	NULL allowed / business validation	Stores consultation scope for this expert_profile record.
8	verification_status	Verification Status	VARCHAR(20)	NULL allowed / business validation	Current status value for verification status.
9	verified_at	Verified At	TIMESTAMP	NULL allowed / business validation	Business timestamp or date value for verified at.
10	verified_by	Verified By	UUID	FK	Foreign key that links this record to the related entity through verified_by.
11	rating_avg	Rating Avg	DECIMAL(3,2)	NULL allowed / business validation	Stores rating avg for this expert_profile record.
12	created_at	Created At	TIMESTAMP	NOT NULL	Timestamp when the record was created.
13	updated_at	Updated At	TIMESTAMP	NULL allowed / business validation	Timestamp when the record was last updated.

Table 297: expert_profiles Data Dictionary

34. expert_credentials

No	Field	Name	Type	Constraint	Description
1	credential_id	Credential ID	UUID	PK, NOT NULL	Primary key that uniquely identifies one credential record.
2	expert_profile_id	Expert Profile ID	UUID	FK	Foreign key that links this record to the related entity through expert_profile_id.

3	credential_type	Credential Type	VARCHAR(50)	NULL allowed / business validation	Classification value for credential type.
4	credential_number	Credential Number	VARCHAR(120)	NULL allowed / business validation	Count or numeric reference used for credential number.
5	issuer	Issuer	VARCHAR(255)	NULL allowed / business validation	Stores issuer for this credential record.
6	issued_date	Issued Date	DATE	NULL allowed / business validation	Business timestamp or date value for issued date.
7	expiry_date	Expiry Date	DATE	NULL allowed / business validation	Business timestamp or date value for expiry date.
8	file_url	File URL	VARCHAR(500)	NULL allowed / business validation	URL or URI value for file URL.
9	review_status	Review Status	VARCHAR(20)	NULL allowed / business validation	Current status value for review status.
10	review_note	Review Note	TEXT	NULL allowed / business validation	Text content recorded for review note.
11	created_at	Created At	TIMESTAMP	NOT NULL	Timestamp when the record was created.
12	updated_at	Updated At	TIMESTAMP	NULL allowed / business validation	Timestamp when the record was last updated.

Table 298: expert_credentials Data Dictionary

35. expert_availability

No	Field	Name	Type	Constraint	Description
1	availability_id	Availability ID	UUID	PK, NOT NULL	Primary key that uniquely identifies one availability record.
2	expert_profile_id	Expert Profile ID	UUID	FK	Foreign key that links this

					record to the related entity through expert_profile_id.
3	start_at	Start At	TIMESTAMP	NULL allowed / business validation	Business timestamp or date value for start at.
4	end_at	End At	TIMESTAMP	NULL allowed / business validation	Business timestamp or date value for end at.
5	channel_type	Channel Type	VARCHAR(20)	NULL allowed / business validation	Classification value for channel type.
6	status	Status	VARCHAR(20)	NOT NULL	Current status value for status.
7	created_at	Created At	TIMESTAMP	NOT NULL	Timestamp when the record was created.
8	updated_at	Updated At	TIMESTAMP	NULL allowed / business validation	Timestamp when the record was last updated.

Table 299: expert_availability Data Dictionary

36. expert_location_shares

No	Field	Name	Type	Constraint	Description
1	location_share_id	Location Share ID	UUID	PK, NOT NULL	Primary key that uniquely identifies one location_share record.
2	expert_profile_id	Expert Profile ID	UUID	FK	Foreign key that links this record to the related entity through expert_profile_id.
3	latitude	Latitude	DECIMAL(10, 7)	NULL allowed / business validation	Location-related value used for latitude.

4	longitude	Longitude	DECIMAL(10, 7)	NULL allowed / business validation	Location-related value used for longitude.
5	accuracy_meters	Accuracy Meters	DECIMAL(10, 2)	NULL allowed / business validation	Stores accuracy meters for this location_share record.
6	availability_status	Availability Status	VARCHAR(20)	NULL allowed / business validation	Current status value for availability status.
7	shared_at	Shared At	TIMESTAMP	NULL allowed / business validation	Business timestamp or date value for shared at.
8	expires_at	Expires At	TIMESTAMP	NULL allowed / business validation	Business timestamp or date value for expires at.
9	consent_reference	Consent Reference	UUID	FK	Foreign key that links this record to the related entity through consent_reference.
10	created_at	Created At	TIMESTAMP	NOT NULL	Timestamp when the record was created.
11	updated_at	Updated At	TIMESTAMP	NULL allowed / business validation	Timestamp when the record was last updated.

Table 300: expert_location_shares Data Dictionary

37. consultation_bookings

No	Field	Name	Type	Constraint	Description
1	booking_id	Booking ID	UUID	PK, NOT NULL	Primary key that uniquely identifies one booking record.
2	requester_user_id	Requester User ID	UUID	FK	Foreign key that links this record to the related entity through

					requester_user_id.
3	expert_profile_id	Expert Profile ID	UUID	FK	Foreign key that links this record to the related entity through expert_profile_id.
4	availability_id	Availability ID	UUID	FK	Foreign key that links this record to the related entity through availability_id.
5	expert_price_id	Expert Price ID	UUID	FK	Foreign key that links this record to the related entity through expert_price_id.
6	shared_summary_id	Shared Summary ID	UUID	FK	Foreign key that links this record to the related entity through shared_summary_id.
7	topic	Topic	VARCHAR(200)	NULL allowed / business validation	Stores topic for this booking record.
8	channel_type	Channel Type	VARCHAR(20)	NULL allowed / business validation	Classification value for channel type.
9	duration_minutes	Duration Minutes	INT	NULL allowed / business validation	Stores duration minutes for this booking record.
10	scheduled_start	Scheduled Start	TIMESTAMP	NULL allowed / business validation	Business timestamp or date value for scheduled start.
11	scheduled_end	Scheduled End	TIMESTAMP	NULL allowed / business validation	Business timestamp or date value for scheduled end.
12	price_snapshot_amount	Price Snapshot Amount	DECIMAL(15,2)	NULL allowed / business validation	Monetary value used for price

					snapshot amount.
13	commission_rate_snapshot	Commission Rate Snapshot	DECIMAL(5,4)	CHECK 0..1	Monetary value used for commission rate snapshot.
14	currency	Currency	CHAR(3)	NULL allowed / business validation	Stores currency for this booking record.
15	cancellation_policy_snapshot	Cancellation Policy Snapshot	TEXT	NULL allowed / business validation	Stores cancellation policy snapshot for this booking record.
16	price_locked_at	Price Locked At	TIMESTAMP	NULL allowed / business validation	Business timestamp or date value for price locked at.
17	status	Status	VARCHAR(25)	NOT NULL	Current status value for status.
18	created_at	Created At	TIMESTAMP	NOT NULL	Timestamp when the record was created.
19	updated_at	Updated At	TIMESTAMP	NULL allowed / business validation	Timestamp when the record was last updated.

Table 301: consultation_bookings Data Dictionary

38. consultation_sessions

No	Field	Name	Type	Constraint	Description
1	session_id	Session ID	UUID	PK, NOT NULL	Primary key that uniquely identifies one session record.
2	booking_id	Booking ID	UUID	FK	Foreign key that links this record to the related entity through booking_id.
3	communication_room_id	Communication Room ID	VARCHAR(255)	NULL allowed / business validation	Stores communication room ID for this session record.

4	started_at	Started At	TIMESTAMP	NULL allowed / business validation	Business timestamp or date value for started at.
5	ended_at	Ended At	TIMESTAMP	NULL allowed / business validation	Business timestamp or date value for ended at.
6	session_status	Session Status	VARCHAR(20)	NULL allowed / business validation	Current status value for session status.
7	expert_summary	Expert Summary	TEXT	NULL allowed / business validation	Stores expert summary for this session record.
8	technical_log_json	Technical Log Json	JSON	NULL allowed / business validation	Stores technical log json for this session record.
9	created_at	Created At	TIMESTAMP	NOT NULL	Timestamp when the record was created.
10	updated_at	Updated At	TIMESTAMP	NULL allowed / business validation	Timestamp when the record was last updated.

Table 302: consultation_sessions Data Dictionary

39. consultation_messages

No	Field	Name	Type	Constraint	Description
1	message_id	Message ID	UUID	PK, NOT NULL	Primary key that uniquely identifies one message record.
2	session_id	Session ID	UUID	FK	Foreign key that links this record to the related entity through session_id.
3	sender_user_id	Sender User ID	UUID	FK	Foreign key that links this record to the related entity through sender_user_id.

4	message_type	Message Type	VARCHAR(20)	NULL allowed / business validation	Text content recorded for message type.
5	message_body	Message Body	TEXT	NULL allowed / business validation	Text content recorded for message body.
6	file_url	File URL	VARCHAR(500)	NULL allowed / business validation	URL or URI value for file URL.
7	sent_at	Sent At	TIMESTAMP	NULL allowed / business validation	Business timestamp or date value for sent at.
8	read_at	Read At	TIMESTAMP	NULL allowed / business validation	Business timestamp or date value for read at.
9	status	Status	VARCHAR(20)	NOT NULL	Current status value for status.

Table 303: consultation_messages Data Dictionary

40. payment_transactions

No	Field	Name	Type	Constraint	Description
1	payment_id	Payment ID	UUID	PK, NOT NULL	Primary key that uniquely identifies one payment record.
2	booking_id	Booking ID	UUID	FK	Foreign key that links this record to the related entity through booking_id.
3	payer_user_id	Payer User ID	UUID	FK	Foreign key that links this record to the related entity through payer_user_id.
4	gateway_name	Gateway Name	VARCHAR(50)	NULL allowed / business validation	Display or reference value for gateway name.
5	gateway_transaction_id	Gateway Transaction ID	VARCHAR(150)	NULL allowed / business validation	Stores gateway transaction ID for this

					payment record.
6	gross_amount	Gross Amount	DECIMAL(15, 2)	CHECK >= 0	Monetary value used for gross amount.
7	gateway_fee	Gateway Fee	DECIMAL(15, 2)	CHECK >= 0	Monetary value used for gateway fee.
8	refund_amount	Refund Amount	DECIMAL(15, 2)	CHECK >= 0	Monetary value used for refund amount.
9	net_paid_amount	Net Paid Amount	DECIMAL(15, 2)	NULL allowed / business validation	Monetary value used for net paid amount.
10	currency	Currency	CHAR(3)	NULL allowed / business validation	Stores currency for this payment record.
11	status	Status	VARCHAR(20)	NOT NULL	Current status value for status.
12	paid_at	Paid At	TIMESTAMP	NULL allowed / business validation	Business timestamp or date value for paid at.
13	refunded_at	Refunded At	TIMESTAMP	NULL allowed / business validation	Business timestamp or date value for refunded at.
14	created_at	Created At	TIMESTAMP	NOT NULL	Timestamp when the record was created.
15	updated_at	Updated At	TIMESTAMP	NULL allowed / business validation	Timestamp when the record was last updated.

Table 304: payment_transactions Data Dictionary

41. commission_records

No	Field	Name	Type	Constraint	Description
1	commission_id	Commission ID	UUID	PK, NOT NULL	Primary key that uniquely identifies one commission record.
2	payment_id	Payment ID	UUID	FK	Foreign key that links this record to the

					related entity through payment_id.
3	expert_profile_id	Expert Profile ID	UUID	FK	Foreign key that links this record to the related entity through expert_profile_id.
4	original_price	Original Price	DECIMAL(15, 2)	CHECK >= 0	Monetary value used for original price.
5	commission_rate	Commission Rate	DECIMAL(5,4)	CHECK 0..1	Monetary value used for commission rate.
6	commission_amount	Commission Amount	DECIMAL(15, 2)	CHECK >= 0	Monetary value used for commission amount.
7	gateway_fee	Gateway Fee	DECIMAL(15, 2)	CHECK >= 0	Monetary value used for gateway fee.
8	refund_amount	Refund Amount	DECIMAL(15, 2)	CHECK >= 0	Monetary value used for refund amount.
9	expert_net_amount	Expert Net Amount	DECIMAL(15, 2)	CHECK >= 0	Monetary value used for expert net amount.
10	eligible_at	Eligible At	TIMESTAMP	NULL allowed / business validation	Business timestamp or date value for eligible at.
11	settlement_status	Settlement Status	VARCHAR(20)	NULL allowed / business validation	Current status value for settlement status.
12	created_at	Created At	TIMESTAMP	NOT NULL	Timestamp when the record was created.
13	updated_at	Updated At	TIMESTAMP	NULL allowed / business validation	Timestamp when the record was last updated.

Table 305: commission_records Data Dictionary

42. expert_reviews

No	Field	Name	Type	Constraint	Description
1	review_id	Review ID	UUID	PK, NOT NULL	Primary key that uniquely identifies one review record.
2	booking_id	Booking ID	UUID	FK	Foreign key that links this record to the related entity through booking_id.
3	reviewer_user_id	Reviewer User ID	UUID	FK	Foreign key that links this record to the related entity through reviewer_user_id.
4	expert_profile_id	Expert Profile ID	UUID	FK	Foreign key that links this record to the related entity through expert_profile_id.
5	rating	Rating	SMALLINT	CHECK	Stores rating for this review record.
6	comment	Comment	TEXT	NULL allowed / business validation	Stores comment for this review record.
7	moderation_status	Moderation Status	VARCHAR(20)	NULL allowed / business validation	Current status value for moderation status.
8	created_at	Created At	TIMESTAMP	NOT NULL	Timestamp when the record was created.
9	updated_at	Updated At	TIMESTAMP	NULL allowed / business validation	Timestamp when the record was last updated.

Table 306: expert_reviews Data Dictionary

43. contribution_points

No	Field	Name	Type	Constraint	Description
----	-------	------	------	------------	-------------

1	point_record_id	Point Record ID	UUID	PK, NOT NULL	Primary key that uniquely identifies one point_record record.
2	user_id	User ID	UUID	FK	Foreign key that links this record to the related entity through user_id.
3	source_type	Source Type	VARCHAR(40)	NULL allowed / business validation	Classification value for source type.
4	source_id	Source ID	UUID	NULL allowed / business validation	Stores source ID for this point_record record.
5	points	Points	INT	NULL allowed / business validation	Stores points for this point_record record.
6	reason	Reason	VARCHAR(255)	NULL allowed / business validation	Text content recorded for reason.
7	recorded_at	Recorded At	TIMESTAMP	NULL allowed / business validation	Business timestamp or date value for recorded at.

Table 307: contribution_points Data Dictionary

44. consultation_price_bands

No	Field	Name	Type	Constraint	Description
1	price_band_id	Price Band ID	UUID	PK, NOT NULL	Primary key that uniquely identifies one price_band record.
2	configured_by	Configured By	UUID	FK	Foreign key that links this record to the related entity through configured_by.
3	channel_type	Channel Type	VARCHAR(20)	NULL allowed / business validation	Classification value for channel type.

4	duration_minutes	Duration Minutes	INT	NULL allowed / business validation	Stores duration minutes for this price_band record.
5	specialty_scope	Specialty Scope	VARCHAR(150)	NULL allowed / business validation	Stores specialty scope for this price_band record.
6	minimum_price	Minimum Price	DECIMAL(15,2)	CHECK >= 0	Monetary value used for minimum price.
7	maximum_price	Maximum Price	DECIMAL(15,2)	CHECK >= 0	Monetary value used for maximum price.
8	commission_rate	Commission Rate	DECIMAL(5,4)	CHECK 0..1	Monetary value used for commission rate.
9	currency	Currency	CHAR(3)	NULL allowed / business validation	Stores currency for this price_band record.
10	effective_from	Effective From	TIMESTAMP	NULL allowed / business validation	Business timestamp or date value for effective from.
11	effective_to	Effective To	TIMESTAMP	NULL allowed / business validation	Business timestamp or date value for effective to.
12	status	Status	VARCHAR(20)	NOT NULL	Current status value for status.
13	created_at	Created At	TIMESTAMP	NOT NULL	Timestamp when the record was created.
14	updated_at	Updated At	TIMESTAMP	NULL allowed / business validation	Timestamp when the record was last updated.

Table 308: consultation_price_bands Data Dictionary

45. expert_consultation_prices

No	Field	Name	Type	Constraint	Description
1	expert_price_id	Expert Price ID	UUID	PK, NOT NULL	Primary key that uniquely identifies one

					expert_price record.
2	expert_profile_id	Expert Profile ID	UUID	FK	Foreign key that links this record to the related entity through expert_profile_id.
3	price_band_id	Price Band ID	UUID	FK	Foreign key that links this record to the related entity through price_band_id.
4	channel_type	Channel Type	VARCHAR(20)	NULL allowed / business validation	Classification value for channel type.
5	duration_minutes	Duration Minutes	INT	NULL allowed / business validation	Stores duration minutes for this expert_price record.
6	price_amount	Price Amount	DECIMAL(15, 2)	CHECK >= 0	Monetary value used for price amount.
7	currency	Currency	CHAR(3)	NULL allowed / business validation	Stores currency for this expert_price record.
8	cancellation_policy	Cancellation Policy	TEXT	NULL allowed / business validation	Stores cancellation policy for this expert_price record.
9	effective_from	Effective From	TIMESTAMP	NULL allowed / business validation	Business timestamp or date value for effective from.
10	effective_to	Effective To	TIMESTAMP	NULL allowed / business validation	Business timestamp or date value for effective to.
11	status	Status	VARCHAR(20)	NOT NULL	Current status value for status.
12	version_no	Version No	INT	NULL allowed / business validation	Stores version no for this expert_price record.

13	created_at	Created At	TIMESTAMP	NOT NULL	Timestamp when the record was created.
14	updated_at	Updated At	TIMESTAMP	NULL allowed / business validation	Timestamp when the record was last updated.

Table 309: expert_consultation_prices Data Dictionary

46. consultation_disputes

No	Field	Name	Type	Constraint	Description
1	dispute_id	Dispute ID	UUID	PK, NOT NULL	Primary key that uniquely identifies one dispute record.
2	booking_id	Booking ID	UUID	FK	Foreign key that links this record to the related entity through booking_id.
3	submitted_by	Submitted By	UUID	FK	Foreign key that links this record to the related entity through submitted_by.
4	resolved_by	Resolved By	UUID	FK	Foreign key that links this record to the related entity through resolved_by.
5	reason_code	Reason Code	VARCHAR(50)	NULL allowed / business validation	Text content recorded for reason code.
6	description	Description	TEXT	NULL allowed / business validation	Descriptive text for the dispute record.
7	evidence_json	Evidence Json	JSON	NULL allowed / business validation	Stores evidence json for this dispute record.
8	status	Status	VARCHAR(20)	NOT NULL	Current status value for status.
9	resolution_type	Resolution Type	VARCHAR(30)	NULL allowed / business validation	Classification value for resolution type.

10	resolution_note	Resolution Note	TEXT	NULL allowed / business validation	Text content recorded for resolution note.
11	submitted_at	Submitted At	TIMESTAMP	NULL allowed / business validation	Business timestamp or date value for submitted at.
12	resolved_at	Resolved At	TIMESTAMP	NULL allowed / business validation	Business timestamp or date value for resolved at.
13	created_at	Created At	TIMESTAMP	NOT NULL	Timestamp when the record was created.
14	updated_at	Updated At	TIMESTAMP	NULL allowed / business validation	Timestamp when the record was last updated.

Table 310: consultation_disputes Data Dictionary

47. refund_records

No	Field	Name	Type	Constraint	Description
1	refund_id	Refund ID	UUID	PK, NOT NULL	Primary key that uniquely identifies one refund record.
2	payment_id	Payment ID	UUID	FK	Foreign key that links this record to the related entity through payment_id.
3	dispute_id	Dispute ID	UUID	FK	Foreign key that links this record to the related entity through dispute_id.
4	approved_by	Approved By	UUID	FK	Foreign key that links this record to the related entity through approved_by.
5	refund_amount	Refund Amount	DECIMAL(15, 2)	CHECK >= 0	Monetary value used for refund amount.

6	currency	Currency	CHAR(3)	NULL allowed / business validation	Stores currency for this refund record.
7	reason	Reason	VARCHAR(255)	NULL allowed / business validation	Text content recorded for reason.
8	gateway_refund_id	Gateway Refund ID	VARCHAR(150)	NULL allowed / business validation	Stores gateway refund ID for this refund record.
9	status	Status	VARCHAR(20)	NOT NULL	Current status value for status.
10	requested_at	Requested At	TIMESTAMP	NULL allowed / business validation	Business timestamp or date value for requested at.
11	processed_at	Processed At	TIMESTAMP	NULL allowed / business validation	Business timestamp or date value for processed at.
12	created_at	Created At	TIMESTAMP	NOT NULL	Timestamp when the record was created.
13	updated_at	Updated At	TIMESTAMP	NULL allowed / business validation	Timestamp when the record was last updated.

Table 311: refund_records Data Dictionary

48. settlement_records

No	Field	Name	Type	Constraint	Description
1	settlement_id	Settlement ID	UUID	PK, NOT NULL	Primary key that uniquely identifies one settlement record.
2	commission_id	Commission ID	UUID	FK	Foreign key that links this record to the related entity through commission_id.
3	expert_profile_id	Expert Profile ID	UUID	FK	Foreign key that links this record to the related entity through

					expert_profile_id.
4	settlement_period_start	Settlement Period Start	DATE	NULL allowed / business validation	Stores settlement period start for this settlement record.
5	settlement_period_end	Settlement Period End	DATE	NULL allowed / business validation	Stores settlement period end for this settlement record.
6	gross_amount	Gross Amount	DECIMAL(15, 2)	CHECK >= 0	Monetary value used for gross amount.
7	commission_amount	Commission Amount	DECIMAL(15, 2)	CHECK >= 0	Monetary value used for commission amount.
8	gateway_fee	Gateway Fee	DECIMAL(15, 2)	CHECK >= 0	Monetary value used for gateway fee.
9	refund_amount	Refund Amount	DECIMAL(15, 2)	CHECK >= 0	Monetary value used for refund amount.
10	expert_net_amount	Expert Net Amount	DECIMAL(15, 2)	CHECK >= 0	Monetary value used for expert net amount.
11	status	Status	VARCHAR(20)	NOT NULL	Current status value for status.
12	settled_at	Settled At	TIMESTAMP	NULL allowed / business validation	Business timestamp or date value for settled at.
13	reference_code	Reference Code	VARCHAR(100)	NULL allowed / business validation	Display or reference value for reference code.
14	created_at	Created At	TIMESTAMP	NOT NULL	Timestamp when the record was created.
15	updated_at	Updated At	TIMESTAMP	NULL allowed / business validation	Timestamp when the record was last updated.

Table 312: settlement_records Data Dictionary

49. triage_assessments

No	Field	Name	Type	Constraint	Description
1	assessment_id	Assessment ID	UUID	PK, NOT NULL	Primary key that uniquely identifies one assessment record.
2	user_id	User ID	UUID	FK	Foreign key that links this record to the related entity through user_id.
3	journey_id	Journey ID	UUID	FK	Foreign key that links this record to the related entity through journey_id.
4	baby_id	Baby ID	UUID	FK	Foreign key that links this record to the related entity through baby_id.
5	symptom_summary	Symptom Summary	TEXT	NULL allowed / business validation	Stores symptom summary for this assessment record.
6	risk_level	Risk Level	VARCHAR(10)	NULL allowed / business validation	Classification value for risk level.
7	recommended_action	Recommended Action	TEXT	NULL allowed / business validation	Stores recommended action for this assessment record.
8	rule_version	Rule Version	VARCHAR(50)	NULL allowed / business validation	Stores rule version for this assessment record.
9	disclaimer_accepted	Disclaimer Accepted	BOOLEAN	NULL allowed / business validation	Boolean flag indicating disclaimer accepted.
10	completed_at	Completed At	TIMESTAMP	NULL allowed / business validation	Business timestamp or date value for completed at.

11	created_at	Created At	TIMESTAMP	NOT NULL	Timestamp when the record was created.
12	updated_at	Updated At	TIMESTAMP	NULL allowed / business validation	Timestamp when the record was last updated.

Table 313: triage_assessments Data Dictionary

50. triage_answers

No	Field	Name	Type	Constraint	Description
1	triage_answer_id	Triage Answer ID	UUID	PK, NOT NULL	Primary key that uniquely identifies one triage_answer record.
2	assessment_id	Assessment ID	UUID	FK	Foreign key that links this record to the related entity through assessment_id.
3	question_code	Question Code	VARCHAR(80)	NULL allowed / business validation	Display or reference value for question code.
4	question_text	Question Text	VARCHAR(500)	NULL allowed / business validation	Stores question text for this triage_answer record.
5	answer_value	Answer Value	TEXT	NULL allowed / business validation	Stores answer value for this triage_answer record.
6	answer_order	Answer Order	INT	NULL allowed / business validation	Stores answer order for this triage_answer record.
7	red_flag_triggered	Red Flag Triggered	BOOLEAN	NULL allowed / business validation	Boolean flag indicating red flag triggered.
8	created_at	Created At	TIMESTAMP	NOT NULL	Timestamp when the record was created.
9	updated_at	Updated At	TIMESTAMP	NULL allowed / business validation	Timestamp when the record

					was last updated.
--	--	--	--	--	-------------------

Table 314: triage_answers Data Dictionary

51. partner_organizations

No	Field	Name	Type	Constraint	Description
1	partner_id	Partner ID	UUID	PK, NOT NULL	Primary key that uniquely identifies one partner record.
2	representative_user_id	Representative User ID	UUID	FK	Foreign key that links this record to the related entity through representative_user_id.
3	partner_type	Partner Type	VARCHAR(30)	NULL allowed / business validation	Classification value for partner type.
4	name	Name	VARCHAR(255)	NULL allowed / business validation	Display or reference value for name.
5	license_number	License Number	VARCHAR(120)	NULL allowed / business validation	Count or numeric reference used for license number.
6	address	Address	VARCHAR(500)	NULL allowed / business validation	Stores address for this partner record.
7	latitude	Latitude	DECIMAL(10,7)	NULL allowed / business validation	Location-related value used for latitude.
8	longitude	Longitude	DECIMAL(10,7)	NULL allowed / business validation	Location-related value used for longitude.
9	verification_status	Verification Status	VARCHAR(20)	NULL allowed / business validation	Current status value for verification status.
10	verified_by	Verified By	UUID	FK	Foreign key that links this record to the related entity

					through verified_by.
11	created_at	Created At	TIMESTAMP	NOT NULL	Timestamp when the record was created.
12	updated_at	Updated At	TIMESTAMP	NULL allowed / business validation	Timestamp when the record was last updated.

Table 315: partner_organizations Data Dictionary

52. partner_expert_links

No	Field	Name	Type	Constraint	Description
1	partner_expert_link_id	Partner Expert Link ID	UUID	PK, NOT NULL	Primary key that uniquely identifies one partner_expert_link record.
2	partner_id	Partner ID	UUID	FK	Foreign key that links this record to the related entity through partner_id.
3	expert_profile_id	Expert Profile ID	UUID	FK	Foreign key that links this record to the related entity through expert_profile_id.
4	relationship_type	Relationship Type	VARCHAR(40)	NULL allowed / business validation	Classification value for relationship type.
5	status	Status	VARCHAR(20)	NOT NULL	Current status value for status.
6	approved_by	Approved By	UUID	FK	Foreign key that links this record to the related entity through approved_by.
7	start_date	Start Date	DATE	NULL allowed / business validation	Business timestamp or date value for start date.

8	end_date	End Date	DATE	NULL allowed / business validation	Business timestamp or date value for end date.
9	created_at	Created At	TIMESTAMP	NOT NULL	Timestamp when the record was created.
10	updated_at	Updated At	TIMESTAMP	NULL allowed / business validation	Timestamp when the record was last updated.

Table 316: partner_expert_links Data Dictionary

53. partner_services

No	Field	Name	Type	Constraint	Description
1	service_id	Service ID	UUID	PK, NOT NULL	Primary key that uniquely identifies one service record.
2	partner_id	Partner ID	UUID	FK	Foreign key that links this record to the related entity through partner_id.
3	service_name	Service Name	VARCHAR(255)	NULL allowed / business validation	Display or reference value for service name.
4	description	Description	TEXT	NULL allowed / business validation	Descriptive text for the service record.
5	price_from	Price From	DECIMAL(15,2)	NULL allowed / business validation	Monetary value used for price from.
6	currency	Currency	CHAR(3)	NULL allowed / business validation	Stores currency for this service record.
7	booking_url	Booking URL	VARCHAR(500)	NULL allowed / business validation	URL or URI value for booking URL.
8	approval_status	Approval Status	VARCHAR(20)	NULL allowed / business validation	Current status value for approval status.
9	created_at	Created At	TIMESTAMP	NOT NULL	Timestamp when the record was created.

10	updated_at	Updated At	TIMESTAMP	NULL allowed / business validation	Timestamp when the record was last updated.
----	------------	------------	-----------	------------------------------------	---

Table 317: partner_services Data Dictionary

54. sponsored_campaigns

No	Field	Name	Type	Constraint	Description
1	campaign_id	Campaign ID	UUID	PK, NOT NULL	Primary key that uniquely identifies one campaign record.
2	partner_id	Partner ID	UUID	FK	Foreign key that links this record to the related entity through partner_id.
3	title	Title	VARCHAR(250)	NULL allowed / business validation	Display or reference value for title.
4	description	Description	TEXT	NULL allowed / business validation	Descriptive text for the campaign record.
5	start_date	Start Date	DATE	NULL allowed / business validation	Business timestamp or date value for start date.
6	end_date	End Date	DATE	NULL allowed / business validation	Business timestamp or date value for end date.
7	sponsor_label	Sponsor Label	VARCHAR(150)	NULL allowed / business validation	Stores sponsor label for this campaign record.
8	approval_status	Approval Status	VARCHAR(20)	NULL allowed / business validation	Current status value for approval status.
9	reviewed_by	Reviewed By	UUID	FK	Foreign key that links this record to the related entity through reviewed_by.

10	created_at	Created At	TIMESTAMP	NOT NULL	Timestamp when the record was created.
11	updated_at	Updated At	TIMESTAMP	NULL allowed / business validation	Timestamp when the record was last updated.

Table 318: sponsored_campaigns Data Dictionary

55. care_facilities

No	Field	Name	Type	Constraint	Description
1	facility_id	Facility ID	UUID	PK, NOT NULL	Primary key that uniquely identifies one facility record.
2	partner_id	Partner ID	UUID	FK	Foreign key that links this record to the related entity through partner_id.
3	name	Name	VARCHAR(255)	NULL allowed / business validation	Display or reference value for name.
4	facility_type	Facility Type	VARCHAR(40)	NULL allowed / business validation	Classification value for facility type.
5	address	Address	VARCHAR(500)	NULL allowed / business validation	Stores address for this facility record.
6	latitude	Latitude	DECIMAL(10,7)	NULL allowed / business validation	Location-related value used for latitude.
7	longitude	Longitude	DECIMAL(10,7)	NULL allowed / business validation	Location-related value used for longitude.
8	phone	Phone	VARCHAR(30)	NULL allowed / business validation	Phone number associated with this record.
9	opening_hours_json	Opening Hours Json	JSON	NULL allowed / business validation	Stores opening hours json for this facility record.

10	source_type	Source Type	VARCHAR(30)	NULL allowed / business validation	Classification value for source type.
11	verification_status	Verification Status	VARCHAR(20)	NULL allowed / business validation	Current status value for verification status.
12	created_at	Created At	TIMESTAMP	NOT NULL	Timestamp when the record was created.
13	updated_at	Updated At	TIMESTAMP	NULL allowed / business validation	Timestamp when the record was last updated.

Table 319: care_facilities Data Dictionary

56. emergency_events

No	Field	Name	Type	Constraint	Description
1	emergency_event_id	Emergency Event ID	UUID	PK, NOT NULL	Primary key that uniquely identifies one emergency_event record.
2	user_id	User ID	UUID	FK	Foreign key that links this record to the related entity through user_id.
3	source_type	Source Type	VARCHAR(30)	NULL allowed / business validation	Classification value for source type.
4	source_reference_id	Source Reference ID	UUID	NULL allowed / business validation	Stores source reference ID for this emergency_event record.
5	risk_level	Risk Level	VARCHAR(10)	NULL allowed / business validation	Classification value for risk level.
6	action_type	Action Type	VARCHAR(30)	NULL allowed / business validation	Classification value for action type.
7	selected_facility_id	Selected Facility ID	UUID	FK	Foreign key that links this record to the related entity

					through selected_facility_id.
8	selected_expert_id	Selected Expert ID	UUID	FK	Foreign key that links this record to the related entity through selected_expert_id.
9	status	Status	VARCHAR(20)	NOT NULL	Current status value for status.
10	opened_at	Opened At	TIMESTAMP	NULL allowed / business validation	Business timestamp or date value for opened at.
11	closed_at	Closed At	TIMESTAMP	NULL allowed / business validation	Business timestamp or date value for closed at.
12	created_at	Created At	TIMESTAMP	NOT NULL	Timestamp when the record was created.
13	updated_at	Updated At	TIMESTAMP	NULL allowed / business validation	Timestamp when the record was last updated.

Table 320: emergency_events Data Dictionary

57. location_snapshots

No	Field	Name	Type	Constraint	Description
1	location_snapshot_id	Location Snapshot ID	UUID	PK, NOT NULL	Primary key that uniquely identifies one location_snapshot record.
2	user_id	User ID	UUID	FK	Foreign key that links this record to the related entity through user_id.
3	context_type	Context Type	VARCHAR(30)	NULL allowed / business validation	Classification value for context type.

4	context_id	Context ID	UUID	NULL allowed / business validation	Stores context ID for this location_snapshot record.
5	latitude	Latitude	DECIMAL(10, 7)	NULL allowed / business validation	Location-related value used for latitude.
6	longitude	Longitude	DECIMAL(10, 7)	NULL allowed / business validation	Location-related value used for longitude.
7	accuracy_meters	Accuracy Meters	DECIMAL(10, 2)	NULL allowed / business validation	Stores accuracy meters for this location_snapshot record.
8	captured_at	Captured At	TIMESTAMP	NULL allowed / business validation	Business timestamp or date value for captured at.
9	expires_at	Expires At	TIMESTAMP	NULL allowed / business validation	Business timestamp or date value for expires at.
10	consent_status	Consent Status	VARCHAR(20)	NULL allowed / business validation	Current status value for consent status.

Table 321: location_snapshots Data Dictionary

58. health_device_connections

No	Field	Name	Type	Constraint	Description
1	connection_id	Connection ID	UUID	PK, NOT NULL	Primary key that uniquely identifies one connection record.
2	user_id	User ID	UUID	FK	Foreign key that links this record to the related entity through user_id.
3	provider_name	Provider Name	VARCHAR(80)	NULL allowed / business validation	Display or reference value for provider name.

4	device_name	Device Name	VARCHAR(150)	NULL allowed / business validation	Display or reference value for device name.
5	scopes_json	Scopes Json	JSON	NULL allowed / business validation	Stores scopes json for this connection record.
6	token_reference	Token Reference	VARCHAR(255)	NULL allowed / business validation	Secure token-related value used by the workflow.
7	consent_granted_at	Consent Granted At	TIMESTAMP	NULL allowed / business validation	Business timestamp or date value for consent granted at.
8	last_synced_at	Last Synced At	TIMESTAMP	NULL allowed / business validation	Business timestamp or date value for last synced at.
9	status	Status	VARCHAR(20)	NOT NULL	Current status value for status.
10	created_at	Created At	TIMESTAMP	NOT NULL	Timestamp when the record was created.
11	updated_at	Updated At	TIMESTAMP	NULL allowed / business validation	Timestamp when the record was last updated.

Table 322: health_device_connections Data Dictionary

59. device_measurements

No	Field	Name	Type	Constraint	Description
1	device_measurement_id	Device Measurement ID	UUID	PK, NOT NULL	Primary key that uniquely identifies one device_measurement record.
2	connection_id	Connection ID	UUID	FK	Foreign key that links this record to the related entity through connection_id.

3	measurement_type	Measurement Type	VARCHAR(40)	NULL allowed / business validation	Classification value for measurement type.
4	value_numeric	Value Numeric	DECIMAL(12, 3)	NULL allowed / business validation	Stores value numeric for this device_measurement record.
5	value_secondary	Value Secondary	DECIMAL(12, 3)	NULL allowed / business validation	Stores value secondary for this device_measurement record.
6	unit	Unit	VARCHAR(30)	NULL allowed / business validation	Stores unit for this device_measurement record.
7	measured_at	Measured At	TIMESTAMP	NULL allowed / business validation	Business timestamp or date value for measured at.
8	source_record_id	Source Record ID	VARCHAR(150)	NULL allowed / business validation	Stores source record ID for this device_measurement record.
9	quality_label	Quality Label	VARCHAR(30)	NULL allowed / business validation	Stores quality label for this device_measurement record.
10	raw_metadata_json	Raw Metadata Json	JSON	NULL allowed / business validation	Stores raw metadata json for this device_measurement record.
11	created_at	Created At	TIMESTAMP	NOT NULL	Timestamp when the record was created.
12	updated_at	Updated At	TIMESTAMP	NULL allowed / business validation	Timestamp when the record was last updated.

Table 323: device_measurements Data Dictionary

60. safety_monitoring_settings

No	Field	Name	Type	Constraint	Description
----	-------	------	------	------------	-------------

1	setting_id	Setting ID	UUID	PK, NOT NULL	Primary key that uniquely identifies one setting record.
2	user_id	User ID	UUID	FK	Foreign key that links this record to the related entity through user_id.
3	is_enabled	Is Enabled	BOOLEAN	NULL allowed / business validation	Boolean flag indicating is enabled.
4	countdown_seconds	Countdown Seconds	INT	NULL allowed / business validation	Count or numeric reference used for countdown seconds.
5	location_sharing_enabled	Location Sharing Enabled	BOOLEAN	NULL allowed / business validation	Boolean flag indicating location sharing enabled.
6	emergency_contact_user_id	Emergency Contact User ID	UUID	FK	Foreign key that links this record to the related entity through emergency_contact_user_id.
7	monitoring_schedule_json	Monitoring Schedule Json	JSON	NULL allowed / business validation	Stores monitoring schedule json for this setting record.
8	sensor_consent_at	Sensor Consent At	TIMESTAMP	NULL allowed / business validation	Business timestamp or date value for sensor consent at.
9	location_consent_at	Location Consent At	TIMESTAMP	NULL allowed / business validation	Business timestamp or date value for location consent at.
10	disclaimer_version	Disclaimer Version	VARCHAR(50)	NULL allowed / business validation	Stores disclaimer version for this setting record.

11	created_at	Created At	TIMESTAMP	NOT NULL	Timestamp when the record was created.
12	updated_at	Updated At	TIMESTAMP	NULL allowed / business validation	Timestamp when the record was last updated.

Table 324: safety_monitoring_settings Data Dictionary

61. safety_events

No	Field	Name	Type	Constraint	Description
1	safety_event_id	Safety Event ID	UUID	PK, NOT NULL	Primary key that uniquely identifies one safety_event record.
2	user_id	User ID	UUID	FK	Foreign key that links this record to the related entity through user_id.
3	setting_id	Setting ID	UUID	FK	Foreign key that links this record to the related entity through setting_id.
4	detected_at	Detected At	TIMESTAMP	NULL allowed / business validation	Business timestamp or date value for detected at.
5	event_type	Event Type	VARCHAR(30)	NULL allowed / business validation	Classification value for event type.
6	confidence_score	Confidence Score	DECIMAL(5,4)	NULL allowed / business validation	Numeric value used for confidence score.
7	peak_acceleration	Peak Acceleration	DECIMAL(10,4)	NULL allowed / business validation	Stores peak acceleration for this safety_event record.
8	angular_velocity	Angular Velocity	DECIMAL(10,4)	NULL allowed / business validation	Stores angular velocity for this

					safety_event record.
9	inactivity_seconds	Inactivity Seconds	INT	NULL allowed / business validation	Stores inactivity seconds for this safety_event record.
10	user_response	User Response	VARCHAR(30)	NULL allowed / business validation	Stores user response for this safety_event record.
11	response_at	Response At	TIMESTAMP	NULL allowed / business validation	Business timestamp or date value for response at.
12	false_positive_reason	False Positive Reason	VARCHAR(80)	NULL allowed / business validation	Text content recorded for false positive reason.
13	status	Status	VARCHAR(20)	NOT NULL	Current status value for status.
14	created_at	Created At	TIMESTAMP	NOT NULL	Timestamp when the record was created.
15	updated_at	Updated At	TIMESTAMP	NULL allowed / business validation	Timestamp when the record was last updated.

Table 325: safety_events Data Dictionary

62. safety_alerts

No	Field	Name	Type	Constraint	Description
1	safety_alert_id	Safety Alert ID	UUID	PK, NOT NULL	Primary key that uniquely identifies one safety_alert record.
2	safety_event_id	Safety Event ID	UUID	FK	Foreign key that links this record to the related entity through safety_event_id.
3	recipient_user_id	Recipient User ID	UUID	FK	Foreign key that links this

					record to the related entity through recipient_user_id.
4	location_snapshot_id	Location Snapshot ID	UUID	FK	Foreign key that links this record to the related entity through location_snapshot_id.
5	alert_reason	Alert Reason	VARCHAR(30)	NULL allowed / business validation	Text content recorded for alert reason.
6	payload_json	Payload Json	JSON	NULL allowed / business validation	Stores payload json for this safety_alert record.
7	delivery_status	Delivery Status	VARCHAR(20)	NULL allowed / business validation	Current status value for delivery status.
8	sent_at	Sent At	TIMESTAMP	NULL allowed / business validation	Business timestamp or date value for sent at.
9	acknowledged_at	Acknowledged At	TIMESTAMP	NULL allowed / business validation	Business timestamp or date value for acknowledged at.
10	created_at	Created At	TIMESTAMP	NOT NULL	Timestamp when the record was created.
11	updated_at	Updated At	TIMESTAMP	NULL allowed / business validation	Timestamp when the record was last updated.

Table 326: safety_alerts Data Dictionary

63. pregnancy_exercises

No	Field	Name	Type	Constraint	Description
1	exercise_id	Exercise ID	UUID	PK, NOT NULL	Primary key that uniquely identifies one exercise record.

2	created_by	Created By	UUID	FK	Foreign key that links this record to the related entity through created_by.
3	title	Title	VARCHAR(200)	NULL allowed / business validation	Display or reference value for title.
4	description	Description	TEXT	NULL allowed / business validation	Descriptive text for the exercise record.
5	trimester_scope	Trimester Scope	VARCHAR(30)	NULL allowed / business validation	Stores trimester scope for this exercise record.
6	difficulty_level	Difficulty Level	VARCHAR(20)	NULL allowed / business validation	Classification value for difficulty level.
7	duration_minutes	Duration Minutes	INT	NULL allowed / business validation	Stores duration minutes for this exercise record.
8	instruction_content	Instruction Content	TEXT	NULL allowed / business validation	Text content recorded for instruction content.
9	media_url	Media URL	VARCHAR(500)	NULL allowed / business validation	URL or URI value for media URL.
10	safety_warning	Safety Warning	TEXT	NULL allowed / business validation	Stores safety warning for this exercise record.
11	supports_posture_analysis	Supports Posture Analysis	BOOLEAN	NULL allowed / business validation	Boolean flag indicating supports posture analysis.
12	status	Status	VARCHAR(20)	NOT NULL	Current status value for status.
13	version_no	Version No	INT	NULL allowed / business validation	Stores version no for this exercise record.
14	created_at	Created At	TIMESTAMP	NOT NULL	Timestamp when the record was created.
15	updated_at	Updated At	TIMESTAMP	NULL allowed / business validation	Timestamp when the record was last updated.

Table 327: pregnancy_exercises Data Dictionary

64. exercise_safety_checks

No	Field	Name	Type	Constraint	Description
1	safety_check_id	Safety Check ID	UUID	PK, NOT NULL	Primary key that uniquely identifies one safety_check record.
2	exercise_id	Exercise ID	UUID	FK	Foreign key that links this record to the related entity through exercise_id.
3	journey_id	Journey ID	UUID	FK	Foreign key that links this record to the related entity through journey_id.
4	user_id	User ID	UUID	FK	Foreign key that links this record to the related entity through user_id.
5	answer_json	Answer Json	JSON	NULL allowed / business validation	Stores answer json for this safety_check record.
6	red_flag_detected	Red Flag Detected	BOOLEAN	NULL allowed / business validation	Boolean flag indicating red flag detected.
7	result_status	Result Status	VARCHAR(20)	NULL allowed / business validation	Current status value for result status.
8	blocked_reason	Blocked Reason	TEXT	NULL allowed / business validation	Text content recorded for blocked reason.
9	completed_at	Completed At	TIMESTAMP	NULL allowed / business validation	Business timestamp or date value for completed at.
10	created_at	Created At	TIMESTAMP	NOT NULL	Timestamp when the record was created.

Table 328: exercise_safety_checks Data Dictionary

65. exercise_sessions

No	Field	Name	Type	Constraint	Description
1	exercise_session_id	Exercise Session ID	UUID	PK, NOT NULL	Primary key that uniquely identifies one exercise_session record.
2	exercise_id	Exercise ID	UUID	FK	Foreign key that links this record to the related entity through exercise_id.
3	journey_id	Journey ID	UUID	FK	Foreign key that links this record to the related entity through journey_id.
4	user_id	User ID	UUID	FK	Foreign key that links this record to the related entity through user_id.
5	safety_check_id	Safety Check ID	UUID	FK	Foreign key that links this record to the related entity through safety_check_id.
6	started_at	Started At	TIMESTAMP	NULL allowed / business validation	Business timestamp or date value for started at.
7	ended_at	Ended At	TIMESTAMP	NULL allowed / business validation	Business timestamp or date value for ended at.
8	paused_seconds	Paused Seconds	INT	NULL allowed / business validation	Stores paused seconds for this exercise_session record.

9	completion_percent	Completion Percent	DECIMAL(5,2)	NULL allowed / business validation	Numeric value used for completion percent.
10	posture_score	Posture Score	DECIMAL(5,2)	NULL allowed / business validation	Numeric value used for posture score.
11	session_status	Session Status	VARCHAR(20)	NULL allowed / business validation	Current status value for session status.
12	warning_count	Warning Count	INT	NULL allowed / business validation	Count or numeric reference used for warning count.
13	summary_json	Summary Json	JSON	NULL allowed / business validation	Stores summary json for this exercise_session record.
14	created_at	Created At	TIMESTAMP	NOT NULL	Timestamp when the record was created.
15	updated_at	Updated At	TIMESTAMP	NULL allowed / business validation	Timestamp when the record was last updated.

Table 329: exercise_sessions Data Dictionary

66. posture_analysis_configs

No	Field	Name	Type	Constraint	Description
1	posture_config_id	Posture Config ID	UUID	PK, NOT NULL	Primary key that uniquely identifies one posture_config record.
2	exercise_id	Exercise ID	UUID	FK	Foreign key that links this record to the related entity through exercise_id.
3	configured_by	Configured By	UUID	FK	Foreign key that links this record to the related entity through configured_by.

4	analysis_mode	Analysis Mode	VARCHAR(20)	NULL allowed / business validation	Classification value for analysis mode.
5	rule_or_model_version	Rule Or Model Version	VARCHAR(10)	NULL allowed / business validation	Classification value for rule or model version.
6	confidence_threshold	Confidence Threshold	DECIMAL(5,4)	CHECK 0..1	Stores confidence threshold for this posture_config record.
7	feedback_level	Feedback Level	VARCHAR(20)	NULL allowed / business validation	Monetary value used for feedback level.
8	config_json	Config Json	JSON	NULL allowed / business validation	Stores config json for this posture_config record.
9	effective_from	Effective From	TIMESTAMP	NULL allowed / business validation	Business timestamp or date value for effective from.
10	effective_to	Effective To	TIMESTAMP	NULL allowed / business validation	Business timestamp or date value for effective to.
11	status	Status	VARCHAR(20)	NOT NULL	Current status value for status.
12	created_at	Created At	TIMESTAMP	NOT NULL	Timestamp when the record was created.
13	updated_at	Updated At	TIMESTAMP	NULL allowed / business validation	Timestamp when the record was last updated.

Table 330: posture_analysis_configs Data Dictionary

67. posture_feedback_events

No	Field	Name	Type	Constraint	Description
1	feedback_event_id	Feedback Event ID	UUID	PK, NOT NULL	Primary key that uniquely identifies one feedback_event record.
2	exercise_session_id	Exercise Session ID	UUID	FK	Foreign key that links this

					record to the related entity through exercise_session_id.
3	posture_config_id	Posture Config ID	UUID	FK	Foreign key that links this record to the related entity through posture_config_id.
4	event_time_ms	Event Time Ms	BIGINT	NULL allowed / business validation	Stores event time ms for this feedback_event record.
5	posture_code	Posture Code	VARCHAR(80)	NULL allowed / business validation	Display or reference value for posture code.
6	confidence_score	Confidence Score	DECIMAL(5,4)	NULL allowed / business validation	Numeric value used for confidence score.
7	severity	Severity	VARCHAR(20)	NULL allowed / business validation	Stores severity for this feedback_event record.
8	feedback_text	Feedback Text	VARCHAR(500)	NULL allowed / business validation	Monetary value used for feedback text.
9	keypoint_summary_json	Keypoint Summary Json	JSON	NULL allowed / business validation	Stores keypoint summary json for this feedback_event record.
10	created_at	Created At	TIMESTAMP	NOT NULL	Timestamp when the record was created.

Table 331: posture_feedback_events Data Dictionary

4.2 Mobile Screen

This section specifies Release 1 mobile screens for Mother, Family Member and Verified Expert roles. Each entry is a detailed, low-fidelity SRS artifact. Existing P0 core screens are ordered first; P1 supporting screens follow. Search, filters, confirmations and protected attachment preview are embedded where they belong rather than being treated as separate screen specifications.

4.2.1 Common Mobile Account Screens

4.2.1.1 Mobile Welcome Screen

Platform: Common Mobile App Feature: MF-01 Account, Trust & Access Control Release: PO
Related UCs: UC-01, UC-03

Purpose: Introduces CareBridge and provides safe entry points for registration or login.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Common Mobile App. Primary role: Guest / Authenticated User. Screen name: "Mobile Welcome". Purpose: Introduces CareBridge and provides safe entry points for registration or login. Related Release 1 use cases: UC-01, UC-03. Use a single portrait mobile frame, 390 x 844 px, with a 24 px safe margin. Do not render signed-in bottom navigation on this guest account screen. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: CareBridge Wordmark: Identifies the application at the top of the entry screen.; Short Purpose Statement: Explains that CareBridge supports maternal and early-childhood care without replacing medical care.; Register Button: Opens the role-appropriate account registration flow.; Login Button: Opens the account login flow.; Terms, Privacy and Safety Links: Opens required terms, data privacy and medical-safety notices.. Interaction and safety constraints: Show only controls authorized for Guest / Authenticated User. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Keep the entry experience short and do not expose private health functions before authentication. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.2.1: Mobile Welcome Screen (to be generated in Stitch)

Field name	Description
CareBridge Wordmark	Identifies the application at the top of the entry screen.
Short Purpose Statement	Explains that CareBridge supports maternal and early-childhood care without replacing medical care.
Register Button	Opens the role-appropriate account registration flow.

Field name	Description
Login Button	Opens the account login flow.
Terms, Privacy and Safety Links	Opens required terms, data privacy and medical-safety notices.

Table 4.2.1: Mobile Welcome Screen

4.2.1.2 Register Account Screen

Platform: Common Mobile App Feature: MF-01 Account, Trust & Access Control Release: PO
Related UCs: UC-01

Purpose: Collects supported account credentials and the initial role needed to create a pending account.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Common Mobile App. Primary role: Guest. Screen name: "Register Account". Purpose: Collects supported account credentials and the initial role needed to create a pending account. Related Release 1 use cases: UC-01. Use a single portrait mobile frame, 390 x 844 px, with a 24 px safe margin. Do not render signed-in bottom navigation on this guest account screen. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Full Name Field: Collects the private account holder name.; Email or Phone Field: Collects the supported unique contact used for OTP and recovery.; Initial Role Selector: Selects Mother, Family Member or Expert Applicant where applicable.; Password and Confirm Password: Collects a policy-compliant password and confirmation.; Terms and Privacy Checkbox: Records consent to account, privacy and safety terms.; Create Account Button: Validates the form and requests OTP delivery.. Interaction and safety constraints: Show only controls authorized for Guest. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Do not expose community identity or health-data settings during registration. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.2.2: Register Account Screen (to be generated in Stitch)

Field name	Description
Full Name Field	Collects the private account holder name.
Email or Phone Field	Collects the supported unique contact used for OTP and recovery.
Initial Role Selector	Selects Mother, Family Member or Expert Applicant where applicable.
Password and Confirm Password	Collects a policy-compliant password and confirmation.
Terms and Privacy Checkbox	Records consent to account, privacy and safety terms.
Create Account Button	Validates the form and requests OTP delivery.

Table 4.2.2: Register Account Screen

4.2.1.3 Verify OTP Screen

Platform: Common Mobile App Feature: MF-01 Account, Trust & Access Control Release: P0
Related UCs: UC-02

Purpose: Validates a one-time password for account activation or another protected account action.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Common Mobile App. Primary role: Guest. Screen name: "Verify OTP". Purpose: Validates a one-time password for account activation or another protected account action. Related Release 1 use cases: UC-02. Use a single portrait mobile frame, 390 x 844 px, with a 24 px safe margin. Do not render signed-in bottom navigation on this guest account screen. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Masked Contact: Shows the destination contact without revealing the full value.; OTP Cells: Accepts the verification code in clearly separated numeric cells.; Expiry and Resend State: Shows the code validity period and whether resend is available.; Verify Button: Submits the OTP and continues only after validation.; Back Control: Returns to the safe previous account screen.. Interaction and safety constraints: Show only controls authorized for Guest. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination

inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Show attempt and resend limits in a compact helper message. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.2.3: Verify OTP Screen (to be generated in Stitch)

Field name	Description
Masked Contact	Shows the destination contact without revealing the full value.
OTP Cells	Accepts the verification code in clearly separated numeric cells.
Expiry and Resend State	Shows the code validity period and whether resend is available.
Verify Button	Submits the OTP and continues only after validation.
Back Control	Returns to the safe previous account screen.

Table 4.2.3: Verify OTP Screen

4.2.1.4 Login Screen

Platform: Common Mobile App Feature: MF-01 Account, Trust & Access Control Release: P0
Related UCs: UC-03

Purpose: Authenticates an approved account and routes the user to the correct role workspace.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Common Mobile App. Primary role: Guest. Screen name: "Login". Purpose: Authenticates an approved account and routes the user to the correct role workspace. Related Release 1 use cases: UC-03. Use a single portrait mobile frame, 390 x 844 px, with a 24 px safe margin. Do not render signed-in bottom navigation on this guest account screen. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Email or Phone Field: Collects the registered contact identifier.; Password Field: Collects the masked password with show or hide control.; Remember Session Toggle: Shows the supported session persistence choice.; Forgot Password Link: Opens the recovery flow.; Login Button: Authenticates credentials and routes by role.; Register Link: Returns a guest to registration.. Interaction and safety constraints: Show only controls authorized for Guest. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive

health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Do not reveal whether an account exists when credentials are invalid. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.2.4: Login Screen (to be generated in Stitch)

Field name	Description
Email or Phone Field	Collects the registered contact identifier.
Password Field	Collects the masked password with show or hide control.
Remember Session Toggle	Shows the supported session persistence choice.
Forgot Password Link	Opens the recovery flow.
Login Button	Authenticates credentials and routes by role.
Register Link	Returns a guest to registration.

Table 4.2.4: Login Screen

4.2.1.5 Request Password Reset Screen

Platform: Common Mobile App Feature: MF-01 Account, Trust & Access Control Release: P0
Related UCs: UC-05

Purpose: Starts password recovery through the registered contact channel.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Common Mobile App. Primary role: Guest. Screen name: "Request Password Reset". Purpose: Starts password recovery through the registered contact channel. Related Release 1 use cases: UC-05. Use a single portrait mobile frame, 390 x 844 px, with a 24 px safe margin. Do not render signed-in bottom navigation on this guest account screen. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Email or Phone Field: Accepts the registered recovery

contact.; Recovery Explanation: Explains that a time-limited recovery proof will be sent if eligible.; Send Recovery Button: Requests reset proof subject to rate limits.; Back to Login Link: Returns to login.. Interaction and safety constraints: Show only controls authorized for Guest. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Use a neutral completion message to avoid account enumeration. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.2.5: Request Password Reset Screen (to be generated in Stitch)

Field name	Description
Email or Phone Field	Accepts the registered recovery contact.
Recovery Explanation	Explains that a time-limited recovery proof will be sent if eligible.
Send Recovery Button	Requests reset proof subject to rate limits.
Back to Login Link	Returns to login.

Table 4.2.5: Request Password Reset Screen

4.2.1.6 Reset Password Screen

Platform: Common Mobile App Feature: MF-01 Account, Trust & Access Control Release: P0
Related UCs: UC-06

Purpose: Sets a new password after a valid recovery proof is confirmed.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Common Mobile App. Primary role: Guest. Screen name: "Reset Password". Purpose: Sets a new password after a valid recovery proof is confirmed. Related Release 1 use cases: UC-06. Use a single portrait mobile frame, 390 x 844 px, with a 24 px safe margin. Do not render signed-in bottom navigation on this guest account screen. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Recovery Proof Status: Shows that the reset request is valid

and time-limited.; New Password Field: Collects the new policy-compliant password.; Confirm Password Field: Confirms the new password.; Password Guidance: Shows compact password policy requirements.; Reset Password Button: Completes the password reset and returns to login..

Interaction and safety constraints: Show only controls authorized for Guest. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Do not display any private profile information on this screen. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.2.6: Reset Password Screen (to be generated in Stitch)

Field name	Description
Recovery Proof Status	Shows that the reset request is valid and time-limited.
New Password Field	Collects the new policy-compliant password.
Confirm Password Field	Confirms the new password.
Password Guidance	Shows compact password policy requirements.
Reset Password Button	Completes the password reset and returns to login.

Table 4.2.6: Reset Password Screen

4.2.1.7 Account Profile and Settings Screen

Platform: Common Mobile App Feature: MF-01 Account, Trust & Access Control Release: P0
Related UCs: UC-04, UC-07, UC-08, UC-09, UC-10, UC-11, UC-14

Purpose: Shows private account details and groups account editing, community identity, notifications, password, logout and lifecycle actions.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Common Mobile App. Primary role: Authenticated User. Screen name: "Account Profile and Settings". Purpose: Shows private account details and groups account editing, community identity, notifications, password, logout and lifecycle actions. Related Release 1 use cases: UC-04, UC-07, UC-08, UC-09, UC-10, UC-11, UC-14. Use a single portrait mobile frame, 390 x 844 px, with a 24 px safe margin. Use a simple top app bar and bottom navigation only for signed-in main-application screens. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white

background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Profile Summary: Shows private display name, role, account status and masked contact information.; Edit Account Action: Opens allowed non-sensitive private account fields.; Community Identity Shortcut: Opens public name, avatar and visibility settings separated from private health data.; Notification and Password Shortcuts: Opens notification preferences and password change controls.; Logout Action: Ends the active device session after an in-screen confirmation state.; Deactivate or Delete Action: Starts account lifecycle action with a retention and care-group impact notice.. Interaction and safety constraints: Show only controls authorized for Authenticated User. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Keep community identity visually separated from account contact information. Destructive actions are in-screen confirmation states, not standalone screens. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.2.7: Account Profile and Settings Screen (to be generated in Stitch)

Field name	Description
Profile Summary	Shows private display name, role, account status and masked contact information.
Edit Account Action	Opens allowed non-sensitive private account fields.
Community Identity Shortcut	Opens public name, avatar and visibility settings separated from private health data.
Notification and Password Shortcuts	Opens notification preferences and password change controls.
Logout Action	Ends the active device session after an in-screen confirmation state.
Deactivate or Delete Action	Starts account lifecycle action with a retention and care-group impact notice.

Table 4.2.7: Account Profile and Settings Screen

4.2.1.8 Notification Center Screen

Platform: Common Mobile App Feature: MF-01 Account, Trust & Access Control Release: P0
Related UCs: UC-12

Purpose: Shows authorized reminders, replies, family alerts and account events with read-state controls.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Common Mobile App. Primary role: Authenticated User. Screen name: "Notification Center". Purpose: Shows authorized reminders, replies, family alerts and account events with read-state controls. Related Release 1 use cases: UC-12. Use a single portrait mobile frame, 390 x 844 px, with a 24 px safe margin. Use a simple top app bar and bottom navigation only for signed-in main-application screens. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Category Tabs or Chips: Filters permitted notification groups such as reminders, community, family and account.; Unread Count: Shows the number of currently unread authorized notices.; Notification List: Shows title, short message, timestamp, type and read state.; Mark All as Read: Marks all eligible notifications as read.; Empty State: Explains when no notifications match the active category.. Interaction and safety constraints: Show only controls authorized for Authenticated User. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Open only destinations the signed-in user is authorized to access. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.2.8: Notification Center Screen (to be generated in Stitch)

Field name	Description
Category Tabs or Chips	Filters permitted notification groups such as reminders, community, family and account.
Unread Count	Shows the number of currently unread authorized notices.
Notification List	Shows title, short message, timestamp, type and read state.

Field name	Description
Mark All as Read	Marks all eligible notifications as read.
Empty State	Explains when no notifications match the active category.

Table 4.2.8: Notification Center Screen

4.2.1.9 Login Sessions Screen

Platform: Common Mobile App Feature: MF-01 Account, Trust & Access Control Release: PO
Related UCs: UC-13

Purpose: Lists active or recent sessions and allows the owner to revoke a selected non-current session.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Common Mobile App. Primary role: Authenticated User. Screen name: "Login Sessions". Purpose: Lists active or recent sessions and allows the owner to revoke a selected non-current session. Related Release 1 use cases: UC-13. Use a single portrait mobile frame, 390 x 844 px, with a 24 px safe margin. Use a simple top app bar and bottom navigation only for signed-in main-application screens. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Current Session Card: Identifies the active device and protects it from accidental removal.; Other Sessions List: Shows device, platform, approximate time and last activity for eligible sessions.; Revoke Session Action: Signs out a selected device after an in-screen confirmation state.; Security Hint: Explains that suspicious sessions should be revoked and the password changed.. Interaction and safety constraints: Show only controls authorized for Authenticated User. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Do not expose precise location or device identifiers beyond what is needed for security. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.2.9: Login Sessions Screen (to be generated in Stitch)

Field name	Description
Current Session Card	Identifies the active device and protects it from accidental removal.
Other Sessions List	Shows device, platform, approximate time and last activity for eligible sessions.
Revoke Session Action	Signs out a selected device after an in-screen confirmation state.
Security Hint	Explains that suspicious sessions should be revoked and the password changed.

Table 4.2.9: Login Sessions Screen

4.2.1.10 Privacy, Consent and Data Sharing Screen

Platform: Common Mobile App Feature: MF-01 Account, Trust & Access Control Release: PO
Related UCs: UC-15, UC-16

Purpose: Creates, reviews and revokes purpose-specific, scoped and time-limited sharing permissions.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Common Mobile App. Primary role: Authenticated User. Screen name: "Privacy, Consent and Data Sharing". Purpose: Creates, reviews and revokes purpose-specific, scoped and time-limited sharing permissions. Related Release 1 use cases: UC-15, UC-16. Use a single portrait mobile frame, 390 x 844 px, with a 24 px safe margin. Use a simple top app bar and bottom navigation only for signed-in main-application screens. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Active Consent List: Shows recipient, selected data scope, purpose, expiry and current status.; Grant New Permission Action: Starts a guided permission form for a family member or verified expert.; Scope Selector: Selects permitted record, summary, log, calendar or alert categories.; Purpose and Expiry Fields: Records why access is granted and when it will expire.; Revoke Action: Revokes an active grant after a clear impact message.. Interaction and safety constraints: Show only controls authorized for Authenticated User. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a

standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Use privacy by default and do not preselect sensitive health, location or device data. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.2.10: Privacy, Consent and Data Sharing Screen (to be generated in Stitch)

Field name	Description
Active Consent List	Shows recipient, selected data scope, purpose, expiry and current status.
Grant New Permission Action	Starts a guided permission form for a family member or verified expert.
Scope Selector	Selects permitted record, summary, log, calendar or alert categories.
Purpose and Expiry Fields	Records why access is granted and when it will expire.
Revoke Action	Revokes an active grant after a clear impact message.

Table 4.2.10: Privacy, Consent and Data Sharing Screen

4.2.2 Mother Mobile App

4.2.2.1 Mother Journey Setup Screen

Platform: Mother Mobile App Feature: MF-02 Mother Care Journey Release: P0 Related UCs: UC-19

Purpose: Initializes the care journey for preconception, pregnancy or postpartum recovery.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Mother Mobile App. Primary role: Mother. Screen name: "Mother Journey Setup". Purpose: Initializes the care journey for preconception, pregnancy or postpartum recovery. Related Release 1 use cases: UC-19. Use a single portrait mobile frame, 390 x 844 px, with a 24 px safe margin. Use a simple top app bar and bottom navigation only for signed-in main-application screens. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Journey Stage Selector: Selects preconception, pregnancy or postpartum recovery.; Key Date Fields: Captures only minimum dates needed for stage-based support.; Privacy Notice: Explains that journey data is private and can be shared only through consent.; Continue Button: Creates the

journey and opens the dashboard.. Interaction and safety constraints: Show only controls authorized for Mother. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Do not imply clinical validation of entered dates. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.2.11: Mother Journey Setup Screen (to be generated in Stitch)

Field name	Description
Journey Stage Selector	Selects preconception, pregnancy or postpartum recovery.
Key Date Fields	Captures only minimum dates needed for stage-based support.
Privacy Notice	Explains that journey data is private and can be shared only through consent.
Continue Button	Creates the journey and opens the dashboard.

Table 4.2.11: Mother Journey Setup Screen

4.2.2.2 Mother Journey Dashboard Screen

Platform: Mother Mobile App Feature: MF-02 Mother Care Journey Release: P0 Related UCs: UC-20, UC-21

Purpose: Shows current stage or week, relevant care priorities, reminders, checklists, reviewed content and safe shortcuts.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Mother Mobile App. Primary role: Mother. Screen name: "Mother Journey Dashboard". Purpose: Shows current stage or week, relevant care priorities, reminders, checklists, reviewed content and safe shortcuts. Related Release 1 use cases: UC-20, UC-21. Use a single portrait mobile frame, 390 x 844 px, with a 24 px safe margin. Use a simple top app bar and bottom navigation only for signed-in main-application screens. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or

personal data. Required components: Stage and Week Summary: Shows the current journey stage, week or postpartum period from user-entered dates.; Update Journey Action: Opens editing for permitted stage and date changes.; Today Priorities: Shows due reminders, checklist items and safe care prompts.; Quick Modules: Links to metrics, postpartum logs, exercises, baby care and community.; Reviewed Content Card: Shows one stage-relevant approved content item with source label.. Interaction and safety constraints: Show only controls authorized for Mother. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Keep clinical language non-diagnostic and show only user-authorized data. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.2.12: Mother Journey Dashboard Screen (to be generated in Stitch)

Field name	Description
Stage and Week Summary	Shows the current journey stage, week or postpartum period from user-entered dates.
Update Journey Action	Opens editing for permitted stage and date changes.
Today Priorities	Shows due reminders, checklist items and safe care prompts.
Quick Modules	Links to metrics, postpartum logs, exercises, baby care and community.
Reviewed Content Card	Shows one stage-relevant approved content item with source label.

Table 4.2.12: Mother Journey Dashboard Screen

4.2.2.3 Maternal Metrics and Trend Screen

Platform: Mother Mobile App Feature: MF-02 Mother Care Journey Release: P0 Related UCs: UC-22, UC-23, UC-24

Purpose: Records, corrects, removes and reviews user-entered maternal observations and non-diagnostic trends.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Mother Mobile App. Primary role: Mother. Screen name: "Maternal Metrics and Trend". Purpose: Records, corrects, removes and reviews user-entered maternal observations and non-diagnostic trends. Related Release 1 use cases: UC-22, UC-23, UC-24. Use a single portrait mobile frame, 390 x 844 px, with a 24 px safe margin. Use a simple top app bar and bottom navigation only for signed-in main-application screens. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Indicator Selector: Selects a supported metric such as weight, blood pressure, glucose or fetal movement note.; Add Metric Form: Collects value, unit, date, source and optional note.; Trend Chart: Displays a selected indicator over time with clear user-entered source labels.; Metric History List: Shows existing entries with edit or delete controls for owner-entered data.; Safety Note: States that trends do not diagnose conditions and may require professional review.. Interaction and safety constraints: Show only controls authorized for Mother. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Embed date and indicator filters within this screen; do not create a separate metric search screen. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.2.13: Maternal Metrics and Trend Screen (to be generated in Stitch)

Field name	Description
Indicator Selector	Selects a supported metric such as weight, blood pressure, glucose or fetal movement note.
Add Metric Form	Collects value, unit, date, source and optional note.
Trend Chart	Displays a selected indicator over time with clear user-entered source labels.
Metric History List	Shows existing entries with edit or delete controls for owner-entered data.

Field name	Description
Safety Note	States that trends do not diagnose conditions and may require professional review.

Table 4.2.13: Maternal Metrics and Trend Screen

4.2.2.4 Postpartum Recovery Logs Screen

Platform: Mother Mobile App Feature: MF-02 Mother Care Journey Release: P0 Related UCs: UC-25, UC-26

Purpose: Records and manages postpartum recovery observations such as sleep, mood, pain note or feeding context.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Mother Mobile App. Primary role: Mother. Screen name: "Postpartum Recovery Logs". Purpose: Records and manages postpartum recovery observations such as sleep, mood, pain note or feeding context. Related Release 1 use cases: UC-25, UC-26. Use a single portrait mobile frame, 390 x 844 px, with a 24 px safe margin. Use a simple top app bar and bottom navigation only for signed-in main-application screens. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Log Type Chips: Selects a supported recovery observation type.; New Log Entry: Collects date, severity or context, note and optional source label.; Recovery Log Timeline: Shows owner-entered log entries with edit and delete controls.; Private-by-Default Notice: Explains that logs are private unless the owner separately shares them.; Safety Escalation Hint: Directs users to professional or emergency care for concerning symptoms.. Interaction and safety constraints: Show only controls authorized for Mother. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Avoid diagnostic language and do not present a medical treatment plan. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.2.14: Postpartum Recovery Logs Screen (to be generated in Stitch)

Field name	Description
Log Type Chips	Selects a supported recovery observation type.
New Log Entry	Collects date, severity or context, note and optional source label.
Recovery Log Timeline	Shows owner-entered log entries with edit and delete controls.
Private-by-Default Notice	Explains that logs are private unless the owner separately shares them.
Safety Escalation Hint	Directs users to professional or emergency care for concerning symptoms.

Table 4.2.14: Postpartum Recovery Logs Screen

4.2.2.5 Preparation Checklist Screen

Platform: Mother Mobile App Feature: MF-02 Mother Care Journey Release: P0 Related UCs: UC-27

Purpose: Creates, edits and completes preparation items for pregnancy, hospital bag, postpartum recovery and early baby care.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Mother Mobile App. Primary role: Mother. Screen name: "Preparation Checklist". Purpose: Creates, edits and completes preparation items for pregnancy, hospital bag, postpartum recovery and early baby care. Related Release 1 use cases: UC-27. Use a single portrait mobile frame, 390 x 844 px, with a 24 px safe margin. Use a simple top app bar and bottom navigation only for signed-in main-application screens. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Checklist Category Tabs: Switches between preparation categories.; Checklist Item List: Shows status, due context and optional note for each item.; Add or Edit Item Action: Creates or updates an owner-managed checklist item.; Completion Toggle: Marks an item complete or incomplete.; Progress Summary: Shows completed and remaining item counts.. Interaction and safety constraints: Show only controls authorized for Mother. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do

not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Do not frame checklist suggestions as medical orders. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.2.15: Preparation Checklist Screen (to be generated in Stitch)

Field name	Description
Checklist Category Tabs	Switches between preparation categories.
Checklist Item List	Shows status, due context and optional note for each item.
Add or Edit Item Action	Creates or updates an owner-managed checklist item.
Completion Toggle	Marks an item complete or incomplete.
Progress Summary	Shows completed and remaining item counts.

Table 4.2.15: Preparation Checklist Screen

4.2.2.6 Pregnancy Exercise Library Screen

Platform: Mother Mobile App Feature: MF-02 Mother Care Journey Release: P0 Related UCs: UC-28

Purpose: Browses reviewed pregnancy exercises with embedded stage, difficulty and duration filters.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Mother Mobile App. Primary role: Mother. Screen name: "Pregnancy Exercise Library". Purpose: Browses reviewed pregnancy exercises with embedded stage, difficulty and duration filters. Related Release 1 use cases: UC-28. Use a single portrait mobile frame, 390 x 844 px, with a 24 px safe margin. Use a simple top app bar and bottom navigation only for signed-in main-application screens. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Search and Filter Bar: Embeds keyword, stage, difficulty and duration controls in the library screen.; Exercise Cards: Show reviewed title, suitable stage, duration, difficulty and compact safety tag.; Exercise Detail Panel: Shows instructions, contraindication note and start action for the selected exercise.; Safety Reminder: States that the library does not assess medical fitness or replace supervision.. Interaction and safety constraints: Show only controls authorized for Mother. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data,

protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Do not add a standalone exercise-search screen. Show only approved exercise content. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.2.16: Pregnancy Exercise Library Screen (to be generated in Stitch)

Field name	Description
Search and Filter Bar	Embeds keyword, stage, difficulty and duration controls in the library screen.
Exercise Cards	Show reviewed title, suitable stage, duration, difficulty and compact safety tag.
Exercise Detail Panel	Shows instructions, contraindication note and start action for the selected exercise.
Safety Reminder	States that the library does not assess medical fitness or replace supervision.

Table 4.2.16: Pregnancy Exercise Library Screen

4.2.2.7 Pre-Exercise Safety Check Screen

Platform: Mother Mobile App Feature: MF-02 Mother Care Journey Release: P0 Related UCs: UC-29

Purpose: Collects required safety answers before an exercise session and blocks entry when configured warning answers are present.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Mother Mobile App. Primary role: Mother. Screen name: "Pre-Exercise Safety Check". Purpose: Collects required safety answers before an exercise session and blocks entry when configured warning answers are present. Related Release 1 use cases: UC-29. Use a single portrait mobile frame, 390 x 844 px, with a 24 px safe margin. Use a simple top app bar and bottom navigation only for signed-in main-application screens. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep

labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Selected Exercise Summary: Shows exercise title, stage suitability and duration.; Safety Questions: Collects required yes or no answers and optional explanation.; Warning Result Area: Explains when the session should not continue and recommends safe next action.; Continue Session Button: Proceeds only when no blocking safety condition is selected.; Cancel Action: Returns to the library without starting a session.. Interaction and safety constraints: Show only controls authorized for Mother. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Use conservative language and never make a fitness or medical clearance decision. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.2.17: Pre-Exercise Safety Check Screen (to be generated in Stitch)

Field name	Description
Selected Exercise Summary	Shows exercise title, stage suitability and duration.
Safety Questions	Collects required yes or no answers and optional explanation.
Warning Result Area	Explains when the session should not continue and recommends safe next action.
Continue Session Button	Proceeds only when no blocking safety condition is selected.
Cancel Action	Returns to the library without starting a session.

Table 4.2.17: Pre-Exercise Safety Check Screen

4.2.2.8 Active Exercise Session and Optional Posture Feedback Screen

Platform: Mother Mobile App Feature: MF-02 Mother Care Journey Release: P0 Related UCs: UC-30

Purpose: Runs an approved exercise with start, pause, resume and complete controls and optional camera-based posture feedback.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Mother Mobile App. Primary role: Mother. Screen name: "Active Exercise Session and Optional Posture Feedback". Purpose: Runs an approved exercise with start, pause, resume and complete controls and optional

camera-based posture feedback. Related Release 1 use cases: UC-30. Use a single portrait mobile frame, 390 x 844 px, with a 24 px safe margin. Use a simple top app bar and bottom navigation only for signed-in main-application screens. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Exercise Progress: Shows current step, elapsed time and remaining duration.; Start, Pause and Complete Controls: Controls the session lifecycle.; Optional Camera Consent: Lets the user enable posture feedback only after clear camera consent.; Posture Feedback Panel: Shows aggregate, non-medical rule or ML feedback without retaining raw video.; Stop for Safety Action: Lets the user end the session immediately.. Interaction and safety constraints: Show only controls authorized for Mother. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. State that raw images or video are not retained and feedback does not replace professional supervision. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.2.18: Active Exercise Session and Optional Posture Feedback Screen (to be generated in Stitch)

Field name	Description
Exercise Progress	Shows current step, elapsed time and remaining duration.
Start, Pause and Complete Controls	Controls the session lifecycle.
Optional Camera Consent	Lets the user enable posture feedback only after clear camera consent.
Posture Feedback Panel	Shows aggregate, non-medical rule or ML feedback without retaining raw video.
Stop for Safety Action	Lets the user end the session immediately.

Table 4.2.18: Active Exercise Session and Optional Posture Feedback Screen

4.2.2.9 Exercise History and Session Result Screen

Platform: Mother Mobile App Feature: MF-02 Mother Care Journey Release: P0 Related UCs: UC-31

Purpose: Shows completed exercise sessions, duration, completion state and aggregate posture feedback.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Mother Mobile App. Primary role: Mother. Screen name: "Exercise History and Session Result". Purpose: Shows completed exercise sessions, duration, completion state and aggregate posture feedback. Related Release 1 use cases: UC-31. Use a single portrait mobile frame, 390 x 844 px, with a 24 px safe margin. Use a simple top app bar and bottom navigation only for signed-in main-application screens. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Session Result Summary: Shows completion status, duration and selected exercise.; Aggregate Feedback: Shows non-medical posture cues or unavailable state when camera feedback was not used.; Safety Note: States that the result is informational and not a fitness assessment.; History List: Shows prior sessions with date, duration and status.; Return to Library Action: Returns to reviewed exercise content.. Interaction and safety constraints: Show only controls authorized for Mother. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Do not show or retain raw camera images or video. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.2.19: Exercise History and Session Result Screen (to be generated in Stitch)

Field name	Description
Session Result Summary	Shows completion status, duration and selected exercise.
Aggregate Feedback	Shows non-medical posture cues or unavailable state when camera feedback was not used.
Safety Note	States that the result is informational and not a fitness assessment.

Field name	Description
History List	Shows prior sessions with date, duration and status.
Return to Library Action	Returns to reviewed exercise content.

Table 4.2.19: Exercise History and Session Result Screen

4.2.2.10 Baby Profiles Screen

Platform: Mother Mobile App Feature: MF-03 Baby Care Journey, Growth and Vaccination

Release: P0 Related UCs: UC-32, UC-33, UC-34

Purpose: Creates, updates, archives and switches the active baby profile.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Mother Mobile App. Primary role: Mother. Screen name: "Baby Profiles". Purpose: Creates, updates, archives and switches the active baby profile. Related Release 1 use cases: UC-32, UC-33, UC-34. Use a single portrait mobile frame, 390 x 844 px, with a 24 px safe margin. Use a simple top app bar and bottom navigation only for signed-in main-application screens. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Baby Profile Cards: Show nickname, birth context, active state and high-level care status.; Add Baby Action: Creates a profile with minimum identity and birth context.; Edit or Archive Action: Updates allowed data or archives a profile without destroying linked history.; Set Active Action: Changes the baby context for journals, growth, vaccination and reminders.; Privacy Notice: Explains profile data is private unless separately shared.. Interaction and safety constraints: Show only controls authorized for Mother. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Use nickname and fictional placeholder data only in the wireframe. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.2.20: Baby Profiles Screen (to be generated in Stitch)

Field name	Description
Baby Profile Cards	Show nickname, birth context, active state and high-level care status.
Add Baby Action	Creates a profile with minimum identity and birth context.
Edit or Archive Action	Updates allowed data or archives a profile without destroying linked history.
Set Active Action	Changes the baby context for journals, growth, vaccination and reminders.
Privacy Notice	Explains profile data is private unless separately shared.

Table 4.2.20: Baby Profiles Screen

4.2.2.11 Baby Care Overview Screen

Platform: Mother Mobile App Feature: MF-03 Baby Care Journey, Growth and Vaccination

Release: P0 Related UCs: UC-35

Purpose: Shows the active baby recent journals, growth, milestones, vaccination state and current prompts.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Mother Mobile App. Primary role: Mother. Screen name: "Baby Care Overview". Purpose: Shows the active baby recent journals, growth, milestones, vaccination state and current prompts. Related Release 1 use cases: UC-35. Use a single portrait mobile frame, 390 x 844 px, with a 24 px safe margin. Use a simple top app bar and bottom navigation only for signed-in main-application screens. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Active Baby Header: Shows selected baby context and switch-profile control.; Recent Daily Care Snapshot: Shows last feeding, sleep, diaper or observation entries.; Growth and Milestone Cards: Show high-level current summaries and shortcut actions.; Vaccination Status Card: Shows recorded and upcoming reference schedule state.; Quick Log Actions: Opens the daily log, growth, milestone or vaccination form.. Interaction and safety constraints: Show only controls authorized for Mother. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-

dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Keep all health information user-entered and clearly source-labelled. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.2.21: Baby Care Overview Screen (to be generated in Stitch)

Field name	Description
Active Baby Header	Shows selected baby context and switch-profile control.
Recent Daily Care Snapshot	Shows last feeding, sleep, diaper or observation entries.
Growth and Milestone Cards	Show high-level current summaries and shortcut actions.
Vaccination Status Card	Shows recorded and upcoming reference schedule state.
Quick Log Actions	Opens the daily log, growth, milestone or vaccination form.

Table 4.2.21: Baby Care Overview Screen

4.2.2.12 Baby Daily Log Screen

Platform: Mother Mobile App Feature: MF-03 Baby Care Journey, Growth and Vaccination

Release: P0 Related UCs: UC-36, UC-37

Purpose: Adds, updates or removes daily feeding, sleep, diaper and supported care observations for the active baby.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Mother Mobile App. Primary role: Mother. Screen name: "Baby Daily Log". Purpose: Adds, updates or removes daily feeding, sleep, diaper and supported care observations for the active baby. Related Release 1 use cases: UC-36, UC-37. Use a single portrait mobile frame, 390 x 844 px, with a 24 px safe margin. Use a simple top app bar and bottom navigation only for signed-in main-application screens. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Log Type Selector: Selects feeding, sleep, diaper, symptom note or another supported category.; Log Entry Form: Collects time, value, unit or duration, source and optional note.; Recent Entries: Shows same-day entries with edit and delete controls for owner-entered data.; Save Action: Creates or updates the selected entry.; Care Safety Hint: Clarifies that the

journal supports observation and is not diagnosis.. Interaction and safety constraints: Show only controls authorized for Mother. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Keep the quick-entry form usable in night mode without creating a separate use case or screen. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.2.22: Baby Daily Log Screen (to be generated in Stitch)

Field name	Description
Log Type Selector	Selects feeding, sleep, diaper, symptom note or another supported category.
Log Entry Form	Collects time, value, unit or duration, source and optional note.
Recent Entries	Shows same-day entries with edit and delete controls for owner-entered data.
Save Action	Creates or updates the selected entry.
Care Safety Hint	Clarifies that the journal supports observation and is not diagnosis.

Table 4.2.22: Baby Daily Log Screen

4.2.2.13 Baby Log Summary Screen

Platform: Mother Mobile App Feature: MF-03 Baby Care Journey, Growth and Vaccination

Release: P0 Related UCs: UC-38

Purpose: Shows recent 24-hour or 7-day patterns for feeding, sleep, diaper and logged observations.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Mother Mobile App. Primary role: Mother. Screen name: "Baby Log Summary". Purpose: Shows recent 24-hour or 7-day patterns for feeding, sleep, diaper and logged observations. Related Release 1 use cases: UC-38. Use a single portrait mobile frame, 390 x 844 px, with a 24 px safe margin. Use a simple top app bar and bottom navigation only for signed-in main-application screens. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6),

medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Time Range Selector: Switches between 24-hour and 7-day summaries.; Pattern Cards: Summarize feeding, sleep, diaper and selected observations.; Simple Trend Visuals: Show user-entered patterns without diagnosing illness.; Source and Data-Gap Labels: Indicate user-entered data and incomplete periods.; Back to Daily Log Action: Returns to entry management.. Interaction and safety constraints: Show only controls authorized for Mother. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Do not infer medical conditions from logged patterns. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.2.23: Baby Log Summary Screen (to be generated in Stitch)

Field name	Description
Time Range Selector	Switches between 24-hour and 7-day summaries.
Pattern Cards	Summarize feeding, sleep, diaper and selected observations.
Simple Trend Visuals	Show user-entered patterns without diagnosing illness.
Source and Data-Gap Labels	Indicate user-entered data and incomplete periods.
Back to Daily Log Action	Returns to entry management.

Table 4.2.23: Baby Log Summary Screen

4.2.2.14 Development Milestones Screen

Platform: Mother Mobile App Feature: MF-03 Baby Care Journey, Growth and Vaccination

Release: P0 Related UCs: UC-39, UC-40

Purpose: Records, updates and removes caregiver-observed development milestones.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Mother Mobile App. Primary role: Mother. Screen name: "Development Milestones". Purpose: Records, updates and removes caregiver-observed development milestones. Related Release 1 use cases: UC-39, UC-40.

Use a single portrait mobile frame, 390 x 844 px, with a 24 px safe margin. Use a simple top app bar and bottom navigation only for signed-in main-application screens. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Milestone Timeline: Shows recorded milestones by age or date.; Add Milestone Action: Records milestone type, observed date and optional note.; Edit or Delete Action: Corrects or removes an owner-entered milestone.; Reference Context: Shows stage-appropriate reference wording without diagnosing developmental disorders.; Professional Consultation Prompt: Suggests professional discussion when the user is concerned.. Interaction and safety constraints: Show only controls authorized for Mother. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Avoid pass or fail wording and do not display diagnostic labels. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.2.24: Development Milestones Screen (to be generated in Stitch)

Field name	Description
Milestone Timeline	Shows recorded milestones by age or date.
Add Milestone Action	Records milestone type, observed date and optional note.
Edit or Delete Action	Corrects or removes an owner-entered milestone.
Reference Context	Shows stage-appropriate reference wording without diagnosing developmental disorders.
Professional Consultation Prompt	Suggests professional discussion when the user is concerned.

Table 4.2.24: Development Milestones Screen

4.2.2.15 Growth Measurements and Trend Screen

Platform: Mother Mobile App Feature: MF-03 Baby Care Journey, Growth and Vaccination
Release: P0 Related UCs: UC-41, UC-42, UC-43

Purpose: Records, updates, removes and reviews growth measurements and charts for the active baby.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Mother Mobile App. Primary role: Mother. Screen name: "Growth Measurements and Trend". Purpose: Records, updates, removes and reviews growth measurements and charts for the active baby. Related Release 1 use cases: UC-41, UC-42, UC-43. Use a single portrait mobile frame, 390 x 844 px, with a 24 px safe margin. Use a simple top app bar and bottom navigation only for signed-in main-application screens. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Measurement Form: Collects date, type, value, unit and source.; Measurement History: Shows owner-entered measurements with edit and delete controls.; Growth Chart: Displays selected measurement trend with clear reference context.; Source Label: Distinguishes home-entered and other supported source types.; Professional Assessment Prompt: Encourages professional assessment for concerns without making conclusions.. Interaction and safety constraints: Show only controls authorized for Mother. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Embed measurement-type and time filters in this screen. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.2.25: Growth Measurements and Trend Screen (to be generated in Stitch)

Field name	Description
Measurement Form	Collects date, type, value, unit and source.
Measurement History	Shows owner-entered measurements with edit and delete controls.
Growth Chart	Displays selected measurement trend with clear reference context.
Source Label	Distinguishes home-entered and other supported source types.

Field name	Description
Professional Assessment Prompt	Encourages professional assessment for concerns without making conclusions.

Table 4.2.25: Growth Measurements and Trend Screen

4.2.2.16 Vaccination Records and Reference Schedule Screen

Platform: Mother Mobile App Feature: MF-03 Baby Care Journey, Growth and Vaccination

Release: P0 Related UCs: UC-44, UC-45

Purpose: Manages user-entered vaccination records and views the reference schedule with linked reminder state.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Mother Mobile App. Primary role: Mother. Screen name: "Vaccination Records and Reference Schedule". Purpose: Manages user-entered vaccination records and views the reference schedule with linked reminder state. Related Release 1 use cases: UC-44, UC-45. Use a single portrait mobile frame, 390 x 844 px, with a 24 px safe margin. Use a simple top app bar and bottom navigation only for signed-in main-application screens. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Reference Schedule Timeline: Shows vaccination names and reference timing for the active baby.; Recorded Status: Shows recorded, pending or unknown status based on user-entered records.; Add or Edit Record Action: Adds or updates a vaccination record with date, source and note.; Reminder Status: Shows linked reminder state and shortcut to create or update a reminder.; Official Record Notice: Clarifies that CareBridge does not replace official vaccination documentation.. Interaction and safety constraints: Show only controls authorized for Mother. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Do not recommend a dose, diagnose suitability or replace pediatric assessment. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.2.26: Vaccination Records and Reference Schedule Screen (to be generated in Stitch)

Field name	Description
Reference Schedule Timeline	Shows vaccination names and reference timing for the active baby.
Recorded Status	Shows recorded, pending or unknown status based on user-entered records.
Add or Edit Record Action	Adds or updates a vaccination record with date, source and note.
Reminder Status	Shows linked reminder state and shortcut to create or update a reminder.
Official Record Notice	Clarifies that CareBridge does not replace official vaccination documentation.

Table 4.2.26: Vaccination Records and Reference Schedule Screen

4.2.2.17 Community Feed Screen

Platform: Mother Mobile App Feature: MF-04 Community Q&A and Moderation Release: P0
Related UCs: UC-46, UC-54

Purpose: Browses moderated community questions and topic collections with embedded search, chips, sort and pagination.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Mother Mobile App. Primary role: User. Screen name: "Community Feed". Purpose: Browses moderated community questions and topic collections with embedded search, chips, sort and pagination. Related Release 1 use cases: UC-46, UC-54. Use a single portrait mobile frame, 390 x 844 px, with a 24 px safe margin. Use a simple top app bar and bottom navigation only for signed-in main-application screens. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Keyword Search: Searches within authorized community content in the feed itself.; Topic Chips and Sort: Filters and orders the feed without leaving the screen.; Question Cards: Show title, topic, anonymous or public label, answer status and moderation or expert labels.; Follow Topic Control: Follows or unfollows a selected topic.; Ask Question Action: Opens the community post composer.; Safety and Source Labels: Explain personal experience, verified expert answer and moderation status.. Interaction and safety constraints: Show only controls authorized for User. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope

and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Do not create a separate community-search or topic-directory screen. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.2.27: Community Feed Screen (to be generated in Stitch)

Field name	Description
Keyword Search	Searches within authorized community content in the feed itself.
Topic Chips and Sort	Filters and orders the feed without leaving the screen.
Question Cards	Show title, topic, anonymous or public label, answer status and moderation or expert labels.
Follow Topic Control	Follows or unfollows a selected topic.
Ask Question Action	Opens the community post composer.
Safety and Source Labels	Explain personal experience, verified expert answer and moderation status.

Table 4.2.27: Community Feed Screen

4.2.2.18 Community Question Detail Screen

Platform: Mother Mobile App Feature: MF-04 Community Q&A and Moderation Release: P0
Related UCs: UC-47, UC-52, UC-53

Purpose: Shows a permitted community question, answer thread, source labels and allowed reactions or bookmarks.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Mother Mobile App. Primary role: User. Screen name: "Community Question Detail". Purpose: Shows a permitted community question, answer thread, source labels and allowed reactions or bookmarks. Related Release 1 use cases: UC-47, UC-52, UC-53. Use a single portrait mobile frame, 390 x 844 px, with a 24 px safe margin. Use a simple top app bar and bottom navigation only for signed-in main-application screens. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow,

photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Question Header: Shows topic, author display mode, time, moderation state and safety labels.; Answer Thread: Shows personal-experience or verified-expert answers with clear source labels.; Helpful Reaction: Adds or removes a permitted helpfulness reaction on an answer.; Bookmark Action: Saves or removes the question from the user private list.; Report Action: Starts the unsafe content or account report flow.; Answer Action: Opens the allowed answer composer.. Interaction and safety constraints: Show only controls authorized for User. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Respect moderation locks and do not expose private identities behind anonymous display. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.2.28: Community Question Detail Screen (to be generated in Stitch)

Field name	Description
Question Header	Shows topic, author display mode, time, moderation state and safety labels.
Answer Thread	Shows personal-experience or verified-expert answers with clear source labels.
Helpful Reaction	Adds or removes a permitted helpfulness reaction on an answer.
Bookmark Action	Saves or removes the question from the user private list.
Report Action	Starts the unsafe content or account report flow.
Answer Action	Opens the allowed answer composer.

Table 4.2.28: Community Question Detail Screen

4.2.2.19 Community Post Composer Screen

Platform: Mother Mobile App Feature: MF-04 Community Q&A and Moderation Release: P0
Related UCs: UC-48, UC-49

Purpose: Creates, edits or deletes the user own topic-based community question with optional anonymous display.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Mother Mobile App. Primary role: User. Screen name: "Community Post Composer". Purpose: Creates, edits or deletes the user own topic-based community question with optional anonymous display. Related Release 1 use cases: UC-48, UC-49. Use a single portrait mobile frame, 390 x 844 px, with a 24 px safe margin. Use a simple top app bar and bottom navigation only for signed-in main-application screens. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Topic Selector: Assigns the post to an enabled community topic.; Question Title and Body: Collects the question content and optional non-sensitive context.; Anonymous Display Toggle: Controls public display while retaining internal accountability.; Safety Notice: Warns against seeking diagnosis, prescription or emergency response from the community.; Submit, Save Draft or Delete Controls: Creates, updates or removes the user own post when permitted.. Interaction and safety constraints: Show only controls authorized for User. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Do not allow posting of sensitive identifiers, protected files or unsafe medical claims. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.2.29: Community Post Composer Screen (to be generated in Stitch)

Field name	Description
Topic Selector	Assigns the post to an enabled community topic.
Question Title and Body	Collects the question content and optional non-sensitive context.
Anonymous Display Toggle	Controls public display while retaining internal accountability.
Safety Notice	Warns against seeking diagnosis, prescription or emergency response from the community.
Submit, Save Draft or Delete Controls	Creates, updates or removes the user own post when permitted.

Table 4.2.29: Community Post Composer Screen

4.2.2.20 Community Answer Composer Screen

Platform: Mother Mobile App Feature: MF-04 Community Q&A and Moderation / MF-05 Expert Contribution Release: P0 Related UCs: UC-50, UC-51, UC-68

Purpose: Creates, edits or deletes a community answer and applies the correct user or verified-expert label.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Mother Mobile App. Primary role: User / Verified Expert. Screen name: "Community Answer Composer". Purpose: Creates, edits or deletes a community answer and applies the correct user or verified-expert label. Related Release 1 use cases: UC-50, UC-51, UC-68. Use a single portrait mobile frame, 390 x 844 px, with a 24 px safe margin. Use a simple top app bar and bottom navigation only for signed-in main-application screens. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Question Context: Shows the target question and topic safely.; Answer Editor: Collects answer content with visible guidance against diagnosis or prescription.; Source Label: Shows personal-experience label or verified-expert identity and safety label.; Submit or Update Action: Posts or revises the answer if it is not locked.; Delete Action: Removes the owner answer where no evidence hold applies.. Interaction and safety constraints: Show only controls authorized for User / Verified Expert. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Verified experts may provide professional guidance within their approved support scope but must not present a diagnosis or prescription. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.2.30: Community Answer Composer Screen (to be generated in Stitch)

Field name	Description
Question Context	Shows the target question and topic safely.
Answer Editor	Collects answer content with visible guidance against diagnosis or prescription.

Field name	Description
Source Label	Shows personal-experience label or verified-expert identity and safety label.
Submit or Update Action	Posts or revises the answer if it is not locked.
Delete Action	Removes the owner answer where no evidence hold applies.

Table 4.2.30: Community Answer Composer Screen

4.2.2.21 Report Unsafe Content or Account Screen

Platform: Mother Mobile App Feature: MF-04 Community Q&A and Moderation Release: PO
Related UCs: UC-55

Purpose: Submits a report about misinformation, dangerous advice, hidden advertising, harassment, privacy breach or suspicious activity.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Mother Mobile App. Primary role: User. Screen name: "Report Unsafe Content or Account". Purpose: Submits a report about misinformation, dangerous advice, hidden advertising, harassment, privacy breach or suspicious activity. Related Release 1 use cases: UC-55. Use a single portrait mobile frame, 390 x 844 px, with a 24 px safe margin. Use a simple top app bar and bottom navigation only for signed-in main-application screens. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Reported Item Summary: Shows the selected post, answer or account context without expanding unnecessary data.; Report Reason Selector: Selects a standardized policy or safety reason.; Optional Evidence Note: Allows a concise explanation without adding sensitive health data.; Urgent Safety Note: Directs immediate danger to emergency services rather than moderation waiting.; Submit Report Button: Submits the report and shows a tracking acknowledgement.. Interaction and safety constraints: Show only controls authorized for User. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Keep the reporter identity confidential from the reported user. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.2.31: Report Unsafe Content or Account Screen (to be generated in Stitch)

Field name	Description
Reported Item Summary	Shows the selected post, answer or account context without expanding unnecessary data.
Report Reason Selector	Selects a standardized policy or safety reason.
Optional Evidence Note	Allows a concise explanation without adding sensitive health data.
Urgent Safety Note	Directs immediate danger to emergency services rather than moderation waiting.
Submit Report Button	Submits the report and shows a tracking acknowledgement.

Table 4.2.31: Report Unsafe Content or Account Screen

4.2.2.22 Verified Expert Directory Screen

Platform: Mother Mobile App Feature: MF-05 Verified Expert Network and Contribution Release: P0 Related UCs: UC-65

Purpose: Browses verified experts with embedded specialty, availability, badge and consented service-area filters.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Mother Mobile App. Primary role: User. Screen name: "Verified Expert Directory". Purpose: Browses verified experts with embedded specialty, availability, badge and consented service-area filters. Related Release 1 use cases: UC-65. Use a single portrait mobile frame, 390 x 844 px, with a 24 px safe margin. Use a simple top app bar and bottom navigation only for signed-in main-application screens. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Keyword and Filter Controls: Embed search, specialty, availability, badge and consented area filters in the directory.; Expert Cards: Show verified name, specialty, approved scope, availability and contribution badges.; Location Consent Hint: Explains that location-aware results require explicit location permission.; Profile Action: Opens the public verified expert profile.; Empty State: Explains when no verified expert matches the active criteria.. Interaction and safety constraints: Show only controls authorized for User. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency,

map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Do not show V2 booking, fee, payment or realtime consultation controls. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.2.32: Verified Expert Directory Screen (to be generated in Stitch)

Field name	Description
Keyword and Filter Controls	Embed search, specialty, availability, badge and consented area filters in the directory.
Expert Cards	Show verified name, specialty, approved scope, availability and contribution badges.
Location Consent Hint	Explains that location-aware results require explicit location permission.
Profile Action	Opens the public verified expert profile.
Empty State	Explains when no verified expert matches the active criteria.

Table 4.2.32: Verified Expert Directory Screen

4.2.2.23 Verified Expert Profile Screen

Platform: Mother Mobile App Feature: MF-05 Verified Expert Network and Contribution Release: P0 Related UCs: UC-66

Purpose: Shows a public verified expert profile, specialties, supported scope, availability and contribution reputation.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Mother Mobile App. Primary role: User. Screen name: "Verified Expert Profile". Purpose: Shows a public verified expert profile, specialties, supported scope, availability and contribution reputation. Related Release 1 use cases: UC-66. Use a single portrait mobile frame, 390 x 844 px, with a 24 px safe margin. Use a simple top app bar and bottom navigation only for signed-in main-application screens. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Professional Header: Shows verified identity, professional title and active trust

status.; Specialty and Service Scope: Shows approved specialties and supported non-emergency scope.; Availability Summary: Shows currently published availability or readiness state.; Contribution Reputation: Shows contribution points and badges as participation recognition, not medical competence.; Safety Boundary: States that profile information is not emergency dispatch or direct clinical treatment.. Interaction and safety constraints: Show only controls authorized for User. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Do not include paid booking, price, payment or V2 session actions. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.2.33: Verified Expert Profile Screen (to be generated in Stitch)

Field name	Description
Professional Header	Shows verified identity, professional title and active trust status.
Specialty and Service Scope	Shows approved specialties and supported non-emergency scope.
Availability Summary	Shows currently published availability or readiness state.
Contribution Reputation	Shows contribution points and badges as participation recognition, not medical competence.
Safety Boundary	States that profile information is not emergency dispatch or direct clinical treatment.

Table 4.2.33: Verified Expert Profile Screen

4.2.2.24 AI Symptom Intake Screen

Platform: Mother Mobile App Feature: MF-06 AI Nurse Assistant and Risk Triage Release: P0
Related UCs: UC-72

Purpose: Collects structured symptom information and permitted context for non-diagnostic risk orientation.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Mother Mobile App. Primary role: Mother. Screen name: "AI Symptom Intake". Purpose: Collects structured symptom information and permitted context for non-diagnostic risk orientation. Related Release 1 use cases: UC-72. Use a single portrait mobile frame, 390 x 844 px, with a 24 px safe margin. Use a

simple top app bar and bottom navigation only for signed-in main-application screens. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Symptom Category Selector: Selects the main concern using structured categories.; Guided Questions: Collect timing, severity, change and limited relevant context.; Context Source Labels: Distinguish user-entered information from device or record context when present.; Emergency Notice: Warns that immediate danger requires emergency care and not waiting for AI.; Continue to Result Button: Runs safe intake processing and opens the triage result.. Interaction and safety constraints: Show only controls authorized for Mother. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Do not show diagnosis, dosage, prescription or deterministic clinical decision claims. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.2.34: AI Symptom Intake Screen (to be generated in Stitch)

Field name	Description
Symptom Category Selector	Selects the main concern using structured categories.
Guided Questions	Collect timing, severity, change and limited relevant context.
Context Source Labels	Distinguish user-entered information from device or record context when present.
Emergency Notice	Warns that immediate danger requires emergency care and not waiting for AI.
Continue to Result Button	Runs safe intake processing and opens the triage result.

Table 4.2.34: AI Symptom Intake Screen

4.2.2.25 Risk Triage Result Screen

Platform: Mother Mobile App Feature: MF-06 AI Nurse Assistant and Risk Triage Release: P0
Related UCs: UC-73, UC-74

Purpose: Shows green, yellow or red non-diagnostic risk orientation and a safe next action.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Mother Mobile App. Primary role: Mother. Screen name: "Risk Triage Result". Purpose: Shows green, yellow or red non-diagnostic risk orientation and a safe next action. Related Release 1 use cases: UC-73, UC-74. Use a single portrait mobile frame, 390 x 844 px, with a 24 px safe margin. Use a simple top app bar and bottom navigation only for signed-in main-application screens. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Risk Orientation Banner: Shows Green, Yellow or Red as a non-diagnostic orientation.; Safe Next Steps: Shows conservative actions appropriate to the result and uncertainty.; Why This Result: Shows a concise summary of the intake factors without medical inference beyond policy.; Emergency Support Action: Opens Emergency Map for red risk or user-selected urgent support.; Disclaimer and Audit Hint: States that AI does not diagnose, prescribe or replace emergency services.. Interaction and safety constraints: Show only controls authorized for Mother. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Use high-contrast but grayscale layout; the label text must be understandable without color. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.2.35: Risk Triage Result Screen (to be generated in Stitch)

Field name	Description
Risk Orientation Banner	Shows Green, Yellow or Red as a non-diagnostic orientation.
Safe Next Steps	Shows conservative actions appropriate to the result and uncertainty.
Why This Result	Shows a concise summary of the intake factors without medical inference beyond policy.
Emergency Support Action	Opens Emergency Map for red risk or user-selected urgent support.

Field name	Description
Disclaimer and Audit Hint	States that AI does not diagnose, prescribe or replace emergency services.

Table 4.2.35: Risk Triage Result Screen

4.2.2.26 Emergency Map and Nearby Facilities Screen

Platform: Mother Mobile App Feature: MF-07 Emergency Map and Nearby Care Support Release: P0 Related UCs: UC-77, UC-78, UC-79

Purpose: Opens nearby facilities, routes, ETA, quick call and navigation after explicit location permission.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Mother Mobile App. Primary role: Mother. Screen name: "Emergency Map and Nearby Facilities". Purpose: Opens nearby facilities, routes, ETA, quick call and navigation after explicit location permission. Related Release 1 use cases: UC-77, UC-78, UC-79. Use a single portrait mobile frame, 390 x 844 px, with a 24 px safe margin. Use a simple top app bar and bottom navigation only for signed-in main-application screens. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Location Permission State: Explains current permission and allows area selection when location is unavailable.; Map Canvas and Facility Markers: Shows eligible nearby care facility markers in a basic grayscale map representation.; Facility List: Shows name, distance, type and availability information where supported.; Route and ETA Panel: Shows route, estimated travel time and selected facility context.; Quick Call and Navigate Actions: Uses supported device capabilities to call or open navigation.; Non-Dispatch Disclaimer: States that CareBridge does not dispatch ambulances or guarantee service availability.. Interaction and safety constraints: Show only controls authorized for Mother. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Keep map search and distance controls inside this screen and show only consented location data. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.2.36: Emergency Map and Nearby Facilities Screen (to be generated in Stitch)

Field name	Description
Location Permission State	Explains current permission and allows area selection when location is unavailable.
Map Canvas and Facility Markers	Shows eligible nearby care facility markers in a basic grayscale map representation.
Facility List	Shows name, distance, type and availability information where supported.
Route and ETA Panel	Shows route, estimated travel time and selected facility context.
Quick Call and Navigate Actions	Uses supported device capabilities to call or open navigation.
Non-Dispatch Disclaimer	States that CareBridge does not dispatch ambulances or guarantee service availability.

Table 4.2.36: Emergency Map and Nearby Facilities Screen

4.2.2.27 Family Emergency Alert and Time-Limited Location Screen

Platform: Mother Mobile App Feature: MF-07 Emergency Map and Nearby Care Support Release: P0 Related UCs: UC-80

Purpose: Creates a time-limited emergency alert for selected family recipients with minimum necessary location and context.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Mother Mobile App. Primary role: Mother. Screen name: "Family Emergency Alert and Time-Limited Location". Purpose: Creates a time-limited emergency alert for selected family recipients with minimum necessary location and context. Related Release 1 use cases: UC-80. Use a single portrait mobile frame, 390 x 844 px, with a 24 px safe margin. Use a simple top app bar and bottom navigation only for signed-in main-application screens. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Selected Recipients: Shows selected eligible family contacts or care-group recipients.; Location Scope and Expiry: Defines accuracy, sharing duration and expiration time.; Alert Context: Allows a short user-entered help message without unnecessary health details.; Send Alert Button: Sends the family alert after explicit consent.; Stop Sharing Action: Stops active location sharing before expiry when the owner chooses.. Interaction and safety constraints: Show only controls authorized for Mother. Maintain privacy by default.

Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Show the sharing scope and expiry prominently. Do not imply emergency dispatch. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.2.37: Family Emergency Alert and Time-Limited Location Screen (to be generated in Stitch)

Field name	Description
Selected Recipients	Shows selected eligible family contacts or care-group recipients.
Location Scope and Expiry	Defines accuracy, sharing duration and expiration time.
Alert Context	Allows a short user-entered help message without unnecessary health details.
Send Alert Button	Sends the family alert after explicit consent.
Stop Sharing Action	Stops active location sharing before expiry when the owner chooses.

Table 4.2.37: Family Emergency Alert and Time-Limited Location Screen

4.2.2.28 Nearby Support Request Screen

Platform: Mother Mobile App Feature: MF-07 Emergency Map and Nearby Care Support Release: P0 Related UCs: UC-81

Purpose: Creates, views or cancels a consented nearby-support request visible only to eligible verified experts.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Mother Mobile App. Primary role: Mother. Screen name: "Nearby Support Request". Purpose: Creates, views or cancels a consented nearby-support request visible only to eligible verified experts. Related Release 1 use cases: UC-81. Use a single portrait mobile frame, 390 x 844 px, with a 24 px safe margin. Use a simple top app bar and bottom navigation only for signed-in main-application screens. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons,

illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Request Purpose and Urgency: Collects minimal request category and a non-diagnostic urgency label.; Consent Summary: Shows what location and context eligible experts may see and for how long.; Request Status: Shows open, accepted, cancelled or expired state.; Cancel Request Action: Cancels an active request and stops new expert responses.; Emergency Disclaimer: States that expert arrival is not guaranteed and emergency services remain separate.. Interaction and safety constraints: Show only controls authorized for Mother. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Do not expose full health records or exact location beyond the configured consent scope. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.2.38: Nearby Support Request Screen (to be generated in Stitch)

Field name	Description
Request Purpose and Urgency	Collects minimal request category and a non-diagnostic urgency label.
Consent Summary	Shows what location and context eligible experts may see and for how long.
Request Status	Shows open, accepted, cancelled or expired state.
Cancel Request Action	Cancels an active request and stops new expert responses.
Emergency Disclaimer	States that expert arrival is not guaranteed and emergency services remain separate.

Table 4.2.38: Nearby Support Request Screen

4.2.2.29 Health Record Timeline and Detail Screen

Platform: Mother Mobile App Feature: MF-08 Personal Health Records and Source Labeling
Release: P1 Related UCs: UC-86

Purpose: Shows authorized maternal and baby records in a timeline with embedded category, date and source filters.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Mother Mobile App. Primary role: Mother. Screen name: "Health Record Timeline and Detail". Purpose: Shows authorized maternal and baby records in a timeline with embedded category, date and source filters. Related Release 1 use cases: UC-86. Use a single portrait mobile frame, 390 x 844 px, with a 24 px safe margin. Use a simple top app bar and bottom navigation only for signed-in main-application screens. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Timeline Filters: Embed date, owner, category and source filters inside the timeline.; Record Cards: Show record title, date, category, source label, attachment state and sharing indicator.; Record Detail Panel: Shows permitted metadata, note and protected attachment preview area.; Add Record Action: Opens the add or edit record screen.; Sharing Badge: Shows active consent scope and expiry when a record is shared.. Interaction and safety constraints: Show only controls authorized for Mother. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Do not create standalone file manager or file viewer screens; attachments remain inside the protected record context. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.2.39: Health Record Timeline and Detail Screen (to be generated in Stitch)

Field name	Description
Timeline Filters	Embed date, owner, category and source filters inside the timeline.
Record Cards	Show record title, date, category, source label, attachment state and sharing indicator.
Record Detail Panel	Shows permitted metadata, note and protected attachment preview area.
Add Record Action	Opens the add or edit record screen.
Sharing Badge	Shows active consent scope and expiry when a record is shared.

Table 4.2.39: Health Record Timeline and Detail Screen

4.2.2.30 Add or Edit Health Record and Attachment Screen

Platform: Mother Mobile App Feature: MF-08 Personal Health Records and Source Labeling

Release: P1 Related UCs: UC-83, UC-84, UC-85

Purpose: Adds a health record and protected attachment, updates permitted metadata, or requests archive or deletion.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Mother Mobile App. Primary role: Mother. Screen name: "Add or Edit Health Record and Attachment". Purpose: Adds a health record and protected attachment, updates permitted metadata, or requests archive or deletion. Related Release 1 use cases: UC-83, UC-84, UC-85. Use a single portrait mobile frame, 390 x 844 px, with a 24 px safe margin. Use a simple top app bar and bottom navigation only for signed-in main-application screens. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Record Owner Selector: Selects mother or active baby within owner permissions.; Title, Category and Date Fields: Collect record metadata needed for safe timeline organization.; Source Label Selector: Labels the origin as user-entered, device-sourced, expert-provided or other supported source.; Protected Attachment Picker: Adds a supported file or image with upload state and size validation.; Save, Archive or Delete Controls: Applies permitted lifecycle action with in-screen confirmation for destructive actions.. Interaction and safety constraints: Show only controls authorized for Mother. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Do not allow attachment access outside record authorization and consent rules. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.2.40: Add or Edit Health Record and Attachment Screen (to be generated in Stitch)

Field name	Description
Record Owner Selector	Selects mother or active baby within owner permissions.
Title, Category and Date Fields	Collect record metadata needed for safe timeline organization.

Field name	Description
Source Label Selector	Labels the origin as user-entered, device-sourced, expert-provided or other supported source.
Protected Attachment Picker	Adds a supported file or image with upload state and size validation.
Save, Archive or Delete Controls	Applies permitted lifecycle action with in-screen confirmation for destructive actions.

Table 4.2.40: Add or Edit Health Record and Attachment Screen

4.2.2.31 Health Summary and Consent Share Screen

Platform: Mother Mobile App Feature: MF-08 Personal Health Records and Source Labeling

Release: P1 Related UCs: UC-87, UC-88

Purpose: Creates a selected-period health summary and shares it or selected records only through an active consent grant.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Mother Mobile App. Primary role: Mother. Screen name: "Health Summary and Consent Share". Purpose: Creates a selected-period health summary and shares it or selected records only through an active consent grant. Related Release 1 use cases: UC-87, UC-88. Use a single portrait mobile frame, 390 x 844 px, with a 24 px safe margin. Use a simple top app bar and bottom navigation only for signed-in main-application screens. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Period and Record Selector: Selects maternal or baby records included in the summary.; Generated Summary Preview: Shows a source-labelled, non-diagnostic summary draft.; Existing Consent Selector: Chooses only an active recipient, purpose, scope and expiry grant.; Share Scope Review: Shows exactly which summary or records will be shared.; Generate or Share Action: Creates the summary or shares selected items after validation.. Interaction and safety constraints: Show only controls authorized for Mother. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Do not create a new implicit consent from this screen; guide users to Privacy, Consent and Data Sharing when no valid

grant exists. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.2.41: Health Summary and Consent Share Screen (to be generated in Stitch)

Field name	Description
Period and Record Selector	Selects maternal or baby records included in the summary.
Generated Summary Preview	Shows a source-labelled, non-diagnostic summary draft.
Existing Consent Selector	Chooses only an active recipient, purpose, scope and expiry grant.
Share Scope Review	Shows exactly which summary or records will be shared.
Generate or Share Action	Creates the summary or shares selected items after validation.

Table 4.2.41: Health Summary and Consent Share Screen

4.2.2.32 Today Tasks and Care Plan Screen

Platform: Mother Mobile App Feature: MF-09 Reminders, Tasks and Care Plan Release: P1

Related UCs: UC-93

Purpose: Shows prioritized due reminders, checklist items and allowed family tasks for today.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Mother Mobile App. Primary role: Mother. Screen name: "Today Tasks and Care Plan". Purpose: Shows prioritized due reminders, checklist items and allowed family tasks for today. Related Release 1 use cases: UC-93. Use a single portrait mobile frame, 390 x 844 px, with a 24 px safe margin. Use a simple top app bar and bottom navigation only for signed-in main-application screens. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Date Header: Shows the selected day and high-level completion progress.; Priority Sections: Groups reminders, checklist items and family tasks by timing or urgency.; Task State Control: Marks a permitted item complete, skipped or needs attention.; Add Reminder Shortcut: Opens the reminder form.; Care Plan Note: States that the plan is organizational and not prescription or treatment planning.. Interaction and safety constraints: Show only controls authorized for Mother. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent

clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Keep sorting and filtering embedded inside the today view. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.2.42: Today Tasks and Care Plan Screen (to be generated in Stitch)

Field name	Description
Date Header	Shows the selected day and high-level completion progress.
Priority Sections	Groups reminders, checklist items and family tasks by timing or urgency.
Task State Control	Marks a permitted item complete, skipped or needs attention.
Add Reminder Shortcut	Opens the reminder form.
Care Plan Note	States that the plan is organizational and not prescription or treatment planning.

Table 4.2.42: Today Tasks and Care Plan Screen

4.2.2.33 Care Reminder Form Screen

Platform: Mother Mobile App Feature: MF-09 Reminders, Tasks and Care Plan Release: P1
Related UCs: UC-89, UC-90, UC-91, UC-92

Purpose: Creates and manages appointment, medicine or vitamin and vaccination reminders.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Mother Mobile App. Primary role: Mother. Screen name: "Care Reminder Form". Purpose: Creates and manages appointment, medicine or vitamin and vaccination reminders. Related Release 1 use cases: UC-89, UC-90, UC-91, UC-92. Use a single portrait mobile frame, 390 x 844 px, with a 24 px safe margin. Use a simple top app bar and bottom navigation only for signed-in main-application screens. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Reminder Type Selector: Selects appointment, medicine or vitamin, vaccination or other supported care reminder.; Title and Schedule Fields: Collect title, date, time, repeat and optional note.; Linked Context: Links to baby vaccination or user-entered appointment context

when selected.; Lifecycle Actions: Save, snooze, complete, skip or delete the reminder where permitted.; No-Prescription Notice: Clarifies that the system does not create medication instruction or dosage.. Interaction and safety constraints: Show only controls authorized for Mother. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Use in-screen confirmation for deletion; do not create a standalone confirmation screen. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.2.43: Care Reminder Form Screen (to be generated in Stitch)

Field name	Description
Reminder Type Selector	Selects appointment, medicine or vitamin, vaccination or other supported care reminder.
Title and Schedule Fields	Collect title, date, time, repeat and optional note.
Linked Context	Links to baby vaccination or user-entered appointment context when selected.
Lifecycle Actions	Save, snooze, complete, skip or delete the reminder where permitted.
No-Prescription Notice	Clarifies that the system does not create medication instruction or dosage.

Table 4.2.43: Care Reminder Form Screen

4.2.2.34 Care Group Management and Family Permissions Screen

Platform: Mother Mobile App Feature: MF-10 Family Sync and Cooperative Care Release: P1
Related UCs: UC-94, UC-95, UC-97, UC-98, UC-99

Purpose: Creates care groups, invites or removes members, defines family permission scope and manages family care tasks.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Mother Mobile App. Primary role: Mother. Screen name: "Care Group Management and Family Permissions". Purpose: Creates care groups, invites or removes members, defines family permission scope and manages

family care tasks. Related Release 1 use cases: UC-94, UC-95, UC-97, UC-98, UC-99. Use a single portrait mobile frame, 390 x 844 px, with a 24 px safe margin. Use a simple top app bar and bottom navigation only for signed-in main-application screens. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Care Group Header: Shows group owner, linked care context and membership count.; Member List: Shows member role, invitation state and permitted status actions.; Invite or Revoke Control: Sends or revokes a family member invitation.; Permission Scope Editor: Defines calendar, tasks, alerts, logs, summaries or record scope per member.; Family Task Panel: Creates, edits or cancels a task assigned to an eligible member.. Interaction and safety constraints: Show only controls authorized for Mother. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Display permission scope and expiry clearly; do not expose data that is not covered by active scope. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.2.44: Care Group Management and Family Permissions Screen (to be generated in Stitch)

Field name	Description
Care Group Header	Shows group owner, linked care context and membership count.
Member List	Shows member role, invitation state and permitted status actions.
Invite or Revoke Control	Sends or revokes a family member invitation.
Permission Scope Editor	Defines calendar, tasks, alerts, logs, summaries or record scope per member.
Family Task Panel	Creates, edits or cancels a task assigned to an eligible member.

Table 4.2.44: Care Group Management and Family Permissions Screen

4.2.2.35 Verified Content Hub Screen

Platform: Mother Mobile App Feature: MF-11 Verified Content and Checklist Hub Release: P1
Related UCs: UC-102

Purpose: Browses approved stage-based articles, FAQs and checklists with embedded keyword, stage and topic controls.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Mother Mobile App. Primary role: User. Screen name: "Verified Content Hub". Purpose: Browses approved stage-based articles, FAQs and checklists with embedded keyword, stage and topic controls. Related Release 1 use cases: UC-102. Use a single portrait mobile frame, 390 x 844 px, with a 24 px safe margin. Use a simple top app bar and bottom navigation only for signed-in main-application screens. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Keyword and Filter Controls: Embed keyword, stage, topic and content-type controls in the hub.; Content Cards: Show title, type, stage, source label, review state and updated date.; Featured Checklist Section: Shows concise action checklists relevant to the selected stage.; Open Content Action: Opens the selected verified content detail.; Safety and Source Labels: Explains reviewed source, version and medical-safety context.. Interaction and safety constraints: Show only controls authorized for User. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Do not create a standalone verified-content search screen; search belongs in this hub. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.2.45: Verified Content Hub Screen (to be generated in Stitch)

Field name	Description
Keyword and Filter Controls	Embed keyword, stage, topic and content-type controls in the hub.
Content Cards	Show title, type, stage, source label, review state and updated date.
Featured Checklist Section	Shows concise action checklists relevant to the selected stage.
Open Content Action	Opens the selected verified content detail.

Field name	Description
Safety and Source Labels	Explains reviewed source, version and medical-safety context.

Table 4.2.45: Verified Content Hub Screen

4.2.2.36 Verified Content Detail Screen

Platform: Mother Mobile App Feature: MF-11 Verified Content and Checklist Hub Release: P1
Related UCs: UC-103

Purpose: Shows a selected approved article, FAQ or checklist with source, review or version state, updated date and safety notes.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Mother Mobile App. Primary role: User. Screen name: "Verified Content Detail". Purpose: Shows a selected approved article, FAQ or checklist with source, review or version state, updated date and safety notes. Related Release 1 use cases: UC-103. Use a single portrait mobile frame, 390 x 844 px, with a 24 px safe margin. Use a simple top app bar and bottom navigation only for signed-in main-application screens. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Content Header: Shows title, type, stage, source label, reviewed date and version state.; Content Body: Shows the approved article, FAQ answer or action checklist.; Source and Safety Note: Explains provenance, review state and relevant caution.; Checklist Completion Controls: Allows local checklist progress where the content type supports it.; Back to Hub Action: Returns to the verified content browse screen.. Interaction and safety constraints: Show only controls authorized for User. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Do not allow direct editing of published content from the user view. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.2.46: Verified Content Detail Screen (to be generated in Stitch)

Field name	Description
Content Header	Shows title, type, stage, source label, reviewed date and version state.
Content Body	Shows the approved article, FAQ answer or action checklist.
Source and Safety Note	Explains provenance, review state and relevant caution.
Checklist Completion Controls	Allows local checklist progress where the content type supports it.
Back to Hub Action	Returns to the verified content browse screen.

Table 4.2.46: Verified Content Detail Screen

4.2.2.37 Expense Planner and Summary Screen

Platform: Mother Mobile App Feature: MF-12 Expense and Preparation Planner Release: P1

Related UCs: UC-109, UC-110, UC-111

Purpose: Records, updates or removes preparation expenses and shows simple household summaries.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Mother Mobile App. Primary role: Mother. Screen name: "Expense Planner and Summary". Purpose: Records, updates or removes preparation expenses and shows simple household summaries. Related Release 1 use cases: UC-109, UC-110, UC-111. Use a single portrait mobile frame, 390 x 844 px, with a 24 px safe margin. Use a simple top app bar and bottom navigation only for signed-in main-application screens. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Expense Entry Form: Collects category, amount, date and optional experience note.; Expense List: Shows entries with edit and delete controls for owner-entered items.; Summary Cards: Show totals by category, month or care period.; Category Filters: Embed category and time filters inside the planner.; Scope Note: Clarifies this is not formal accounting, billing, insurance or financial advice.. Interaction and safety constraints: Show only controls authorized for Mother. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or

list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Use fictional monetary values only in the wireframe. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.2.47: Expense Planner and Summary Screen (to be generated in Stitch)

Field name	Description
Expense Entry Form	Collects category, amount, date and optional experience note.
Expense List	Shows entries with edit and delete controls for owner-entered items.
Summary Cards	Show totals by category, month or care period.
Category Filters	Embed category and time filters inside the planner.
Scope Note	Clarifies this is not formal accounting, billing, insurance or financial advice.

Table 4.2.47: Expense Planner and Summary Screen

4.2.2.38 Connected Devices and Synchronisation Screen

Platform: Mother Mobile App Feature: MF-13 Connected Device and Health Data Integration

Release: P1 Related UCs: UC-112, UC-113, UC-115

Purpose: Connects supported devices or platforms, imports selected observations and disconnects or deletes imported data.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Mother Mobile App. Primary role: Mother. Screen name: "Connected Devices and Synchronisation". Purpose: Connects supported devices or platforms, imports selected observations and disconnects or deletes imported data. Related Release 1 use cases: UC-112, UC-113, UC-115. Use a single portrait mobile frame, 390 x 844 px, with a 24 px safe margin. Use a simple top app bar and bottom navigation only for signed-in main-application screens. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Supported Connection Cards: Show supported device or platform options and current connection status.; Permission Scope Review: Shows indicators selected for import such as sleep, steps, heart rate or blood pressure where supported.; Sync Status: Shows last sync, data availability and error state.; Disconnect and Delete Action: Stops

future sync and requests deletion of imported data under policy.; Source and Quality Notice: Explains that device observations are not clinical measurements by default.. Interaction and safety constraints: Show only controls authorized for Mother. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Do not access a device or import data until the user explicitly accepts the scope. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.2.48: Connected Devices and Synchronisation Screen (to be generated in Stitch)

Field name	Description
Supported Connection Cards	Show supported device or platform options and current connection status.
Permission Scope Review	Shows indicators selected for import such as sleep, steps, heart rate or blood pressure where supported.
Sync Status	Shows last sync, data availability and error state.
Disconnect and Delete Action	Stops future sync and requests deletion of imported data under policy.
Source and Quality Notice	Explains that device observations are not clinical measurements by default.

Table 4.2.48: Connected Devices and Synchronisation Screen

4.2.2.39 Device Data Trend and Quality Screen

Platform: Mother Mobile App Feature: MF-13 Connected Device and Health Data Integration
Release: P1 Related UCs: UC-114

Purpose: Shows trends, source labels, data gaps and quality cautions for imported device observations.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Mother Mobile App. Primary role: Mother. Screen name: "Device Data Trend and Quality". Purpose: Shows trends, source labels, data gaps and quality cautions for imported device observations. Related Release 1 use cases: UC-114. Use a single portrait mobile frame, 390 x 844 px, with a 24 px safe margin. Use

a simple top app bar and bottom navigation only for signed-in main-application screens. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Indicator Selector: Switches the displayed imported observation type.; Trend Chart: Shows available data trend with time range control.; Source and Quality Labels: Show originating device or platform, last sync, gaps and warning state.; Context Note: Explains that the view supports observation, not diagnosis.; Connection Shortcut: Returns to connection and sync settings.. Interaction and safety constraints: Show only controls authorized for Mother. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Do not generate an abnormal medical diagnosis from device data. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.2.49: Device Data Trend and Quality Screen (to be generated in Stitch)

Field name	Description
Indicator Selector	Switches the displayed imported observation type.
Trend Chart	Shows available data trend with time range control.
Source and Quality Labels	Show originating device or platform, last sync, gaps and warning state.
Context Note	Explains that the view supports observation, not diagnosis.
Connection Shortcut	Returns to connection and sync settings.

Table 4.2.49: Device Data Trend and Quality Screen

4.2.2.40 Emergency Contacts and Smart Activity Monitoring Screen

Platform: Mother Mobile App Feature: MF-14 Smart Activity Monitoring and Safety Support

Release: P1 Related UCs: UC-116, UC-117, UC-118, UC-119

Purpose: Manages emergency contacts and configures, enables or disables optional phone IMU safety monitoring.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Mother Mobile App. Primary role: Mother. Screen name: "Emergency Contacts and Smart Activity Monitoring". Purpose: Manages emergency contacts and configures, enables or disables optional phone IMU safety monitoring. Related Release 1 use cases: UC-116, UC-117, UC-118, UC-119. Use a single portrait mobile frame, 390 x 844 px, with a 24 px safe margin. Use a simple top app bar and bottom navigation only for signed-in main-application screens. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Emergency Contact List: Shows alert recipients, order and verification state.; Add or Edit Contact Action: Adds, verifies, reprioritizes or removes a recipient.; Monitoring Configuration: Sets sensor consent, countdown duration, alert recipients and optional location-sharing behavior.; Monitoring Status Toggle: Enables or disables monitoring after permissions and configuration are valid.; Limitations Notice: States that monitoring is optional, device-dependent and not certified emergency detection.. Interaction and safety constraints: Show only controls authorized for Mother. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Do not promise continuous background detection or dispatch service. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.2.50: Emergency Contacts and Smart Activity Monitoring Screen (to be generated in Stitch)

Field name	Description
Emergency Contact List	Shows alert recipients, order and verification state.
Add or Edit Contact Action	Adds, verifies, reprioritizes or removes a recipient.
Monitoring Configuration	Sets sensor consent, countdown duration, alert recipients and optional location-sharing behavior.
Monitoring Status Toggle	Enables or disables monitoring after permissions and configuration are valid.

Field name	Description
Limitations Notice	States that monitoring is optional, device-dependent and not certified emergency detection.

Table 4.2.50: Emergency Contacts and Smart Activity Monitoring Screen

4.2.2.41 Safety Check and Incident Alert Screen

Platform: Mother Mobile App Feature: MF-14 Smart Activity Monitoring and Safety Support
Release: P1 Related UCs: UC-120

Purpose: Handles a suspected fall or strong impact using a safety countdown and I am OK, Need Help or emergency-support actions.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Mother Mobile App. Primary role: Mother / Phone IMU Sensor. Screen name: "Safety Check and Incident Alert". Purpose: Handles a suspected fall or strong impact using a safety countdown and I am OK, Need Help or emergency-support actions. Related Release 1 use cases: UC-120. Use a single portrait mobile frame, 390 x 844 px, with a 24 px safe margin. Use a simple top app bar and bottom navigation only for signed-in main-application screens. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Suspected Event Summary: Shows the detected time and general event type without claiming certainty.; Safety Countdown: Shows the time remaining before configured alert escalation.; I Am OK Action: Cancels escalation and records a safety confirmation.; Need Help Action: Triggers selected-family alert and emergency support handoff according to consent.; Emergency Support Shortcut: Opens Emergency Map and nearby support actions.. Interaction and safety constraints: Show only controls authorized for Mother / Phone IMU Sensor. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. This must be a prominent emergency-safe screen with a non-dispatch and non-certified detection disclaimer. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.2.51: Safety Check and Incident Alert Screen (to be generated in Stitch)

Field name	Description
Suspected Event Summary	Shows the detected time and general event type without claiming certainty.
Safety Countdown	Shows the time remaining before configured alert escalation.
I Am OK Action	Cancels escalation and records a safety confirmation.
Need Help Action	Triggers selected-family alert and emergency support handoff according to consent.
Emergency Support Shortcut	Opens Emergency Map and nearby support actions.

Table 4.2.51: Safety Check and Incident Alert Screen

4.2.2.42 Safety Event History and False Positive Feedback Screen

Platform: Mother Mobile App Feature: MF-14 Smart Activity Monitoring and Safety Support

Release: P1 Related UCs: UC-121

Purpose: Shows prior safety events and lets the user report a false positive with a reason.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Mother Mobile App. Primary role: Mother. Screen name: "Safety Event History and False Positive Feedback". Purpose: Shows prior safety events and lets the user report a false positive with a reason. Related Release 1 use cases: UC-121. Use a single portrait mobile frame, 390 x 844 px, with a 24 px safe margin. Use a simple top app bar and bottom navigation only for signed-in main-application screens. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Safety Event Timeline: Shows event time, response, alert state and final status.; Event Detail Panel: Shows permitted context and location-sharing state for the selected event.; False Positive Feedback: Lets the user mark an event as false detection and select or enter a reason.; Monitoring Shortcut: Returns to monitoring configuration.; Privacy Notice: Explains the retention and access boundary for safety event data.. Interaction and safety constraints: Show only controls authorized for Mother. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-

dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Do not reveal family recipient data beyond authorized context. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.2.52: Safety Event History and False Positive Feedback Screen (to be generated in Stitch)

Field name	Description
Safety Event Timeline	Shows event time, response, alert state and final status.
Event Detail Panel	Shows permitted context and location-sharing state for the selected event.
False Positive Feedback	Lets the user mark an event as false detection and select or enter a reason.
Monitoring Shortcut	Returns to monitoring configuration.
Privacy Notice	Explains the retention and access boundary for safety event data.

Table 4.2.52: Safety Event History and False Positive Feedback Screen

4.2.3 Family Member Mobile App

4.2.3.1 Care Group Invitation Screen

Platform: Family Member Mobile App Feature: MF-10 Family Sync and Cooperative Care Release: P1 Related UCs: UC-96

Purpose: Lets a recipient review and accept or reject an active care-group invitation after authentication.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Family Member Mobile App. Primary role: Family Member. Screen name: "Care Group Invitation". Purpose: Lets a recipient review and accept or reject an active care-group invitation after authentication. Related Release 1 use cases: UC-96. Use a single portrait mobile frame, 390 x 844 px, with a 24 px safe margin. Use a simple top app bar and bottom navigation only for signed-in main-application screens. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Invitation Summary: Shows group name, inviter, care context and invitation expiry.; Permission Preview: Shows the data and actions the recipient may receive after acceptance.; Accept Button: Joins the active care group after confirmation.; Reject Button: Declines the

invitation without joining.; Privacy Note: Explains that accepting does not grant access beyond the shown scope.. Interaction and safety constraints: Show only controls authorized for Family Member. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Do not show private maternal or baby data before the invitation is accepted. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.2.53: Care Group Invitation Screen (to be generated in Stitch)

Field name	Description
Invitation Summary	Shows group name, inviter, care context and invitation expiry.
Permission Preview	Shows the data and actions the recipient may receive after acceptance.
Accept Button	Joins the active care group after confirmation.
Reject Button	Declines the invitation without joining.
Privacy Note	Explains that accepting does not grant access beyond the shown scope.

Table 4.2.53: Care Group Invitation Screen

4.2.3.2 Family Shared Care Home Screen

Platform: Family Member Mobile App Feature: MF-10 Family Sync and Cooperative Care Release: P1 Related UCs: UC-100, UC-101

Purpose: Shows assigned task updates plus shared calendar, data and alerts within the family member active permission scope.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Family Member Mobile App. Primary role: Family Member. Screen name: "Family Shared Care Home". Purpose: Shows assigned task updates plus shared calendar, data and alerts within the family member active permission scope. Related Release 1 use cases: UC-100, UC-101. Use a single portrait mobile frame, 390 x 844 px, with a 24 px safe margin. Use a simple top app bar and bottom navigation only for signed-in main-application screens. Render exactly one complete screen, not a user-flow

diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Care Group Selector: Selects one joined care group when the member belongs to several.; Assigned Tasks: Shows tasks assigned to the signed-in member with status update controls.; Shared Calendar: Shows authorized reminders and calendar items.; Shared Data and Alerts: Shows only permitted summaries, records or emergency alerts.; Permission Scope Badge: Displays active scope and expiry for sensitive information.. Interaction and safety constraints: Show only controls authorized for Family Member. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Embed shared data and alert viewing inside the care-group home rather than creating multiple standalone browsing screens. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.2.54: Family Shared Care Home Screen (to be generated in Stitch)

Field name	Description
Care Group Selector	Selects one joined care group when the member belongs to several.
Assigned Tasks	Shows tasks assigned to the signed-in member with status update controls.
Shared Calendar	Shows authorized reminders and calendar items.
Shared Data and Alerts	Shows only permitted summaries, records or emergency alerts.
Permission Scope Badge	Displays active scope and expiry for sensitive information.

Table 4.2.54: Family Shared Care Home Screen

4.2.4 Verified Expert App

4.2.4.1 Expert Profile and Verification Documents Screen

Platform: Verified Expert App Feature: MF-05 Verified Expert Network and Contribution Release: P0 Related UCs: UC-60, UC-61, UC-62

Purpose: Creates or updates professional information and submits or replaces required verification evidence.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Verified Expert App. Primary role: Expert Applicant / Verified Expert. Screen name: "Expert Profile and Verification Documents". Purpose: Creates or updates professional information and submits or replaces required verification evidence. Related Release 1 use cases: UC-60, UC-61, UC-62. Use a single portrait mobile frame, 390 x 844 px, with a 24 px safe margin. Use a simple top app bar and bottom navigation only for signed-in main-application screens. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Professional Profile Form: Collects specialty, experience, approved support scope and public profile information.; Credential Upload Area: Adds required certificates and supporting documents with upload status.; Public Preview: Shows a compact preview of the data that may appear in the directory after approval.; Save Draft or Submit Action: Saves a draft or submits evidence for review.; Privacy and Review Notice: Explains that documents are protected and profile changes may require review.. Interaction and safety constraints: Show only controls authorized for Expert Applicant / Verified Expert. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Do not show revenue, paid consultation, commission or partner actions. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.2.55: Expert Profile and Verification Documents Screen (to be generated in Stitch)

Field name	Description
Professional Profile Form	Collects specialty, experience, approved support scope and public profile information.
Credential Upload Area	Adds required certificates and supporting documents with upload status.
Public Preview	Shows a compact preview of the data that may appear in the directory after approval.

Field name	Description
Save Draft or Submit Action	Saves a draft or submits evidence for review.
Privacy and Review Notice	Explains that documents are protected and profile changes may require review.

Table 4.2.55: Expert Profile and Verification Documents Screen

4.2.4.2 Expert Verification Status and Renewal Screen

Platform: Verified Expert App Feature: MF-05 Verified Expert Network and Contribution Release: P0 Related UCs: UC-63

Purpose: Shows verification result, missing information, expiry status and renewal submission action.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Verified Expert App. Primary role: Expert Applicant / Verified Expert. Screen name: "Expert Verification Status and Renewal". Purpose: Shows verification result, missing information, expiry status and renewal submission action. Related Release 1 use cases: UC-63. Use a single portrait mobile frame, 390 x 844 px, with a 24 px safe margin. Use a simple top app bar and bottom navigation only for signed-in main-application screens. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Verification Status Banner: Shows pending, approved, needs correction, expired or restricted state.; Credential Checklist: Shows required documents and review feedback.; Expiry Information: Shows expiration date when applicable.; Renew or Resubmit Action: Opens the eligible submission path.; Support Scope Notice: Explains that only approved status can provide verified-expert actions.. Interaction and safety constraints: Show only controls authorized for Expert Applicant / Verified Expert. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Keep review feedback factual and do not expose internal reviewer-only data. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.2.56: Expert Verification Status and Renewal Screen (to be generated in Stitch)

Field name	Description
Verification Status Banner	Shows pending, approved, needs correction, expired or restricted state.
Credential Checklist	Shows required documents and review feedback.
Expiry Information	Shows expiration date when applicable.
Renew or Resubmit Action	Opens the eligible submission path.
Support Scope Notice	Explains that only approved status can provide verified-expert actions.

Table 4.2.56: Expert Verification Status and Renewal Screen

4.2.4.3 Expert Availability and Service Scope Screen

Platform: Verified Expert App Feature: MF-05 Verified Expert Network and Contribution Release: P0 Related UCs: UC-64

Purpose: Configures availability, support methods, service area and public readiness state for directory and nearby support eligibility.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Verified Expert App. Primary role: Verified Expert. Screen name: "Expert Availability and Service Scope". Purpose: Configures availability, support methods, service area and public readiness state for directory and nearby support eligibility. Related Release 1 use cases: UC-64. Use a single portrait mobile frame, 390 x 844 px, with a 24 px safe margin. Use a simple top app bar and bottom navigation only for signed-in main-application screens. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Availability Status: Sets online, offline or scheduled availability.; Support Method Selector: Shows approved non-emergency support methods.; Service Area and Location Scope: Defines broad service area and optional time-limited location visibility.; Availability Schedule: Sets available periods and expiry of temporary readiness.; Save Availability Button: Publishes permitted availability data.. Interaction and safety constraints: Show only controls authorized for Verified Expert. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination

inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Do not include V2 booking, price, payments or realtime channel creation. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.2.57: Expert Availability and Service Scope Screen (to be generated in Stitch)

Field name	Description
Availability Status	Sets online, offline or scheduled availability.
Support Method Selector	Shows approved non-emergency support methods.
Service Area and Location Scope	Defines broad service area and optional time-limited location visibility.
Availability Schedule	Sets available periods and expiry of temporary readiness.
Save Availability Button	Publishes permitted availability data.

Table 4.2.57: Expert Availability and Service Scope Screen

4.2.4.4 Expert Question Queue Screen

Platform: Verified Expert App Feature: MF-05 Verified Expert Network and Contribution Release: P0 Related UCs: UC-67

Purpose: Shows community questions matched to verified specialties and approved public support scope.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Verified Expert App. Primary role: Verified Expert. Screen name: "Expert Question Queue". Purpose: Shows community questions matched to verified specialties and approved public support scope. Related Release 1 use cases: UC-67. Use a single portrait mobile frame, 390 x 844 px, with a 24 px safe margin. Use a simple top app bar and bottom navigation only for signed-in main-application screens. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Question Queue Filters: Embed topic, specialty match and status filters in the queue.; Question Cards: Show topic, safety flag, anonymized display and answer status.; Question Detail Preview: Shows the permitted context needed before answering.; Answer Action: Opens the verified expert answer composer.; Queue Notice: Explains that community answers do not create

a private consultation relationship.. Interaction and safety constraints: Show only controls authorized for Verified Expert. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Show only questions within verified specialty and moderation rules. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.2.58: Expert Question Queue Screen (to be generated in Stitch)

Field name	Description
Question Queue Filters	Embed topic, specialty match and status filters in the queue.
Question Cards	Show topic, safety flag, anonymized display and answer status.
Question Detail Preview	Shows the permitted context needed before answering.
Answer Action	Opens the verified expert answer composer.
Queue Notice	Explains that community answers do not create a private consultation relationship.

Table 4.2.58: Expert Question Queue Screen

4.2.4.5 Expert Answer Composer Screen

Platform: Verified Expert App Feature: MF-05 Verified Expert Network and Contribution Release: P0 Related UCs: UC-68

Purpose: Posts a verified expert answer with badge, source label, safety boundary and contribution attribution.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Verified Expert App. Primary role: Verified Expert. Screen name: "Expert Answer Composer". Purpose: Posts a verified expert answer with badge, source label, safety boundary and contribution attribution. Related Release 1 use cases: UC-68. Use a single portrait mobile frame, 390 x 844 px, with a 24 px safe margin. Use a simple top app bar and bottom navigation only for signed-in main-application screens. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat

rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Question Context: Shows the selected community question and topic.; Expert Answer Editor: Collects the answer under verified expert identity.; Verified and Source Labels: Shows professional verification and approved safety/source labels.; Safety Guidance: Reminds the expert not to diagnose, prescribe or misuse free answers for promotion.; Post or Update Button: Submits the answer for allowed publication state.. Interaction and safety constraints: Show only controls authorized for Verified Expert. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Do not include private contact, paid-consultation or realtime-session actions. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.2.59: Expert Answer Composer Screen (to be generated in Stitch)

Field name	Description
Question Context	Shows the selected community question and topic.
Expert Answer Editor	Collects the answer under verified expert identity.
Verified and Source Labels	Shows professional verification and approved safety/source labels.
Safety Guidance	Reminds the expert not to diagnose, prescribe or misuse free answers for promotion.
Post or Update Button	Submits the answer for allowed publication state.

Table 4.2.59: Expert Answer Composer Screen

4.2.4.6 Contribution Points and Badges Screen

Platform: Verified Expert App Feature: MF-05 Verified Expert Network and Contribution Release: P0 Related UCs: UC-69

Purpose: Shows contribution points, badges and qualifying community activities as participation recognition.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Verified Expert App. Primary role: Verified Expert. Screen name: "Contribution Points and Badges". Purpose: Shows contribution points, badges and qualifying community activities as participation recognition. Related Release 1 use cases: UC-69. Use a single portrait mobile frame, 390 x 844 px, with a 24 px safe margin. Use a simple top app bar and bottom navigation only for signed-in main-application screens. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Contribution Summary: Shows total points and current badges.; Qualifying Activity List: Shows approved public-answer and contribution activities.; Badge Criteria Hint: Explains how badges are earned or retained.; Trust Boundary Note: States that points and badges do not prove clinical competence or replace verification.. Interaction and safety constraints: Show only controls authorized for Verified Expert. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. No leaderboard should encourage unsafe volume-based advice. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.2.60: Contribution Points and Badges Screen (to be generated in Stitch)

Field name	Description
Contribution Summary	Shows total points and current badges.
Qualifying Activity List	Shows approved public-answer and contribution activities.
Badge Criteria Hint	Explains how badges are earned or retained.
Trust Boundary Note	States that points and badges do not prove clinical competence or replace verification.

Table 4.2.60: Contribution Points and Badges Screen

4.2.4.7 Nearby Availability and Support Requests Screen

Platform: Verified Expert App Feature: MF-07 Emergency Map and Nearby Care Support Release: P0 Related UCs: UC-82

Purpose: Lets an eligible expert opt in to nearby availability and respond to minimal-context nearby support requests.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Verified Expert App. Primary role: Verified Expert. Screen name: "Nearby Availability and Support Requests". Purpose: Lets an eligible expert opt in to nearby availability and respond to minimal-context nearby support requests. Related Release 1 use cases: UC-82. Use a single portrait mobile frame, 390 x 844 px, with a 24 px safe margin. Use a simple top app bar and bottom navigation only for signed-in main-application screens. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Nearby Availability Toggle: Opt in or out of nearby support eligibility with expiry and location scope.; Request List and Map Preview: Shows eligible requests with minimum necessary context and distance.; Request Detail: Shows permitted urgency label, time, consent scope and route option.; Accept, Decline or Stop Responding: Records the expert response without guaranteeing arrival.; Safety Boundary: States the feature is not emergency dispatch or a substitute for emergency services.. Interaction and safety constraints: Show only controls authorized for Verified Expert. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Do not disclose precise location or health records beyond the request consent scope. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.2.61: Nearby Availability and Support Requests Screen (to be generated in Stitch)

Field name	Description
Nearby Availability Toggle	Opt in or out of nearby support eligibility with expiry and location scope.
Request List and Map Preview	Shows eligible requests with minimum necessary context and distance.
Request Detail	Shows permitted urgency label, time, consent scope and route option.

Field name	Description
Accept, Decline or Stop Responding	Records the expert response without guaranteeing arrival.
Safety Boundary	States the feature is not emergency dispatch or a substitute for emergency services.

Table 4.2.61: Nearby Availability and Support Requests Screen

4.3 Web Screen

This section specifies Release 1 web screens for Expert, Moderator, Content Admin and System Admin roles. Each entry uses a consistent desktop portal shell for signed-in screens and excludes Version 2 consultation, payment, realtime, partner and operation-impact dashboard workflows.

4.3.1 Common Web Account Screens

4.3.1.1 Web Portal Login Screen

Platform: Common Web Portal Feature: MF-01 Account, Trust & Access Control Release: P0
Related UCs: UC-01, UC-02, UC-03, UC-05, UC-06

Purpose: Authenticates supported Expert, Moderator, Content Admin and System Admin accounts, with optional registration or recovery entry where applicable.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Common Web Portal. Primary role: Web Portal User. Screen name: "Web Portal Login". Purpose: Authenticates supported Expert, Moderator, Content Admin and System Admin accounts, with optional registration or recovery entry where applicable. Related Release 1 use cases: UC-01, UC-02, UC-03, UC-05, UC-06. Use a single desktop web frame, 1440 x 900 px. Use a centered authentication card; do not render the signed-in portal shell. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: CareBridge Wordmark: Identifies the portal.; Email or Phone Field: Collects the registered account identifier.; Password Field: Collects the masked password.; Login Button: Authenticates and routes by role.; Recovery and Registration Links: Open the supported account flows.. Interaction and safety constraints: Show only controls authorized for Web Portal User. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime

session or partner workflow. Use a centered authentication card and avoid displaying sensitive operational data before login. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.3.1: Web Portal Login Screen (to be generated in Stitch)

Field name	Description
CareBridge Wordmark	Identifies the portal.
Email or Phone Field	Collects the registered account identifier.
Password Field	Collects the masked password.
Login Button	Authenticates and routes by role.
Recovery and Registration Links	Open the supported account flows.

Table 4.3.1: Web Portal Login Screen

4.3.1.2 Web Account and Session Settings Screen

Platform: Common Web Portal Feature: MF-01 Account, Trust & Access Control Release: P0
Related UCs: UC-04, UC-07, UC-08, UC-09, UC-10, UC-11, UC-12, UC-13, UC-14, UC-15, UC-16

Purpose: Groups account profile, community identity, notification preferences, sessions, consent and account lifecycle settings.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Common Web Portal. Primary role: Web Portal User. Screen name: "Web Account and Session Settings". Purpose: Groups account profile, community identity, notification preferences, sessions, consent and account lifecycle settings. Related Release 1 use cases: UC-04, UC-07, UC-08, UC-09, UC-10, UC-11, UC-12, UC-13, UC-14, UC-15, UC-16. Use a single desktop web frame, 1440 x 900 px. Use a simple top header and left sidebar for the signed-in portal shell. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Profile and Role Summary: Shows permitted private account fields and role status.; Settings Navigation: Switches among profile, community identity, notifications, password, sessions, privacy and lifecycle sections.; Consent Summary: Shows active permissions with recipient, scope, purpose and expiry.; Session Management: Shows active sessions and revoke controls.; Lifecycle Controls: Provides logout, deactivate or delete actions with in-screen confirmation.. Interaction and safety constraints: Show

only controls authorized for Web Portal User. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Use a desktop shell with header and sidebar. Enforce role scope for all settings. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.3.2: Web Account and Session Settings Screen (to be generated in Stitch)

Field name	Description
Profile and Role Summary	Shows permitted private account fields and role status.
Settings Navigation	Switches among profile, community identity, notifications, password, sessions, privacy and lifecycle sections.
Consent Summary	Shows active permissions with recipient, scope, purpose and expiry.
Session Management	Shows active sessions and revoke controls.
Lifecycle Controls	Provides logout, deactivate or delete actions with in-screen confirmation.

Table 4.3.2: Web Account and Session Settings Screen

4.3.2 Verified Expert Web Portal

4.3.2.1 Expert Portal Dashboard Screen

Platform: Expert Web Portal Feature: MF-05 Verified Expert Network and Contribution Release: P0 Related UCs: UC-64, UC-67, UC-69

Purpose: Summarizes professional profile status, availability, question queue activity, contribution recognition and nearby-support readiness.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Expert Web Portal. Primary role: Verified Expert. Screen name: "Expert Portal Dashboard". Purpose: Summarizes professional profile status, availability, question queue activity, contribution recognition and nearby-support readiness. Related Release 1 use cases: UC-64, UC-67, UC-69. Use a single desktop web frame, 1440 x 900 px. Use a simple top header and left sidebar for the signed-in portal shell. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS

wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Profile and Verification Summary: Shows current expert status and any pending action.; Availability Card: Shows public readiness and shortcut to edit service scope.; Question Queue Summary: Shows matched open questions and answer activity.; Contribution Summary: Shows points, badges and qualifying activities.; Nearby Support Readiness: Shows opt-in state without exposing user data.. Interaction and safety constraints: Show only controls authorized for Verified Expert. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Do not include V2 paid consultation, revenue, commission or realtime modules. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.3.3: Expert Portal Dashboard Screen (to be generated in Stitch)

Field name	Description
Profile and Verification Summary	Shows current expert status and any pending action.
Availability Card	Shows public readiness and shortcut to edit service scope.
Question Queue Summary	Shows matched open questions and answer activity.
Contribution Summary	Shows points, badges and qualifying activities.
Nearby Support Readiness	Shows opt-in state without exposing user data.

Table 4.3.3: Expert Portal Dashboard Screen

4.3.2.2 Expert Professional Profile and Documents Screen

Platform: Expert Web Portal Feature: MF-05 Verified Expert Network and Contribution Release: P0 Related UCs: UC-60, UC-61, UC-62, UC-63

Purpose: Creates or updates expert professional profile and manages required verification document submissions.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Expert Web Portal. Primary role: Expert Applicant / Verified Expert. Screen name: "Expert Professional Profile and Documents". Purpose: Creates or updates expert professional profile and manages required verification document submissions. Related Release 1 use cases: UC-60, UC-61, UC-62, UC-63. Use a single desktop web frame, 1440 x 900 px. Use a simple top header and left sidebar for the signed-in portal shell. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Professional Profile Form: Collects specialty, experience, service scope and public professional description.; Credential Document Area: Shows submitted documents, expiry and replacement controls.; Verification Status Panel: Shows pending, approved, correction or renewal state.; Public Directory Preview: Previews approved public professional information.; Submit or Update Action: Submits changes under review policy.. Interaction and safety constraints: Show only controls authorized for Expert Applicant / Verified Expert. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Protect uploaded evidence and do not expose reviewer-only information. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.3.4: Expert Professional Profile and Documents Screen (to be generated in Stitch)

Field name	Description
Professional Profile Form	Collects specialty, experience, service scope and public professional description.
Credential Document Area	Shows submitted documents, expiry and replacement controls.
Verification Status Panel	Shows pending, approved, correction or renewal state.
Public Directory Preview	Previews approved public professional information.
Submit or Update Action	Submits changes under review policy.

Table 4.3.4: Expert Professional Profile and Documents Screen

4.3.2.3 Expert Availability and Nearby Support Screen

Platform: Expert Web Portal Feature: MF-05 Verified Expert Network and Contribution / MF-07 Emergency Map Release: P0 Related UCs: UC-64, UC-82

Purpose: Configures availability and service scope, and manages nearby-support opt-in and eligible request responses.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Expert Web Portal. Primary role: Verified Expert. Screen name: "Expert Availability and Nearby Support". Purpose: Configures availability and service scope, and manages nearby-support opt-in and eligible request responses. Related Release 1 use cases: UC-64, UC-82. Use a single desktop web frame, 1440 x 900 px. Use a simple top header and left sidebar for the signed-in portal shell. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Availability Calendar: Sets current and scheduled readiness states.; Service Scope Controls: Defines allowed support methods and public service area.; Nearby Availability Scope: Controls location visibility, accuracy and expiry for nearby support.; Eligible Support Request List: Shows minimal request context with accept or decline controls.; Safety Notice: States no dispatch or arrival guarantee.. Interaction and safety constraints: Show only controls authorized for Verified Expert. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Do not reveal exact user location or private health records beyond consented request context. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.3.5: Expert Availability and Nearby Support Screen (to be generated in Stitch)

Field name	Description
Availability Calendar	Sets current and scheduled readiness states.
Service Scope Controls	Defines allowed support methods and public service area.
Nearby Availability Scope	Controls location visibility, accuracy and expiry for nearby support.

Field name	Description
Eligible Support Request List	Shows minimal request context with accept or decline controls.
Safety Notice	States no dispatch or arrival guarantee.

Table 4.3.5: Expert Availability and Nearby Support Screen

4.3.2.4 Expert Question Queue and Answer Editor Screen

Platform: Expert Web Portal Feature: MF-05 Verified Expert Network and Contribution Release: P0 Related UCs: UC-67, UC-68

Purpose: Browses specialty-matched community questions and composes verified expert answers.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Expert Web Portal. Primary role: Verified Expert. Screen name: "Expert Question Queue and Answer Editor". Purpose: Browses specialty-matched community questions and composes verified expert answers. Related Release 1 use cases: UC-67, UC-68. Use a single desktop web frame, 1440 x 900 px. Use a simple top header and left sidebar for the signed-in portal shell. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Queue Filters: Embed specialty, topic, safety and answer-status filters.; Question List: Shows matched questions with moderation and anonymous-display state.; Question Detail Pane: Shows permitted context and reply thread.; Answer Editor: Composes a verified expert response with safety/source label.; Post Action: Submits the answer when allowed.. Interaction and safety constraints: Show only controls authorized for Verified Expert. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Use a two-pane desktop layout and keep all content inside the community and expert-scope boundary. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.3.6: Expert Question Queue and Answer Editor Screen (to be generated in Stitch)

Field name	Description
Queue Filters	Embed specialty, topic, safety and answer-status filters.
Question List	Shows matched questions with moderation and anonymous-display state.
Question Detail Pane	Shows permitted context and reply thread.
Answer Editor	Composes a verified expert response with safety/source label.
Post Action	Submits the answer when allowed.

Table 4.3.6: Expert Question Queue and Answer Editor Screen

4.3.3 Moderator Web Workspace

4.3.3.1 Moderator Dashboard and Moderation Queue Screen

Platform: Moderator Web Portal Feature: MF-04 Community Q&A and Moderation Release: P0
Related UCs: UC-56

Purpose: Shows reported, automatically flagged and pre-publication community content awaiting review.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Moderator Web Portal. Primary role: Moderator. Screen name: "Moderator Dashboard and Moderation Queue". Purpose: Shows reported, automatically flagged and pre-publication community content awaiting review. Related Release 1 use cases: UC-56. Use a single desktop web frame, 1440 x 900 px. Use a simple top header and left sidebar for the signed-in portal shell. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Queue Summary Cards: Show pending, high-risk, reported and pre-publication counts.; Moderation Filters: Embed content type, safety level, topic, queue state and date filters.; Queue Table: Shows item summary, report reason, flag source, time and current state.; Priority Safety Marker: Highlights content requiring urgent review using text and grayscale emphasis.; Open Item Action: Opens the selected moderation item.. Interaction and safety constraints: Show only controls authorized for Moderator. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the

owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Do not include operations-impact dashboard metrics or partner workflows. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.3.7: Moderator Dashboard and Moderation Queue Screen (to be generated in Stitch)

Field name	Description
Queue Summary Cards	Show pending, high-risk, reported and pre-publication counts.
Moderation Filters	Embed content type, safety level, topic, queue state and date filters.
Queue Table	Shows item summary, report reason, flag source, time and current state.
Priority Safety Marker	Highlights content requiring urgent review using text and grayscale emphasis.
Open Item Action	Opens the selected moderation item.

Table 4.3.7: Moderator Dashboard and Moderation Queue Screen

4.3.3.2 Moderation Item and Report Resolution Screen

Platform: Moderator Web Portal Feature: MF-04 Community Q&A and Moderation Release: P0
Related UCs: UC-57, UC-58

Purpose: Reviews a selected item or report, applies moderation action and resolves a content or account report under policy.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Moderator Web Portal. Primary role: Moderator / System Admin. Screen name: "Moderation Item and Report Resolution". Purpose: Reviews a selected item or report, applies moderation action and resolves a content or account report under policy. Related Release 1 use cases: UC-57, UC-58. Use a single desktop web frame, 1440 x 900 px. Use a simple top header and left sidebar for the signed-in portal shell. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Reported Content Context: Shows post or answer content, moderation history and report summary.; Policy and Safety Evidence: Shows relevant policy categories, flag reason and

permitted evidence.; Moderation Actions: Approve, hide, lock, label, request revision or escalate the item.; Enforcement Panel: Applies warning, posting restriction, suspension escalation or no action when authorized.; Resolution Note and Audit Notice: Records the decision reason and states it will be audited.. Interaction and safety constraints: Show only controls authorized for Moderator / System Admin. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Protect reporter identity and avoid exposing unrelated user health data. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.3.8: Moderation Item and Report Resolution Screen (to be generated in Stitch)

Field name	Description
Reported Content Context	Shows post or answer content, moderation history and report summary.
Policy and Safety Evidence	Shows relevant policy categories, flag reason and permitted evidence.
Moderation Actions	Approve, hide, lock, label, request revision or escalate the item.
Enforcement Panel	Applies warning, posting restriction, suspension escalation or no action when authorized.
Resolution Note and Audit Notice	Records the decision reason and states it will be audited.

Table 4.3.8: Moderation Item and Report Resolution Screen

4.3.3.3 Community Topics and Visibility Screen

Platform: Moderator Web Portal Feature: MF-04 Community Q&A and Moderation Release: P0
Related UCs: UC-59

Purpose: Creates, updates, reorders or disables community topics and visibility rules.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Moderator Web Portal. Primary role: Moderator / Content Admin. Screen name: "Community Topics and Visibility". Purpose: Creates, updates, reorders or disables community topics and visibility rules. Related

Release 1 use cases: UC-59. Use a single desktop web frame, 1440 x 900 px. Use a simple top header and left sidebar for the signed-in portal shell. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Topic List: Shows topic name, description, order and active visibility state.; Topic Editor: Creates or updates permitted topic metadata.; Display Order Control: Reorders topics used by the community feed.; Visibility Toggle: Enables, disables or archives a topic under policy.; Save Changes Button: Applies the updated topic configuration.. Interaction and safety constraints: Show only controls authorized for Moderator / Content Admin. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Topic changes should not silently alter existing post history or user consent. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.3.9: Community Topics and Visibility Screen (to be generated in Stitch)

Field name	Description
Topic List	Shows topic name, description, order and active visibility state.
Topic Editor	Creates or updates permitted topic metadata.
Display Order Control	Reorders topics used by the community feed.
Visibility Toggle	Enables, disables or archives a topic under policy.
Save Changes Button	Applies the updated topic configuration.

Table 4.3.9: Community Topics and Visibility Screen

4.3.4 Content Admin Web Workspace

4.3.4.1 Content Dashboard and Library Screen

Platform: Content Admin Web Portal Feature: MF-11 Verified Content and Checklist Hub Release: P1 Related UCs: UC-104, UC-105

Purpose: Shows verified-content drafts, published items, review state and library navigation.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Content Admin Web Portal. Primary role: Content Admin. Screen name: "Content Dashboard and Library". Purpose: Shows verified-content drafts, published items, review state and library navigation. Related Release 1 use cases: UC-104, UC-105. Use a single desktop web frame, 1440 x 900 px. Use a simple top header and left sidebar for the signed-in portal shell. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Content Summary Cards: Show drafts, pending review, published and archived item counts.; Content Filters: Embed type, stage, topic, status, source and version filters.; Content List: Shows article, FAQ and checklist title, status, last update and owner.; Create Content Action: Opens the content editor for a new draft.; Open Detail Action: Opens the selected draft or published version.. Interaction and safety constraints: Show only controls authorized for Content Admin. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Do not expose V2 partner sponsored-content workflows. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.3.10: Content Dashboard and Library Screen (to be generated in Stitch)

Field name	Description
Content Summary Cards	Show drafts, pending review, published and archived item counts.
Content Filters	Embed type, stage, topic, status, source and version filters.
Content List	Shows article, FAQ and checklist title, status, last update and owner.
Create Content Action	Opens the content editor for a new draft.
Open Detail Action	Opens the selected draft or published version.

Table 4.3.10: Content Dashboard and Library Screen

4.3.4.2 Content Editor and Source Management Screen

Platform: Content Admin Web Portal Feature: MF-11 Verified Content and Checklist Hub Release: P1 Related UCs: UC-104, UC-105

Purpose: Creates or updates verified articles, FAQs and checklists with source, stage, topic and safety metadata.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Content Admin Web Portal. Primary role: Content Admin. Screen name: "Content Editor and Source Management". Purpose: Creates or updates verified articles, FAQs and checklists with source, stage, topic and safety metadata. Related Release 1 use cases: UC-104, UC-105. Use a single desktop web frame, 1440 x 900 px. Use a simple top header and left sidebar for the signed-in portal shell. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Content Type Selector: Selects article, FAQ or checklist.; Title and Body Editor: Collects the content text and action guidance.; Source and Review Metadata: Captures source labels, reviewed date, safety note and version context.; Stage and Topic Mapping: Maps the draft to applicable care stages and topics.; Save Draft or Submit Review Action: Saves a draft or sends it for publication review.. Interaction and safety constraints: Show only controls authorized for Content Admin. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Do not create unlabelled medical advice or unsafe clinical claims. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.3.11: Content Editor and Source Management Screen (to be generated in Stitch)

Field name	Description
Content Type Selector	Selects article, FAQ or checklist.
Title and Body Editor	Collects the content text and action guidance.
Source and Review Metadata	Captures source labels, reviewed date, safety note and version context.

Field name	Description
Stage and Topic Mapping	Maps the draft to applicable care stages and topics.
Save Draft or Submit Review Action	Saves a draft or sends it for publication review.

Table 4.3.11: Content Editor and Source Management Screen

4.3.4.3 Content Review, Publication and Archive Screen

Platform: Content Admin Web Portal Feature: MF-11 Verified Content and Checklist Hub Release: P1 Related UCs: UC-106, UC-107

Purpose: Reviews a content version and approves, rejects, requests revision, unpublishes or archives it.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Content Admin Web Portal. Primary role: Content Admin / System Admin. Screen name: "Content Review, Publication and Archive". Purpose: Reviews a content version and approves, rejects, requests revision, unpublishes or archives it. Related Release 1 use cases: UC-106, UC-107. Use a single desktop web frame, 1440 x 900 px. Use a simple top header and left sidebar for the signed-in portal shell. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Version Comparison: Shows current draft against the previous published version when available.; Review Checklist: Shows source, safety, stage mapping and policy review checks.; Decision Controls: Approve, reject, request revision, unpublish or archive the version.; Decision Note: Records a concise decision reason for author and audit use.; Publication State: Shows draft, pending, published, withdrawn or archived state.. Interaction and safety constraints: Show only controls authorized for Content Admin / System Admin. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Published content must retain version and review context; do not silently overwrite it. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.3.12: Content Review, Publication and Archive Screen (to be generated in Stitch)

Field name	Description
Version Comparison	Shows current draft against the previous published version when available.
Review Checklist	Shows source, safety, stage mapping and policy review checks.
Decision Controls	Approve, reject, request revision, unpublish or archive the version.
Decision Note	Records a concise decision reason for author and audit use.
Publication State	Shows draft, pending, published, withdrawn or archived state.

Table 4.3.12: Content Review, Publication and Archive Screen

4.3.4.4 Content Categories and Stage Topic Mapping Screen

Platform: Content Admin Web Portal Feature: MF-11 Verified Content and Checklist Hub Release: P1 Related UCs: UC-108

Purpose: Maintains categories, stage mappings, topic tags and display order used by the content hub.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Content Admin Web Portal. Primary role: Content Admin. Screen name: "Content Categories and Stage Topic Mapping". Purpose: Maintains categories, stage mappings, topic tags and display order used by the content hub. Related Release 1 use cases: UC-108. Use a single desktop web frame, 1440 x 900 px. Use a simple top header and left sidebar for the signed-in portal shell. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Category List: Shows category name, type, active state and display order.; Stage Mapping Grid: Maps content categories to preconception, pregnancy, postpartum or baby-care stages.; Topic Tag Management: Creates, updates or deactivates topic tags.; Display Order Controls: Reorders categories and tags shown in the content hub.; Save Configuration Button: Applies approved mapping changes.. Interaction and safety constraints: Show only controls authorized for Content Admin. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or

list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Avoid changes that break existing content links or hide safety information. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.3.13: Content Categories and Stage Topic Mapping Screen (to be generated in Stitch)

Field name	Description
Category List	Shows category name, type, active state and display order.
Stage Mapping Grid	Maps content categories to preconception, pregnancy, postpartum or baby-care stages.
Topic Tag Management	Creates, updates or deactivates topic tags.
Display Order Controls	Reorders categories and tags shown in the content hub.
Save Configuration Button	Applies approved mapping changes.

Table 4.3.13: Content Categories and Stage Topic Mapping Screen

4.3.4.5 Approved AI Knowledge Sources Screen

Platform: Content Admin Web Portal Feature: MF-06 AI Nurse Assistant and Risk Triage Release: P0 Related UCs: UC-75

Purpose: Registers, reviews, versions or disables approved knowledge sources available to AI retrieval and safe answer composition.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: Content Admin Web Portal. Primary role: Content Admin / System Admin. Screen name: "Approved AI Knowledge Sources". Purpose: Registers, reviews, versions or disables approved knowledge sources available to AI retrieval and safe answer composition. Related Release 1 use cases: UC-75. Use a single desktop web frame, 1440 x 900 px. Use a simple top header and left sidebar for the signed-in portal shell. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Knowledge Source List: Shows source name, owner, review state, version and active status.; Source Detail and Evidence: Shows reference metadata, coverage and safety notes.; Review Controls: Approve, disable, version or request updates to a source.; Usage Scope: Shows which AI topics or safe-response contexts the source supports.; Audit Note: Records that changes

are traceable for safety review.. Interaction and safety constraints: Show only controls authorized for Content Admin / System Admin. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Do not allow unreviewed community content to become an approved AI knowledge source. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.3.14: Approved AI Knowledge Sources Screen (to be generated in Stitch)

Field name	Description
Knowledge Source List	Shows source name, owner, review state, version and active status.
Source Detail and Evidence	Shows reference metadata, coverage and safety notes.
Review Controls	Approve, disable, version or request updates to a source.
Usage Scope	Shows which AI topics or safe-response contexts the source supports.
Audit Note	Records that changes are traceable for safety review.

Table 4.3.14: Approved AI Knowledge Sources Screen

4.3.5 System Admin Web Workspace

4.3.5.1 System Administration Dashboard and User Accounts Screen

Platform: System Admin Web Portal Feature: MF-01 Account, Trust & Access Control Release: P0
Related UCs: UC-17

Purpose: Reviews account status and applies authorized role, access and misuse restrictions under separation-of-duties controls.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: System Admin Web Portal. Primary role: System Admin. Screen name: "System Administration Dashboard and User Accounts". Purpose: Reviews account status and applies authorized role, access and misuse restrictions under separation-of-duties controls. Related Release 1 use cases: UC-17. Use a single desktop web frame, 1440 x 900 px. Use a simple top header and left sidebar for the signed-in portal shell. Render exactly one complete screen, not a user-flow diagram. Visual style: basic

academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Administration Summary: Shows account status, verification workload and high-level access-risk counts.; User List Filters: Embed role, status, verification and risk indicator filters.; User Account Table: Shows account, role, status and permitted action shortcuts.; Role and Access Panel: Applies authorized role or access changes with reason capture.; Restriction Controls: Applies account restriction under policy and audit notice.. Interaction and safety constraints: Show only controls authorized for System Admin. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Do not include an operation-impact dashboard, partner management, payment or consultation modules. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.3.15: System Administration Dashboard and User Accounts Screen (to be generated in Stitch)

Field name	Description
Administration Summary	Shows account status, verification workload and high-level access-risk counts.
User List Filters	Embed role, status, verification and risk indicator filters.
User Account Table	Shows account, role, status and permitted action shortcuts.
Role and Access Panel	Applies authorized role or access changes with reason capture.
Restriction Controls	Applies account restriction under policy and audit notice.

Table 4.3.15: System Administration Dashboard and User Accounts Screen

4.3.5.2 Expert Verification Queue and Trust Status Screen

Platform: System Admin Web Portal Feature: MF-05 Verified Expert Network and Contribution
Release: P0 Related UCs: UC-70, UC-71

Purpose: Reviews expert verification submissions and restricts, suspends or reinstates expert trust status when policy requires it.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: System Admin Web Portal. Primary role: System Admin. Screen name: "Expert Verification Queue and Trust Status". Purpose: Reviews expert verification submissions and restricts, suspends or reinstates expert trust status when policy requires it. Related Release 1 use cases: UC-70, UC-71. Use a single desktop web frame, 1440 x 900 px. Use a simple top header and left sidebar for the signed-in portal shell. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Verification Queue: Shows pending, renewal, correction and expired verification submissions.; Credential Review Pane: Shows protected evidence, expiry dates and required review checks.; Decision Controls: Approve, request supplementation, reject or renew the submission.; Trust Status Panel: Applies, changes or lifts restriction and suspension status with reasons.; Audit and Separation Notice: Shows that decisions are tracked and must follow permission boundaries.. Interaction and safety constraints: Show only controls authorized for System Admin. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Do not expose expert private documentation to users outside administrator scope. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.3.16: Expert Verification Queue and Trust Status Screen (to be generated in Stitch)

Field name	Description
Verification Queue	Shows pending, renewal, correction and expired verification submissions.
Credential Review Pane	Shows protected evidence, expiry dates and required review checks.
Decision Controls	Approve, request supplementation, reject or renew the submission.
Trust Status Panel	Applies, changes or lifts restriction and suspension status with reasons.

Field name	Description
Audit and Separation Notice	Shows that decisions are tracked and must follow permission boundaries.

Table 4.3.16: Expert Verification Queue and Trust Status Screen

4.3.5.3 Sensitive Access and Security Events Screen

Platform: System Admin Web Portal Feature: MF-01 Account, Trust & Access Control Release: P0
Related UCs: UC-18

Purpose: Reviews sensitive record access, abnormal login, permission changes and protected-file events for follow-up or investigation.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: System Admin Web Portal. Primary role: System Admin. Screen name: "Sensitive Access and Security Events". Purpose: Reviews sensitive record access, abnormal login, permission changes and protected-file events for follow-up or investigation. Related Release 1 use cases: UC-18. Use a single desktop web frame, 1440 x 900 px. Use a simple top header and left sidebar for the signed-in portal shell. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Event Filters: Embed event type, time range, actor, target type and severity filters.; Security Event Table: Shows event time, actor, action, target category, result and risk indicator.; Event Detail Panel: Shows permitted audit context and relevant policy flags.; Investigation Note: Records follow-up decision and action status.; Privacy Boundary: States that audit review does not authorize unnecessary health-data browsing.. Interaction and safety constraints: Show only controls authorized for System Admin. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Do not turn each audit event into a separate screen or expose raw protected data by default. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.3.17: Sensitive Access and Security Events Screen (to be generated in Stitch)

Field name	Description
Event Filters	Embed event type, time range, actor, target type and severity filters.
Security Event Table	Shows event time, actor, action, target category, result and risk indicator.
Event Detail Panel	Shows permitted audit context and relevant policy flags.
Investigation Note	Records follow-up decision and action status.
Privacy Boundary	States that audit review does not authorize unnecessary health-data browsing.

Table 4.3.17: Sensitive Access and Security Events Screen

4.3.5.4 AI Risk and Red-Flag Rules Screen

Platform: System Admin Web Portal Feature: MF-06 AI Nurse Assistant and Risk Triage Release: P0 Related UCs: UC-76

Purpose: Maintains conservative risk thresholds, red-flag patterns, safe fallback actions and rule version history.

Stitch Wireframe Prompt

Create a low-fidelity grayscale wireframe mockup for CareBridge. Platform: System Admin Web Portal. Primary role: System Admin. Screen name: "AI Risk and Red-Flag Rules". Purpose: Maintains conservative risk thresholds, red-flag patterns, safe fallback actions and rule version history. Related Release 1 use cases: UC-76. Use a single desktop web frame, 1440 x 900 px. Use a simple top header and left sidebar for the signed-in portal shell. Render exactly one complete screen, not a user-flow diagram. Visual style: basic academic SRS wireframe with white background (#FFFFFF), light gray panels (#F2F2F2), light gray input and table cells (#E6E6E6), medium gray 1 px borders (#777777), black text (#111111), Arial or Inter, flat rectangles, 4 px maximum corner radius, no gradients, shadows, glow, photographs, colorful icons, illustrations, decorative patterns, 3D or brand-color accents. Keep labels in English. Use realistic but fictional placeholder data only; do not use real medical or personal data. Required components: Rule List: Shows rule name, active version, risk level, last review and status.; Risk Threshold Editor: Defines conservative input conditions and green, yellow or red orientation.; Red-Flag Pattern Editor: Maintains patterns that require emergency or offline-care guidance.; Safe Fallback Actions: Defines non-diagnostic next actions when certainty is low or service fails.; Version and Review History: Shows version, reviewer and effective date.. Interaction and safety constraints: Show only controls authorized for System Admin. Maintain privacy by default. Display source labels for user-entered, device-sourced, expert-provided, or system-generated information when relevant. When sensitive health data, protected files, location, device observations or shared-care context is visible, show a concise consent scope and expiry badge. Never present AI output as a diagnosis, prescription, dosage or independent clinical decision. On risk, emergency, map and safety states, display a clear non-dispatch

emergency-care disclaimer. Embed search, filters, sort and pagination inside the owning browse or list screen. Do not add a standalone search screen, confirmation-dialog screen, file-manager screen, Version 2 paid consultation, payment, commission, realtime session or partner workflow. Do not provide fields for diagnosis, prescription, dosage or independent clinical decision-making. Produce exactly one grayscale wireframe screen suitable for insertion into a Software Requirement Specification document.

Figure 4.3.18: AI Risk and Red-Flag Rules Screen (to be generated in Stitch)

Field name	Description
Rule List	Shows rule name, active version, risk level, last review and status.
Risk Threshold Editor	Defines conservative input conditions and green, yellow or red orientation.
Red-Flag Pattern Editor	Maintains patterns that require emergency or offline-care guidance.
Safe Fallback Actions	Defines non-diagnostic next actions when certainty is low or service fails.
Version and Review History	Shows version, reviewer and effective date.

Table 4.3.18: AI Risk and Red-Flag Rules Screen

5. Non-Functional Requirements

This section specifies the non-functional requirements that govern how CareBridge must operate across the Mobile App, Expert Portal, Admin & Moderation Portal, Partner Portal, backend services and external integrations. These requirements apply to all functional requirements, workflows, screens and business rules unless a stricter requirement is defined in a specific use case.

5.1 External Interfaces

CareBridge communicates with human users, mobile device capabilities, web browsers and external software services. All interfaces must enforce role-based access control, explicit consent, auditability and healthcare-safety boundaries defined in the functional specification and appendix.

5.1.1 User Interfaces

- **NFR-UI-01:** The mobile app must support Mother, Family Member/Caregiver and Verified Expert workflows through responsive screens for account setup, care journey tracking, community Q&A, AI triage, emergency map, consultation, family sync, connected devices, pregnancy exercise and safety monitoring.
- **NFR-UI-02:** The web portals must support Verified Expert, System Admin, Moderator, Content Admin and Partner Representative workflows through role-specific dashboards, queues, detail screens, configuration screens and audit views.

- **NFR-UI-03:** All screens that collect sensitive health, location, payment or identity data must display clear validation, confirmation or safety messages before the user submits or shares the data.
- **NFR-UI-04:** Safety-critical outputs, including AI triage, emergency alerts, fall-detection events, nearby-care guidance and consultation boundaries, must be visually distinguishable from normal informational messages.

5.1.2 Software and Service Interfaces

- **NFR-SI-01:** The backend must expose secure APIs for mobile and web clients using authenticated requests, role checks, consent checks and request validation before processing protected data.
- **NFR-SI-02:** VNPay payment integration must support consultation fee payment, payment status callback, commission calculation, dispute/refund tracing and reconciliation without exposing raw card or bank credentials to CareBridge.
- **NFR-SI-03:** Firebase Cloud Messaging or an equivalent notification service must support reminder, community, consultation, moderation, emergency and family-alert notifications with retry and delivery-state logging.
- **NFR-SI-04:** TrackAsia Map Service must support nearby facility search, distance calculation, ETA, route display and location-based emergency or nearby-expert workflows only after location permission is granted.
- **NFR-SI-05:** Gemini AI Service and the RAG knowledge service must support symptom-intake structuring, safe guidance generation and verified-content retrieval while preserving disclaimers that the system does not diagnose, prescribe or replace emergency care.
- **NFR-SI-06:** ZegoCloud and Firebase Realtime services must support chat, voice, video or realtime consultation sessions for authorized participants and must expose connection status, session errors and termination events to the application.
- **NFR-SI-07:** The audit/security log interface must record consent changes, expert access, payment/refund events, admin actions, moderation actions, AI/safety events and security incidents in an immutable or tamper-evident manner.

5.1.3 Hardware, Device and Sensor Interfaces

- **NFR-DI-01:** Health device and smartwatch integrations must collect only user-approved data types, such as heart rate, steps, sleep, activity or related health trend data, and must allow users to disconnect the device and stop future sync.
- **NFR-DI-02:** Phone IMU sensor access for smart activity monitoring and fall-detection support must be permission-based, configurable and disableable by the user.
- **NFR-DI-03:** Camera access for MediaPipe posture analysis must be requested only during pregnancy exercise posture-check workflows and must not be used for unrelated background monitoring.
- **NFR-DI-04:** When sensor, wearable, camera or location permission is denied, unavailable or interrupted, the system must provide a graceful fallback such as manual input, limited-mode guidance or an explanatory message.

5.1.4 Communication Interfaces

- **NFR-CI-01:** All client-server and server-service communication that carries account, health, consultation, payment, location or audit data must use encrypted transport such as HTTPS/TLS.
- **NFR-CI-02:** External-service requests must include timeout handling, retry where appropriate and safe degradation so that third-party failures do not corrupt CareBridge data or bypass safety rules.
- **NFR-CI-03:** The system must record integration errors with correlation information sufficient for troubleshooting while avoiding exposure of sensitive health, payment or authentication secrets in logs.

5.2 Quality Attributes

The following quality attributes define measurable expectations for usability, reliability, performance, security, safety, maintainability, scalability, compatibility and data quality. They are mandatory acceptance criteria for CareBridge implementation and testing.

5.2.1 Usability

- **NFR-US-01:** A first-time Mother user must be able to register, verify OTP, set up a basic journey profile and reach the home dashboard within 10 minutes without staff assistance under normal network conditions.
- **NFR-US-02:** A returning user must be able to log in, view today tasks/reminders, view notifications and open a care record within 3 minutes.
- **NFR-US-03:** A Family Member/Caregiver must be able to accept an invitation, view shared care information and complete an assigned task within 5 minutes after receiving the invitation.
- **NFR-US-04:** A Verified Expert must be able to update availability, view a consultation request and open the shared health summary within 5 minutes after login.
- **NFR-US-05:** Admin, moderator and content users must be able to process normal queue items using consistent list, filter, detail, decision and confirmation patterns across all web portals.
- **NFR-US-06:** All forms must provide field-level validation messages, preserve user-entered data after recoverable errors and clearly distinguish required fields from optional fields.
- **NFR-US-07:** The system must provide Vietnamese UI text for end users and may provide English labels where required by academic documentation, admin operation or external-service terminology.
- **NFR-US-08:** Destructive or sensitive actions such as deleting accounts, revoking permissions, disconnecting devices, cancelling consultation, refund decision, suspending users or sending emergency alerts must require confirmation before execution.

5.2.2 Reliability

- **NFR-RL-01:** Core authenticated functions, including login, profile access, care records, reminders, community feed, expert directory and notifications, must target at least 99.5% monthly availability during the pilot period, excluding scheduled maintenance.

- **NFR-RL-02:** Safety-related functions, including emergency alert creation, emergency map opening, family alert delivery attempt and safety-monitoring event recording, must target at least 99.9% monthly backend availability during the pilot period.
- **NFR-RL-03:** Planned maintenance should be scheduled during low-traffic hours and announced to affected users at least 24 hours in advance when user-facing downtime is expected.
- **NFR-RL-04:** Mean Time To Repair for critical production incidents affecting login, health record access, emergency alert processing or consultation payment should not exceed 4 hours during the pilot support window.
- **NFR-RL-05:** The system must prevent loss of confirmed user-created records after successful save, including health logs, baby logs, vaccination records, growth measurements, community posts, expert answers, consultation summaries and payment/refund records.
- **NFR-RL-06:** Notification, AI, map, payment, realtime and device-sync failures must be captured as recoverable errors with user-safe messages and retry/reconciliation paths where applicable.
- **NFR-RL-07:** Critical defects are defined as defects that cause unauthorized data disclosure, complete loss of confirmed data, inability to access emergency/safety functions, incorrect payment settlement, or inability for major roles to perform their primary workflow.
- **NFR-RL-08:** No release candidate may contain known unresolved critical defects; significant defects must have approved workaround or deferral records before pilot release.

5.2.3 Performance

- **NFR-PF-01:** Under normal pilot load, 95% of authenticated API requests for standard create, read, update and delete operations must complete within 2 seconds, excluding third-party service latency.
- **NFR-PF-02:** Login, OTP verification and session refresh must complete within 3 seconds for 95% of successful requests under normal network conditions.
- **NFR-PF-03:** List screens such as community feed, notifications, consultation list, expert directory, moderation queue and admin user list must render the first page of results within 3 seconds for 95% of requests.
- **NFR-PF-04:** AI symptom-intake structuring and RAG answer generation should return an initial result within 10 seconds for 95% of requests; if processing takes longer, the UI must show progress or retry guidance.
- **NFR-PF-05:** Emergency map, nearby-care search, nearby-expert search and route/ETA calculation should return usable results within 5 seconds for 95% of requests when the map service is available.
- **NFR-PF-06:** Push notification delivery must be attempted within 30 seconds after the triggering backend event for reminders, consultation updates, moderation results and family alerts.
- **NFR-PF-07:** Realtime consultation connection setup for chat, voice or video should complete within 10 seconds for 95% of attempts when the realtime service and user network are available.

- **NFR-PF-08:** Health device sync should process newly available device data within 5 minutes of a scheduled sync cycle or within 30 seconds for a user-triggered manual sync, excluding external device-platform delay.
- **NFR-PF-09:** The pilot deployment must support at least 1,000 registered users, 200 daily active users, 50 verified experts, 20 partner representatives and 20 admin/moderation/content staff without architectural redesign.

5.2.4 Security and Privacy

- **NFR-SP-01:** The system must enforce authentication, role-based authorization and consent checks before allowing access to personal, health, family, consultation, payment, location or moderation data.
- **NFR-SP-02:** Sensitive data must be protected in transit with TLS and protected at rest using database, storage or platform-level encryption appropriate to the deployment environment.
- **NFR-SP-03:** Passwords must never be stored in plain text and must be protected using secure hashing and salting mechanisms.
- **NFR-SP-04:** OTP, access token and refresh session mechanisms must enforce expiration, retry limits, session revocation and suspicious-activity logging.
- **NFR-SP-05:** Users must be able to grant, view, limit and revoke data-sharing permissions for family members and experts where applicable; revocation must block future access except for legally or operationally retained records.
- **NFR-SP-06:** Location, camera, motion-sensor and device-health data must be collected only after explicit permission and only for the declared feature purpose.
- **NFR-SP-07:** Audit logs must record security-relevant events, admin actions, consent events, expert access, payment/refund actions, safety events and incident handling without exposing unnecessary sensitive payloads.
- **NFR-SP-08:** Account deletion and deactivation must follow retention notices for audit, financial, medical-record, investigation or legal-safety constraints defined by the appendix and business rules.

5.2.5 Safety and Compliance

- **NFR-SF-01:** CareBridge must consistently present itself as a healthcare support platform, not a hospital, diagnosis system, prescription system or emergency medical service.
- **NFR-SF-02:** AI-generated and RAG-generated outputs must include safety boundaries, red-flag escalation guidance and instructions to seek professional or emergency care when risk indicators are detected.
- **NFR-SF-03:** Community answers and expert answers must be distinguishable, and user-generated community content must not be presented as verified medical advice.
- **NFR-SF-04:** Pregnancy exercise, posture analysis, growth tracking, vaccination reminders, AI triage and device trends must be framed as supportive guidance and tracking, not treatment decisions.

- **NFR-SF-05:** Emergency, fall-detection and nearby-support workflows must prioritize fast access to call, navigation, family alert and emergency instruction actions over non-critical content.
- **NFR-SF-06:** Moderation and abuse-prevention workflows must allow staff to hide unsafe content, resolve reports, warn or suspend accounts and escalate safety cases according to role permission.

5.2.6 Maintainability and Supportability

- **NFR-MT-01:** Functional modules must be organized around the documented domains: account/access, mother and baby care, community, expert consultation, AI/RAG, reminders, family sync, partner governance, safety monitoring, content, moderation, payment and audit.
- **NFR-MT-02:** External integrations must be isolated behind service adapters so that mock, sandbox or replacement services can be used during academic testing or third-party outage.
- **NFR-MT-03:** Configuration for AI safety rules, consultation price bands, posture analysis, content categories, notification templates and moderation rules must be manageable without code changes where the functional specification defines admin screens.
- **NFR-MT-04:** Logs and monitoring must provide enough information to trace a failed request across client, backend and third-party service boundaries using correlation identifiers.
- **NFR-MT-05:** The system must maintain traceability between major features, use cases, screens, business rules and test cases so requirement changes can be reviewed safely.

5.2.7 Scalability and Capacity

- **NFR-SC-01:** The system architecture must allow horizontal scaling of backend APIs, notification workers, AI/RAG processing workers and realtime-session coordination components when user volume increases.
- **NFR-SC-02:** Data storage must support growth in health records, baby logs, community content, notifications, audit logs, uploaded files, consultation summaries and device-sync data without requiring changes to user workflows.
- **NFR-SC-03:** Search, filter and queue screens must use pagination or equivalent lazy loading to avoid loading unbounded records into the client.
- **NFR-SC-04:** The MVP must support Hanoi-first pilot operation and leave room for later national scaling without depending on nationwide healthcare database integration in the pilot scope.

5.2.8 Compatibility and Portability

- **NFR-CP-01:** Mobile functions must support commonly used modern Android and iOS devices capable of running the selected development framework, camera access, push notification and standard permission prompts.
- **NFR-CP-02:** Web portals must support current versions of Chrome, Edge and Safari on desktop or laptop screen sizes used by experts, administrators, moderators, content admins and partner representatives.

- **NFR-CP-03:** The UI must remain usable on common mobile screen sizes and must not require landscape orientation for normal mother, family or expert mobile workflows except where camera/posture workflow explicitly benefits from it.
- **NFR-CP-04:** If a third-party service is unavailable in a local development or demo environment, the system must support sandbox, mock or manual fallback behavior sufficient to demonstrate the related workflow safely.

5.2.9 Data Quality, Integrity and Auditability

- **NFR-DQ-01:** The system must validate required fields, data type, range and ownership before saving account, health, baby, vaccination, growth, payment, consultation, community, partner or safety records.
- **NFR-DQ-02:** Health trends, growth charts, vaccination schedules, device data and AI/RAG context must identify data sources where relevant, such as user-entered, expert-created, device-synced, imported or system-generated.
- **NFR-DQ-03:** The system must prevent duplicate processing of payment callbacks, notification jobs, device-sync events and emergency/safety events where duplicate processing could mislead users or corrupt records.
- **NFR-DQ-04:** Important records must keep created time, updated time, actor identity and status history where needed for audit, dispute, moderation, expert verification, consent or safety investigation.
- **NFR-DQ-05:** Reports and dashboards must use consistent definitions for active users, verified experts, contribution points, consultation revenue, commission, partner performance, safety events and impact metrics.

6. Requirement Appendix

This appendix defines the cross-module business rules, common requirements, reusable application messages and additional constraints that govern the CareBridge system. These items are derived from the CareBridge use cases, context diagrams and main workflow diagrams. They are not independent use cases; they constrain how use cases, screens, services and background processes must behave.

The rules in this appendix apply unless a more specific use case or functional requirement defines a stricter behavior.

6.1 Business Rules

Business rules define mandatory conditions and constraints that must be enforced by the CareBridge application, backend services and supporting integrations. Each rule includes traceability to the relevant use case group or workflow

6.1.1 Account, Identity, Role, and Access Control

ID	Business Rule	Rule Definition
BR-01	Role-Based Access Control	Before running any protected function, the backend must verify identity, active session, account status, role and permission scope. Hiding a UI control is not sufficient authorization.

BR-02	Actor Generalization	Mother, Family Member, Verified Expert, Moderator, Content Admin, Admin and Partner Representative inherit common User functions, but each role may access only role-specific functions assigned to it.
BR-03	Protected Function Authentication	The system must block protected functions when the user is not logged in, the session has expired, or the account is blocked, suspended, deactivated or deleted.
BR-04	OTP and Identity Re-verification	The system must re-verify identity before sensitive actions such as registration verification, password reset, password change, account deletion, role-sensitive changes and payment-sensitive decisions.
BR-05	Session Lifecycle	A successful login creates an active session. Logout, account deactivation, password reset and session revocation must invalidate affected sessions or tokens immediately.
BR-06	Account State Separation	Account deactivation and deletion are different states. Deactivation blocks normal use while preserving the account; deletion starts controlled removal subject to verification, retention, audit and dispute rules.
BR-07	Administrative Role Management	Only authorized Admin users may create staff accounts, assign roles, change permissions or restrict accounts. The system must prevent privilege escalation outside the approved role model.
BR-08	Sensitive Action Audit	The system must record audit events for password changes, session revocation, permission changes, health-data access, moderation, expert verification, payment/refund decisions, partner approval, safety alerts and security investigation.

Table 555: Account, Identity, Role, and Access Control Business Rules

6.1.2 Privacy, Consent, Health Data, and

ID	Business Rule	Rule Definition
BR-09	Data Ownership and Source Label	Each mother record, baby record, health summary, uploaded file, device metric, expert note, verified content item and partner item must retain the owner and source label where relevant.
BR-10	Consent Attributes	Every sharing permission must store owner, recipient, purpose, allowed data scope, start time, expiry time and status.
BR-11	Consent Expiry and Revocation	When permission expires or is revoked, the recipient must immediately lose future access. Required audit evidence may still be retained according to retention rules.
BR-12	Sharing History	The system must keep a reviewable history of grant, update, expiry and revoke events so the data owner can see who had access, what was shared and when.
BR-13	Privacy by Default	Health information, baby information, consultation data, exact location, online status, safety events and device metrics must be private by default.
BR-14	Family Permission Granularity	Joining a care group does not provide full access. Family Members may view or update only data, tasks, reminders or alerts explicitly shared by the Mother.
BR-15	Expert Access to Shared Data	Verified Expert status does not automatically grant access to personal health data. Experts may view only shared summaries or records within approved purpose, scope and validity period.

BR-16	File Ownership and Authorized Access	Before upload, preview, download, archive or deletion, the system must check file ownership, linked record, current consent and access scope.
BR-17	Archive and Soft Delete	The system should archive or soft-delete data instead of physical deletion when the data is tied to consultation, audit, investigation, payment, dispute or retention requirements.
BR-18	Minimum Necessary Disclosure	Family alerts, nearby-support requests, dashboards, partner reports and impact reports must disclose only the minimum necessary data for the current purpose

Table 556: Privacy, Consent, Health Data, and Uploaded Medical Documents

6.1.3 Community, Moderation, Content, and Partner Governance

ID	Business Rule	Rule Definition
BR-19	Community Content Ownership	Posts and answers from normal users must be treated as personal experience or emotional support. The system must not present them as diagnosis, treatment, prescription, dosage instruction, or emergency advice.
BR-20	Anonymous Display Boundary	Anonymous mode hides the display identity from normal community viewers but does not hide the author from authorized moderation, audit and safety processes.
BR-21	Community Safety Pre-check	Before publishing or routing community content, the system must check spam, prohibited content, unsafe keywords and red-flag health risk.
BR-22	Red-Flag Community Routing	Potential emergency, self-harm, severe symptom or unsafe-care content must be routed to AI triage or emergency guidance before normal community discussion or paid consultation promotion.
BR-23	Moderation Decision Governance	Moderators may approve, hide, reject, warn, restrict or suspend according to reason codes and audit policy. Users must receive clear status or rejection messages when applicable.
BR-24	Verified Expert Answer Identification	Answers must be labeled as Verified Expert Answers only when the author is a verified expert and the expert status is valid for the answer context. The system must not display ordinary user answers, AI responses or partner content as verified expert answers.
BR-25	Accepted Answer and Correctness Boundary	A question owner or authorized moderation process may mark an expert answer as accepted, helpful or selected. This marker means the answer was selected for the question context; it must not be presented as guaranteed medical correctness, diagnosis or prescription.
BR-26	Verified Content Versioning	Articles, FAQs, checklists and exercise content must have category, version, status, reviewer and publication state. Only approved versions may be published or indexed for RAG.
BR-27	Outdated or Unsafe Content Handling	Reported outdated, misleading or unsafe verified content must be reviewed and may be revised, unpublished and de-indexed from AI/RAG knowledge.
BR-28	Partner and Sponsored Content Labeling	Partner clinic/service listings and sponsored content must be reviewed before publication and clearly labeled as partner or sponsored content.

Table 557: Community, Moderation, Content, and Partner Governance

6.1.4 Expert Verification, Booking, Consultation, Payment, Pricing, Dispute, Refund, and Realtime Communication

ID	Business Rule	Rule Definition
BR-29	Expert Verification Before Practice	An expert must not appear as verified or provide official consultation until profile, credentials and required documents are approved by authorized admin users.
BR-30	Expert Status Enforcement	Expired, rejected, suspended, revoked or unverified experts must be hidden from booking, nearby expert display and expert-answer routing.
BR-31	Availability and Booking Consistency	Booking must use current expert availability and prevent double booking, invalid time slots and booking with ineligible experts.
BR-32	Consultation Data Sharing Consent	Before consultation, the system must check whether the user has shared the selected health summary or records with the expert and must not expose unshared data.
BR-33	Payment-Gated Booking State	A paid consultation must not be confirmed until payment status is successful or otherwise explicitly approved by a controlled rule. Failed, expired or cancelled payments must not create active consultations.
BR-34	Refund, Dispute and No-show Governance	Cancellation, no-show, dispute and refund decisions must follow time, status, evidence and payment rules and must be auditable.
BR-35	Commission Calculation	Expert revenue and platform commission must be calculated from completed eligible consultations according to the active commission configuration and refund state.
BR-36	Consultation Pricing Control	Experts may set or update prices only within admin-configured price bands. Mothers must see the current consultation price before booking or payment.
BR-37	Realtime Communication Eligibility	Chat, voice call or video call must be opened only for eligible consultation, support or approved contact contexts and must respect account, role, status and privacy rules.
BR-38	Consultation Summary Boundary	Consultation summaries must be created by authorized experts and shown only to eligible users within the consultation context. They are support notes, not diagnosis or prescription.

Table 558: Expert Verification, Booking, Consultation, Payment, Pricing, Dispute, Refund, and Realtime Communication

6.1.5 AI, Health Guidance, Exercise, Growth, and Reminder Safety

ID	Business Rule	Rule Definition
BR-39	AI Safety Boundary	AI/RAG/Triage output must be presented as support guidance only. It must not diagnose, prescribe, replace doctors, replace hospitals or delay emergency/offline care.

BR-40	Red-Flag Overrides Normal Flow	When red-flag indicators are detected, emergency guidance, nearby care and safe action must take priority over community posting, self-care tips or paid consultation promotion.
BR-41	AI Input Scope Control	AI intake may use only the symptom text, selected records, device data or RAG sources that the user has approved for that interaction.
BR-42	AI Fallback and Escalation	If AI service is unavailable, incomplete or unsafe, the system must show fallback guidance and safe next actions without fabricating medical advice.
BR-43	Exercise Safety Pre-check	Pregnancy exercise sessions must require safety disclaimer and contraindication checklist before starting. Failed safety checks must block the exercise session and advise offline care when needed.
BR-44	Posture Feedback Boundary	MediaPipe posture feedback must be treated as exercise-form assistance only, not medical assessment. Low camera quality must trigger retry, reposition or continue-without-posture options.
BR-45	Growth and Health Trend Disclaimer	Growth charts, health metrics and device trends provide tracking and visualization only. Abnormal values should trigger guidance to consult experts or offline care, not diagnosis.
BR-46	Reminder and Checklist Safety	Reminders and checklists must not override urgent care guidance. Missed reminders should be visible but must not imply medical failure or emergency status unless connected to a safety workflow.
BR-47	Health Content Source Control	Health guidance, checklist, FAQ and exercise instructions must use approved content sources and show appropriate safety wording when used in app, AI/RAG or expert workflows.

Table 559: AI, Health Guidance, Exercise, Growth, and Reminder Safety

6.1.6 Map, Location, Nearby Support, and Connected Devices

ID	Business Rule	Rule Definition
BR-48	Location Permission Requirement	Exact user or expert location may be collected or shared only after explicit location permission and purpose-specific consent.
BR-49	Manual Location Fallback	If location permission is denied or unavailable, emergency and nearby-care workflows must provide manual location entry, quick call or default safety guidance.
BR-50	Map Service Boundary	Map, route, distance and ETA results are guidance only. CareBridge must not guarantee ambulance dispatch, facility availability or medical response.
BR-51	Nearby Expert Opt-in	Experts appear in nearby-support results only if verified, eligible, currently opted in and within configured radius, duration and visibility scope.
BR-52	Nearby Request Minimum Data	Nearby support requests must share only minimum necessary request data until an expert accepts, rejects, expires or the requester cancels the request.
BR-53	No Sponsored Blocking Before Safe Action	Emergency and red-flag workflows must not place ads, sponsored content or partner promotion before safe action options such as emergency guidance, nearby care or quick call.

BR-54	Connected Device Consent and Scope	Smartwatch or wearable data import must explain data types, sync frequency, source labels, retention and limitations before connection or manual import.
BR-55	Device Data Quality Rules	Missing, stale, duplicated or conflicting device data must be labeled and must not be fabricated, silently overwritten or used as confirmed medical evidence.

Table 560: Map, Location, Nearby Support, and Connected Devices

6.1.7 Smart Activity Monitoring and Safety Support (MF21)

ID	Business Rule	Rule Definition
BR-56	Voluntary Activation	MF-21 must be off by default. It may start only after the Mother reads the limitations, gives the required permissions, completes the safety settings, and confirms consent. Sharing location remains optional.
BR-57	Sensor and Background Limitation	The system must not promise that MF-21 will run continuously on every phone. Where practical, raw motion data should be processed on the device and should not be uploaded continuously unless there is a clear need and valid consent.
BR-58	Suspected Event Wording	Motion analysis may create only a “suspected fall/impact” event. The system must not say that the user is injured, unconscious, or definitely in an emergency.
BR-59	Mandatory Safety Check	After detecting a suspected event, the system must open the safety-check screen with sound, vibration, and a countdown. The user must be able to select “I am safe”, “Need help”, or “Emergency support”.
BR-60	Family Alert Conditions	The system may send a family alert only when the Mother selects “Need help” or does not respond before the countdown ends. The alert may go only to configured contacts and may include location only when permission exists.
BR-61	Safety History and False Positives	The system must keep the minimum event details and the user's response in safety history. The user may mark an event as a false positive, but the original sensor-event information and any alert already sent must not be rewritten.
BR-62	MF-21 to MF-19 Escalation	MF-21 may open the MF-19 emergency-support flow so the user can call, search for nearby care, or start navigation. This action does not dispatch an ambulance or guarantee that help will arrive.

Table 561: Smart Activity Monitoring and Safety Support

6.2 Common Requirements

This section specifies common requirements shared across multiple modules. They support consistent behavior for authentication, authorization, validation, privacy, external-service interaction, safety handling, auditing and traceability.

ID	Common Request	Rule Definition
CR-01	Authentication Check	Every protected UI action, API endpoint and background action must verify the current session and account state before execution.

CR-02	Authorization Check	Authorization must be enforced server-side using role, ownership, consent and feature permission checks.
CR-03	Input Validation	All user input, upload metadata, payment status, location values and external-service responses must be validated before storage or use.
CR-04	Error Handling	The system must show safe, user-understandable errors without exposing internal stack traces, tokens or sensitive data.
CR-05	Audit Event Creation	Sensitive actions must create append-only audit events with actor, action, target, timestamp, result and reason where applicable.
CR-06	Consent Gate	Any access to shared health data, location, device metrics, consultation summary or safety alert must pass active consent and scope checks.
CR-07	Notification Preference Handling	Noncritical notifications must respect user preferences. Safety-critical notifications may follow special mandatory-delivery rules.
CR-08	External Service Fallback	If payment, AI, map, notification, realtime or storage services fail, the system must show a safe fallback state and preserve data consistency.
CR-09	Traceability	Each rule, message and common requirement should be traceable to at least one UC group or workflow section.
CR-10	Data Minimization	Only collect, display and transmit the minimum data needed for the active use case and purpose.
CR-11	Source and Timestamp Labeling	Records imported from device, uploaded by user, generated by expert or created by system must include source and timestamp when relevant.
CR-12	Status Consistency	State changes such as pending, approved, rejected, cancelled, expired, completed, refunded, suspended and revoked must be consistent across UI and backend.
CR-13	Idempotency for Sensitive Actions	Payment confirmation, refund, notification retry, device sync retry and alert delivery retry must avoid duplicate side effects.
CR-14	Retention and Deletion	Deletion, archive, anonymization and retention must preserve required audit, dispute and safety evidence while honoring privacy requirements.
CR-15	Localization	User-facing messages should support approved Vietnamese wording while preserving message codes and safety meaning.
CR-16	Accessibility and Readability	Critical safety, consent, payment and error messages must be visible, readable and not hidden behind decorative content.
CR-17	Versioned Content	Verified content, exercise configuration, AI/RAG knowledge and safety wording must be version-controlled and auditable.
CR-18	Role-specific Navigation	Users must be routed to the correct dashboard and feature set after login based on role and onboarding state.
CR-19	Sensitive Data Masking	OTP targets, contact details, health data previews and dashboard indicators must be masked or summarized when full display is not authorized.
CR-20	Safe Default State	When the system cannot verify permission, service response, expert status or safety state, it must default to the safer restrictive behavior.
CR-21	No Fabricated Results	The system must not fabricate payment confirmation, AI output, device metrics, location ETA or notification delivery status when an external service fails.

CR-22	Report Consistency	Section 6 business rules, messages and common requirements must remain consistent with context diagrams, workflows, use case list and functional requirements.
-------	--------------------	--

Table 562: Common Requirements

6.3 Application Messages List

This section defines reusable, safety-critical, and cross-module messages. The UI may provide approved Vietnamese translations, but the message code, severity and safety meaning must remain unchanged. Function-specific success messages may be added in screen specifications.

Message type and display color convention used in this section:

1. Inline message (neutral/dark gray or red when blocking): text displayed inside the screen area near the related content or action.
2. Inline red under field: validation text displayed in red under or near an input field; the related input may also use a red border until corrected.
3. Toast message (green/blue/red): temporary feedback displayed after an action; green indicates success, blue indicates neutral information and red indicates failure.
4. Confirmation dialog (neutral or red destructive action): a dialog that asks the user to confirm a sensitive or irreversible action before continuing.
5. Modal dialog (neutral/amber/red): a blocking dialog used when the user must read or resolve the message before continuing; amber/red may be used for warning or safety-critical cases.
6. Banner / safety notice (amber or red): a persistent and visible notice; amber is used for caution/limitation and red is used for urgent, emergency or safety-critical action.
7. Badge / label (blue/green/amber/gray): a short visual label; blue or green identifies verified content/expert answers, amber identifies accepted or selected answers, and gray identifies partner or sponsored content.

6.3.1 Common Validation and Access Messages

Message Code	Message Type	Context	Content
MSG-COM-01	Toast message (red)	Authentication	Your session has expired. Please log in again to continue.
MSG-COM-02	Inline message (red)	Authorization	You do not have permission to perform this action or view this data.
MSG-COM-03	Modal dialog (neutral)	Sensitive action	Please verify your identity before continuing this sensitive action.
MSG-COM-04	Inline message (red)	Account state	This account is currently restricted, deactivated, suspended or deleted.
MSG-COM-05	Toast message (green)	Logout/session	The selected session has been signed out successfully.

Table 563: Common Validation and Access Messages

6.3.2 Privacy and Consent Messages

Message Code	Message Type	Context	Content
MSG-PRI-01	Inline red under field	Consent	This data cannot be shared until scope, purpose and expiry are selected.
MSG-PRI-02	Toast message (green)	Consent revoked	Sharing permission has been revoked. Future access is no longer allowed.
MSG-PRI-03	Inline message (red)	Shared data	The selected recipient no longer has access to this record.
MSG-PRI-04	Inline message (red)	Uploaded medical document access	This uploaded medical document is private and can be opened only by authorized users.

Table 564: Privacy and Consent Messages

6.3.3 Community, Content and Partner Governance Messages

Message Code	Message Type	Context	Content
MSG-COMM-01	Banner / safety notice (red)	Community red flag	This content may describe an urgent health concern. Please review emergency guidance before posting.
MSG-COMM-03	Inline message (red)	Moderation rejected	This content cannot be published because it violates CareBridge community or safety rules.
MSG-COMM-04	Badge / label (blue/green)	Verified expert answer	This answer was written by a verified expert. It is support information and does not replace professional medical care.
MSG-COMM-02	Toast message (blue)	Moderation pending	Your content has been submitted and is waiting for moderation review.
MSG-COMM-05	Badge / label (amber)	Accepted answer	This answer has been marked as accepted or helpful for this question. It must not be treated as diagnosis or prescription.
MSG-CONT-01	Inline message (red)	Verified content	This content version is not approved and cannot be published or indexed for AI/RAG.
MSG-PART-01	Badge / label (gray)	Sponsored content	This is sponsored or partner-provided content and is labeled for transparency.

Table 565: Community, Content and Partner Governance Messages

6.3.4 Expert, Consultation, Payment and Realtime Messages

Message Code	Message Type	Context	Content
MSG-EXP-01	Inline message (red)	Expert verification	This expert profile is not verified or is currently restricted. Booking is unavailable.
MSG-EXP-02	Banner / safety notice (amber)	Credential expiry	Expert verification requires renewal before this profile can be shown as verified.
MSG-CON-01	Inline message (red)	Booking	The selected time slot is no longer available. Please choose another slot.
MSG-CON-02	Toast message (red)	Payment	Payment was not completed. The consultation has not been confirmed.
MSG-CON-03	Toast message (blue)	Refund	Your refund or dispute request has been submitted for review.
MSG-CON-04	Banner / safety notice (amber)	Shared summary	No health summary is shared with this expert. You may continue without sharing or grant access.
MSG-CON-05	Toast message (red)	Realtime	The call or chat cannot start because the consultation is not eligible or the service is unavailable.
MSG-PRICE-01	Inline red under field	Pricing	The expert price must stay within the admin-configured price band.

Table 566: Expert, Consultation, Payment and Realtime Messages

6.3.5 Community, Content and Partner Governance Messages

Message Code	Message Type	Context	Content
MSG-EXP-01	Inline message (red)	Expert verification	This expert profile is not verified or is currently restricted. Booking is unavailable.
MSG-EXP-02	Banner / safety notice (amber)	Credential expiry	Expert verification requires renewal before this profile can be shown as verified.
MSG-CON-01	Inline message (red)	Booking	The selected time slot is no longer available. Please choose another slot.
MSG-CON-02	Toast message (red)	Payment	Payment was not completed. The consultation has not been confirmed.
MSG-CON-03	Toast message (blue)	Refund	Your refund or dispute request has been submitted for review.
MSG-CON-04	Banner / safety notice (amber)	Shared summary	No health summary is shared with this expert. You may continue without sharing or grant access.

MSG-CON-05	Toast message (red)	Realtime	The call or chat cannot start because the consultation is not eligible or the service is unavailable.
MSG-PRICE-01	Inline red under field	Pricing	The expert price must stay within the admin-configured price band.

Table 567: Expert, Consultation, Payment and Realtime Messages

6.4 Other Requirements

Other cross-cutting constraints capture rules that must be stated in the SRS but should not be modeled as standalone use cases.

ID	Requirement	Descrip
OR-01	Safety Boundary	CareBridge is a maternal and early childhood health support platform, not a hospital information system, medical diagnosis tool, prescription system or emergency dispatch service.
OR-02	No Diagnosis / No Prescription	No AI, expert summary, community answer, device trend, posture feedback or growth chart may be worded as official diagnosis or prescription.
OR-03	Audit Immutability	Audit records must be append-only for normal operators. Corrections must be made by additional records, not by editing or deleting existing evidence.
OR-04	Raw IMU Retention Limit	Continuous raw phone IMU streams must not be stored long term. Only suspected event metadata and safety confirmation results should be retained as required.
OR-05	External Service Dependency	Payment, map, AI, notification, realtime, storage and health-device services are external dependencies and require fallback, retry and reconciliation handling.
OR-06	Partner Transparency	Sponsored or partner-provided items must be labeled and must not be prioritized above safety actions in emergency or red-flag flows.

OR-07	Data Protection	Encryption, access control, consent checks, source labels, retention rules and audit logging are security/data requirements, not separate end-user use cases.
OR-08	Localization	Approved Vietnamese UI wording may be used, but message codes, severity and safety meaning must remain consistent.
OR-09	Consistency with Diagrams	Business rules must stay consistent with the Web/Mobile context diagrams and the 14 main workflow diagrams.
OR-10	Consistency with Use Cases	Rules that UC guidance marks as cross-cutting must remain in Section 6 rather than being drawn as independent use cases.

Table 568: Other Requirements